



TELEPHONE: + 34 93 699 50 00



FAX:

+ 34 93 699 45 12



E-MAIL:

General: info@deltalab.es

Marketing: marketing@deltalab.es



MAIL:

Plaza de la Verneda, 1 (Pol. Ind. La Llana) P.O. BOX 195 08191 RUBÍ (BARCELONA) SPAIN



WEB:

www.deltalab.es





Dear customers, distributors, suppliers and friends,

We present the new edition of the DELTALAB General Catalogue, the result of two years of an intense work, which contains a wide range of products consisting of thoroughly organized items by specialty.

The catalogue that you have in your hands is the result in the shape of products, of DELTALAB's vocation of service and adaptation to the needs that are generated constantly in our industry. In response to market demands, we have incorporated new interesting items, always meeting the highest quality standards.

And in our desire for overall continuous development, I would also remark another important aspect. We redesigned the catalog image making it more minimalist and conceptual. We have also added new images that give a great aesthetic value. Also we grouped and rearranged some of the chapters in order to improve usability and content balance.

DELTALAB is a leading company in Spain in the manufacturing of single use labware and it is established as a reference company in the sector at an international level. It is certified by business volume data: over 60% of the turnover came from the sales in export markets.

Our constant ambition for self-improvement, our efforts and our investment policy in R&D allow us to improve year by year our production processes, and work intensively on product development. This is the DELTALAB commitment, working on innovation as value added of the company in an increasingly competitive market. And so, it is reflected in our slogan: "your needs, our innovations".

José Sáez Mateos Managing Director





Iso Quality Standard



UNE-EN ISO 9001 Quality management systems. Requirements.



UNE-EN ISO 13485 Medical devices. Quality management systems. Requirements for regulatory purposes.



CE mark in medical devices (MDD) and medical devices for in vitro diagnostics (IVD)

The CE Mark concerns those products that, affecting human beings, may involve a risk to their health. For that reason they must comply with the following requirements:

- Do not compromise health or clinical status.
- Do not compromise user's health or safety.
- Do not compromise the health or safety of third persons.
- · Minimize the risk of errors of use thanks to its ergonomic features and the environment in which it is intended to be used.
- It takes into account the technical knowledge, expertise, experience, training, and medical and physical conditions of the intended users.

The CE Mark implies installations, and work organisation, and systems are submitted to rigorous both internal and external quality controls. It guarantees the highest level of safety when using those products intended for the sample collection and handling of human samples, both for the patient and for the user.

CE Mark guarantees that those products will be freely commercialised within the European Union. Since December the 7th 2003, the only sanitary products (MDD and IVD) that can be commercialised are those ones with CE Mark.

Medical devices clasification

1. MEDICAL DEVICES (MDD)

According to Directive 93/42/EEC (modified by Directive 2007/47/EC), a MDD Medical Device is "any instrument, device, equipment, material or any other article, including software, developed to be used with human beings, alone or in association with others, with the aim of: diagnosis, prevention, tracing, treatment or relief of a disease; diagnosis, tracing, treatment, relief or compensation in case of an injury or disability; investigation, substitution or modification of anatomy or physiological processes; fertility control; and not designed to perform the main action expected to be obtained by immunizing, pharmaceutical or metabolical means inside or over the human body".

This Directive establishes 4 categories of medical devices; DELTALAB manufactures and markets temporary use product, so: Class I (sterile and non sterile) and Class IIa:

• Class I (Non sterile): Low risk, invasive medical devices

It includes: non sterile tongue depressors, speculums... Inspection gloves are considered both as MDD Class I and PPE (Personal Protective Equipment). All of them are self certification products.

Class I (sterile): Low risk, sterile, invasive medical devices

Among others: sterile tongue depressors, sterile Ayre spatulas, sterile cytology brushes... They are certified by a Notified Body.

• Class IIa (sterile): Moderated risk medical devices

Swabs with and without medium are included in this Class. They are certified by a Notified Body.

CE Mark allows swabs to be used in surgical operations as an invasive product.

2. MEDICAL DEVICES FOR IN VITRO DIAGNOSTIC (IVD)

According to Directive 98/79/EC, a medical device for in vitro diagnostic is "any medical device consisting on a reagent, reactive product, calibrator, monitoring equipment, instrumentation and equipment cases, sample containers, instrument, apparatus or system designed by the manufacturer to be used in vitro to examine human samples, including blood and tissues donations; with the aim of collecting information: about a physiologic or pathologic status, or a congenital anomaly; or to determine the security and compatibility with potential receptors; or in order to supervise therapeutical measures".

They are considered low risk devices as they are not in touch with human beings and are self certification products.

As an example, sample containers, collection tubes, capillary tubes, etc. are included.

PPE PRODUCTS (Personal Protective Equipment)

According to Directive 89/686/EEC, a PPE product is "any device or system a human being is wearing or using, which protects him from one or various risks that may endanger his health and/or security". Security glasses and protection gloves are included. Inspection gloves are considered both as MDD and PPE as they protect the user while being a low risk, invasive medical device.



Sales conditions

- Goods are shipped at the purchaser's risk and on the purchaser's account.
- Goods are considered delivered, and at the responsibility of the purchaser, from their shipment from our warehouses. Any problem of transport or other is therefore the entire responsibility of the purchaser from that moment.
- MINIMUM ORDERING QUANTITY (per total order) according to our price list.
- MINIMUM ORDER QUANTITY (per product): In those products indicated, it is possible to acquire minor quantities from the standard selling unit. In those cases, a surcharge is applied in order to cover the manipulation costs. Please see specific conditions in the Sales Conditions attached to the Price List or contact our Sales Department. Orders over the standard selling unit must be multiple of this standard selling unit.
- Orders must be in complete case quantities unless indicated otherwise.
- Goods delivered are the property of DELTALAB S.L. until paid in full.
- ORDERS should be confirmed in writing, stating our code number, description, total quantity and referring to any quotation number that may apply. In case of doubt the code number will have preference over its description.
- For any **CLAIM** please quote our delivery note reference, batch number, product code, total quantity and any other relevant information. Any claim should be made within 60 days of receipt of goods.
- **RETURN OF GOODS**: We only accept the return of goods due to a problem caused by DELTALAB S.L. *A prior authorisation* must be obtained from our commercial department, in which we will indicate which freight forwarder should be used.
- Weights and measures in this catalogue are in Kgs and cubic meters.

For reasons which are beyond our control we reserve the right to modify the product specifications.

- Please contact us for any SPECIAL PRODUCT. In addition to the products listed we can offer a wide range of colour caps, different concentrations of anticoagulants, different configurations of our stainless steel racks and private labelled products.
- **SALES CONDITIONS:** You can see our complete sales conditions at DELTALABs Price list.
- Please contact our commercial department for more information.

DELTALAB S.L.

All rights reserved. No part of this catalogue may be reproduced without the prior permission in writing of DELTALAB S.L.



Nomenclator

Stackable	
	- <u>1</u>] =
Graduated	-2 -3
Alphanumerical Identification	A-Z O-1
Aprilation to a footblood of	
Floats	
	DNAse RNAse
DNAse and RNAse free	FREE
	S¶″
Autoclavable (see more information about autoclaving in page 11)	[2]
	X
Suitable for PCR	PCR
	X
Suitable for QPCR (Real Time PCR)	QPCR
	₹ ∃
Innovative product	INNOVATIVE
	NEW
New	****
	Ç€(IVD)
Directive 98/79/EC (in Vitro Diagnostic medical device)	
Directive 93/42/EEC (low and moderate risk medical device)	C€(MDD)
Directive 93/42/EEC (Moderate risk, sterile medical device)	(€ class IIa
	(€ (EPI)
Directive 89/686/EEC (Personal protective equipment)	(C (EPI)

The technical data (dimensions, capacities and resistances) specified in this catalogue are approximated.

The resistances (to chemical products, temperature, centrifugation...) may suffer variations depending on external conditions. We recommend to test the product under the forecasted conditions of use. Please do not hesitate to ask for samples if necessary.







_0	Technical information	8	i
_1	Microbiology	14	<u>@</u>
_2	Molecular biology	48	٠٠.
3	Haematology and biochemistry	66	
_4	Histology, microscopy and staining	84	٤
_5	Sampling and transport	108	
_6	Tubes, microtubes and cryovials	118	
	Sample containers	146	
8	Liquid handling	166	•
9	Sample storage	198	
10	Wilmut system	230	W
11	Hygiene, safety and general labware	238	
12	Laboratory and industrial packaging	278	
13	Index	298	A-Z

	PE	LD	PE	HD	Р	P	Р	s	TI	PX	Al	3S	PM	IMA	Р	С	P۱	/C	PT	FE
°C	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50
1.4 - dioxane							×	×					×	×			×	×	•	•
Acetaldehyde	•	×	•		•	×	×	×	•	×	×	×	×	×		×		×	•	
Acetic acid	•	•	•			•			•		×	×	×	×	•		•		•	
Acetone		×	•			•	×	×					×	×	×	×	×	×	•	
Acrylonitrile	•	•	•			×	×	×	×	×			×	×	×	×	×	×	•	•
Adipic acid	•	•	•	•	•	•	•	•	•	•			•	•	•	•	•			•
Aluminium chloride			•						•				•		×	×	•			•
Aluminium hydroxide									•							×		•		
Allyl alcohol	•	•	•			•	•		•				×	×	•			×		•
Amino acids	•	•	•	•	•	•	•	•	•						•	•	•	•	•	•
Ammoniac	•		•			•		×	•						×	×				
Ammonic hydroxide (30 %)		•	•			•		×	•		•				×	×	•		•	•
Ammonium chloride		•					•		•			•					•		•	
Amyl acetate		×	•			×	×	×					•		×	×	×	×	•	•
Amyl alcohol	•	•		-	•	•				-						•			•	
Amyl chloride	*	*	×	×	×	×	×	*	×	×			×	×	×	*	×	*	•	
Aniline	•	•		•	•	•	*	*	•		×	×	*	*	ï	*	*	*		
Aqua regia	*	×	×	×		×	-	*			×	×	×	×	×	×		-		
Benzaldehyde	•	•	•	-		-	×	×			×	×	×	×		×	×	×		
Benzene		*					×	*			×	×	×	×	×	×	*	*		
Benzine		*					×	×	•		-		•	•	•	×	•	•		
Benzoid acid	*	*	×	×	*	*	*	×	×	×			×	×		•				
Boric acid (10 %)	•	•	•	-	•	•	•	•	•	•		•	•	•				-		
Bromine	×	*	×	×	×	×	×	×	×	×			×	×	×	×	×	×		
Bromoform	*	*	×	×	*	*	×	×	×	×			×	×	×		×	×		
									•				×	-		*	×	×	-	
Butyl acetate			•				×	*			*	*		*	×	*			•	
Butyl alcohol	•	•	•	•	•	•		*	•		×	*		×				•		
Calcium chloride	•	•	•			•		•		<u> </u>			•			•		×	•	
Calcium hydroxide	•	•			•	•	•			•			•	-	*	*	-	•	•	
Calcium hypoclorit	•	•	•	•	•	•	•	•	•				•	•		*	•	*	•	
Carbon sulphate	*	*	*	*	*	*	*	*	*	*			×	*	*	*	*	*	•	
Carbon tetrachloride	*	×		*	*	×	*	*	×	*				*	×	*	*	*	•	
Citric acid			•						•		•	•		-	•					
Cupric sulphate	-	•	-	•	-	•	•	•	_	•			-	•	-	•	-			
Chlorine (10 %)	-	*		*	-	*	*	*	_	*	_			*				*	•	
Chlorine water	*	*	×	*	×	*	×	×		×	•	•	*	×	*	*	*	×	•	-
Chloroform	×	×	•		×	×	×	×		×	×	×	×	×	×	×	×	×		
Chlorydic acid (35 %)	•	-	•	•	•	•			•	•				*	×	×		*	•	•
Chromic acid (10 %)	•	•	•	•	•	•			•	•				×	•		•		•	
Chromic acid (50 %)	•		•				×	×					×	×		×	•	×	•	
Chromic sulphate blend	•	×	•	×	×	×				×			×	×	×	×	•		•	•
Decahydronaphtalene					×	×	×	×		×			×	×	×	×	•		•	•
Dichlorobenzene		×		×		×	×	×	×	×	×	×	×	×	×	×	×	×	•	•
Diethylene glycol	•	•	•	•	•	•		×	•	•			×	×	•		×	×	•	•
Dimethyl formamide	•	•	•	•	•	•	×	×	•	•			×	×	×	×		×	•	•
Dimethyl sulphoxide	•	•	•	•	•	•	×	×	•	•			×	×	×	×	×	×	•	•
Ether	×	×		×		×	×	×	×	×	×	×	×	×	×	×	×	×	•	•
Ethyl acetate	•	•	•	•	•	•	×	×	•		×	×			×	×	×	×	•	•
Ethyl alcohol (100%)	•	•	•	•	•	•	×	×	•				×	×	•		•		•	•
Ethylene chloride	×	×	×	×		*	×	×	×	×	×	*	×	×	×	×	×	×	•	•
Ethylene oxide						×	×	×		×			×	×		×		×	•	•
Fluorine	×	*	×	×	×	*	×	×		×			×	×	•		•	•	•	•
Formaldehyde (40 %)	•	•	•	•	•	•	×	×	•	•		×	×	×	•			×	•	•
Formic acid (98-100 %)	•	•	•	•	•	•	•		•	•			×	×	•		×	×	•	•
Fuel oil		×	•		•	•	×	×						×	•		×	×	•	•
Glacial acetic acid	•		•	•	•		×	×	•		×	×			×	×	•		•	•
Glycerine	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•
Glycol	•	•	•	•	•	•	•	•	•	•			•	•	•	•	•	•	•	•
				*																



	PE	LD	PE	HD	F	P	P	s	TE	PΧ	A	BS	РМ	MA	P	C	P	VC	PT	FE
°C	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50	20	50
Hexane	×	×	•					×		×	×	×	•	•	×	×		×	•	•
Hydrofluoric acid (40 %)	•	•	•	•	•	•	•	•	•	•			×	×	×	×		×	•	•
Hydrofluoric acid (70 %)		×	•				×	×	•				×	×	×	×	×	×	•	
Hydrogen peroxide (35%)	•	•	•			•	•	•	•	•	×	×	×	×	•	•	•		•	
lodine	×	×	×	×	•	•		×	•				×	×		×	×	×	•	•
Isobutyl alcohol	•	•	•	•	•	•			•	•				×	•	•	•		•	•
Isopropyl alcohol	•	•	•	•	•	•			•	•				×	•	•	•		•	•
Isopropyl benzene		×	•			×	×	×	×	×			×	×	×	×	×	×	•	•
Lactic acid	•	•	•	•	•	•	•	•	•	•	•	•		×	•	•			•	•
Mercurous chloride	•	•	•	•	•	•	•		•	•			•	•	•	•	×	×	•	•
Mercury	•	•	•	•	•	•	•	•	•	•			•	•	•	•	•	•	•	•
Methyl acetate		×			•		×	×	•	•	×	×			×	×	*	×	•	•
Methyl alcohol	•		•	•	•	•		×	•	•		×	×	×	•		•		•	•
Methyl propyl ketone	•		•	•	•		×	×					×	×	×	×	*	×	•	•
Methylene chloride	×	×		×	×	×	×	×	*	×	*	×	×	×	×	×	*	×	•	•
Mineral oil	•		•	•	•	•	•	•	•	•			•	•	•	•	•	•	•	•
Monochloroacetic acid	•	•	•	•	•	•	×	×	•	•				×		×	•	•	•	•
Nitric acid (10 %)	•	•	•	•	•	•	×	×	•	•	•	•	•		•		•		•	•
Nitric acid (50 %)				×		×	×	×		×					•			×	•	•
Nitric acid (70 %)	×	×	×	×	×	×	×	×		×	×	×		×	×	×	×	×	•	•
Nitrobenxene	×	×		×	×	×	×	×	×	×	×	×	×	×	×	×	×	×	•	•
Oxalic acid	•	•	•	•	•	•	•	•	•	•	•		•	•	•	•	•	•	•	•
Ozone	•	•	•	•	•	•			•	•			•	•	•	•	•		•	•
Perchloric acid	•	×	•	×	•	×	×	×		×	×	×	×	×	×	×		×	•	•
Perchloroethylene	×	×	×	×	×	×	×	×	×	×				×	×	×	×	×	•	•
Phenol (100 %)	•		•	•	•	•	×	×					×	×	×	×	×	×	•	•
Phosphoric acid (85 %)	•	•	•	•	•	•	•		•	•	•	•	×	×	•	•	•		•	•
Phtalat dibutylic		×		×	•		×	×	•				×	×	×	×	*	×	•	•
Potassium chloride	•	•	•	•	•	•			•	•			•	•	•	•	•		•	•
Potassium hydroxide	•	•	•	•	•	•			•	•	•	•	•	•	×	×			•	•
Potassium permanganate	•	•	•	•	•	•	•	•	•	•			•	•	•	•	•	•	•	•
Propylene glycol	•	•	•	•	•	•	•	•	•	•	•	•			•		×	×	•	•
Propylene oxide	•	•	•	•	•	•	×	×	•	•					×	×	×	×	•	•
Pyridine	-		•				*	×	-		*	×	×	×	×	×		×	-	-
Salicylic acid	•	•	•	•	•	•	•	•	•	•								*	•	•
Salicylic aldehyde	•	•	•	•	•	•	*	*	•	•							×	×	•	•
Silver acetate	•	•	•	•	•	•			•	•					•	•			•	•
Silver nitrate	•	•	•	•	•	•			•	•			•	•	•	•			•	•
Sodium acetate	•	•	•	•	•	•	•	•	•	•			*	×	•	•			•	•
Sodium dichromate	•	•	•	•	•	•	•	•	•	•				•			•	•	•	•
Sodium hydroxide	•		•	•	•	•	•	•	•	•		•		•	×	*	•	•	•	•
Sulfuric acid (60 %)	-	•	•	•	•	•	*	*	•	•	44		*	*	•	•	•	*	•	•
Sulfuric acid (98 %) Tartaric acid		×		×	*	×	×	×	•	•	×	×	×	×	×	×	*	×	•	•
	•	•	•	•	•	•	•	•	•	•			•	•	•	•	•	•	•	•
Tetrahydrofuran Toluene	*	*		×	*	×	×	*		*	×	*	*	×	*	×	×	×	•	•
Trichloroethane	•	×		•	*	×	×	×	*	×	×	×	×	×	×	×	*	×	•	•
Trichloroethylene	×	×		×		×	×	×		×	×	×	×		×	×	×	×	•	
		×		*	*		*	*	*		*	*	*	×		*	*	*	•	
Triethylene glycol Tripropylene glycol	•									•								*	•	
Turpentine	-	×	_	_		-	_	-	_		•	•							•	
Urea				*	*	*	×	×			×	×	•		×	×	•	•	•	_
Vinylidene chloride	*	-		-		-	•	•	_					_	×	×	•	×	_	•
Xylol	*	×		×	×	×	×	×	×	×	•	•	×	×	×	×	×	×	•	•
Zinc chloride (10 %)		×		×	*	×	×	×		×	×	×	×	×	*	×	×	*	•	
			•	•	•	•		•	•	•	×	•	×	*	•	•			•	•
Zinc sulphate (10 %)											*	×								

- = Excellent / very good chemical resistance = Good / limited chemical resistance
- = Poor chemical resistance



	POLYSTYRENE	HIGH DENSITY POLYETHYLENE	LOW DENSITY POLYETHYLENE	POLYPROPYLENE	
GENERAL PROPERTIES	PS	PEHD	LDPE	PP STANDAR	
Structure	Amorphous structure with stable dimensions $\left[-\left(CH_{2}-CH_{n}\right)-CH_{2}-CH-\right]m$	The catalytic controlled polymerisation of ethylene minimises the production of compact, lateral structures resulting in a pliable but rigid material.	Ethylene is polymerised under super pressure, which produces a higher number of less compact lateral structures resulting in a material which is more flexible than HDPE.	Polypropylene properties are similar to polyethylene. A methyl group is bonded to each carbon pair, which greatly improves the thermal resistance. Polyethylene becomes brittle when in contact with oxidants for long periods of time.	
	$\left[- CH - CH_{2} - (CH - CH_{2}) - \right] m$	$\left[-\operatorname{CH_2-CH_2-} \right] \operatorname{n}$	$\left[\!$	[- CH(CH ₃) - CH ₂ -]	
Transparency	Optically clear	Opaque	Opaque	Translucent	
Mechanical resistance	Rigid / Brittle	Rigid	Flexible	Resists Deformation	
Density	1,05 gr/cm ²	0,955 gr/cm ²	0,924 gr/cm ²	0,91 gr/cm ²	
Max. temperature	70° C	80° C	75° C	121° C	
Min. temperature	-10° C	-50° C	-50° C	0° C*	
Autoclavable	NO	NO	NO	YES	
Gas sterilization	DOSES LIMITED	YES	YES	YES	
Gamma irradiation Sterilization	YES	YES	YES	DOSES LIMITED	
Beta irradiation Sterilization	YES	YES	YES	YES	
Water absorption	< 0,05 %	< 0,01 %	< 0,01 %	< 0,02 %	
Properties	-Bright finish		Light weighted.Slight electrostatic changeImpact resistant	-Virtually inert	

 $^{^{\}star}$ There are special blends of PP like those from our cryovials, which are able to withstand up to -196 $^{\circ}$ C.

	POLYSTYRENE	HIGH DENSITY POLYETHYLENE	LOW DENSITY POLYETHYLENE	POLYPROPYLENE	
GENERAL CHEMICAL RESISTANCE	PS	HDPE	LDPE	PP STANDARD	
Oils	GOOD	GOOD	GOOD	GOOD	
Acids	LIMITED* GOOD*		GOOD*	GOOD*	
Alcohol	GOOD	GOOD	GOOD	GOOD	
Bases	GOOD	GOOD	GOOD	GOOD	
Ketones	POOR	LIMITED	LIMITED	LIMITED	
Esters	POOR	LIMITED LIMITED		LIMITED	
Fats	GOOD	GOOD	GOOD	GOOD	
Hydrocarbons -aromatics	POOR	GOOD	LIMITED	GOOD	
Hydrocarbons Chloride	See chemical resistance table**	LIMITED	LIMITED	See chemical resistant table**	
Hydrocarbons Halogenated	POOR	LIMITED	LIMITED	LIMITED	
Metals (Cu, Mn, Co)	See chemical resistance table**	See chemical resistant table**	See chemical resistant table**	LIMITED	
Oxidants	See chemical resistance table**	LIMITED/POOR	LIMITED/POOR	LIMITED/POOR	

^{*} Poor resistance to strong oxidant acids




 $^{^{\}star\star}$ Chemical resistances of plastics are detailed in a specific table

Sterilisation methods

METHOD	PROCEDURE	RECOMMENDED FOR	PRECAUTIONS	LIMITATIONS
Dry heat	Direct action of dry heat 190 °C for 80' or 160 °C (130-160').	Glassware, metals and liquids.	High temperatures may damage fragile metals.	The own material's limitations. High heat soaking may unacceptabbly affect material properties.
Autoclave,Super Heated Steam under pressure	Action of 3 elements: temperature, water steam and pressure 121 °C (15'), 126 °C (10')	Glassware, fabrics, liquids etc. All materials resistant to heat above 121 °C and moisture.	Not recommended for most ordinary plastics. When autoclaving capped containers, total closure shall be avoided	Generally used for the sterilisation of small articles.
Ethylene Oxide ("EO Gas")	Used alone, associated with freon or carbon dioxyde. 55 to 60 °C (2-3h), 27 °C to 33 °C (5h30)	Generally all materials, some exceptions.	Requires a subsequent ventilation to ensure sterilised products are free of residual gases that may be toxic.	Ethylene oxide is toxic and explosive.
Gamma irradiation (Radioactive cobalt source)	By radiations emitted by a radioactive source	Widely used in industries to sterilise single use material. Sterilisation doses are calculated from bioburden counts.	Limitations in some applications, as some material properties may be unacceptably altered by this method.	Effects on material are cumulative, therefore gamma sterilised products are not resterilisable, by most conventional methods, after first use.
Beta irradiation (non radioactive particle accelerator)	High energy electron beam	Widely used in industries to sterilise single use material. Sterilisation doses are calculated from bioburden counts.	Limitation in some applications, as some material properties may be unacceptably altered by this method.	Affects on material are cumulative, therefore beta sterilised products are not resterilisable, by most conventional methods, after first use.
Chemical (Formol- Aldehydes)	Chemical products are steam and heat	Special applications, minimal industrial use.	Toxic	Requires specialised sterilisation equipment and procedures.

In an **ASEPTIC PRODUCTION**, manufacturing processes, packaging methods and environment are designed to prevent bacterial contamination (bioburden). Our **EUROTUBO®** containers are injected at a temperature of 250 °C and are automatically connected to the assembling and packaging processes, without manual intervention, so as to guarantee an aseptic production.

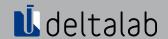

Glove properties chart

	LATEX	NITRILE	VINYL
Mechanical resistance			
- Perforation			
- Tearing			
Adequate for foodstuff handling	×		
Comfort		0	

Glove chemical resistance chart

Diluted mineral acids			
Hydrochloric acid			
Chromic acid	×		
Nitric acid			
Perchloric acid			
Phosphoric acid and Sulphuric acid			
Concentrated acids			
Chlorhidric acid			
Chromic acid	×		
Nitric acid			
Sulphuric acid	×	*	
Hydrocarbon and oil by products			
Aniline			×
White spirit	×		
Styrene	×		•
Gasoline, Hexane, Paraffine, Kerosene	×		



Maximum resistance of nitrile gloves 125 °C.

: Excelent; : Limited; *

: Good;: Not recommended

These indications are solely for information. It is strongly recommended that the user carries out preliminary testings.

CENTRIFUGATION:

Conversion G - R.P.M.

The relative centrifugal force (RCF) can be determined using the formula:

FCR= 1,118 x 10-6 x r x n²


r= radius of the rotor (mm); distance between the axis of the rotor and the farest wall of the tube. n= rotating speed (revolutions per minute).

The result is expressed in terms of acceleration (g); 1 g is equal to 9,807 m/s².

It is recommended that the centrifuge caps fit in size and shape to the tubes to be centrifuged.

Conversion chart xg - R.P.M.

r xg	1,000 xg	1,500 xg	2,000 xg	2,500 xg	3,000 xg	3,500 xg	4,000 xg	4,500 xg	5,000 xg	10,000 xg	15,000 xg
50 mm	rpm 4,227	5,177	5,978	6,683	7,321	7,908	8,454	8,967	9,452	13,367	16,371
75 mm	3,451	rpm 4,227	4,881	5,457	5,978	6,457	6,903	7,321	7,717	10,914	13,367
100 mm	2,989	3,661	rpm 4,227	4,726	5,177	5,592	5,978	6,340	6,683	9,452	11,576
125 mm	2,673	3,274	3,781	rpm 4,227	4,630	5,001	5,347	5,671	5,978	8,454	10,354
150 mm	2,440	2,989	3,451	3,859	rpm 4,227	4,566	4,881	5,177	5,457	7,717	9,452
175 mm	2,259	2,767	3,195	3,572	3,913	rpm 4,227	4,519	4,793	5,052	7,145	8,751
200 mm	2,113	2,588	2,989	3,342	3,661	3,954	rpm 4,227	4,483	4,725	6,683	8,185
225 mm	1,993	2,440	2,818	3,151	3,451	3,728	3,985	rpm 4,227	4,456	6,301	7,717
250 mm	1,890	2,315	2,673	2,989	3,274	3,537	3,781	4,010	rpm 4,227	5,978	7,321

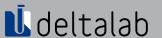


METRIC EQUIVALENCES

ANGLO-SAXON MEASURES	METRIC SYSTEM						
Length: 1 inch (in) = 2.54 cm 12 in = 1 foot (ft) = 30.48 cm 3 ft = 1 yard (yd) = 91.44 m 1,760 yd = 1 mile = 5,280 ft = ±1.6093 km	1 centimeter (cm) = 0.01 metro (m) 1 millimeter (mm) = 0.1 cm 1 micron (µm) = 0.001 mm 1 nanometer (nm) = 0.001 µm						
Volume: 1 cubic inch (in³) = 16.387 cm³ 1,728 in³ = 1 cubic foot (ft³) = 0.0283 m³ 27 ft³ = 1 cubic yard (yd³) = 0.7646 m³	$1,000 \text{ mm}^3 = 1 \text{ cubic centimeter (cm}^3)$ $1,000 \text{ cm}^3 = 1 \text{ cubic decimeter (dm}^3)$ $1,000 \text{ dm}^3 = 1 \text{ cubic meter (m}^3)$						
Capacity: 1 ounce EUA (fl.oz) = 29.57 ml 1 gallon EUA = 3.7854 l 1 gill = 0.142 l 4 gills = 1 pint = 0.5683 l 2 pints (pt) = 1 cuarto (qt) = 1.1365 l 4 qt = 1 brit gallon = 4.5461 l	1 microliter (µI)= 0.001 mI 1 milliliter (mI) = 0.001 I 1 centiliter (cI) = 0.01 I 1 deciliter (dI) = 0.1 I 1 cubic decimeter (dm³) = 1 I 1 hectoliter (hI) = 100 I						
Weight: 1 ounce (oz) = 28,349 g (1 g = 0.035 oz) 16 oz = 1 pound (lb) = 0.4536 kg 112 lb = 1 hundredweight (cwt) = 50.80 kg	1 kilogramme (kg) = 1,000 g 1 milligramme(mg) = 0.001 g 1 microgramme (µg) = 0.001 mg						

PRESSURE: CONVERSION

IS EQUAL TO:							
1 at (technical atmosphere)	0.980665 bar						
1 Atm (physical atmosphere)	1.01325 bar						
1 Atm	1.033 at						
1 bar	105 Pa (Pascal)						
1 Torr	1.3332 mbar (milibar)						

HOW TO FIGURE OUT VOLUMES

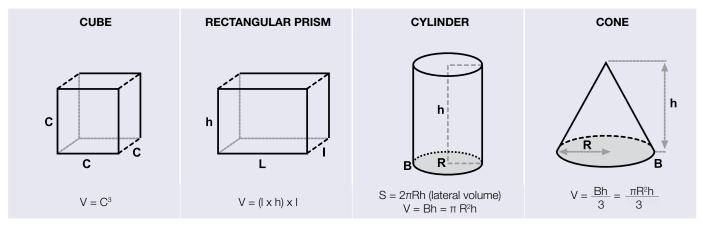

CONVERSION (°F - °C - °K)

$^{\circ}$ F (Farenheit) = ($^{\circ}$ C x 1,8) + 32	°C (Celsius) = (°F – 32) x 0,556	$^{\circ}$ K (Kelvin) = $^{\circ}$ C + 273,15

USEFUL INFORMATION: TRANSPORT

Containers: capacity:	20 feet: 28 m3 maximum
	40 feet: 58 m3 maximum
	40 HC feet: 68 m3 maximum
Pallets:	maximum height land transport: 2.40 m
	maximum height sea transport: 2.10 m
	width x depth european pallet: 80 x 120 cm
	width x depth american pallet: 100 x 120 cm





Swabs for microbiological sampling

Swabs are used for biological sampling.

Specially used for processing samples which, after being coloured, will be analysed by microscopy.

Also suitable for isolations in culture medium.

Another important use is its capacity to spread dishes by dissemination techniques (for example for susceptibility testing by Kirby-Bauer or by E-test).

Swabs class IIa

The swabs for the collection and transport of microbiological samples are considered invasive surgical-type product and with a use, in compliance with the essential requirements of the "Directive 93/42/EEC. Medical devices". According to the classification rules contained in this Directive, the swabs are part of Class IIa.

The **Agencia Española del Medicamento y Productos Santarios** (AEMPS – ON 0318) is the institution who certifies these swabs.

Quality standards - Requirements fulfilled by swabs

- 1. UNE-EN ISO 556-1 Sanitary products sterilisation. Requirements to be designed as "STERILE". Part 1: sterilised sanitary products requirements.
- 2. UNE-EN ISO 868-2 Packaging materials for medical devices sterilized in its final phase. Part 2: sterilization wrap. Requirements and test methods.
- 3. UNE-EN ISO 11737-1 Sterilization of medical devices. Microbiological methods. Part 1: determination of a population of microorganisms.
- 4. UNE-EN ISO 11737-2 Sterilization of medical devices. Microbiological methods. Part 2: sterility tests conducted to validate a sterilization process.
- 5. UNE-EN ISO 11135-1 Sterilization of medical devices. Ethylene oxide. Part 1: requirements for development, validation and control of the routine of a sterilization process of medical devices.
- 6. UNE-EN ISO 11137-1 Sterilization of health care products. Radiation. Part 1: requirements for the development, validation, and control of the routine of a sterilization process of PS.
- 7. UNE-EN ISO 11137-2 Sterilization of health care products. Radiation. Part 2: setting the sterilization dose.
- 8. UNE-EN ISO 15223-1 Sanitary products. Symbols to use on labels, labelling and information to be supplied. Part 1: general requirements.
- 9. UNE-EN ISO 14971 Sanitary products. Application of risk on management.
- 10. UNE-EN ISO 13485 Sanitary products. Systems of quality management. Requirements for regulatory purposes.

Presentations










- 1. **In bulk**: the non-sterile swabs presented in bags of 100 units.
- 2. **Peel-pack package**: single bag. Paper (fiber free) which can be peeled. It consists in one side of medical paper and another side of polyethylene.
- 3. **Flow-pack package**: single bag. Polypropylene bioriented bag.
- 4. Package in polypropylene tube: The product is identified on the label.

All sort of package have the following parameters printed on: product code, product description, lot number, expiry date, CE mark, manufacturer name and address, sterilisation method and single use mark ②.



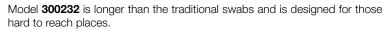

For use in gynaecology when taking endocervical cells while using the speculum.

Dimensions:

Wood

300232 shaft 200 x 2.5 mm. Tip 5 mm Ø.

Other models dimensions:


Wood: 150×2.2 mm (tip $\varnothing \pm 5$ mm). Plastic: 150×2.5 mm (tip $\varnothing \pm 5$ mm). Aluminium: 150×0.9 mm (tip $\varnothing \pm 1.5$ mm).


They are supplied in bags of 100 units, excepting code **300205**, supplied bulk.

code	description	case quantity	case weight	case volume
300232	extra large wood (200 mm) + cotton	80 x 100	5.85	0.030
300230	standard wood (150 mm) + cotton	100 x 100	4.40	0.032
300260	polypropylene + cotton	40 x 100	2.20	0.010
300270	snappable polystyrene + cotton	40 x 100	2.50	0.016
300268	snappable polystyrene + polyester	40 x 100	5.40	0.030
300243	aluminium + cotton	10 x 100	3.00	0.005
300205*	wooden shaft only	1 x 1,000	0.41	0.001

* Non **CE** product Expiry date: 60 months

All swabs have a round, compact and shred resistant head for a proper sampling avoiding fiber loose.

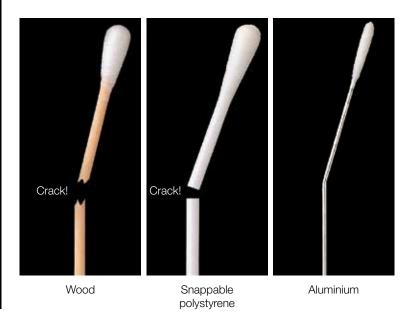

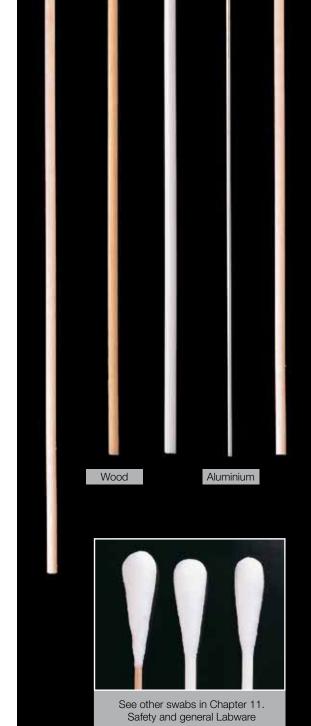

Wood

Polystyrene



Sterile swabs individually wrapped in peel-pack (sterilized by ethylene oxide) or flow pack (sterilized by radiation), depending on the model.

For use when samples do not need to be transported.

We recommend our transport swabs in tubes when sample transport is required (see next page).

code	description	package	case quantity	case weight	case volume
300200	wood + cotton	peel-pack	2 x 1,000	2.6	0.029
310200	wood + cotton	flow-pack	2 x 1,000	2.6	0.029
300201	snappable PS + cotton	peel-pack	2 x 1,000	3.5	0.029
300202	snappable PS + viscose	peel-pack	2 x 1,000	3.5	0.029
310202	snappable PS + viscose	flow-pack	2 x 1,000	2.02	0.029
300203	aluminium + cotton	peel-pack	2 x 1,000	2.4	0.029
310253.1	aluminium + viscose	flow-pack	2 x 1,000	2.3	0.029

Cases per pallet: 54.

1. Flow-pack A: Wood

2. Peel-pack B: Snappable polystyrene

C: Aluminium

Expiry date: 48 month from sterilization date.



Sterile swabs (2 units)

(€ class IIa (MDD)

Two Sterile swabs wrapped in peel-pack (sterilized by ethylene oxide) or flow pack (sterilized by radiation), depending on the model.

One swab is designed for cleaning the sampling area.


The other swab is designed for sample collection. For use when samples do not need to be transported.

We recommend our transport swabs in tubes when sample transport is required

1,000 peel-packs/flow-packs with 2 units each one per case.

code	description	package	case quantity	case weight	case volume	cases per pallet
300210	wood + cotton	peel-pack	1,000	2.54	0.025	54
300211.1	snappable PS + cotton	peel-pack	1,000	2.72	0.025	54
310211.1	snappable PS + cotton	flow-pack	1,000	2.90	0.029	54

Expiry date: 48 month from sterilization date.

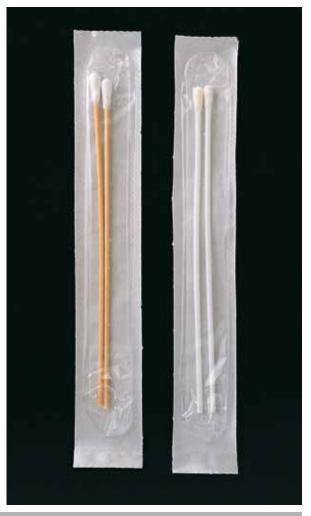







Sterile swab in round tube

(€ class IIa (MDD)


Sterile dry swabs supplied in shockproof round bottom polypropylene tube, with a label sealing the cap. Dimensions of tube: \emptyset 13 x 165 mm. Sterilised by ethylene oxide.

code	description	case quantity	case weight	case volume	cases per pallet
300250	wood + cotton	4 x 500	15.00	0.065	24
300259	wood + viscose	4 x 500	14.00	0.070	24
300261	snappable polystyrene + cotton	4 x 500	14.20	0.070	24
300252	snappable polystyrene + viscose	4 x 500	14.20	0.070	24
300251	aluminium + cotton	4 x 500	14.00	0.070	24
300253	aluminium + viscose	4 x 500	13.90	0.070	24
300263	snappable polystyrene + polyester	4 x 500	14.20	0.070	24
300250.1	wood + pure cotton	4 x 500	14.00	0.070	24
300254	twisted aluminium + viscose	4 x 500	13.40	0.066	24
300265	polystyrene + flocked polyester	4 x 500	14.20	0.065	24

Expiry date: 48 month from sterilization date.



Twisted aluminium: This is a new material excepcionally flexible to use it in sensitive zones as nasopharyngeal, ocular or for pediatric uses.

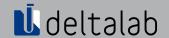





Human DNA free certified swabs, steriles

Human DNA free Certified. Sterile.

The swab is supplied in a polypropylene tube, which protects the sample up to the laboratory prior to its analysis. The stick of the swab is made of polystyrene while the head is produced with viscose, polyester or cotton according to the code. The tube is labeled indicating code, description, lot, expiry date and providing an identifying area to note down collection details (site, date, etc.). Moreover, the label seals the tube with the cap of the swab, acting like a tamper evident system.

	code	description	selling unit	case quantity	case weight	case volume
Nav	300252DNA	polystyrene + viscose	500	4 x 500	14.20	0.066
NEW	300263DNA	polystyrene + polyester (dacron® type)	500	4 x 500	14.20	0.066
WEW	300261DNA	polystyrene + cotton	500	4 x 500	14.20	0.066

Expiry date: 48 months from sterilisation date







Collection and transport system for microbiological samples

- This innovative systems brings together a **flocked swab** with a **liquid transport medium**
- The sample is suspended on the medium getting higher feasibility and reliability
- Compatible with automatic inoculation and streaking, as well as with the manual ones
- All products of δswab® range consist in:

Flocked swab

- Guarantees the maximum absorption and elution of the sample in the medium
- Offers a great versatility, adaptability and comfort for the user



Standard

Urethral

Nasopharyngeal

Pediatric

Liquid media


- Concept change in the collection and transport of microbiological samples
- The liquid medium makes possible the obtention of different homogeneous innoculum of one sample

Transport media

 Allows the transport of multiple microorganisms assuring its feasibility until the analysis

Enrichment media

• Designed for the isolation and culture of bacteria

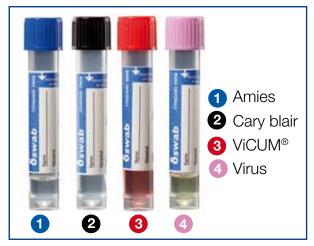



The **\delta**swab® liquid media range has been designed following the methodology of the CLSI Standard, M40-A (Quality Control of Microbiology Transport Systems: approved standard)

Instructions













Transport swabs with media



Safety

 Paper less, plastic peel-pack is an extraordinary safe system to maintain pack integrity and sterility

Comfort

- Space saving
- Collapsible internal case
- Easy wharehouse identification
- Selling unit easy to commercialize
- New standardised pack quantity for all our swabs
- The peel-pack has a small slit at the top to allow easy opening

Information

- Full range of transport medias listed on the package
 Clear dual CE Marking IVD + MDD printed in both
- Clear dual CE Marking IVD + MDD printed in both, peel-pack and case
- User's guide provided in both, peel-pack and case

Design

- Attractive peel-pack and friendly users guide
- Eye-catching case designed exclusively for swabs
- Drawings of a swab and a tube on both the peelpack and the case for first sight identificacion

Traceability

- Each single peel-pack is printed with the catalogue number, bar code, lot mumber, expiry date and product description
- Both external and internal cases are also identified with catalogue number, bar code, lot number, expiry date and product description.

Environmentally friendly

 Peel pack made of one material only: easy recycling!

(€ Swabs

As the two main components of a transport swab have different purposes and different ways of interaction with the patient, both have different CE treatment:

- 1. **Swab**: for microbiological sample collection. As a surgical invasive product of temporary use, is classified as CLASS IIa product in compliance with Guideline 93/42/EEC (Directive of Sanitary Products) amended by Directive 2007/47/EC.
- 2. **Tube with transport media**: container facilitator of survival and transport of biological samples prior to analysis.

Classified as in vitro diagnostic IVD (according to Directive 98/79/EC).

Transport media

The media used to transport the sample, has the appropriate chemical composition that maintains the bioburden levels with the minimum physiological activity, achieving thus the minimal quantitative and qualitative variation of the sample from the collection to its analysis.

Components:

Each set consists in a tube and a swab, packed in a peel-pack bag (A) and sterilised by radiation.

- (B) **Tube with transport media:** Made of rigid polypropylene with round bottom.
- (C) Polyethylene **cap** specially designed to close hermetically; non slippery surface when using gloves.
- (D) Coloured **label** in each tube with the following parameters printed: transport media type, sterilisation method, manufacturer name, CE mark, "single use" and "see the instructions" marks, lot number and expiry date, as well as spaces to write patient name, date, time, medical report number, source and kind of sample. The label seals the cap and the tube.
- (E) **Swab** fitted in a second cap with a shaft which ends in a tip of synthetic or organic material hardly adhered to the shaft. The cap provides a leakproof seal after sample collection.

Both the peel-pack package and the label of the tube are printed with the following information: transport media type, sterilization method, manufacturer name, CE mark, "single use" and "see the instructions" marks, lot number, and expiry date. Peel pack also includes bar code.





Amies. Sterile

Swab: (class IIa (MDD)
Tube with media: ((IVD)

AMIES media is a modification of Cary Blair media, which is at the same time a modification of Stuart media.

Basically the glycerophosphate is replaced by an inorganic phosphate, and the methylene blue by pharmacological charcoal.

Calcium and magnesium are also added, thus maintaining the permeability of the bacterial cell.

This media assures the viability of organisms such as:

Neisseria sp.

Haemophilus sp.

Corynebacterium sp.

Trichomonas vaginalis

Streptococcus pyogenes

Streptococcus pneumoniae

Shigella flexneri

Salmonella typhi

Brucella abortus

Enterobacterias, etc.

The presence of charcoal on the media neutralises inhibitors and bacterial toxins, and improves the recovery ratio of Neisseria gonorrhoeae and Vibrio cholerae

Some microorganisms can resist on the media for up to 3 days, although the recovery of microorganisms is better if cultured in the first 24 hours. Amies transport swabs are available with or without charcoal.

Swabs are sterilised by radiation. Peel-pack dimensions: 38 x 210 mm. Expiry date: 30 months after sterilisation date.

	code	description	case quantity	case weight	case volume	cases per pallet
	300287	snappable polystyrene + viscose	6 x 100	8.5	0.055	32
	300285	snappable polystyrene + viscose (with charcoal)	6 x 100	9.18	0.052	32
	300281	aluminum + viscose	6 x 100	9.15	0.057	32
	300281/1	aluminum + viscose (with charcoal)	6 x 100	9.15	0.052	32
NEW	300281T	twisted aluminum + viscose	6 x 100	8.24	0.057	32
NEW	300281TC	twisted aluminum + viscose (with charcoal)	6 x 100	8.24	0.057	32



Liquid Amies. Sterile

Swab: (€ class IIa (MDD) Tube with media: (€ (IVD)

This is a variation of the traditional Amies media swab where the media is presented on liquid form without agar.

Suitable for conservation of pathogens such as:

Haemophilus sp.

Corynebacterium sp.

Trichomonas vaginalis

Streptococcus pyogenes

Streptococcus pneumoniae

Shigella flexneri

Salmonella typhi

Brucella abortus

Staphylococcus epidermidis

Escherichia coli, etc.

Media is supplied in its liquid form, in a sponge, suitable for direct extensions on slides. **Swabs are sterilised by radiation.**

Peel-pack dimensions: 38 x 210 mm.

Expiry date: 30 months after sterilisation date.

code	description	case quantity	case weight	case volume	cases per pallet
300284	snappable polystyrene + viscose	6 x 100	7.5	0.057	32

Stuart swabs. Sterile

Swab: (€ class IIa (MDD) Tube with media: (€ (IVD)

The modified Stuart media allows the conservation and transportation of a large number of pathological microorganisms, such as:

Neisseria gonorrhoeae

Heamophilus influenzae

Neisseria meningiditis

Bordetella pertusis

Corynebacterium diphteriae

Trichomonas vaginalis

Staphylococcus aeureus

Streptococcus sp.

Salmonella sp.

Shigellas sp.

etc.

The most unstable organisms will remain viable for up to 24 hours and other for several days. The media is reduced due to the presence of thioglycolate, which difficults the enzymatic reactions of the bacteria.

The multiplication of the bacteria is prevented due to the lack of nitrogen in the media.

Swabs are sterilised by radiation.

Dimensions of the peel-pack: 38 x 210 mm. Expiry date: 30 months from sterilisation date.

	code	description	quantity	weight	volume	pallet
	300290	wood + cotton	6 x 100	8.5	0.052	32
	300291	aluminium + cotton	6 x 100	8.5	0.052	32
	300295	snappable polystyrene + viscose	6 x 100	8.5	0.052	32
į	300295C	snappable polystyrene + viscose (with charcoal)	6 x 100	8.9	0.057	32
į	300292	twisted aluminium + viscose	6 x 90	7.6	0.057	32



Cary Blair. Sterile

Swab: (€ class IIa (MDD) Tube with media: (€ (IVD)

Cary Blair is another modification of Stuart media.

The glycerophosphate has been substituted by inorganic phosphate since glycerophosphate is a metabolite for some bacteria, so they could grow and disgrise the pathogens amount.

The methylene blue has also been removed and the pH increased to 8.4.

This is a media originally developed for faecal samples but is also used successfully for anaerobic transport, such as:

Neisseria gonorrhoeae

Vibrio cholerae

Vibrio parahaemolyticus

Heamophilus influenzae

Neisseria meningiditis

Bordetella pertusis

Streptococcus pneumoniae

Shigella flexneri

Pasteurella pestis

Campylobacter Spp., etc.

Swabs are **sterilised by radiation**.

Dimensions of the peel-pack: 38 x 210 mm. Expiry date: 30 months from sterilisation date.

code	description	case quantity	case weight	case volume	cases per pallet
300280	wood + cotton	6 x 100	8.5	0.052	32
300280.2	snappable polystyrene + viscose	6 x 100	8.5	0.052	32









Media for virus. Sterile

Swab: (€ class IIa (MDD) Tube with media: (€ (IVD)

This viral transport media preserves the specimen during transportation to the laboratory.

Suitable, among many others:

Papiloma

Pseudorrabia

Influenza aviar (H7N1)

Influenza A (H1N1) o H1N1/09 pandemic

Suis herpesuirus, etc.

The media maintains the sample viable for up to 72 hours and in some cases for longer time.

The addition of antimicrobial substances inhibit the growth of bacteria and fungi. Ideal for nasal, pharyngeal, ocular, skin and mucus samples.

Media is supplied in its liquid form, in a sponge, suitable for direct extensions on slides.

Swabs are sterilized by radiaton.

Expiry date: 18 months from sterilisation date.

code	description	case quantity	case weight	case volume
300297	snappable polystyrene + polyester	6 x 100	7.8	0.052
300294	aluminium + polyster	6 x 100	7.8	0.052

Minimum order quantity: 100 units.

Cases per pallet: 32.



Media for Chlamydia. Sterile

Swab: (€ class IIa (MDD) Tube with media: (€ (IVD)

Liquid media, suitable for Chlamydia.

For cervical samples it is recommended to firstly use a dry swab to clear the cervical channel.

Media is supplied in its liquid form, in a sponge, suitable for direct extensions on slides. **Sterilised by radiaton**.

Expiry date: 12 months from sterilisation.

code	description	case quantity	case weight	case volume
300299	snappable polystyrene + polyester	6 x 100	7.8	0.052

Minimum order quantity: 100 units.

Cases per pallet: 32.



90 x 14 mm Petri Dish

Made in polystyrene.


Supplied in groups of 20 units, packaged in heat sealed bags.

Code 200200 is aseptic.

Code 200209 is sterile by radiation.

Suitable for automatic filling machines.

code	description	sterile	case quantity	case weight	case volume	cases per pallet		
With three vents								
200200	Ø 90 x 14 mm	no	25 x 20	7.71	0.071	28		
200209	Ø 90 x 14 mm	RADIATION	25 x 20	7.50	0.071	28		
Not vented dishes, for the cultivation of anaerobic								
200200.4	Ø 90 x 14 mm	no	25 x 20	6.90	0.071	28		
200200.4	Ø 90 x 14 mm	no	25 x 20	6.90	0.071	28		

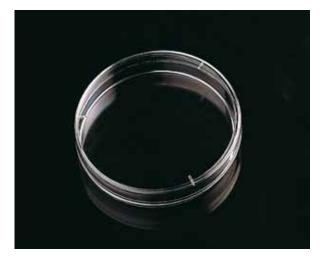

90 x 14 mm Petri Dish, two compartments

Made in polystyrene. **Aseptic**.

With three vents.

Supplied in groups of 20 units, packaged in heat sealed bags.

Suitable for automatic filling machines.

code	description	case quantity	case weight	case volume	cases per pallet
200202	Ø 90 x 14 mm 2 compartments	25 x 20	7.46	0.072	28



90 x 14 mm Petri Dish, three compartments

Made in polystyrene. **Aseptic**.

With three vents.

Supplied in groups of 20 units, packaged in heat sealed bags. Suitable for automatic filling machines.

code	description	case quantity	case weight	case volume	cases per pallet
200203	Ø 90 x 14 mm 3 compartments	25 x 20	7.82	0.070	28



90 x 15 mm Petri Dish, four compartments

Made in polystyrene. **Aseptic**.

Vented.

Supplied in groups of 24 units, packaged in heat sealed bags. Suitable for automatic filling machines.

code	description	case quantity	case weight	case volume	cases per pallet
200210	Ø 90 x 15 mm 4 compartments	24 x 25	9.20	0.077	28







Petri Dish 90 x 15 mm with internal cross

Made in polystyrene. Aseptic.

With four vents. Stackable.

With a mark on a side that helps orientation and a writing area on the bottom. Supplied in groups of 25 units, packaged in heat sealed bags.

Suitable for automatic filling machines.

Internal cross guarantees the totally flat base.

code	description	case quantity	case weight	case volume	cases per pallet
200200.5	Ø 90 x 15 mm internal cross	24 x 25	7.50	0.076	20



140 x 20 mm Petri Dish

Made in polystyrene. With three vents.

Supplied in groups of 10 units, packaged in heat sealed bags.

Suitable for automatic filling machines.

Code 200214 is aseptic.

Code 200219 is sterile by radiation.

code	description	sterile	case quantity	case weight	case volume	cases per pallet
200214	Ø 140 x 20 mm	no	12 x 10	5.84	0.068	30
200219	Ø 140 x 20 mm	RADIATION	12 x 10	5.18	0.066	30



120 x 120 mm squared Petri Dish

Made in polystyrene. **Aseptic**.

With four vents.

Supplied in groups of 10 units, packaged in heat sealed bags.

Suitable for automatic filling machines.

code	description	case quantity	case weight	case volume	cases per pallet
200204	120 x 120 mm	24 x 10	11.18	0.070	28



55 x 14 mm Petri Dish

Made in polystyrene. With four vents.

Supplied in groups of 15 units, packaged in heat sealed bags.

Suitable for automatic filling machines.

Code 200201 is aseptic.

Code 200201.B is sterile by radiation.

code	description	sterile	case quantity	case weight	case volume	cases per pallet
200201	Ø 55 x 14 mm	no	80 x 15	8.36	0.068	30
200201.B	Ø 55 x 14 mm	RADIATION	80 x 15	8.40	0.068	30



Ø 100 mm Petri dish baskets

These stainless steel baskets are suitable for dishes of diameter up to 10 cm. Two models available to hold 10 or 20 Petri dishes.



^{*}Lenght x width x height. Heigh with handle: 360 mm.

We can supply customer-designed models.



Moulded 90 or 100 mm Petri dish rack

Useful for safe transport and incubation 90 or 100 mm Petri dishes and helpful to streamline bench top sample processing.

Each rack holds up to forty-two Petri dishes with lids, configured in six columns of seven plates.

Up to 4 racks can be stacked together, and all interlock to prevent accidental knock over.

Access slots under the bottom Petri dishes allows a stack of dishes to be safely removed from the rack using a secure ring that places a finger under the bottom dish and thumb on top of the lid if the top dish.


A large centre divider provides space for labelling and incorporates a handle. All parts are moulded with white polypropylene and are steam **autoclavable** at 121 °C.

Racks are supplied individually bagged, flat, and can be easily assembled (instructions included).

code	dimensions	case	case	case
	mm	quantity	weight	volume
H-611	330 x 210 x 178	1	0.350	0.004





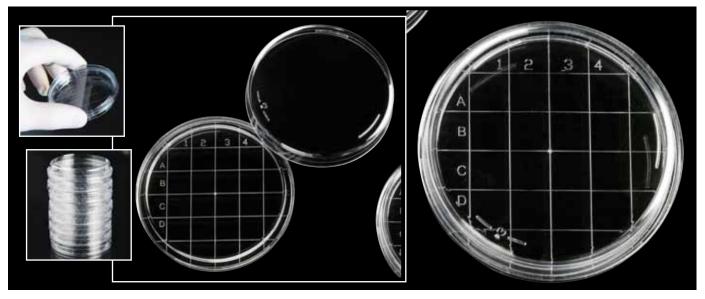

Contact plates

Aseptic production.

Used in the pharmaceutical or food-processing industry, hospital environments for determining bacterial contamination of surfaces such as the skin, operating tables, refrigerated gondolas, and work surfaces.

They are manufactured in transparent polystyrene.

The moulded grid at the bottom makes it possible to identify the detected contamination per cm² and it facilitates the counting of colonies. **Stable stacking:** The shape of the lid makes stacking perfectly stable during transport and incubation, and saves space on work surfaces. Dimensions of the dish: 65,7 mm x 14,7 mm height.

code	description	sterile	case quantity	case weight	case volume	cases per pallet
200208	Rodac contact plate	no	25 x 20	5.00	0.034	70
200218	Rodac contact plate	RADIATION	25 x 20	5.00	0.034	70



See chapter 5. Sampling and transport



See chapter 4. Histology, microscopy and staining



Moulded contact plate rack

Useful for safe transport and incubation of Rodac contact plates and helpful to streamline bench top sample processing.

Each rack holds up to 60 contact plates with lids (minimum diameter: 65 mm, maximum diameter: 72 mm) configured in six columns of ten plates. Up to 4 racks can be stacked together, and all interlock to prevent accidental knock over.

Access slots under each stack of contact plates allows a stack of plates to be safely removed from the rack using a secure ring that places a finger under the bottom plate and thumb on top of the lid if the top plate. A large centre divider provides space for labelling and incorporates a handle.

All parts are moulded with white polypropylene and are steam ${\bf autoclavable}$ at 121 °C.

Racks are supplied individually bagged, flat, and can be easily assembled (instructions included).

code	dimensions	case	case	case
	mm	quantity	weight	volume
H-610	266 x 165 x 178	1	0.33	0.002



Stackable Petri dish incubation tray

This polypropylene tray increases overall capacity of incubator.

Designed to hold five 100 mm diameter or fifteen 60 mm diameter Petri dishes, it's provided with two large label areas for critical sample identification by marker or label.

With 16 vents.

Not autoclavable.

code	dimensions	bag	bag	bag	
	mm	quantity	weight	volume	
H-615	251 x 237 x 35	3	0.75	0.0105	





Microtitre plates, sterile

Standard 96 well-plate made of crystal-clear polystyrene.

Three models with different well shapes are available.

In general terms, flat bottom plates are used when optical reading is required (for example, using a spectrophotometer); rounded and conical bottom plates are very useful when total recovery of the sample is required; they also offer better performance than flat bottom plates for sample separation through centrifugation.

Some applications:

- Diluting
- Serial sample preparation
- Protein precipitation
- Hemagglutination
- Conical bottom plates may be used for complement fixation
- Flat bottom plates are also suitable for cell growth in suspension, and for ELISA applications, when there is no need of a pre-treated plate (which enhances the antigen-antibody adhesion to the surface of the plate)

Polystyrene lids are supplied separately, see code 900015.

Those plates can also be covered using the sealing mat or the adhesive sealing film presented on page 62.

Available in two versions: non sterile in bags of ten units, or sterilized and individually wrapped.

Plate dimensions: 127.7 x 85.8 mm.


Well plate estimate volume:

900010: 0.281 μl. **900011**: 0.404 μl. **900012**: 0.219 μl.

code	description	sterile	case quantity	case weight	case volume
900010	microtitre plate "U" form	RADIATION	50	2.33	0.013
900010.1	microtitre plate "U" form	non	5 x 10	2.33	0.014
900011	microtitre plate "∐" form	RADIATION	50	2.33	0.013
900011.1	microtitre plate "∐" form	non	5 x 10	2.33	0.013
900012	12 microtitre plate "V" form		50	2.33	0.013
900012.1	microtitre plate "V" form	non	5 x 10	2.33	0.013
900015	lid for microtitre plate	RADIATION	50	1.16	0.013
900015.1	lid for microtitre plate	non	5 x 10	1.16	0.013






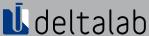









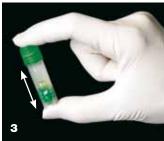



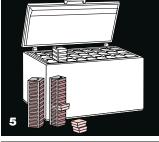

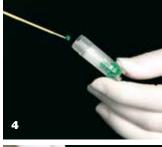





CRYOINSTANT: Porous beads for microbiological culture preservation

CryoInstant is an excellent preservation system for storing microbiological culture (for example fungus on sporulation phase), consisting on **2 ml cryovials** containing 25 porous beads and cryopreservative-added broth.

With these vials we can:

- Get a perfect freezing storage system.
- Obtain up to 25 replicates from the original microorganism, to use them during many years
- Ease the recovery of the stored strain streaking an inoculated bead directly onto solid media
- Prevent changes and repetitive thawing / freezing of the rest of beads
- Avoid the growth of ice when recovering,
- Improve recovery rates of frozen cultures by adding cryopreservative to broth
- Minimise the risk of cross contamination,
- Save space in the freezers when storing.

The external thread cryovials are made of polypropilene (page 145). Their smooth inner surface prevents from contamination. Their caps embody a silicone washer to ensure a positive leakproof seal. Tubes and caps withstand up to **–190 °C**.

In order to ease the classification of the samples, we offer the possibility to get cryovials with caps and beads in five different colours: every code displays caps and beads from a single and same colour (excepting code **409113/6**, which is an assortment), allowing a fast and easy identification of the samples, and providing a system whereby users are able to code different bacterial species, different samples, different laboratories, etc. And even more: with colour coded inserts (described on page 144) a further classification can be made.

Cryovials are presented in carton boxes (resistant to **-100 °C**; dimensions: 150x150x55 mm) with 100 units; boxes are furtherly described in page 206. Each box is supplied plastic-wrapped and labelled with code reference, batch, expiration date and specifying the colour of the caps and beads.

code	beads and cap colour	case quantity	case weight	case volume
409113/1	white	100	0.53	0.002
409113/2	blue	100	0.53	0.002
409113/3	yellow	100	0.53	0.002
409113/4	red	100	0.53	0.002
409113/5	green	100	0.53	0.002
409113/6	assorted: 5 colours x 20 cryovials	100	0.53	0.002

Way of use:

- 1 Take the sample from the strain with a loop (please see our loops on page 32 and 33)
- 2 Inoculate the cryovial by putting the loop in contact with the broth
- 3 Close the cryovial and shake it gently in order to permeate the strain into the beads
- **4** Extract the cryopreservative-added broth with a Pasteur pipette (see our plastic pipettes on pages 189-193)
- 5 Close the cryovial and take it to the freezer
- 6 Every time we need to reproduce the strain, we will take a bead with a loop or a forceps
- 7 We will put it down on a Petri dish with culture medium, attempting to put in contact all the surface of the bead with the medium



Nichrome loops microincinerator

Desktop microincinerator for sterilization of Nichrome loops.

The device reaches a burning temperature of about 1000 °C inside, but their external temperature does not exceed 30 °C.

In this way there is no danger of burns and can be manipulated with gloves without any problem.

After reaching the sterilization temperature in 7 minutes, the loop is fully sterilized in only 2 seconds.

The material of the tube where the sterilization is performed is quartz so that it prevents the formation of aerosols.

Dimensions: 150x250x215 mm.

	code	description	case quantity	case weight	case volume
Waw	DEL001	loops microincinerator	1	5.0	0.015





Nichrome wire loops

Nichrome wire loops for microbiology.

Comparing them to the traditional loops made of platinum and rhodium, the nickel and chrome alloy results much more competitive in terms of durability and price. Moreover, it offers a totally smooth surface, and rapid cooling after heat sterilisation.

Non-calibrated loops (mod. A and B):

Non-calibrated nichrome loops, without holder, saving money and storage space. Supplied in groups of 25 units, in one tube, specifying code, description, and lot.

Calibrated loops (mod. C and D):

Calibrated nichrome loops. Supplied individually in a capped tube, labeled with the code, description and lot.

Calibrated and non-calibrated loops are easily joinable to the PVC and aluminium holder (mod. E) with an easy and fast screw movement.


Dimensions of the aluminium holder 150 mm.


Non-calibrated loops with holder (mod. F and G):

Non-calibrated loop made of nichrome alloy.

Light aluminium holder (**mod. G**, golden colour, **mod. F**, green) partially protected by a transparent, insulating, non-skid plastic. Dimensions, 275 mm.

mod.	code	description	quantity
А	302762	Non-calibrated loop 2 mm Ø	25
В	302764	Non-calibrated loop 4 mm Ø	25
С	302771	Calibrated loop 1 µl	1
D	302772	Calibrated loop 10 μl	1
Е	302780	Loop holder, PVC and aluminium	1
F	302792	Non-calibrated loop (2 mm Ø) with holder	5
G	302794	Non-calibrated loop (4 mm Ø) with holder	5

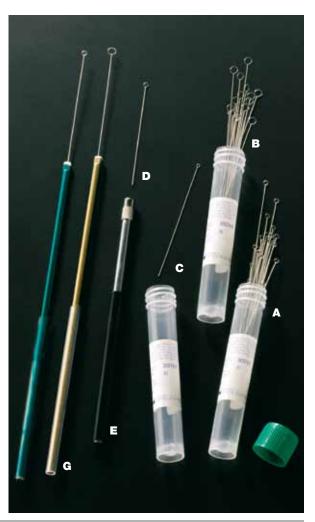

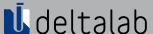





Sterile calibrated loops

Manufactured in flexible polystyrene.

Sterilized by radiation.

Two sizes available: 1 μ l and 10 μ l. Very flexible loop, with rounded tip.

Not recommended for colony counting.

Sterile inoculating loops.

Batch and expiry date printed on the bag.

code	μΙ	presentation	colour	case quantity	case weight	case volume
302703	1	10 ud./peel-pack	green	4 x 500	5.00	0.012
302723	1	20 ud./peel-pack	green	4 x 1,000	7.40	0.039
302704	10	10 ud./peel-pack	blue	4 x 500	5.00	0.037
302724	10	20 ud./peel-pack	blue	4 x 1,000	7.80	0.039



Digralsky spreader. Sterile

Designed for surface spreading. Made of white polystyrene.

Radiation sterilized.

Total length: 149 mm, base length: 40 mm.

Base length allows spreading liquid samples over any type of Petri dish and avoids contact with the dish walls.

The curved Spreader extremity and rounded corners avoid damaging media surface while spreading.

The tip bent and rounded edges minimize the possibility of breakage of the solid medium during extension.

It allows their use in places where sterilization before use is not possible (for instance when samples are collected far from the lab).

Available in peel-pack containing 1 or 5 units, marked with code, expiry date and lot.

code	presentation	case quantity	case weight	case volume
200500	5 units peel-pack	200 x 5	3.6	0.025
200510	1 unit peel-pack	500	3.1	0.025





Sterile culture tubes in polystyrene

Tubes supplied with either a two position ribbed polyethylene cap, which can be left loose for aerobic work or sealed for anaerobic cultures.

They are biologically inert, exempt from mold release agents, and with stand up to ${\bf 1,400~xg}$ and ${\bf 70~^{\circ}C}.$

Packaged in self-standing resealable zip-lock bags of 125 units.

code	dimensions mm	case quantity	case weight	case volume
300807	12 x 75	8 x 125	4.24	0.033
300808	17 x 100	8 x 125	7.14	0.060







PosSealed position for anaerobic cultures





EUROTUBO® 12 ml screw cap tubes, round bottom

C€ (IVD)

Made of **autoclavable** (121 °C) transparent polypropylene or polystyrene. Green cap made of high density polyethylene. Dimensions: 15 x 102 mm. The external skirt allows the tubes to remain **free-standing.** Supplied screwed. Recommended volume: 12 ml. The sterile model (ethylene oxide) is supplied individually in flow-pack bag, with indication of batch, expiration date, etc. According to the guidelines for sterile products. Resistance to centrifugation: **PS: 7,500 xg. PP: 15,000 xg.**

Attention: For autoclaving, the cap should be loose on the thread and not tightly fitted.

	mod.	code	description	maximum volume ml	volume recom. ml	case quantity	case weight	case volume	cases pallet
	1	301403	polystyrene	14.4	12	6 x 250	12.66	0.082	20
	2	301402	STERILE EO polystyrene individually wrapped	14.4	12	6 x 250	13.00	0.096	16
121°C	3	401402	STERILE EO polypropylene individually wrapped	14.2	12	6 x 250	11.80	0.096	16
121°C	4	401403	high transparency polypropylene	14.2	12	6 x 250	11.00	0.082	20

Dimensions (±0,09):

code	external cap Ø mm A	external tube Ø mm B	internal tube Ø mm C	length with cap mm D	length without cap mm E	D————	E
301402, 301403	20.9	16.3	14.4	103.9	102.5	A	\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\
401402, 401403	20.9	16.2	14.3	102.9	101.5		₩±↓



15 ml and 50 ml centrifugal tubes

C€ (IVD)

 $\label{eq:made_of_polypropylene} \mbox{Made of polypropylene, suitable for both clinical and research applications.}$

DNAsa, RNAsa and pyrogen free. They are also free from natural rubber and heavy metals.

High transparency of the material for a clear visualization during experiments specially for molecular biology and animal tissue culture.

Tube and cap designed with the system of flat threads for a complete leakproof.

Highly smooth hydrophobic surface for minimum disturbance during centrifugation.

Silk-screen blue graduation in the tube and large white frosted portion for easy writing.

Autoclavable at 121 °C.

Centrifugation resistance: 14.000 xg.

Available models: 15 ml non-skirted and 50 ml skirted and non-skirted.

code	description	presentation	sterile	case weight	case volume	cases per pallet
15 ml tubes						
429940	non-skirted tube	500 tubes in bulk	NO	500	4.5	0.0264
429945	non-skirted tube	20 bag to 25 tubes.	NO	500	4.6	0.0264
429942	non-skirted tube	20 bag to 25 tubes.	STERILE R	500	4.6	0.0264
50 ml tubes						
429930	non-skirted tube	20 bag to 25 tubes.	NO	500	8	0.0663
429931	non-skirted tube	20 bag to 25 tubes.	STERILE R	500	8	0.0663
429950	skirted tube	20 bag to 25 tubes.	NO	500	9	0.0663
429951	skirted tube	20 bag to 25 tubes.	STERILE R	500	9	0.0663





EUROTUBO® 15 ml conical tubes

C€ (IVD)

Tubes made of transparent polypropylene, conical bottom, suitable for centrifugation tests in immunology, microbiology, etc. Continuous thread, external moulded graduations in 0,5 ml increments.

Frosted area (55 x 10 mm).

Blue cap in polyethylene with hermetical closure thanks to its internal elastic obturation.

Code 429910 is autoclavable (121°C) with the cap not being closed, just placed on the thread. Autoclave not recommended for the codes 429920 and 429946, as they have been sterilised by radiation.

Code 429920 is supplied sterile in 25 microns polyethylene bags (100 units).

The bag specifies IVD, lot number, expiration date, etc.

Code **429946** is supplied sterile in individual bags. Resistance to centrifugation: **7,000 xg** (non sterile model) and **5,000 xg** (sterile models). Length (with cap): 120 mm. External mouth diameter: 17 mm. External cap diameter: 20,9 mm.

	code	description	case quantity	case weight	case volume	cases per pallet
21 25 ₹	429910	non sterile	500	3,8	0,029	54
	429920	STERILE	5 x 100	3,9	0,029	54
	429946	STERILE R individually bagged	500	3,9	0,038	40







EUROTUBO® 50 ml conical tubes

C€ (IVD)

Tubes made of transparent polypropylene, conical bottom, suitable for centrifugation tests in immunology, microbiology, etc. Perfect to analyse Koch bacillus. Continuous thread, external moulded graduations in 5 ml increments. Blue cap in high density polyethylene with hermetical closure thanks to its internal elastic obturation. Skirted and non skirted versions, sterilised by radiation and non sterilised versions. Sterile codes are either supplied in individual bags (codes 429926, 429927), in bags of 100 units (codes 429926.10, 429927.10), or in bags of 25 units (codes 429926.25, 429927.25) and resist up to 7,000 xg. Autoclave not recommended as they have been ionised. Non sterile codes withstand up to 12,000 xg and are autoclavable (121 °C) with the cap loose on the thread and not tightly fitted.

Every version is supplied capped, excepting code 429900SP, which is supplied uncapped, with the cap in a separated bag.

code		description		internal tube Ø mm	external tube Ø mm	length with cap mm	case quantity	case weight	case volume	cases per pallet
429900		non skirted		27.2	29.5	117.5	500	7.9	0.075	20
429900SP	no	non skirted, unscrewed cap		27.2	29.5	117.5	500	7.9	0.082	20
429901		skirted		27.2	29.5	117.5	500	8.5	0.082	20
429926	STERILE	non skirted, individually wrapped	34.4	27.2	29.5	117.5	500	8.2	0.082	20
429926.25	STERILE	non skirted, bag 25 units	34.4	27.2	29.5	117.5	20x25	8.2	0.082	20
429926.10	STERILE	non skirted, bag 100 units	34.4	27.2	29.5	117.5	5x100	7.9	0.082	20
429927	STERILE	skirted, individually wrapped	34.4	27.2	29.5	117.5	500	8.8	0.082	20
429927.25	STERILE	skirted, bag 25 units	34.4	27.2	29.5	117.5	20x25	8.8	0.082	20
429927.10	STERILE	skirted, bag 100 units	34.4	27.2	29.5	117.5	5x100	8.5	0.082	20



Please find racks for these tubes in Chapter 9. Sample Storage



POLYPROPYLENE RACKS:

- W-018 • 19568
- 19570


METAL RACKS:

- R-292 R-293
- R-281 R-282
 - R-283







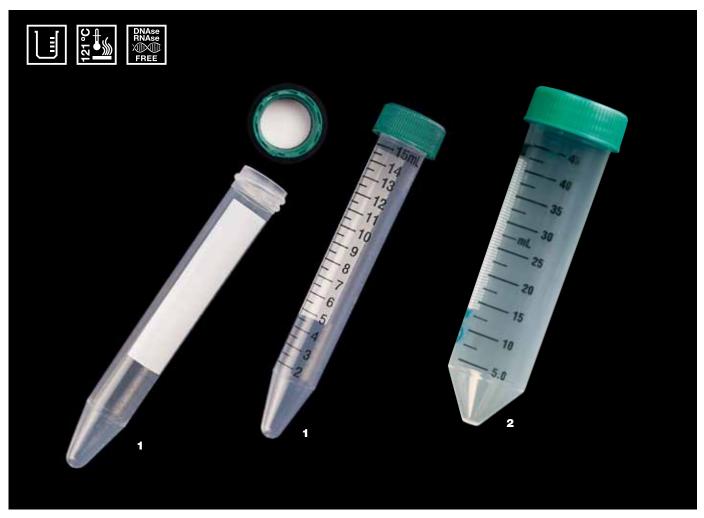

15 ml and 50 ml conical high resistance tubes. Sterile

C€ (IVD)

Tubes made of transparent, copolymer polypropylene.

Green caps made of polyethylene with an internal liner which ensures leakproofness.

Tubes feature a solvent resistant white panel and black graduations for use both for clear or dark samples.

DNAse, RNAse, endotoxins and metal free.

They are sterile by radiation and withstand temperatures down to -90 °C (15 ml) and -80 °C (50 ml)...

Resistance to centrifugation: 17,000 xg (15 ml) and 20.000 xg (50 ml).


Autoclavable tube (the liner of the cap is not autoclavable)

Supplied in bags of 50 units.

Bags are printed with instructions for use and feature a double closure: a first tamper evident seal that helps guaranteeing sterility, and a secondary zip-lock, resealable closure.

	code	presentation	capacity ml	dimensions mm	case quantity	case weight	case volume
1	409920	bag 50 units	15	17 x 118	10 x 50	4.12	0.035
2	409922	bag 50 units	50	29.6 x 114.6	10 x 50	8.16	0.076



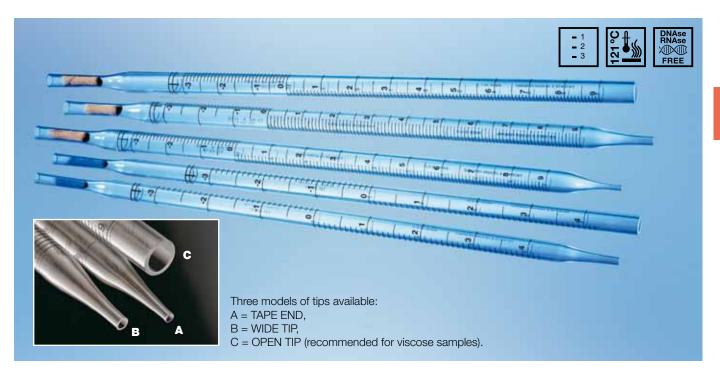

Sterile serological pipettes

Single use, made of crystal polystyrene. Accuracy > 98 %.

Sterilised by radiation. Made in one piece for an optimum precision, except code 90043 (50 ml).

No liquid retention: as tip is not sealed, liquid is not held at the joint. Total dispensing guaranteed.

Non pyrogenic (ISO 10993-11), non cytotoxic (ISO 10993-5) and non haemolytic (ISO 10993-4). **DNase and RNase free.** Volumes range from 1 to 50 ml. Pipettes include a coloured cotton plug for better identification. The peel-pack have a **fiber-free** packaging including a strip for easy opening. Black graduations are bright, black and **permanent**. 5, 10 and 25 ml models have negative graduations for extra capacity, and a second reverse scale (ascending and descending).

	code	capacity ml	presentation	cotton colour	case quantity	tip	graduation ml	negative graduation ml	total capacity ml	case weight	case volume	cases palet
	900030	1	1 peel-pack	•	4 x 500	Α	0/0.9	until –0.3	1.3	10.00	0.068	24
	900031	1	bag of 25	•	4 x 1,000	А	0/0.9	until –0.3	1.3	14.52	0.071	30
	900130	1	1 peel-pack	•	4 x 500	С	0/0.9	until –0.3	1.3	10.10	0.071	30
	900131	1	bag of 25	•	4 x 1,000	С	0/0.9	until –0.3	1.3	14.00	0.071	25
_	900032	2	1 peel-pack	•	4 x 500	А	0/1.8	until –0.4	2.4	11.60	0.068	30
	900033	2	bag of 25	•	4 x 1,000	А	0/1.8	until –0.4	2.4	17.20	0.071	30
	900034	5	1 peel-pack	•	4 x 250	А	0/4	until –3	8	12.20	0.071	24
	900038	5	bag of 25	•	4 x 375	А	0/4	until –3	8	15.00	0.071	30
NEW	900144N	5	1 peel-pack	•	4 x 250	В	0/4	until –3	8	12.20	0.072	30
NEW	900134N	5	1 peel-pack	•	4 x 250	С	0/4	until –3	8	12.20	0.072	30
	900036	10	1 peel-pack	•	4 x 200	А	0/9	until –3	13	11.07	0.071	24
	900037	10	bag of 25	•	4 x 250	А	0/9	until –3	13	13.08	0.071	30
NEW	900146N	10	1 peel-pack	•	4 x 200	В	0/9	until –3	13	11.20	0.072	30
NaW	900136N	10	1 peel-pack	•	4 x 200	С	0/9	until –3	13	11.20	0.072	30
	900041	25	1 peel-pack	•	4 x 100	Α	0/23	until –10	35	11.20	0.072	24
	900043	50	1 peel-pack	•	4 x 25	А	1/2	until –10	60	3.60	0.072	60




Find trays, rotary stands and other products for pipettes in Chapter 11. **Safety and General Labware**



Please find racks for these tubes in Chapter 9. Sample Storage



See our plastic Pasteur pipettes in chapter 8. Liquid Handling

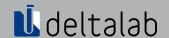

EUROTUBO® Pipetting bulb

This one-handed design is the simplest safety pipette filler to use available. Manufactured in natural orange rubber.

Approx. drawn capacity: 25 ml.

Single hand use, only two operating points. Evacuate via the automatic valve. Standard model, accommodates all pipettes. Ability to clean inside of bulb by removing patented valve and rinsing out.

code	code description		bag weight	bag volume
19200	pippeting bulb	1	0.05	0.004



1. Evacuate the air by 2. Intake by pressing pressing the bulb, as the drawing indicates.




on point $B \bigcirc \Lambda$.



3. Drain the liquid by pressing on point



4. Blow out the pipette by pressing as the drawing indicates (point C).

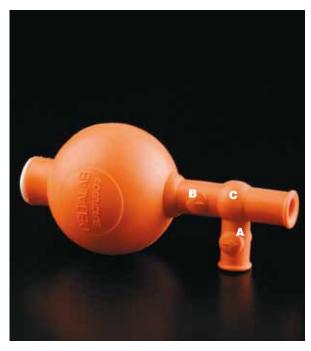

Pipetting bulb

Made of rubber.

Used to avoid mouth pipetting and contamination risk.

Can be opened, cleaned and autoclaved.

Ideal for Wintrobe and Westergren pipettes.

code	description	bag quantity	bag weight	bag volume
19201	red pippeting bulb	1	0.03	0.0001



Pipette pumps

Several models for various pipette volumes.

Designed for fast and efficient pipetting with simple, one handed operation. Pipettes fit smoothly into collar.

Rotate the knurled thumb wheel on the side for precision filling or dispensing, and press the fast release lever for quick emptying.

Easy to use and easily disassembled for cleaning.

Sizes are colour coded.


Pipette pumps resist acids and alkalies.

code	description	colour	bag quantity	bag weight	bag volume
W-100	up to 2 ml	blue	1	0.05	0.0002
W-110	up to 10 ml	green	1	0.05	0.0002
W-120	up to 25 ml	red	1	0.05	0.0002







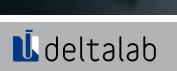



Sterile Whirl-Pak® blender bags with and without filter

Bags made of a low density polyethylene blend, resulting in an extra-resistant and transparent bag.

Leakproof closure with several metallic rounded sticks (see how to use it on figures on page 110).

Codes 200373, 200374 and 200376 feature a filter layer of finely perforated polyethylene, to separate the liquid and solids. This allows the easy pipetting of the sample.


There are 6.45 holes per cm², each one measuring 330 microns diameter.

Bags with a **one-piece seam**, avoiding the possible risks of the loss of the corners bags.

They feature a write-on strip (excluding codes 200342 y 200343).

Ethylene Oxyde sterilised. Made with materials suitables for alimentary use.

mod.	code	description	capacity ml	dimensions cm	tickness microns	case quantity	case weight	case volume
1	200342	standard bag	390	13 x 19	76	500	2.50	0.0170
1	200343	standard bag	720	15 x 23	102	500	3.88	0.0168
2	200351	bag with write-on strip	1,650	19 x 30	102	500	5.66	0.0182
3	200373	bag with filter and write-on strip	720	15 x 23	102	250	1.98	0.0166
3	200376	bag with filter and write-on strip	1,650	19 x 30	102	250	3.88	0.0180
3	200374	bag with filter and write-on strip	2,041	19 x 38	102	250	4.28	0.0170

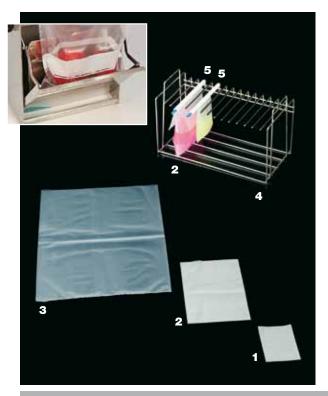

Sterile bags for STOMACHER® blenders

Made of polyethylene. Thickness 7-8/100 mm. Ideal for mixing and blending bacterial samples or when collecting toxic substances from food, fabrics, etc. **Irradiation sterilised**. Supplied in bags of 25 units. Code **15006** is a stainless steel rack for 14 bags. Side handles incorporated.

To close bags. May be used with the rack **15006**. Compounded by a blue tubular piece (230 mm) and a white clip (197 mm) which tightens the bag around the blue tube.

Code **15007** is a clip to adjust bags onto the rack.

Also available non sterile versions. Please ask for the minimum order quantity and delivery time. Add an "S" at the end of the code.

mod.	code	description	capacity ml	case quantity	case weight	case volume
1	15001	bag 100 x 155 mm	80	100 x 25	5.00	0.014
2	15003	bag 180 x 300 mm	400	20 x 25	4.06	0.010
3	15004	bag 380 x 580 mm	3,500	400	12.60	0.030
4	15006	rack 390 x 200 x 240 mm	-	1	1.62	0.039
5	15007	clips	-	10	0.15	0.001



Sterile homogeniser lateral filter bag

Filter bag suitable for PCR and small volumes, it can be used with very short pipette.

The unique and patented Pull-Up system allows to pinch the bag and pull up the filtered liquid for easy pipetting.

Nonwoven side filter having a porosity of 50 microns, rigid and transparent. The filtration is performed during homogenization, it is not necessary to wait for the sedimentation of particles.

Compatible with any lab blender.

Made with materials suitables for alimentary use.

	code	packaging	case quantity	case weight	case volume
NEW	15008	homogeniser lateral filter bag 400 ml	15x20	5.3	0.009



Sterile homogeniser full size filter bag

Bags incorporate a full size filter dividing them into two identical pouches. Insert the sample in one compartment. After sample process, solid particles remain in their original compartment while liquid has filtered into the second compartment, thus allowing a convenient liquid pipetting. Bag capacity: 2 000 ml

Sterilised by radiation. Graduated.

code	code packaging		case weight	case volume	
15005	195 x 300 mm in bags of 10 units	500	7.2	0.04	



Tissue grinders-homogenisers

Borosilicate glass vessels and serrated plunger. Tips: head made of TEFLON, shaft made of stainless steel. Fully **autoclavable**.

The distance between the glass and the plunger is $\pm~200~\mu m.$

Glass vessels:

code	volume ml	body Ø mm	height mm	quantity	weight	volume
196102	2	8	120	1	0.01	0.001
196105	5	12	132	1	0.01	0.001
196110	10	16	150	1	0.02	0.001
196115	15	19	155	1	0.03	0.001

Plunger-serrated tips:

•	-				
code	for tube of mm	height mm	quantity	weight	volume
196302	2	230	1	0.06	0.001
196305	5	235	1	0.06	0.001
196310	10	270	1	0.06	0.001
196315	15	270	1	0.02	0.001







Surface sampling kits

Designed according to the international standard **ISO 18593**. which is taken as reference on European Directives and Rules about microbiology of food and animal feeding stuff, horizontal methods for sampling techniques from surfaces.


The kit includes:

- A bottle containing 40 ml of sterile De-Neutralising solution
- A disposable sterile template, supplied in a flow pack bag
- A sponge or a swab, depending on the code

code	packaging	case quantity	case weight	case volume
200396.V	surface sampling kit with sponge	20	3.00	0.025
200397.V	surface sampling kit with a swab	20	3.00	0.025





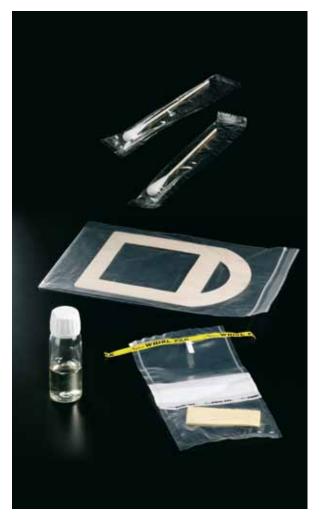

Meat or turkey carcass sampling kit

Designed according to the international standard ISO 17604, which is taken as reference on European Directives and Rules about microbiology of food and animal feeding stuff – carcass sampling for microbiological analysis.

The kit includes:

- A bottle containing 25ml of sterile buffered peptone water
- A disposable sterile template, supplied in a flow pack bag
- A sponge or 2 swabs, depending on the code

code	packaging	case quantity	case weight	case volume
200393.V	carcass sampling kit with sponge	20	3.00	0.025
200394.V	carcass sampling kit with 2 swabs	20	3.00	0.025
200393P	template	20	0.75	0.0051









Sterile jars for water sampling

Transparent squared jars made of PET.

Wide mouth and grips at two sides, easing the sample collection. Red cap made of polypropylene with internal liner made of polexan with a label that performs as tamper evident. **Sterilised by radiation**.

Each bottle is labelled specifying description, code, lot, and expiry date.

Internal mouth diameter mm: 55

Jar weight (gr):

500 ml jar: 44 g - 1,000 ml jar: 65 g - 1,500 ml jar: 72 g

Available with or without thiosulfate sodium

With liquid sodium thiosulfate (24 mg/l): Ideal to analyze water meant for human consumption, swimming pools, and any other water where chlorine may modify the composition of the sample while being transported.

Empty, without thiosulfate: Ideal to sample water for physico-chemical analysis, as well as for other sampling requiring sterile bottles.

mod.	code	description	sterile	capacity ml	case quantity	case weight	case volume	cases per pallet
1	282340	jar with thiosulfate	RADIATION	500	44	2.68	0.036	40
2	282341	jar with thiosulfate	RADIATION	1,000	48	3.72	0.076	20
3	282350	jar without thiosulfate	RADIATION	500	44	2.68	0.036	40
4	282351	jar without thiosulfate	RADIATION	1,000	48	3.72	0.076	20

^{*}Minimum ordering quantity: 1 pallet

We can supply any other concentration according to each country's standards.

We can also dose thiosulfate in other jars from our catalogue.

Optionally presented in individual bags if requested.



See industrial funnels on chapter 12. Laboratory and industrial packaging



Please ask for customized labes







Sterile bottles for water sampling

Rectangular body and cap made of polyethylene. The minor sides are gripped for a better handling.

The larger sides are flat for labelling. Tamper evident red cap with an internal joint. Hermetical closure.

Sterilised by radiation.

Each bottle is labelled specifying description, code, lot, and expiry date.

Dimensions:

500 ml bottle: $83 \times 65 \times 135$ mm - 1,000 ml bottle: $83 \times 65 \times 235$ mm.

Mouth diameter intern: 28 mm. Unitary weight (empty bottle):

500 ml bottle: 40 g - 1,000 ml bottle: 61 g.

Available with or without thiosulfate sodium

With liquid sodium thiosulfate (24 mg/l): Ideal to analyse water meant for human consumption, swimming pools, and any other water where chlorine may modify the composition of the sample while being transported.

Empty, without thiosulfate: Ideal to sample water for physico-chemical analysis, as well as for other sampling requiring sterile bottles.

code	description	sterile	capacity ml	case quantity	case weight	case volume	cases per pallet
282320	bottle with thiosulfate	RADIATION	500	24	1.10	0.025	48
282321	bottle with thiosulfate	RADIATION	1,000	20	1.37	0.033	42
282323	bottle with thiosulfate	RADIATION	500	111	6.15	0.110	16
282323.BU	bottle with thiosulfate , individually wrapped	RADIATION	500	111	6.15	0.110	16
282324	bottle with thiosulfate	RADIATION	1,000	74	5.63	0.110	16
282324.BU	bottle with thiosulfate , individually wrapped	RADIATION	1,000	74	5.63	0.110	16
282330	bottle without thiosulfate	RADIATION	500	24	1.10	0.025	48
282331	bottle without thiosulfate	RADIATION	1,000	20	1.37	0.033	42
282333	bottle without thiosulfate	RADIATION	500	111	6.10	0.110	16
282334	bottle without thiosulfate	RADIATION	1,000	74	5.58	0.110	16

We can supply any other concentration according to each country's standards.

In 250 ml or 125 ml bottles if requested.

Optionally presented in individual bags if requested.



Please ask for customized labels



Sodium thiosulfate sterile bags Whirl-Pak

Specially designed for sampling chlorinated water: waste water, surface and irrigation water supplies, recreational waters, potable water.

Contains 3 tablets of 10 mg sodium thiosulfate.

As an added feature, these bags have a convenient. Write-on strip and markings for 100 ml and 300 ml measurements.

These markings help the technician to collect the correct amount of water as required by the testing procedure.

The write-on strip will accept virtually any writing instrument and minimizes ink smearing on the surface.

code	capacity	dimensions	case	case	case
	ml	cm	quantity	weight	volume
292604	for 300	11.5 x 23	100	0.54	0.0035



Long-handled dippers

Polypropylene cup and polyethylene handle.

For convenient sampling from tanks, vats, streams, ponds, lakes, etc. These dippers are light, easy to use and portable. The polypropylene cups have threaded fittings that screw onto the handle.


Two pouring spouts on the 500 ml dipper make it useful for left or right handed people to pour from either side.

The end of the handle has a hook for hanging.

Autoclavable container (up to 121°C).

code	description	quantity	weight	volume
19575	91 cm handle, 500 ml container	1	0.30	0.009
19576	183 cm handle, 500 ml container	1	0.75	0.020
19577	container 500 ml (fits either handle)	1	0.30	0.005

Please ask the commercial department for other lengths.

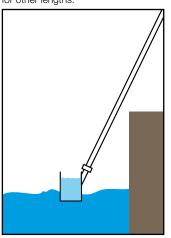











Water control: absence/presence of Coliforms (E. Coli)

125 ml bottle, made of PET, with white ribbed tamper evident cap made of PP. It contains dehydrated culture medium in order to detect the presence of coliforms in human consumption, industrial or leisure water (pools, etc.).

Sterile by radiation.

It allows the water testing directly on a bottle.


The culture medium is based on chromogenic substrates, which make possible the coliform reading; and complemented with fluorogenic substrates, which presumptively identify the presence of E. Coli.

Process

- 1. A sample of water shall be taken in a sterile bottle, without additives, or with a neutralising (like thiosulfate) in case of chlorinated water.
- 2. A part of this sample (100 ml) is poured into the bottle containing the dehydrated culture medium; then shake gently until the medium is dissolved.
- 3. Incubate for 24 hours at 37 °C.
- 4. After incubation, the colour will reveal the results:
 - A Yellow: the sample is negative for non enteric bacteria, total coliforms and E. Coli.
 - B Muddy yellow: the sample is positive for enteric bacteria (not coliforms).
 - C Turquoise blue: the presence of total coliforms is confirmed.
- 5. To discard escherichia place the sample under a UV light in a dark environment; if the sample gets fluorescent it contains E.coli; if there is not fluorescency, then it is another coliform, not E. Coli.

Every case includes two boxes with ten bottles each.

code	description	case quantity	case weight	case volume
282300	absence/presence kit	2 x 10	0.84	0.012

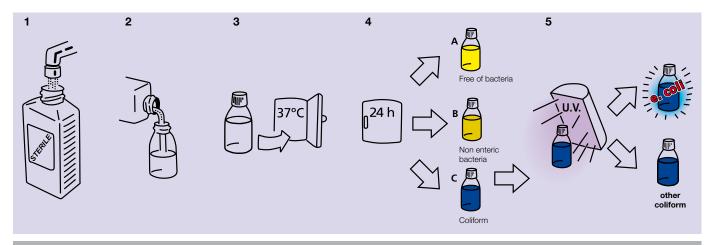







Screw thread microtubes

Made of **autoclavable** polypropylene, they can be used at extreme temperatures from **-190 °C to +121 °C**. Two versions available: in transparent polypropylene, or opaque brown **(UV resistant**, designed to be used with light sensitive samples).

Certified RNAse, DNAse and pyrogen free. Tubes and caps can be centrifuged at 17,000 xg. Caps are supplied separately, see below. Dimensions: $11 \times 44 \text{ mm.}$

mod.	code	volume ml	skirt	case quantity	case weight	case volume				
Transpare	Transparent polypropylene									
1	409110.1	0.5	yes	1,000	1.46	0.010				
2	409110.2	1.5	yes	1,000	1.45	0.010				
3	409110.3	1.5	no	1,000	1.13	0.010				
4	409110.4	2.0	yes	1,000	1.34	0.010				
Brown po	lypropylene									
5	409113.1	0.5	yes	1,000	1.53	0.010				
6	409113.2	1.5	yes	1,000	1.34	0.010				
7	409113.3	1.5	no	1,000	1.34	0.009				
8	409113.4	2	yes	1,000	1.34	0.009				



Caps for screw thread microtubes

Made of polypropylene.

Caps are available either with an attachment loop or without loop.

Both models have a sealing O-ring (red) of silicone to ensure a positive leakproof seal.

For sample identification, colour coding inserts can be placed upon caps (made of polypropylene).

Cap dimensions: 13 x 8 mm.

mod.	code	description	case quantity	case weight	case volume
Transpare	ent polypropyler	ne			
1	409007.N	cap with loop	1,000	0.40	0.003
2	409008.N	cap without loop	1,000	0.40	0.003
Brown po	olypropylene				
3	409007.M	brown cap with loop	1,000	0.40	0.003
4	409008.M	brown cap without loop	1,000	0.40	0.003
Inserts					
5	409111A	blue insert	500	0.58	0.003
6	409111AM	yellow insert	500	0.58	0.003
7	409111B	white insert	500	0.58	0.003
8	409111R	red insert	500	0.58	0.003



Screw cap microtubes, with cap. Sterile

Tubes and caps in medical grade, transparent polypropylene.

The cap embodies a non-reactive ethylene-propylene o-ring.

Suitable for autoclave, liquid nitrogen (gaz) and boiling processes.

Perfect for long term storage.

Withstand temperatures from -190 °C. Centrifugation Resistance: 20,000 xg.


DNAse, RNAse, DNA, and PCR inhibitors free.


Graduated models feature a frosted area for writing.


The non-graduated model incorporates an external grip for an easy handling. Microtubes are supplied capped, in bags of 50 units.

code	description		graduation	case quantity	case weight	case volume
409115/4	0.5 ml skirted	STERILE	×	50 x 50	5.90	0.030
409115/2	1.5 ml non skirted	STERILE	\checkmark	50 x 50	5.90	0.030
409115/6	2 ml skirted	STERILE	\checkmark	50 x 50	5.90	0.030
409115/3	2 ml non skirted	STERILE	\checkmark	50 x 50	5.90	0.030



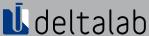

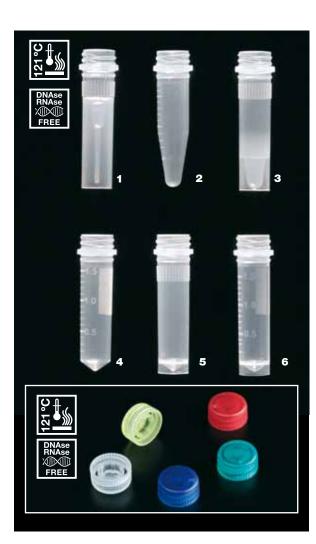

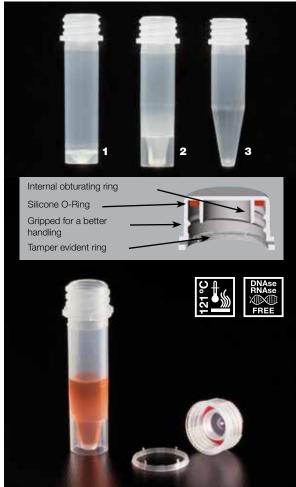

Screw thread microtubes

Made of transparent polypropylene.

Suitable for use in liquid nitrogen, **autoclave** and for boiling applications, and can be used at temperatures down to -190 °C.

Certified RNAse, DNAse and PCR inhibitors free.

Products ideal for long term sample storage.

Withstand centrifugation at **20,000 xg**. The codes **409111/4**, **409111/5** and **409111/6**, with external moulded slots for better handling with gloves.

Dimensions: 10.3 x 44.5 mm (except code **409111/2**: 10.3 x 43.6 mm).

Caps are supplied separately, see below.

mod.	code	description	graduation	case quantity	case weight	case volume
1	409111/4	0.5 ml skirted	×	500	0.70	0.002
2	409111/2	1.5 ml	\checkmark	500	0.60	0.002
3	409111/5	1.5 ml skirted	×	500	0.70	0.002
4	409111/3	2.0 ml	\checkmark	500	0.70	0.002
5	409111/6	2.0 ml skirted	×	500	0.70	0.002
6	409111/7	2.0 ml skirted	\checkmark	500	0.70	0.002

Caps for screw thread microtubes

Made of medical grade polypropylene.

Feature an internal O-ring to ensure leakproof seal.

Dimensions 13.0 x 6.0 mm.

Certified RNAse, DNAse and PCR inhibitors free.

code colour		case quantity	case weight	case volume
409112/0	natural	500	0.024	0.0001
409112/1 blue		500	0.024	0.0001
409112/2	409112/2 green		0.024	0.0001
409112/4	red	500	0.024	0.0001
409112/6	yellow	500	0.024	0.0001







Screw cap until locking ring clicks over serrated tube neck.

Contents are now protected until cap is removed and tamper-evident ring is detached.

Screw thread tamper evident microtubes

Microtubes and caps are made of autoclavable ultra clear polypropylene. Ribbed cap with internal silicone O-ring for a positive **leakproof seal**. Super fast ¼ turn thread design. **Tamper evident** seal which notices if microtube has been opened.

Rnase, Dnase and Pyrogen free.

They can be used at extreme temperatures from -196 °C to 121 °C.

Resistance to centrifugation: 17,000 xg.

Microtube dimensions: 11 x 44 mm.

Cap dimensions (with tamper-evident ring): 15 x 9 mm.

Used on:

- Test of fertility and DNA testing
- · Packaging of diagnostic kits and reagents
- Forensic laboratories

	mod.	code	volume ml	case quantity	case weight	case volume
	1	409110.4T	2	1,000	2.06	0.013
	2	409110.2T	1.5	1,000	2.06	0.013
N=W	3	409110.5T	1.5	1,000	0.65	0.013





Human DNA free certified swabs, steriles

Human DNA free Certified. Sterile.

The swab is supplied in a polypropylene tube, which protects the sample up to the laboratory prior to its analysis. The stick of the swab is made of polystyrene while the head is produced with viscose, polyester or cotton according to the code. The tube is labeled indicating code, description, lot, expiry date and providing an identifying area to note down collection details (site, date, etc.). Moreover, the label seals the tube with the cap of the swab, acting like a tamper evident system.

	code	code description		case quantity	case weight	case volume
Nav	300252DNA polystyrene + viscose		500	4 x 500	14.20	0.066
NEW	300263DNA polystyrene + polyester (dacron® type)		500	4 x 500	14.20	0.066
WEW	300261DNA polystyrene + cotton		500	4 x 500	14.20	0.066

Expiry date: 48 months from sterilisation date





0.2 ml Real Time PCR tubes

Tubes made of polypropylene, featuring attached hinged caps.

Flat caps are easily pierceable and offer optical quality, thus allowing their application in **Real Time PCR.**

Available in strips of 8 tubes (see code **4095.1NP** in the following page). Certified **RNAse, DNAse and PCR inhibitors free.**

code description		case quantity	case weight	case volume
4094.5N	PCR 0.2 ml QPCR tube	1,000	0.25	0.003





0.2 ml PCR tubes

Tubes made of polypropylene, featuring attached hinged caps. Caps are flat and easily pierceable.

See strips of these tubes on the following page (codes **4094.3N** and **4094.4N**). Certified **RNAse, DNAse and PCR inhibitors free.**

code	colour	case quantity	case weight	case volume
4094.1N	.1N natural		0.25	0.003
4094.1A	4094.1A blue		0.25	0.003
4094.1R red		1,000	0.25	0.003
4094.1V green		1,000	0.25	0.003
4094.1AM	yellow	1,000	0.25	0.003

See minimum order and delivery time for other colors.



0.2 ml PCR tubes

Made of polypropylene.

Tubes with attached domed cap.

Certified RNAse, DNAse and PCR inhibitors free.

code		description	colour	case quantity	case weight	case volume
4095.9	4095.9N individual tub		natural	1,000	0.25	0.003



0.2 ml PCR tubes in strips

Made of natural colour polypropylene.

Certified RNAse, DNAse and pyrogen free.

Different models available:


- Strip of 8 or 12 tubes with its corresponding strip of 8 or 12 domed caps.
- Strip of 8 tubes with its attached hinged strip of domed caps.
- Strip of 8 tubes with its attached hinged strip of flat caps, suitable for Real Time PCR.
- Strip of 8 tubes.
- Strip of 8 flat caps, suitable for Real Time PCR.

Cap and tube strips allow easy handling and inventory of tubes and caps. Cap strips prevent cross-contamination from tube to tube.

Code **4095.7N** is also suitable for plates.

	code	description	case quantity	case weight	case volume
1	4095.2N	8 tube strip + 8 domed cap strip	125 strips	0.30	0.0035
2	4095.7N	strip of 8 flat caps	125 strips	0.20	0.0028
3	4095.6N	8 tube strip	125 strips	0.20	0.0036
4	4095.1NP	8 tube strip attached to 8 flat cap strips	125 strips	0.27	0.0036
5	4095.1N 8 tube strip attached to 8 domed cap strips		125 strips	0.27	0.0036
6	4095.4N 12 tube strip attached to 12 domed cap strips		80 strips	0.20	0.0028

Ask for another colours.

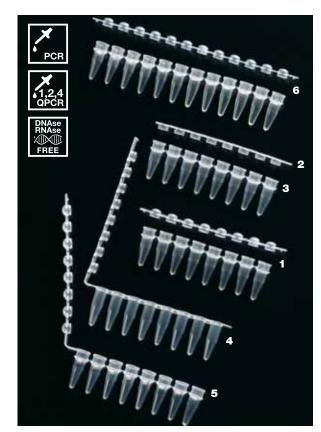

0.2 ml PCR tubes in strips

Tubes made of polypropylene, featuring attached hinged caps.


Suitable for Real Time PCR.

See individual tubes on page 51 (4094.1N). Caps are flat and easily pierceable. Certified RNAse, DNAse and pyrogen free.

In bags of 10 strips.

Available in standard height (code **4094.3N**) and low profile (code **4094.4N**) that minimizes the effects of condensation and it allows working with small samples even less than 20 μ l.

	code	description	case quantity	case weight	case volume
1	4094.3N	strip of 8 tubes and 8 attached caps	120 strips	0.27	0.0036
2	2 4094.4N strip of 8 tubes and 8 attached caps		120 strips	0.20	0.0028

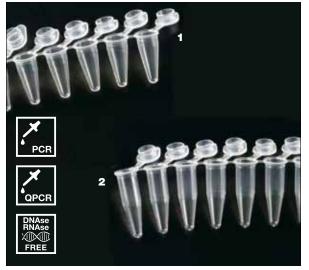

0.2 ml PCR tubes in strips

Made of autoclavable polypropylene.

Certified RNAse, DNAse and pyrogen free.


The strip includes 8 integral 0.2 tubes with ultra thin sidewalls and bottoms with individually attached caps.

Available with either **flat (needle pierceable) or dome-topped individually** attached hinged caps.

Every cap embodies a shield in order to prevent contamination when being opened.

While easily opened and closed with one hand, their positive sealing will fully protect the contents from evaporation during the whole thermal cycle. These strips can be cut to any length while each tube has its own attached cap.

	code	сар	colour	case quantity	case weight	case volume
1	4096.2N	flat cap	natural	125 strips	0,40	0,055
2	4096.3N	domed cap	natural	125 strips	0,40	0,055

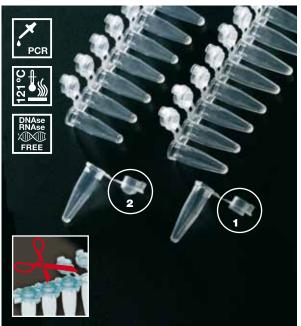



Real Time PCR



Polimerase Chain Reaction is one of the most common techniques in every lab where Molecular Biology is being used. Its applications comprise Diagnostics, Genetic and Prenatal testing, Tissue Typing, Forensics, Pharmacological evaluation, among others.

One of the most useful PCR techniques is the Real Time PCR, also known as Quantitative PCR (QPCR). Its main advantage in front of standard PCR is that it saves time as the quantification is being held during DNA amplification. While standard PCR needs a quantifying process after DNA amplification, QPCR gives results in real time.

PCR standard consumables are mainly made of transparent polypropylene. After some years of experience, laboratory technicians have found that **crosstalking among transparent wells may affect the real time quantification. Opaque white wells and tubes** have proved to avoid this well-to-well crosstalking, thus **ensuring an exact** and **reliable quantification**.

Opace wells also allows luminosity absorption.

The following pages include the newest innovations in this brand new area; like **bi-mould technologies**, which allow the manufacture of an opaque tube attached to an optically clear cap.

White well technology is available in both strips of tubes or also 96 well plates.



0.2 ml, Real Time PCR tubes in strips

Tubes made of polypropylene.

Strip of 8 tubes, each one featuring an attached flat cap.


Strips are manufactured by biomaterial molding, so the tubes are made of opaque white PP, while caps are compounded of transparent, optically clear

Opaque white tubes perform the highest performance in **Real Time PCR**, as they avoid crosstalking between wells.

Specially conceived for **Real Time PCR**. Caps are easily pierceable. **Certified RNAse, DNAse and PCR inhibitors free**.

In bags of 10 strips.

	code	description	case quantity	case weight	case volume
4094.5BP		strip of 8 white 0.2 ml QPCR tubes and 8 attached caps	120	0.26	0.004

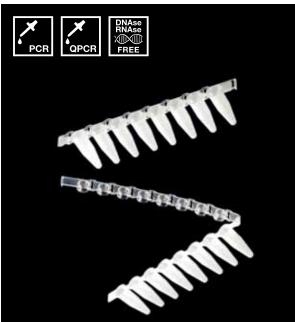

0.2 ml, Real Time PCR tubes in strips

Tubes made of polypropylene.

Strip of 8 tubes, featuring an attached strip of 8 flat caps. Strips are manufactured by biomaterial molding, so the tubes are made of opaque white PP, while caps are compounded of transparent, optically clear PP. Opaque white tubes perform the highest performance in Real Time PCR, as they avoid crosstalking between wells. Specially conceived for **Real Time PCR**. Caps are easily pierceable.

Certified RNAse, DNAse and pyrogen free.

code	description	case quantity	case weight	case volume
4095.1BP	strip of 8 white 0.2 ml QPCR tubes and 8 caps	125	0.23	0.004



Graduated 0.5 ml PCR tubes

Made of autoclavable polypropylene.

Certified RNAse, DNAse and pyrogen free.

Attached hinged caps are flat and easily pierceable.

Tubes are easily opened and closed with one hand.

Tubes have moulded-in graduations in 0.1 ml increments from 0.1 up to 0.6 ml, and a frosted panel on their side for writing or labelling.

code	colour	case quantity	case weight	case volume
4094.2N	natural	1,000	0.38	0.0046
4094.2V	green	1,000	0.38	0.0046
4094.2AM	yellow	1,000	0.38	0.0046
4094.2A	blue	1,000	0.38	0.0046
4094.2R	red	1,000	0.38	0.0046



Graduated 0.5 ml PCR tubes

Same features as the above tubes, but with a domed cap.

code	colour	case quantity	case weight	case volume
4095.5N	natural	1,000	0.46	0.0050



384 wells plate skirted

Made of transparent polypropylene. 384 wells plate.

Each well has a capacity of 50 µl.

A low rim around the top of each well helps to prevent accidental cross-contamination

All wells are thin walled for an excellent thermal transfer.

Orientation cut is at position A 24 (upper right).

A black printed alphanumeric grid helps sample identification.

RNAse, DNAse, DNA and PCR inhibitors free.

Suitable for PCR and real time PCR (QPCR).

Dimensions according to the **SBS** standard.

code	description	case quantity	case weight	case volume
900384	384 PCR well plate	10 x 10	3.20	0.026


Minimun order quantity: 10





PCR Plates – COMPATIB	ILITY CHART	900111 and 900111B	900110 unskirted	900122 semi skirt	900112 semi-skirt	900113B semi skirt	900123 skirted	900384 skirted	9000098 unskirted	900095 white, skirted 900093 black, skirted
Well Shape	Well Shape		9001	9001	9001	9001	9001	3006	0006	0006
	rox. 21 mm, max. capacity approx. 350 µl)	•		•					•	
	pprox. 16 mm, max. capacity approx. 200 µl)		•		•	•	•	50 µl		•
Manufacturer	Model									
Apollo Brand	100									
A 11 1 1 1 1	ATC401	•							•	•
Applied Biosystems	0700									
Thermal Cyclers	2700 9600	•		•				•	•	
	9700	•						•		
	9800 "Fast"				•			_		
	Veriti 0.1 ml									
	Veriti 384							•		
"Real Time" Cyclers	vond de i									
	5700	•		•						
	PRISM 7000	•		•					•	
	7300	•							•	
	7500	•							•	
	7500 "Fast"				•					
	7700	•		•					•	
	Step One plus				•					
Sequencers	PRISM 310	•		•						
	PRISM 3100	•		•						
	3130 (XL)	•		•						
	3700	•		•				•		
	PRISM 3730 (XL)	•		•						
Amersham										
Sequencers	MegaBace 500						•			•
	MegaBace 1000						•			
n 1	MegaBace 4000							•		
Beckman	050									
Sequencers Biometra	CEQ	•								
Thermal Cyclers	Uno	•	•				•		•	
memai Oydieis	Uno II									
	T1 Thermal Cycler						•			
	Tgradient	•								
	Trobot	•	•	•				•	•	
Bio-Rad/MJ Research										
Thermal Cyclers	Gene Cycler	•								
·	PTC-100	•	•	•			•	•	•	•
	PTC-200	•	•	•			•	•	•	•
	PTC-220/221	•	•				•		•	•
	PTC-225 Tetrad	•	•	•			•	•	•	•
	Dyad/Dyad Disciple	•	•	•			•	•		
	iCycler	•		•			•		•	
	MyCycler	•		•					•	
	Mini Gradient	•	•							
	Personal	•		•						
	DNA Engine family	•	•	•			•	•		
<u>"" </u>	C1000/S1000	•	•	•			•	•	•	•
"Real Time" Cyclers	Opticon/Opticon2	•	•				•			-
	Chromo-4		•				•			-
	iCycler	•		•			•	•		
	MyiQ	•		•			•			
	iQ5	•		•			•			
	CFX96						•			
0	CFX384						_	•		+-
Sequencers	BaseStation						•			



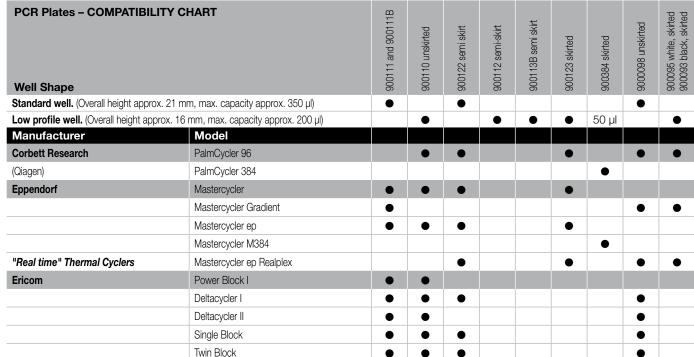







96-well flexible PCR plate

Made of **autoclavable** polypropylene.

It is suitable for both real time and standard PCR thermocyclers.

This 96-well PCR plate is thin-walled for rapid thermal transfer. Flexible, the plate could be easily cut into sections of 24, 32 or 48 tubes.

An alphanumeric grid helps sample identification, and to facilitate orientation, the bottom right corner of the plate is cut away.

Dimensions: 120 x 73 x 20 mm. **Certified RNAse and DNAse free.**

Dimensions according to the SBS standard

code	colour	case quantity	case weight	case volume
900098	natural	10	0.21	0.0024





96-well low profile skirted PCR plate

Made of polypropylene. Each "low profile" well has a capacity of 200 μ l. A low rim around the top of each well helps to prevent accidental cross-contamination and makes easy the sealing with film.

All wells are thin walled for an excellent thermal transfer.

This plate has a skirt approximately 15 mm high that can be filled using automatic fluid handling systems or standard multichannel pipettors.

A black printed alphanumeric grid helps sample identification.


Orientation cut is at position A 12 (upper right).

RNAse, DNAse, DNA and PCR inhibitors free.

Suitable for PCR and real time PCR (QPCR).

Dimensions according to the SBS standard.

code	colour	case quantity	case weight	case volume
900123	natural	10 x 10	2.96	0.030

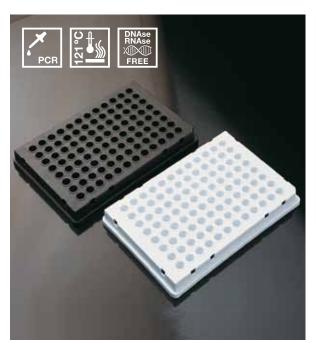

Opaque skirted 96 well PCR plate

Made of **autoclavable** polypropylene.

Alphanumeric identification.

The skirt around the plate provides a labelling area.

Conical-bottomed wells.

The plate can be handled by robotic equipment and is ideal with automated pipetting systems.

To facilitate orientation, the bottom right corner of the plate is cut away.

Well dimensions: 5.5 mm. Well depth: 14 mm.

Code 900093: Black plate for fluorescent PCR.

Code 900095: White plate for luminescence (optical absorption) PCR.

Dimensions: 126 x 84 x 15 mm.

Certified RNAse and DNAse free.

code	description	case quantity	case weight	case volume
900093	96 well black plate	5 x10	1.34	0.010
900095	96 well white plate	5 x 10	1.34	0.010



96 well standard plate

Made of transparent polypropylene.

96 wells standard plate with a capacity of 350 µl each well.

A low rim around the top of each well helps to prevent accidental cross-contamination and makes easy the sealing with film.

All wells are thin walled for an excellent thermal transfer.

A black printed alphanumeric grid helps sample identification.

Orientation cut is at position A 12 (upper right).

RNAse, DNAse, DNA and PCR inhibitors free.

Suitable for PCR and real time PCR (QPCR).

Dimensions according to the **SBS** standard.

code	colour	case quantity	case weight	case volume
900111	natural	10 x 10	2.72	0.030



96-well low profile PCR plate

Made of transparent polypropylene.

Low profile well (15 mm high).

Perfect to work with samples of 100 µl or even less.

A low rim around the top of each well helps to prevent accidental cross-contamination.

The orientation cut is at position H 12 (lower right).


A printed alphanumeric grid helps sample identification.

DNAse, RNAse, DNA and PCR inhibitors free.

It is suitable for both standard or **Real Time PCR**.

Dimensions according to the SBS standard.

code	description	case quantity	case weight	case volume
900110	low profile PCR plate	5 x 20	1.9	0.018

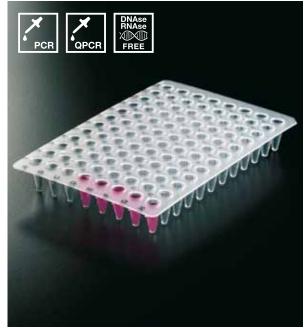

96 well semi-skirt plate

Made of transparent polypropylene.

This plate has a semi skirt ±7.5 mm high.

Each well has a capacity of 200 µl and embodies a low rim around its top which prevents accidental cross-contamination and it makes easy the sealing with foils


Orientation cut is at position A 12 (upper right).

Alphanumeric identification printed in black.

DNAse, RNAse, DNA and PCR inhibitors free.

Dimensions according to the SBS standard.

code	description	case quantity	case weight	case volume
900122	semi-skirt PCR plate	10	0.32	0.0027

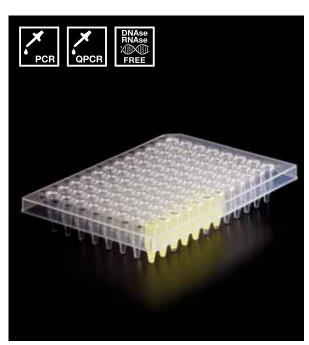





96 wells white standard plate

Made of white polypropylene. 96 wells standard plate.

Each well has a capacity of 350 µl.

A low rim around the top of each well helps to prevent accidental cross-contamination and makes easy the sealing with film.

All wells are thin walled for an excellent thermal transfer.

Its flexible design allows the plate to be easily cut into sections of 24, 32 or 48 tubes.

Orientation cut is at position A 12 (upper right).

A black printed alphanumeric grid helps sample identification.

RNAse, DNAse, DNA and PCR inhibitors free.

Suitable for PCR and real time PCR (QPCR).

Dimensions according to the SBS standard.

code	description	case quantity	case weight	case volume
900111B	white PCR plate	10 x10	2.96	0.030



96 wells low profile, semi-skirted white plate

Specially designed for Roche thermocycler.

Made of white polypropylene. Semi-skirt "low profile" 96 well. Each well has a capacity of 200 μ l.

A low rim around the top of each well helps to prevent accidental cross-contamination and makes easy the sealing with film.

Orientation cut is at position H 12 (bottom right).


A black printed alphanumeric grid helps sample identification.

RNAse, DNAse, DNA and PCR inhibitors free.

Suitable for PCR and real time PCR.

Dimensions according to the SBS standard.

code	description	case quantity	case weight	case volume
900113B	semi-skirt PCR plate	10 x 10	2.96	0.030

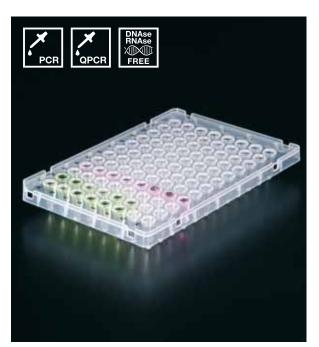

96 well, low profile, semi-skirted plate

Specially designed for ABI "Fast" thermocycler.

Made of transparent polypropylene.

Semi-skirt "low profile" 96 well.

Each well has a capacity of 200 μ l.

A low rim around the top of each well helps to prevent accidental cross-contamination.

Orientation cut is at position A 1 (upper left).

A black printed alphanumeric grid helps sample identification.

RNAse, DNAse, DNA and PCR inhibitors free.

Suitable for PCR and real time PCR (QPCR).

Dimensions according to the SBS standard.

ı	code	description	case quantity	case weight	case volume
	900112	semi-skirt PCR plate	10 x 10	2.96	0.030



PCR sealing mat

Manufactured from a non reactive rubber, this PCR sealing mat ensures a secure, leak free seal during the PCR cycling process.

It is designed to seal **96 well PCR plates**, but may also be cut to fit 24/32/48 well plates.

The mat is marked "this side up" for a proper positioning into the plate (sharp cone side up).

This sealing mat may be sterilised in an **autoclave** or cleaned by immersion in a bleach solution.

The PCR sealing mat can be re-used as often as required.

DNAse, RNAse, DNA and PCR inhibitors free.

code	description	case quantity	case weight	case volume
900305	PCR sealing mat	10 x 5	1.20	0.012



4095.7N



Adhesive sealing film

Sealing film for use with microplates, multiwell plates and microtiter plates. Advantages:

- Minimises the risk of contamination or reagent spillage during ELISA or PCR processes.
- 2. Minimises the risk of contamination from tube to tube and from plate to plate.
- 3. Prevents sample evaporation.

A 5 mm wide strip (opaque white) at the lateral edges of the film helps pull the film from its protective paper and prevents it sticking onto fingers. The film is thermostable and functional from **-70 °C to 95 °C** at 75% humidity.

RNAse and DNAse free. DMSO resistant.

We recomend to use the "roller" 900330 to ensure a perfect seal.

code	description	case quantity	case weight	case volume
900300	adhesive film	1 x 100	0.23	0.001
900330	sealing roller	1	0.014	0.0001

DNASE RNASE FREE 900330

Aluminium sealing foil

This type of material is ideal for manual sealing during PCR work, microtiter plates or manipulation and file plates.

For high throughput applications from **-86 °C to 95 °C**.

Adhesive backing.


Pierceable with a pipette tip for easy access to sample.

DMSO resistant.

It is recommended to use the "roller" code **900330** for ensuring a perfect sealing, eliminating the danger of evaporation.


RNAsa, DNAsa and DNA free.

code	description	case quantity	case weight	case volume
900320	aluminium foil	1 x 100	0.28	0.001
900330	sealing roller	1	0.014	0.0001

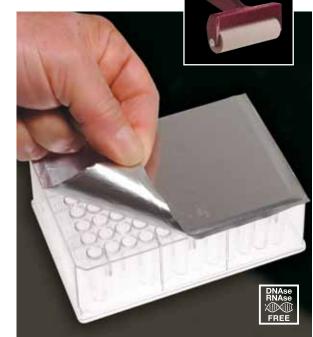





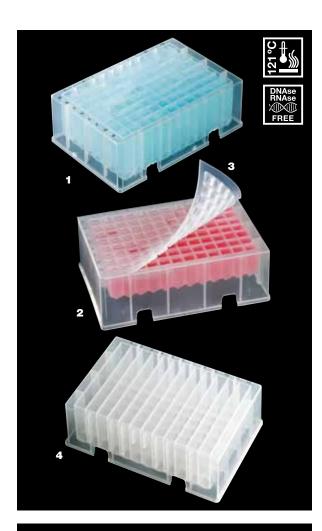

Deep well plate (96 square wells or 12 rectangular channels)

Made of medical grade polypropylene.

A printed alphanumeric grid helps sample identification, and to facilitate orientation the bottom left corner of the plate is cut away. Conical bottom wells.

Widely used for liquid handling and storage, whether with manual pipetting or robotic handling.

Withstand temperatures from -90 °C to 121 °C (Autoclave).

Dimensions according to SBS.

Certified free from detectable RNase, DNase, DNA and PCR inhibitors.

Codes 900195 and 900196 are supplied in 2 bags of 5 units each.

Code **900309** is a sealing mat manufactured from autoclavable silicone. Exclusive use with the square well plates (**900195** and **900196**).

Pierceable so the sample can be aspirated without removing the mat. Supplied in a bag of 5 units.

Code **900197** proves to be extremely convenient as a reservoir thanks to its 12 high volume channels. Supplied in a bag of 5 units.

mod.	code	description	case quantity	case weight	case volume
1	900195	96 well plate 2.2 ml	2 x 5	1.18	0.009
2	900196	96 well plate 1.1 ml	2 x 5	0.95	0.007
3	900309	sealing mat for 96 squared wells	1 x 5	0.16	0.001
4	900197	12 channel plate 22 ml	1 x 5	0.45	0.005



96 deep well plate (96 round wells)

Made of medical grade polypropylene, it is comprised of 96 round bottom wells, each one of 1.2 ml capacity. Printed alphanumeric grid helps sample identification, and to facilitate orientation the bottom right corner of the plate is cut away.

Widely used for liquid handling and storage, whether with manual pipetting or robotic handling.

It with stands temperatures from **–90 °C to 121 °C** (autoclave).

Dimensions according to SBS.

Certified free from detectable RNase, DNase, DNA and PCR inhibitors. It can be sealed with the plastic sealing film **900300** or the pierceable aluminum foil **900320**.

code	description	case quantity	case weight	case volume
900155	96 well plate 1.2 ml	10	1.41	0.009

96 deep well plate (square well, round bottom)

Made of polypropylene. The 2.1 ml well capacity (2.0 when capped) plate is used mainly for compound storage and enzyme assays. An alphanumeric grid helps sample identification, and to facilitate orientation the bottom right corner of the plate is cut away.

DMSO resistant. Can be centrifuged up to 6000 xg.

Withstand temperatures up to -150 °C. Autoclavable.

Each well measures 8.3 x 8.3 mm.

Plate dimensions: 41.6 mm high x 127.8 mm long x 85.5 mm wide.

We recommend to protect the samples with the sealing mat code **900306** (DMSO resistant).

cod	le	description	case quantity	case weight	case volume
9001	70	96 deep well plate, natural colour	4	0.47	0.002
9003	306	sealing mat	24	0.15	0.002





Sample storage system

96 round bottom tubes, supplied in twelve strips of eight tubes each, held in a rack with lid. Rack, lid, tubes and caps are made of **autoclavable** polypropylene, manufactured with the SBS standard footprint. Compatible with all robotic applications and multichannel pipetting systems. Alphanumeric numbered wells, indelibly printed in black, allowing identification in short light conditions. **DNAse, RNAse, DNA and PCR inhibitors free**.


Used for serial dilutions, mixing, storage and harvesting of cells, cell growth for cell culture assays and DNA screening, and as an ideal long-term storage system.

Caps are sold separately in strips of eight units; every cap has a tag to make an easy opening and closing.

Tube dimensions (height): 8 x 44 mm

Rack dimensions (with lid): 128 x 86 x 48 mm

code	description	case quantity	case weight	case volume
409009	rack with 96 tubes (8 x 12)	10	1.63	0.0140
409010	strips of 8 tubes each	125 strips	0.27	0.0073
409011	strips of 8 caps each	125 strips	0.11	0.0012
409012	loose tubes	1,000	0.63	0.0073

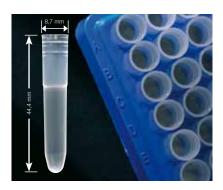

Sample storage system

Compact sample storing system. It consists on a blue rack with a translucent lid, holding 96 loose round bottom tubes (1,2 ml) arranged in 8 x 12.

Autoclavable and stackable, it resists up to **-100 °C**, and embodies a moulded alphanumeric identification.


Rack, lid and tubes made of polypropylene. Caps made of low density polyethylene.

Caps are sold apart in strips of eight units; every cap has a tag to make an easy opening and closing.

Tube dimensions (height, closed): 8,7 x 45 mm. Rack dimensions (with lid): 118 x 82 x 50 mm.

DNAse, RNAse free.

code	description	case quantity	case weight	case volume
409008	rack with 96 tubes (8 x 12)	10	1.24	0.0100
409005	strips of 8 caps each	120 strips	0.01	0.0001
409008.1	tubes in bulk	1,000	1.00	0.0500

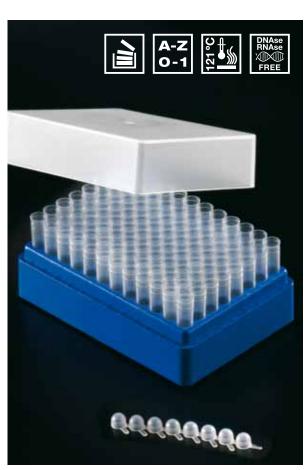







It consists on a white rack and a transparent lid, with 96 tubes (12 x 8) 1.2 ml (capped, 1.1 ml).

Tubes and rack are manufactured in polypropylene, being autoclavable. Caps in non autoclavable polyethylene.

Caps are acquired apart, in strips on eight units.

It is ideal to work with robots and multichannel pipetting systems, as well as for sample transport, storage, or freezing (it can withstand up to **-80 °C**).

It stands up the majority of chemical agents.

Both lid and rack embody an alphanumeric identification.


The top left edge of the lid is cut for an exact orientation.


Tube dimensions (capped): 9 x 48 mm.

Rack dimensions (with lid): 126 x 81 x 53 mm.

According to SBS standard.

code	description	case quantity	case weight	case volume
409004	rack with 96 tubes (8 x 12)	10	1.81	0.0080
408002	tubes of 8 caps each	120 strips	0.60	0.0040
408005	strips of 8 caps each	120 strips	0.01	0.0001
408003	loose tubes	5 x 960	0.80	0.0200

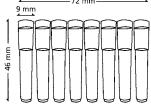

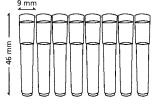

Storage rack with 2 ml tubes

Compatible with most robotic Workstation, this polypropylene storage rack can be used with most cell harvesters and leading 8 – & 12 – chanel pipettors. It contains 96 removable polypropylene square tubes in a 8 x 12 configuration, each having a 2.1 ml capacity.

Although the tubes are square, the bottom is round to facilitate emptying. Tubes and rack are **autoclavable** and they are ideal for storage of blood and other biological samples at temperatures from -30°C up to 70°C.

Racks are stackable to save storage space. According to SBS standard.

code	description	case quantity	case weight	case volume
418003	96 well storage rack with tubes	10	1.81	0.008









100% Traceability

Our blood collection tubes are individually identified with its code, LOT number and expiry date, so the traceability from its manufacture to the final consumer is guaranteed.

Raw materials traceability.

Process traceability.

Final product traceability.

We manufacture our blood collection tubes by following:

Directive 98/79/EC "In vitro" diagnostic medical devices.

- ISO 6710 Single use containers for blood collection
- UNE-EN ISO 14971 for the Risks Management of Sanitary Products
- **UNE-EN 980** for Labelling and Marking
- UNE-EN ISO 13485, medical devices quality management systems





Expiry date our tubes, anticoagulants and other

•	Serum tubes (granules and gel):	24 months
•	Serotub with granules and preservative:	12 months
•	Lithium heparin:	24 months
•	lodoacetate heparin:	13 months
•	Edta:	24 months
•	Citrate for coagulation:	15 months
•	Anticoagulants in containers:	24 months
•	Thrombocyte count and osmotic brittleness:	24 months
•	Reticulocyte vital staining:	24 months

Procedures and controls

During the manufacturing process a lot of controls and guides are used, insisting on:

- Reproducibility dosing control
- Quantity and quality dosing control (flame photometer)
- Watertightness control (vacuum, centrifugation and fluorescein method)

Each case includes a brochure with recommendations of use.



0.8 thick and 19 mm wide. Non toxic, USP, Class VI, silicone. Autoclavable.


code	description	quantity	weight	volume
GS-01 roll of 50 meters		1 x 50	0.85	0.0026
GS-02	box with 10 bands 0.5 m long each	1	0.10	0.001









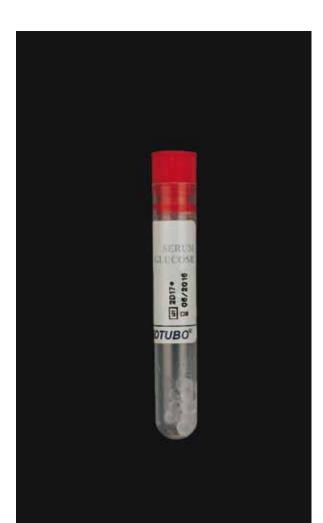

Serum glucose serotub

Made of a clear polypropylene, supplied capped and labeled with an inert polyethylene cap.

To permit a fast serum separation of the blood clot, the tubes have an inert additive inside which speeds up the coagulation. This accelerator allows a mixture with the blood almost immediate and avoids the unnecessary sample

The coagulation speed is much higher than the obtained with the glass tubes or the other existing tubes on the market.

The inert granules are rounded to avoid erythrocytes lesions during centrifugation and minimize the risk of haemolysis. These granules are located between the clot and the separated serum working as a retaining wall. It assures a comfortable pipette action or serum decanting. (It isn't a watertight barrier). Inside the tube there is also a little amount of a glucose preservative, lithium iodoacetat, which maintains the sample stable up to 6 days and allows to make most of the biochemical measurements with an only tube (the exceptions are CPK and Lithium); so it's possible to save one of the two tubes used normally to determine routine tests and glucose respectively.

With the ergonomic features of the cap the design achieves a high reliability without avoiding the main characteristic of the cap: it's very easy to use.

On the label of each tube the lot number, expiry date and volume are printed.

Presentation: packed in black bags in order to maintain them out of the light because the iodum is photosensitive.

Keep these tubes at a temperature below 35 °C.

Supplied in plastic racks.

	code	description	case quantity	case weight	case volume	pallet quantity
	600602	13 x 75 for 4 ml blood	12 x 100	4,92	0.033	42

Ask for minimum quantity and delivery time for tube 16 x 100, 10 ml blood.



Citrate for coagulation

(€ (IVD)

Made of clear polypropylene, supplied capped and labeled.

Our sodium citrate, 3.8% or 3.2% concentration, buffered and sterile has a ratio citrate: blood 1:9 and is highly recommended for coagulation tests. According to the prevailing rules, this liquid anticoagulant permits to determine the prothrombin time (Quick) up to 12 hours after sample collection.

The cap is not only easy-to-use but also assures a watertight closing; allowing a comfortable and reliable capping.

A volume indication mark, lot number and expiry date on the label of each tube, allows to assure the traceability of the product.

Tubes are supplied in plastic racks.

code	type of coagulant	description	case quantity	case weight	case volume	pallet quantity
601102	buffered to 3.8%	13 x 75 round for 4 ml blood	12 x 100	5.30	0.033	42
601103	buffered to 3.8%	13 x 75 round for 2.5 ml blood	12 x 100	5.30	0.033	42
601203	buffered to 3.2%	13 x 75 round for 2.5 ml blood	12 x 100	5.30	0.033	42





Serum tubes with clot accelerator and gel serum separator

C€ (IVD)

Made of a clear polypropylene, supplied capped and labelled with an inert polyethylene cap.

The tubes are designed and processed to permit a fast serum and blood clot separation.

Each tube contains a special inert additive which accelerates the coagulation, and the result is a **fast clot retraction**.

The coagulation speed is, doubtlessly, much higher than the obtained with the glass tubes or the other existing tubes on the market.

The inert gel is located, after centrifugation, between the clot and the obtained serum and it works as a totally watertight barrier.


The ergonomic features of the cap design results in a very reliable and easy-to-use cap.

On the label of each tube the lot number, expiry date and volume are printed.

Supplied in plastic racks.

We strongly recommend its usage for biochemistry, routine tests, special biochemistry, markers, hormones, immunology (tube without anticoagulants).

code	description	case quantity	case weight	case volume	pallet quantity
600801	13 x 75 round for 4 ml	12 x 100	5.80	0.033	48
600800	15 x 100 round for 9 ml	6 x 120	6.00	0.040	36

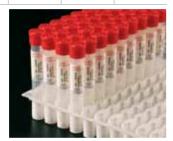



Serum tubes with clot accelerator and gel serum separation


Clotting speed increased: with our the tubes processed with gel serum separator and accelerator we improve the clotting speed achieving a 20-25% advantage over glass tubes.

The serum is obtained after 12 minutes, depending on the particular working conditions.


Serum volume increased: The volume obtained is higher than the obtained with granules tubes or glass tubes, achieving an advantage around 20%. The gel serum separator tube is the preferred option when analysis requires maximum serum yield from the blood sample.

Recommended in case the objective is to obtain larger serum volume.

Note: Before using, the gel serum separator remains at the bottom, even when the tube is bent.

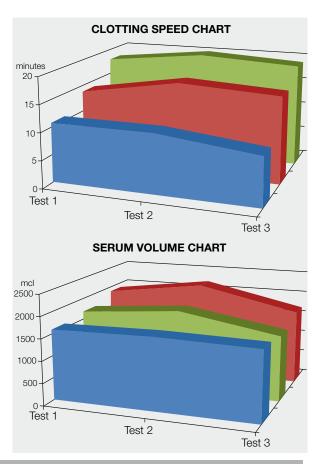



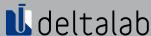

(€ (IVD)



Serum tubes with clot accelerator and granule serum separator

Made of a clear polypropylene, supplied capped and labeled with an inert polyethylene cap.

The tubes are designed and processed to permit a **fast serum and blood clot separation.**

Each tube contains a special inert additive which accelerates the coagulation, and the result is a fast clot retraction. The coagulation speed is much higher than the obtained with the glass tubes or the other existing tubes on the market.


Allows the obtention of serum in 12 min.

The inert **granules are rounded** to avoid cellular lesions during centrifugation and minimize the risk of haemolysis. These granules are located between the clot and the separated serum working as a retaining wall. It assures a comfortable pipette action or serum decanting.(It isn't a watertight barrier). The ergonomic features of the cap design results in a very reliable and easy-to-use cap.

On the label of each tube the lot number, expiry date and volume are printed.

code	description	case quantity	case weight	case volume	pallet quantity
600400	13 x 75 round for 4 ml	2 x 1,000	6.90	0.045	32
600300	16 x 100 round for 9 ml	2 x 500	6.10	0.045	32
707094*	bottle of 750 g with separator granules	20	17.50	0.045	40

^{*} Not **C€** Product

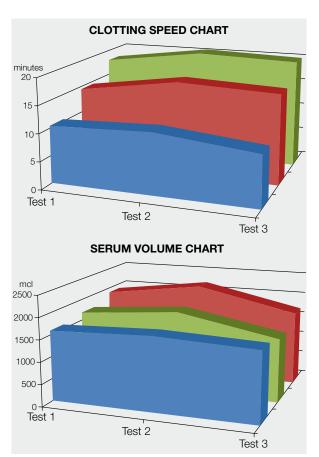


Serum tubes with clot accelerator and granules

Clotting speed increased: with our the tubes processed with granules and accelerator we achieve a clotting speed higher than the obtained not only with glass tubes but also with gel tubes. The advantage achieved is around 40-50% over glass tubes.

Serum can be obtained starting from 12 minutes, depending on working conditions.

Serum volume increased: The volume obtained is a little higher than the obtained with glass tubes but fewer than the obtained with gel serum separator tubes.

Due to its effectiveness and low price this product is highly recommendable for all routine serum tests.

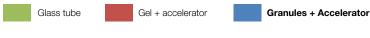

C€ (IVD)

Made of a clear polypropylene, supplied capped and labeled, indicating **filling line**, **lot number**, **and expiry date**.

The anticoagulant pulverization inside the tube optimizes the mixture and avoids the unnecessary blood dilution.

The anticoagulant mechanism is the inhibition of the thrombin action.

Supplied in plastic racks.

Type of coagulant: **spray**.

code	description	case quantity	case weight	case volume	pallet quantity
601802	13 x 75 round for 4 ml blood	12 x 100	4.10	0.033	42
601803	13 x 75 round for 2.5 ml blood	12 x 100	4.10	0.033	42
601810	16 x 100 round or 10 ml blood	6 x 120	4.46	0.038	30



Iodoacetate lithium + heparin lithium tubes

C€ (IVD)

Made of a clear polypropylene, supplied capped and labeled. On the label of each tube a fill line indicates the level of blood required, as well is printed the lot number and expiry date.

The anticoagulant and preservative pulverization inside the tube optimizes the mixture and avoids the unnecessary blood dilution.

The blend anticoagulant-glucose preservative is ideal for biochemical tests, so it's possible to determine most of the biochemical parameters with only one tube and preserve the product stable for 4 days. It is recommended to maintain the tubes in a dark and cool place (at room temperature).

The tubes are packed in black bags in order to maintain them out of the light because the iodum is photosensitive.

Supplied in plastic racks.

Type of coagulant: liquid.

code	description	case quantity	case weight	case volume	pallet quantity
602002	13 x 75 round for 4 ml blood	12 x 100	4.40	0.033	42
602003	13 x 75 round for	12 x 100	4.30	0.033	42







Edta: tripotassium

Made of clear polypropylene, supplied capped and labeled.

Tetraceticetilendiamin acid, a tripotassium salt, works as an anticoagulant thanks to its capacity to fix the blood calcium. Because the anticoagulant is pulverized, it allows a mixture with the blood almost immediate.

The quantity of additive is very small so there are no dilution mistakes (this may occur in the tubes with big additive liquid solution volumes). There is no risk of anticoagulant loss when uncapping because it is adhered to the tube walls.

The cap shape, both internal and external design, assures a comfortable and reliable capping.

A volume indication mark, lot number and expiry date on the label of each tube, allows to assure the traceability of the product.

Supplied in plastic racks.

Type of coagulant: spray.

code	description	case quantity	case weight	case volume	pallet quantity
601612	15 x 50 fat for 4 ml blood	10 x 120	5.10	0.035	36
601613	15 x 50 flat for 2.5 ml blood	10 x 120	4.50	0.035	36
601603	13 x 75 round for 2.5 ml blood	12 x 100	4.10	0.033	42
601702	13 x 75 round for 4 ml blood	12 x 100	4.10	0.033	42



Edta: tripotassium. Rubber cap

(€ (IVD)

Tubes made of clear polypropylene, supplied capped and labeled. Mauve, pierceable (but not pierced) and plugged cap made of thermoplastic rubber. Suitable for hematological automatic machines.

The label indicates code, volume, lot number and expiry date, assuring total traceability.

Tubes supplied in trays of 100 units.

Type of coagulant: liquid.

mod.	code	description	case quantity	case weight	case volume	pallet quantity
1	611603	13 x 75 for 3 ml blood	10 x 100	5.30	0.034	48
2	611604	13 x 80 for 3 ml blood	8 x 100	3.00	0.030	50





Blood-plasma. Edta tubes: dipotassium

(€ (IVD)

Made of clear polypropylene, supplied capped and labelled.

Tetraaceticetilendiamin acid, a dipotassium salt, works as an anticoagulant thanks to its capacity to fix the blood calcium.


Because the anticoagulant is pulverized, it allows a mixture with the blood almost immediate. The quantity of additive is very small so there are no dilution mistakes (this may occur in the tubes with big additive liquid solution volumes). There is no risk of anticoagulant loss when uncapping because it is adhered to the tube walls. The cap shape, both internal and external shape, assures a comfortable and reliable capping.

A volume indication mark, lot number and expiry date on the label of each tube, allows to assure the traceability of the product.

Supplied in plastic racks.

type of coagulant: spray .					
code	description	case quantity	case weight	case volume	pallet quantity
601402	13 x 75 round for 4 ml blood	12 x 100	4.30	0.033	42
601403	13 x 75 round for 2.5 ml blood	12 x 100	3.90	0.033	42
601412	15 x 50 flat for 4 ml blood	10 x 120	4.98	0.033	36
601413	15 x 50 flat for 2.5 ml blood	10 x 120	4.26	0.033	36



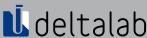

(€ (IVD)

Capped and labeled transparent polypropylene tubes.

Filled with 3.8% sodium citrate (anticoagulant).

The citrate: blood ratio (according to the Westergren method) is 1:4, so the tubes contain 0.4 ml of stable 3.8% sodium citrate solution.

The printed fill line indicates 2.0 ml so 1.6 ml of blood will be added.

The cap design results in a comfortable and reliable handling, specially because it is possible working with gloves without slipping.

Volume fill line, expiry date and batch number are printed on the label; so the product traceability is guaranteed.

Supplied in 100 units plastic racks.


code	description	case quantity	case weight	case volume
601006	round tube 13 x 75 mm for 2 ml (1.6 ml of blood)	12 x 100	4.90	0.033

Pallet quantity: 42 cases.





Comparation of the ESR reading obtained with the Micro system and standard type

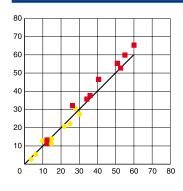

ESR reading with MICRO system (mm)

 Micro system reading at 50 minutes and with standard at 60 minutes.


Micro system reading at 100 minutes and with standard at 120 minutes.

Sedirate micro system


(€ (IVD)

System for the determination of the erythrocyte sedimentation rate. Specially recommended for **Paediatrics**. The system consists of a tube and a pipette. The tube includes a pierceable and re-sealable stopper (same model as the stoppers presented on all our tubes on next page) filled with 0.08 ml of trisodic citrate 0.106 M for 0.32 ml of blood according to standards of Westergren standard method. The pipette of pressure filling has an inner diameter of 1.25 mm and is graduated. Once blood and citrate are mixed together, introduce the pipette into the tube (without removing the stopper). The blood will automatically reach the 0 level. The results obtained are comparable to those obtained with the standard method (macro).

code	description	case quantity	case weight	case volume
27	set pipette + tube	400	3.28	0.029

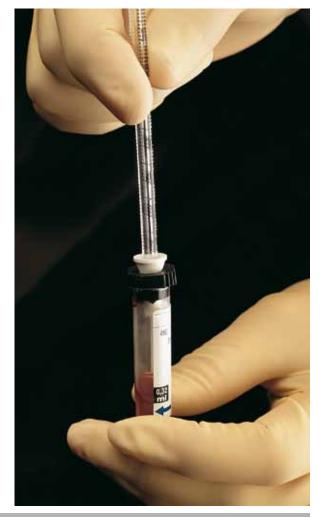



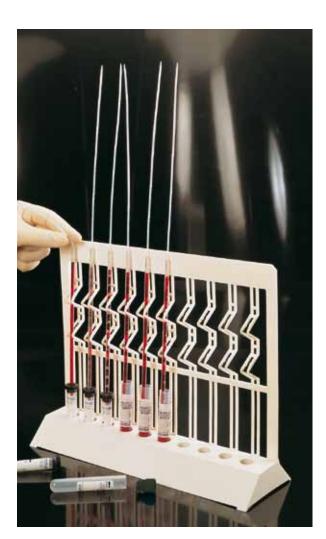

E.S.R. semi micro system TAKIVES with self-levelling system

C€ (IVD)

Pipette graduated from 0 to 160 mm, with a total length of 200 mm and an inner diameter of 2.5 mm, according to the Westergren Method.

The system accepts a total volume of 1 ml.

The plunger must be drawn up manually up to a limit inside the pipette to ensure an adequate volume of blood-citrate blend.

Two tubes available, both made of high transparency polypropylene:

Code **1164**: for 1 ml of total volume (0.75 ml of blood); with a pierceable rubber cap that allows inserting the pipette without taking off the cap.

Code **601006**: for 2 ml of total volume (1.6 ml of blood); with a polyethylene cap that shall be taken off before inserting the pipette.

code	description	case quantity	case weight	case volume
1360	graduated pipette	2 x 500	3.35	0.008
1361*	stand for 10 pipettes	5	3.20	0.023
601006	tube 13 x 75 mm with sterile citrate for 2 ml (1.6 ml of blood)	12 x 100	4.90	0.033
1164	tube 12 x 55 mm with citrate with pierceable cap	10 x 100	3.00	0.020

^{*} Non **(€** product



E.S.R. with self-levelling system

C€ (IVD)

Consists of a polystyrene pipette with self-filling system using a plunger suitable for 12 mm or 13 mm \varnothing tubes.

Graduated from 0 to 180 mm. 1.25 ml blood-citrate mixture is enough for determination.

Westergren method.

code	description	case quantity	case weight	case volume
29	E.S.R. pipette, 230 mm long graduated up to 180 mm	3 x 200	3.00	0.028
601006	tube 13 x 75 mm with sterile citrate for 2 ml (1.6 ml of blood)	12 x 100	4.90	0.033

The E.S.R. system consists of a holder and a set of pipettes with a polyethylene filter

Pipettes may be filled using an automatic suction pump, a pippetting bulb or pump, or a hand hold pippeter (both manual or electronic).

The holder accommodates up to 10 pipettes. Plastic pipettes are graduated up to 180 mm. Tubes are filled with sterile citrate (0.106 M). Westergren method.

code	description	case quantity	case weight	case volume
28	pipette graduated up to 180 mm, length 230 mm	600	1.72	0.0039
1361*	holder to 10 places	5	3.20	0.023
601006	tube 13 x 75 mm with sterile citrate for 2 ml (1.6 ml of blood)	12 x 100	4.90	0.033

^{*} Non **(€** product



Mailing container

Mailing container with safety screw blue cap, both manufactured in polyethylene.

Leak proof.

Dimensions uncapped 128 x 30 mm. Ideal for 10 ml tubes.

The tube includes a piece of absorbent paper to prevent any leakage. Container and cap are sold separately.

code	description	case quantity	case weight	case volume
401301	mailing container	500	7.70	0.090
401302	screw cap	500	2.00	0.095



Blood group test plate

Made of ultra clear polystyrene. Features 10 numbered wells.


Stackable.

Dimensions: 160 x 40 x 6 mm.

Each case contains 50 bags of 10 units each.

code	description	case quantity	case weight	case volume
128030	blood group plate	25 x 10	3.70	0.014



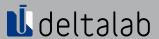



Blood collection tubes with pierceable and re-sealable rubber caps

Manufactured from transparent polypropylene. Robust and resistant to breakage. Tubes feature a unique flexible rubber cap which has two cross cuts and is specially designed to reseal automatically after introduction or withdrawal of the sample. This system eliminates the need to remove the stopper when introducing or withdrawing blood samples. Tubes can be used with most manual, semi-automatic and automatic blood sampling machines, and are suitable for transportation by intra hospital pneumatic systems.



Blood collection tubes

(€ (IVD)

code	type	description	expiration months	case quantity	case weight	case volume
620300	serum separator	16 x 100 with granules for 10 ml	36	6 x 100	4.20	0.034
620400	serum separator	13 x 80 with granules for 4 ml	36	8 x 100	3.20	0.030
621611	edta tripotassium	16 x 55 skirted for 2.0 ml	24	10 x 100	3.50	0.034
621613	edta tripotassium	13 x 80 for 2.5 ml	24	8 x 100	2.90	0.030
621102	citrate for coagulation	12 x 70 for 2 ml	15	8 x 100	2.40	0.030

Pallet quantity: 50 cases.



Special pediatrics

(€ (IVD)

code	type	description	expiration months	case quantity	case weight	case volume
620200	serum separator	12 x 55 with granules for 2 ml	36	10 x 100	3.10	0.025
621610	edta tripotassium	12 x 55 for 1 ml	18	10 x 100	2.75	0.025
621101	citrate for coagulation	12 x 55 for 1 ml	15	10 x 100	2.80	0.025

Pallet quantity: 50 cases.

See speed sedimentation citrate for pediatrics in the page 74





PROBLEM		POSSIBLE CAUSES	REASON FOR MISUSE	SOLUTION
		Incomplete or inhibited coagulation. Non-homogeneous blend of blood and coagulant	The tube was not inverted for 5 times	After tube is filled with blood and capped, gently invert the tube for 5 times
	Poor or incomplete serum separation.	Early centrifugation, insufficie nt allowance coagulation time prior to centrifugation	The recommended 30 min. standing time for coagulation was not observed	Wait for 30 min before centrifugation
The barrier is not clearly defined	Barrier is intact and well defined	Centrifugation xg-setting under appropriate value	The tube was not centrifuged at or above the minimum, appropriate xg-setting	a) Set the centrifuge between 1,000-1,500 xgb) Confirm centrifuge is correctly calibrated
20		Less centrifugation time than recommended	The tube was not centrifuged for the required time	Centrifuge the tube at the appropriate xg-force for at least 10 min
	Barrier is skewed and poorly formed, separation is incomplete	Refrigerated centrifuge	The temperature of the centrifuge was lower than the recommended	 a) Keep the centrifuge at a temperature between 24 °C and 26 °C (around 77 °F). b) Thermally insulate the tube and avoid contact with the centrifuge and centrifuge rotor
		xg-force centrifugation exceeded the tube construction	The tube was subjected to xg-forces exceeding 1,500 xg	Maintain centrifuge at or below 1,500 xg
	fracture during fugation	Centrifuge not positioned level and/or foreign matter or debris inside the rotor of the centrifuge	a) Shock absorbers or tube stand-off cushions absent. b) The inside of the rotor is not clean	Install shock absorbers or tube stand-off cushions as required. Remove any foreign agent or debris that may be inside the centrifuge
		Incomplete or inhibited coagulation. Non-homogeneous blend of blood and coagulant	The tube was not inverted for 5 times	Invert the tube for 5 times
	served within the ion barrier	Insufficient time before centrifugation	The recommended 30 min standing time for coagulation was not observed	Wait for 30 minutes before centrifugation
		Excessive xg-force centrifugation	The tube was centrifuged at above 1,500 xg	Set the centrifuge at a speed below 1,500 xg
Processes of	fibrin in serum	Incomplete or inhibited coagulation. Non-homogeneous blend of blood and coagulant	The tube was not inverted for 5 times	Invert the tube for 5 times
Presence of	iibiii iii seruiii	Early centrifugation, insufficient allowance coagulation time prior to centrifugation	The recommended 30 min standing time for coagulation was not observed	Wait for 30 minutes before centrifugation
		Poor serum quality due to haemolysis	Blood subject to aggressive agitation or rough handling	Handle tube with moderation. Invert gently. Never shake severely or agitate the tube
AbilOffilal	assay results	Excessive xg-force centrifugation	The tube was centrifuged above 1,500 xg	Set the centrifuge to a maximum of 1,500 xg xg-force







Anticoagulants and preservatives in containers

Available in bottles of 15 ml.

The dosage of one drop (15 ml = 300 drops) is enough for 5 ml of blood. Both products have preservatives for stability.

code	description			
705000	705000 lithium heparin 15 ml			
705010	tripotassium edta 15 ml	10		





Special techniques

Liquid for thrombocyte count determination:

Due to its optical characteristics it identifies the thrombocyte avoiding confusion with other cells. This reagent also prevents adhesion and aggregation of the thrombocytes.

Full instructions are included with the kit.

Test of osmotic brittleness of the erythrocytes:

The test for the erythrocyte osmotic fragility detects the resistance of these cells to haemolysis, in hypotonical solutions with decreasing concentration of sodium chloride.

This set contains 2 complete tests and each one is composed of 12 tubes with stable and buffered solutions. Full instructions are included in the kit.

code	description	case quantity	case weight	case volume
800000	thrombocyte counting kit 1 x 50 tubes	30 kits	5.50	0.046
802000	osmotic brittleness of erythrocytes kit 2 x 12 plus 2 lithium heparine tubes	30 kits	7.50	0.046

Minimum order quantity: 1 kit.



Reticulocyte staining kit

(€ (IVD)

This simple to use kit consists of a tube containing 100 μ I of stable buffered bright cresil blue stain which allows the determination of the erythrocyte count. Two to three drops of blood are added directly to the tube and incubated for 10 minutes at room temperature.

The erythrocytes become a pale blue colour making them easy to identify. Full instructions are included with the kit.

Tube made of polypropylene and cap of mauve polyethylene.

code	description	case quantity	case weight	case volume
801000	reticulocyte count kit (1 x 50 tubes)	30 kits	5.50	0.046

Minimum order quantity: 1 kit.



Haematological stains

(€ (IVD)

code	description	case quantity	case weight	case volume
808000	Eosin, solution, May Grunwald, 250 ml bottle Danger symbols: R11-23/24/25- 39/23/24/25. S 7-16-36/37-45	16	4.30	0.015
808001	Eosin, solution, May Grunwald, 1000 ml bottle Danger symbols: R11-23/24/25- 39/23/24/25. S 7-16-36/37-45	12	13.00	0.046
808100	Eosin, solution modified Giemsa, 250 ml bottle Danger symbols: R11-23/24/25- 39/23/24/25. S 7-16-36/37-45	16	5.15	0.015
808101	Eosin, solution modified Giemsa,1000 ml bottle Danger symbols: R11-23/24/25- 39/23/24/25. S 7-16-36/37-45	12	13.00	0.046
808200	Eosin, solution Wright, 250 ml bottle Danger symbols: R11-23/24/25- 39/23/24/25. S 7-16-36/37-45	16	4.40	0.015
805013	fast stain for blood extensions, 250 ml containers, Kit with 2 containers A and 2 containers stain B	12 Kits	20.00	0.046



Please see more information (page 99 about these stains at Chapter 4. **Histology, Microscopy and Staining**.



Microhaematocrit capillary tubes

(€ (IVD)

Soda neutral glass tubes with a colour-coded print for an easier identification, with sodium heparin (red) or without heparin (blue).

Supplied inside a glass tube with a plastic cap, specifying code, lot and expiry date.

Two lengths available: 70 or 75 mm.

Internal diameter: 1.15 mm, external diameter: 1.55 mm.

They meet ISO 12772.

code	description	volume µl	case quantity	case weight	case volume
7301	75 mm with heparin	75	10 x 100	0.36	0.001
7311	75 mm without heparin	75	10 x 100	0.36	0.001
7401	70 mm with heparin	70	10 x 100	0.36	0.001
7411	70 mm without heparin	70	10 x 100	0.36	0.001




Capillary tubes sealing wax

Vinyl plastic wax, in a numbered (from 1 to 24) plastic holder.

Suitable for all capillary tubes, specially for those that have to be centrifuged, like microhaematocrit capillary tubes.

Must be preserved in temperatures from 8 °C to 30 °C.

1	code	description	case quantity	case weight	case volume
	7600	sealing wax	3	0.050	0.001

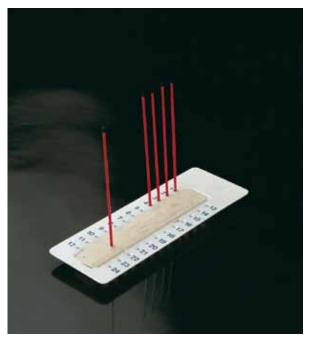









4 5 6 7

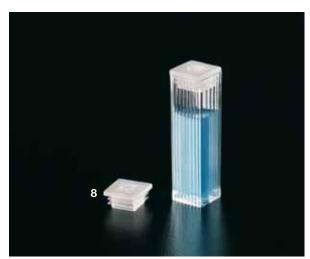

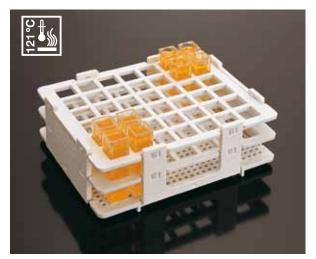

Spectrophotometer cuvettes

Disposable cuvettes suitable for most of the open spectrophotometers. Homogeneous measures, specially of the surface crossed by the light beam, assuring an optimum transmission level on the whole visible spectral.

The material used avoids any possible measurement interference. Due to the strict quality control during the manufacture process, a high reliability is assured. The maximum absorption variations are $\pm\,1\%$.

The two sides not crossed by the light beam are ribbed to an easy identification of the cuvette position inside the spectrophotometer measurement chamber, resulting in an easy positioning and removal.

They are supplied in a dust proof, expandable polystyrene box (100 units per box) with lid.

Dimensions: 12.55 x 12.65 x 44.55 mm (± 0,1 mm).

Light path: 10 mm.

Standard cuvettes

Made of polystyrene for assays in the visible spectral range (340 to 800).

	code	description	case quantity	case weight	case volume
1	302000	4.5 ml macro	5 x 100	1.60	0.018
2	302100	1.5 ml micro	5 x 100	1.60	0.018
3	302400	2.5 ml semimicro	5 x 100	1.60	0.018

Special UV grade cuvettes

Made of UV grade PMMA for accuracy throughout UV/Vis range (280 to 800).

Model 303100 with four clear sides is ideal for fluorometry and nephefelometry.

code		description	case quantity	case weight	case volume
4	303100	4.5 ml 4 sides transparent	5 x 100	1.60	0.018
5	303102	4.5 ml macro	5 x 100	1.60	0.018
6	303101	2.5 ml semimicro	5 x 100	1.60	0.018
7	303103	1.5 ml micro	5 x 100	1.60	0.018

Cap

In polyethylene. Suitable for the macro cuvettes, 4.5 ml (302000, 303100, and 303102)

	code	description	case quantity	case weight	case volume
8	304000	polyethylene cap	1,000	0.20	0.0024



Please see our rack for the spectrophotometer cuvettes, code M-100, on Chapter 9. Sample Storage.

- Polypropylene
- Autoclavable
- Assembling rack that saves space
- Capability: 42 cuvettes





Sample cups

Polystyrene cups, with polyethylene caps.

mod.	code	description	type	Ø ext. mm	h mm	case quantity	case weight	case volume
1	900024	0.50 ml	Gemsaec, Kone Lab 20	13.55	24.50	10 x 1,000	12.70	0.069
2	900023	0.25 ml	Centrifichem	13.70	16.40	14 x 1,000	15.00	0.069
3	900022	1.50 ml	Technicon	13.80	22.60	10 x 1,000	10.60	0.069
4	192503	1.50 ml	Technicon RA 1000	13.90	24.70	12 x 1,000	10.50	0.060
5	910022	2.00 ml	Technicon	13.75	24.90	10 x 1,000	10.00	0.070
6	900026	push on cap suitable, for all models	-	16.00	8.50	20 x 1,000	3.50	0.069

Cobas Mira cuvette segments

Multicuvettes made of PMMA.

Rack made of polypropylene (colour: red).


For use on Cobas Mira.

Light path: 6 mm.

Complete rack with 15 segments of 12 cuvettes.

code	description	case quantity	case weight	case volume
900019	complete rack: 15 segments of 12 cuvettes	30	9.00	0.040



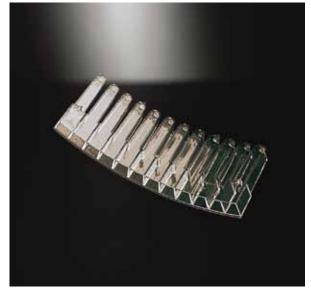

Coulter counter cups

Coulter counter cups: single use polystyrene cups suitable for any Coulter counter for cell counting.

Two models available, with or without lid.


Volume: 20 ml.

Dimensions: 30 x 56 mm (Ø x h).

code	description	case quantity	case weight	case volume	
200103	coulter counter cup with lid	1,000	7.80	0.087	
200102	coulter counter cup without lid	1,400	8.70	0.120	





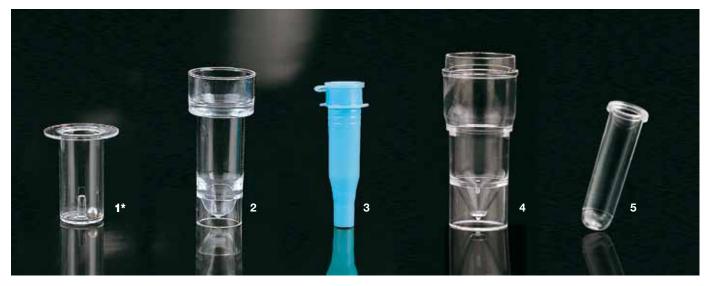

Sample cups

Cups are made of polystyrene, except code 900008 which is made of HDPE.

mod.	code	code description		Ø ext. mm	h mm	case quantity	case weight	case volume
1*	900005.1	1.5 ml, cup with steel	Amelung	18.65 / 11.25	22.8	12 x 1,000	9.40	0.092
2	900020	2.50 ml	Hitachi	16.75	38.0	6 x 1,000	10.75	0.069
3	900008	0.70 ml	Cobas bio	7.65	35.5	12 x 1,000	8.35	0.067
4	910023	4.00 ml	Technicon	16	37.9	6 x 1,000	12.20	0.069
5	910026	0.5 ml (0.8 ml total volume)	Sysmex	9.95 / 7.94	29.8	10 x 2,000	13.50	0.056

^{* 12} bags of 1,000 sample cups and 12 dispenser containers of 1,000 balls.



Scintillation vials

Single use vials made of high density polyethylene to minimize solvent losses. Screw caps ensure a leakproof seal.

Compatible with most liquid scintillation counters available on the market. 2 models available: code **900100**, standard volume 20 ml, and code **900101**, volume 4 ml, designed for insertion in the 20 ml vial so as to minimize the volume of scintillation liquid.

Dimensions: Vial 20 ml: $26.5 \times 58.5 \text{ mm}$ (Ø x h) Minivial 4 ml: $13.71 \times 53.15 \text{ mm}$ (Ø x h).

code	description	case quantity	case weight	case volume
900100	scintillation vial (20 ml)	1,000	7.90	0.069
900101	minivial (4 ml) for the above vial	2,000	5.90	0.042



Beakers for cell counters


Single use beakers ideal for dilutions and cell counting in Toa and Royco-Hycel cell counters.

Both models are made of polyethylene.

code	baker type	Ø mouth mm	Ø base mm	h mm	vol. ml.	case quantity	case weight	case volume
1988314	Toa	30.5	25	68.6	30 ml	3,000	6.54	0.039
1988315	Royco-Hycel	30.5	26	57.5	25 ml	3,000	6.50	0.039





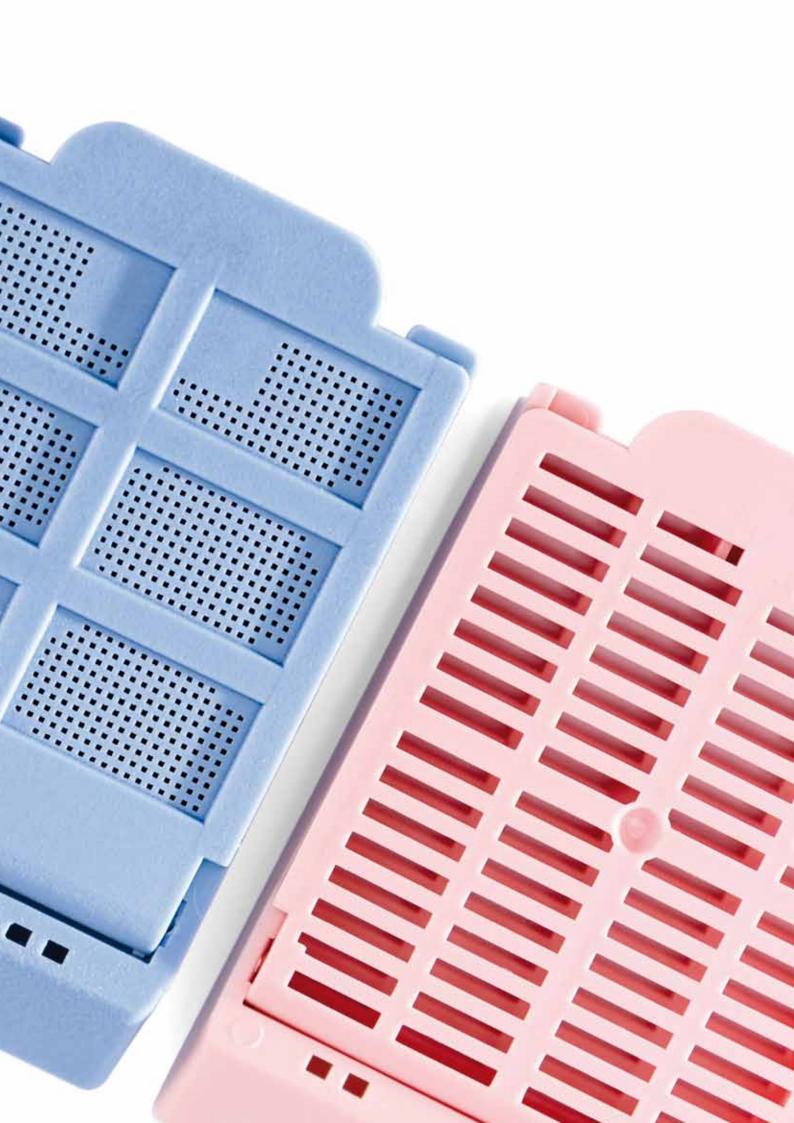

Polyethylene wide mouth jars with lid and insert plug

High density Polyethylene jars with screw cap and insert plug. Jars and plugs are natural transluscent; lids are black.

Upon request, jars can be supplied individually wrapped for a minimum order of 3 cases for each code, and sterile for a minimum order of 6 cases, for each code.

code	capacity ml	Ø neck mm	Ø base mm	heigth mm	case quantity	case weight	case volume
202811	30	25	32	53	100	1.00	0.011
202810	60	38	52	48	200	3.30	0.045
202808	90	51	60	51	150	3.90	0.045
202809	125	51	60	62	200	6.00	0.082
202801	170	51	60	80	160	6.46	0.082
202802	250	54	67	100	125	5.24	0.078
202814	400	60	74	124	130	8.00	0.140
202803	500	67	80	131	120	8.40	0.140
202821	500	86	103	93	95	9.20	0.140
202813	750	86	103	127	75	7.70	0.140
202818	1,000	86	103	157	50	7.00	0.140



See more capacities on pages 280-281

Jars with lid

Body and lid in white **autoclavable** polypropylene.

Self sealing security lid.

Made with materials suitables for alimentary use.

Code 241014 embodies a plastic handle.

Jars are supplied uncapped.

code	nominal capacity ml	Ø int. mouth x Ø base x height mm	weight g	case quantity	case weight	case volume	casses per pallet
241010	250	84.5 x 79.5 x 62	18	150	3.03	0.042	32
241011	550	107 x 98 x 80	32	250	5.01	0.103	16
241013	1,000	112.5 x 112 x 137.5	46	150	7.68	0.078	20
241014	1,560	133.5 x 134 x 140	58	150	4.42	0.130	16



Tamper evident, buckets with lid

Disposable containers made of **autoclavable** polypropylene. Made with materials suitables for alimentary use. Include tight fitting and leakproof lids. All models include a sturdy and convenient white plastic handle for easier carrying. Buckets are supplied with lid apart.

code	capacity I	Ø int. mouth mm	height mm	body weight g	case quantity	case weight	case volume	casses per pallet
222802	3	184.00	137.50	105.00	80	9.00	0.140	16
222803	4.5	210.00	156.00	134.00	48	10.46	0.130	16
222804	5.6	211.00	194.00	152.00	45	11.58	0.130	16
222805	10.6	251.00	263.00	314.00	20	10.50	0.130	16









EUROTUBO® Leak proof screw cap containers

(€ (IVD)

Specially designed for histology, they are provided with double internal security closure. Manufactured in translucent polypropylene with polyethylene ribbed yellow cap. They have moulded graduations in increments of 50 ml, and a ring in the upper part of the container that prevents dripping when decanting liquids. Its leakproof closure makes them ideal for liquid (including alcohols) or solid sample storage. They meet the standard UNE-EN 14401:2005 Rigid plastic containers. Methods to test the effectiveness of closures. Supplied uncapped. Other cap colours available upon request.

mod.	code	capacity ml	external mouth Ø mm	height mm	case quantity	case weight	case volume	cases per pallet
1	202845	250	89	57	100	4.68	0.046	32
2	202846	500	90	110	200	10.96	0.084	20
3	202847	1,000	111	139	100	11.40	0.097	16



Security screw cap containers

(€ (IVD)


Ultra clear polypropylene containers with yellow polyethylene cap.

Watertight closure: the cap, with an innovative shape, has a double closure. Graduated Container.

The base and the cap are ribbed so the handling with gloves is comfortable and reliable.


There is a ring in the upper part of the body, in order to avoid the dripping in case of liquid decanting.

Six different capacities from 20 to 120 ml.

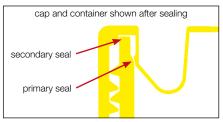



Find more information about this product on Chapter 7. **Sample containers**

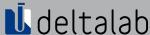

Paraffin (€ (IVD)

It has a low dimethyl sulphoxyde content and a quick penetration of tissues for $2\ \mbox{micron}$ cuts.

Melting point: 56/58 °C. Can be used with usual solvents. Ideal for histology.

	mod. code		description	case quantity	case weight	case volume
	1	440139	Paraplast Plus® - 1 kg bag	8 bags	8.6	0.021
Nav	2	440139P	paraffin - 1 kg bag	6 bags	6.5	0.021







Histology mould

Single use clear polypropylene mould. Can be re-used.

The smooth interior ensures the paraffin does not adhere to the mould. The special material used ensures excellent thermal exchange. Available in various sizes

External dimensions: 50 x 37.2 x 12 mm. Suitable for the most of cassettes.

mod.	code	well dimensions mm	case quantity	case weight	case volume
1	440131	15 x 15 x 5	1,000	2.20	0.013
2	440132	24 x 24 x 5	1,000	2.20	0.013
3	440133	30 x 24 x 5	1,000	2.20	0.013



Metal tray histology

Material: stainless steel.

External dimensions: $52 \times 35 \times 11 \text{ mm}$.

mod.	code	well dimensions mm	bag quantity	bag weight	bag volume
1	192931	7 x 7 x 7	10	0.08	0.0001
2	192932	15 x 15 x 7	10	0.08	0.0001
-	192933	24 x 24 x 7	10	0.09	0.0001
3	192934	30 x 24 x 7	10	0.09	0.0001
4	192935	37 x 24 x 7	10	0.09	0.0001




Tissue embedding sponge

Foam sponge, very useful for small specimens in biopsy tests. Fits standard cassettes.


Avoids sample loss though drain holes when sample is processed. Does not affect paraffin drainage.

code	dimensions mm	colour	case quantity	case weight	case volume
192922.04	32 x 26 x 3	blue	500	0.06	0.002

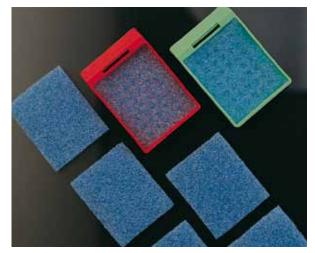





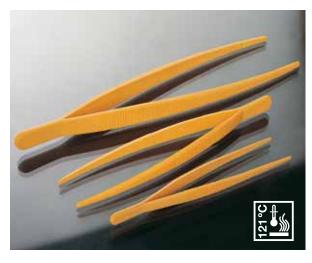

Forceps

Made of P.O.M.


Colour: orange.

Self-sprung forceps with rounded ends.

Code **19503**, the inner part of the tip is ribbed to enhance sampling. Ribbed for a better handling.

Autoclavable. Flexible.

code	length mm	case quantity	case weight	case volume
19500	115	5	0.23	0.0006
19501	145	5	0.05	0.00005
19503	250	5	0.17	0.0043

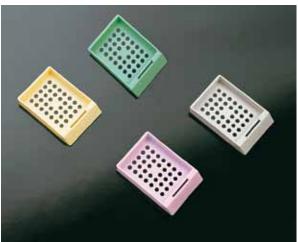

Cassettes without lid

Made of polyacetyle.

With 34 round holes of 2 mm Ø.

Resistant to the solvents used in histological processes. It is possible to write on them.

External dimensions: 40 x 28 x 6 mm.

code	colour	case quantity	case weight	case volume
440141	white	3 x 1,000	7.00	0.047
440147	green	3 x 1,000	7.00	0.047
440144	pink	3 x 1,000	7.00	0.047
440143	yellow	3 x 1,000	7.00	0.047



Lids for cassettes

Stainless steel lids.

Fit all standard cassettes.

Suitable for **440141** and series.

code	description	case quantity	case weight	case volume
440149	lid for cassettes	25	0.07	0.0005



Cassette for biopsy 6 compartments

Material: polyacetyle.

Designed for small biopsies, each cassette has six individual cells to hold six samples.

Each cell is 7 x 12 mm. The lid allows multiple opening and closing. A frosted area allows writing on the cassette.

Each cassette has near 3,000 micro holes of 0.26 mm diameter.

code	colour	case quantity	case weight	case volume
44240	pink	4 x 250	5.20	0.026
44241	white	4 x 250	5.20	0.026
44245	blue	4 x 250	5.20	0.026

Ask for minimum quantity and delivery time of other colours.



CASSETTES HISTOSETTE I

- Made of polyacetal.
- Air vents allowing more efficient filling with paraffin.
- Safety closure.
- External dimensions: 40 x 28 x 6 mm (with closed lid)
- With one-piece integral lid.
- To lock the cassette, just snap apart the lid, turn it round and lock it into the base of the cassette.
- They may be opened and closed as many times as required.

HOW TO USE




Cassettes pour tissus. Histosette I

Writing area at 30° or 45° angle, depending on the code. 64 cuts (1 x 5 mm) in base and lid, distributed 16 x 4.

code 30°	code 45°	colour	case quantity	case weight	case volume
454141	454111	white	3 x 500	5.40	0.039
454147	454117	green	3 x 500	5.40	0.039
454140	454110	pink	3 x 500	5.40	0.039
454143	454113	yellow	3 x 500	5.40	0.039
454148	454118	blue	3 x 500	5.40	0.039

Ask for minimum quantity and delivery time for other colours.

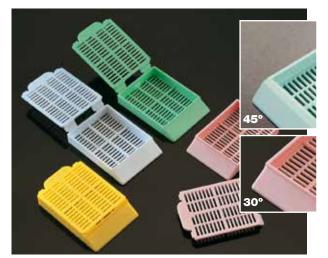

Cassettes for biopsy. Histosette I

Writing area at 45° angle. 224 draining holes (1 x 1 mm).

code	colour	case quantity	case weight	case volume
554111	white	3 x 500	6.35	0.034
554117	green	3 x 500	6.35	0.034
554113	yellow	3 x 500	6.35	0.034

Ask for minimum quantity and delivery time for other colours.



Cassettes for tissue

Made of acetal polymer. With an integral lid. Small tab on the lid for easy opening and locking.

code	colour	case	case	case
2340	23.041	quantity	weight	volume
554141	white	4 x 500	7.58	0.064
554147	green	4 x 500	7.58	0.064
554140	pink	4 x 500	7.58	0.064
554143	yellow	4 x 500	7.58	0.064
554148	blue	4 x 500	7.58	0.064

Ask for minimum quantity and delivery time for other colours.






CASSETTES HISTOSETTE II

- Made of polyacetal.
- Hinged lid: can be opened and closed as often as necessary.
- Supplied assembled, saving time and making storage easier.
- External dimensions: 41 x 28 x 6 mm.
- Air vents allowing more efficient filling with paraffin.
- Easy handling with one hand.
- Safety closure (see figure).
- Writing area slanted at a 45° angle, suitable to be used with most labelling instruments.



HOW TO USE

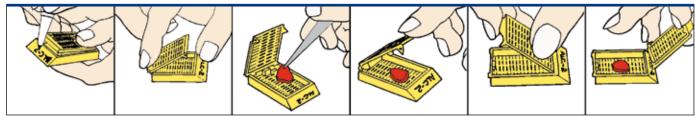











Cassettes for tissue. Histosette II

Material: acetal polymer. 64 cuts in base and lid.

Cassettes and lids are available: pre-mounted, or with the lid separated in another case.

code	colour	case quantity	case weight	case volume					
Cassettes w	Cassettes with the lid pre-mounted:								
44140	pink	3 x 500	5.00	0.026					
44141	white	3 x 500	5.00	0.026					
44143	yellow	3 x 500	5.00	0.026					
44145	blue	3 x 500	5.00	0.026					
44147	green	3 x 500	5.00	0.026					
Cassettes w	ith the lid separated:								
441140	pink	2 x 500*	3.56	0.025					
441141	white	2 x 500*	3.56	0.025					
441143	yellow	2 x 500*	3.56	0.025					
441145	blue	2 x 500*	3.56	0.025					

^{*2} dispensing boxes with 500 cassettes each plus one box containing 1,000 lids.

Cassettes for biopsy. Histosette II

Made of acetal polymer. 224 holes. 1 mm square openings to maximize fluid exchange and to ensure proper drainage.

Cassettes and lids are available: pre-mounted, or with the lid separated in another case.

code	colour	case quantity	case weight	case volume			
Cassettes w	ith the lid pre-mounted:						
445140	pink	3 x 500	5.20	0.026			
445141	white	3 x 500	5.20	0.026			
445143	yellow	3 x 500	5.20	0.026			
445145	blue	3 x 500	5.20	0.026			
445147	green	3 x 500	5.20	0.026			
Cassettes w	Cassettes with the lid separated:						
451241	white	2 x 500*	3.90	0.026			
451243	yellow	2 x 500*	3.90	0.026			
451245	blue	2 x 500*	3.90	0.026			

^{*2} dispensing boxes with 500 cassettes each plus one box containing 1,000 lids.



CASSETTES SLIMSETTE

- Made of polyacetal.
- Hinged lid: can be opened and closed as often as necessary.
- Supplied assembled, saving time and making storage easier.
- Hinged lid may be separated if required.
- Dimensions: 41 x 28.5 x 6 mm.


- Easy handling with one hand.
- Safety closure.
- Air vents allowing more efficient filling with paraffin.
- Writing area slanted at a 45° angle, suitable with most labeling instruments.
- Only 6 mm high, allowing more cassettes to be stacked in automated labeling machines and in storage cabinets.
- Two versions available: single cassettes and in stacks for Leica® and Sakura® labeling instruments.

Cassettes for tissue. Slimsette

Quick drainage thanks to 114 vents each measuring 5 x 1 mm.


Ask for minimum quantity and delivery time for other colours.

code	colour	case quantity	case weight	case volume		
Single cassettes (The main case contains three dispenser boxes of 500 cassettes each)						
555141	white	3 x 500	5.01	0.025		

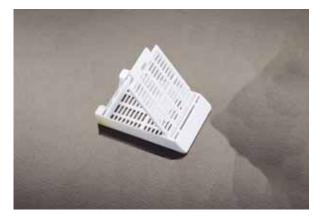

code	colour	case quantity	case weight	case volume		
Stacks for Leica® and Sakura® (each case contains 50 stacks of 40 cassettes)						
555251 white		50 x 40	7.12	0.023		
555258 blue		50 x 40	7.12	0.023		

Ask for minimum quantity and delivery time for other colours.

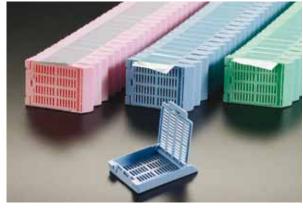

Cassettes for biopsy. Slimsette

Quick drainage thanks to its 224 vents (0.8 x 0.8 mm).

adion drainage thanks to its 22 i volta (c.s x c.e iiiii).						
code	colour case case quantity weight v					
Single cassettes (The main case contains three dispenser boxes of 500 cassettes each)						
44151	white	3 x 500	5.01	0.025		

Ask for minimum quantity and delivery time for other colours.



code	colour	case quantity	case weight	case volume	
Stacks for Leica® and Sakura® (each case contains 50 stacks of 40 cassettes)					
44351 white		50 x 40	6.84	0.022	
44355	blue	50 x 40	6.84	0.022	

Ask for minimum quantity and delivery time for other colours.





- Made of polyacetal.
- Cassettes type Richard-Allan® and Surgipath® Multi-Cassette (Leica)
- Hinged lid: can be opened and closed as often as necessary.
- Supplied assembled, saving time and making storage easier.
- Dimensions: 42 x 28 x 6 mm.

- Large tab for convenient an easy opening
- Easy handling with one hand.
- Safety closure.
- Writing area at 45°, suitable for most automatic labelers.
- Air vents allowing more efficient filling with paraffin.



Cassettes for tissues. Swingsette

Designed to hold tissue specimens. 64 holes of approximately 0.76 x 5.06 mm; distributed 16 x 4.

code	colour	case quantity	case weight	case volume
415141	white	3 x 500	5.44	0.025
415143	yellow	3 x 500	5.44	0.025
415147	green	3 x 500	5.44	0.025

Ask for minimum quantity and delivery time for other colours.





Cassettes for biopsy. Swingsette

Designed to hold biopsy specimens. 224 holes of approximately 0.76 x 0.76 mm.

code	colour	case quantity	case weight	case volume
414141	white	3 x 500	5.44	0.025
414143	yellow	3 x 500	5.44	0.025
414147	green	3 x 500	5.44	0.025
	9			

Ask for minimum quantity and delivery time for other colours.



Mega cassettes

Made from polyacetal.

With an integral lid.


Designed for larger tissue samples.

Holes are 1 x 5 mm.

Cassette dimensions: 26 x 40 x 13 mm.

code	colour	case quantity	case weight	case volume
44125	white	500	2.57	0.025

Ask for minimum quantity and delivery time for other colours.

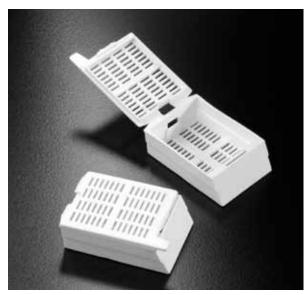



Storage cabinet

Made of ABS. Marble color.

This storage cabinet with lid includes 5 drawers and is designed for the storage of small units such as slides or embedding rings.

Drawers are designed to hold ABS racks code 19921.

The storage cabinet can hold a total of 5,000 slides, as each drawer can hold 10 racks and each rack can hold 100 slides.

Each slot on the rack is numbered.

Polystyrene covers for the racks are available (code 19920).

A four wheel support (code 19931) is also available for the cabinet.

code	description	case quantity	case weight	case volume
19930	420 x 420 x 730 mm	1	13.64	0.164
19931	4 wheels	1	0.02	0.00005
19921	19921 cabinet rack		0.37	0.003
19920	cover	2	0.19	0.001





Economical storage box

This box, complete with a lid, is specially designed for the storage of histology cassettes.

Made of heavy-duty cardboard for long term storage.

Each box stores up to 250 cassettes or 160 embedding rings.

Dimensions: 455 x 235 x 50 mm.

code	description	case quantity	case weight	case volume
DS42	waterproof cardboard	18	5.3	0.130



Modular storage drawers


Made of high impact polystyrene, this set consists of 6 stackable drawers (7 slots per drawer) providing permanent storage and identification of up to 165 embedding rings or 250 cassettes per drawer (total capacity: 1,500 cassettes).

Dimensions: 400 x 230 x 340 mm. Front: black colour. Sides: grey.

code	description	case quantity	case weight	case volume
DS40	set of 6 drawers	1	6.14	0.037





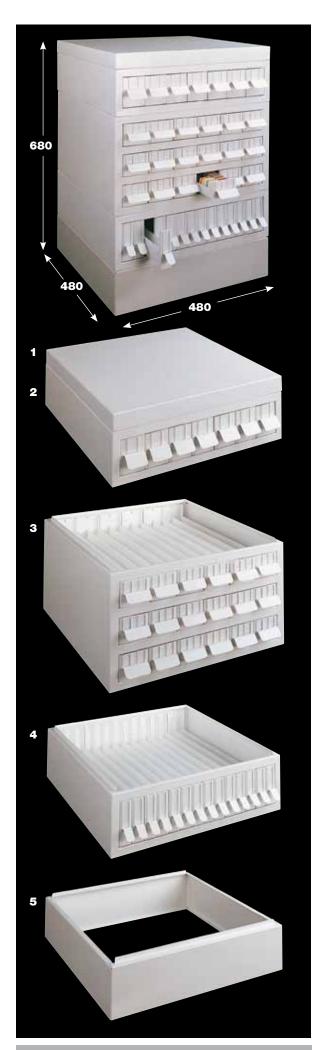

Laboratory storage cabinet

Designed especially for the storage of microscope slides, photo slides, histology cassettes and embedding rings.

Universal models interchangeable with other cabinets on the market. All parts of the cabinet are stackable.

Manufactured from enameled metal. Colour grey.

You can make up your own composition from one or several of the following elements:

- Modular photo slide drawers
- Modular cassette drawers
- Modular microscope slide drawers

All modules are stackable.



Modular top (1)

code	external dimensions mm		unit	unit	
	height	length	width	weight	volume
DS61	50	480	480	2	0.056

Modular photo slide drawers (2)

Cabinet consisting of 7 independent drawers to hold a total of 1,700 transparencies of 50×60 mm.

Dimensions of a drawer: 58 x 475 mm.

code	external dimensions mm		unit	unit	
	height	length	width	weight	volume
DS12	120	480	480	9.54	0.056

Modular cassette drawers (3)

Cabinet consisting of 6 double drawers laid on 3 levels.

Each drawer has 2 compartments.

The cabinet can hold approximately 1800 embedding rings.

Internal dimensions of each drawer: 68 x 475 mm.

code	exte	external dimensions mm			unit	
	height	length	width	weight	volume	
DS35	255	480	480	24.32	0.091	

Modular microscope slide drawers (4)

Consisting of 14 sliding drawers to hold a total of approximately 6,500 microscope slides of 26×76 mm.


Internal drawer dimensions: 27x475 mm.

code	exte	external dimensions mm			unit	
code	height	length	width	weight	volume	
DS11	140	480	480	15	0.056	

Modular base (5)

To fit all modular cabinets and drawers.

code	exte	external dimensions mm			unit	
	height	length	width	weight	volume	
	DS62	120	480	480	3.72	0.090

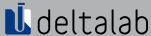


STAR FROST® ADHESIVE slides

Standard size: 26 x 76 mm. Thickness: 1 mm. Made in Germany. High definition, pre-washed glass. Hydrofophobic. No stain residues, no proteins in the surface.

Feature a white, 20 mm. marking area. This special printed area is resistant to all common chemicals used in the laboratories, and may be marked with pencil, ballpoint pen, marker and different laboratory-printers (e.g. Leica XPS). Thanks to the single side printing the user always knows the best surface to be used.

Star Frost® Adhesive slides are treated with silane-coating, which ensures cell-and-tissue-adhesion and supports the forming of covalent bonds between the tissue and the glass surface through electrostatic attraction. So tissue and cell substrates stick rapidly on the surface without using further adhesives or protein coating media, saving costs and time. This adhesive coating is restricted to the top surface and therefore to the usable area only. As the underside is untreated, the slides do not tend to stick when stacked or transported as with other competitive products.

StarFrost® Adhesive slides do not show striae and disturbing background staining such as the blue or red background obtained with Haematoxylin and Eosin-staining of tissues on Albumin coated glass, or the brown background obtained during the Immunoperoxidase or In-situ-DNS-Hybridisation process. They are ideal for molecular-hybridization, as being DNASe and RNAse free there is no need of washing them with chromic-acid.

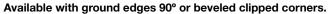

Supplied in boxes of 50 units.

Final package (20 boxes) by tropical pack.

Conforms with ISO 8037/1 (european standard).

Expiry date: 14 months after manufacturing date.

code	edges	corners	case quantity	case weight	case volume
A100018	grounded 90°	no clipped	20 x 50	4.98	0.003
A100019	grounded bevelled	clipped	20 x 50	4.98	0.003



TROPICAL PACK:

Aluminium film designed to protect slides from humidity.

STAR FROST® slides

Standard size: 26 x 76 mm. Thickness: 1 mm.


Made in Germany of high definition, pre-washed glass. Hydrofophilic (epoxy). Feature a colour-coded, 20 mm marking area.

This special colour-printed area is resistant to all common chemicals used in the laboratories, and may be marked with pencil, ballpoint pen, marker and different laboratory-printers (e.g. Leica XPS). Thanks to the single side printing the user always knows the best surface to be used.

Ground edges at 45° (angle) and clipped corners. Supplied in boxes of 50 units. Final package (20 boxes) by tropical pack.

Expiry date: 14 months after manufacturing date.

code	strip colour	case quantity	case weight	case volume
A100010	white	20 x 50	4.98	0.003
A100011	green	20 x 50	4.98	0.003
A100013	blue	20 x 50	4.98	0.003
A100014	yellow	20 x 50	4.98	0.003

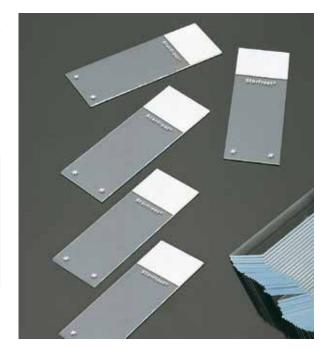












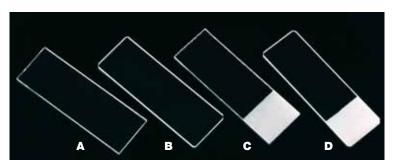

EUROTUBO® slides

Standard dimensions: 26 x 76 mm (\pm 0.2). Clear, pre-washed glass, 1.1 mm (\pm 0.1) thick.

Ideal for routine analysis.

Displayed in cases of 50 units.

Each case is vacuum packed, easy to open, with no need for scissors.

Each case includes a bag of silica salt.

Expiry date: 15 months after manufacturing date.

	code	edges	corners	matt area	case quantity	case weight	case volume
Α	D100001	cut 90°	no clipped	no	50 x 50	14.25	0.013
В	D100002	grounded 45°	clipped	no	50 x 50	14.20	0.010
С	D100003	cut 90°	no clipped	yes	50 x 50	14.20	0.010
D	D100004	grounded 45°	clipped	yes	50 x 50	14.25	0.013



EUROTUBO® cover slides

Specific use: protection of the sample and lens.

Clear, pre-washed glass. Thickness: No 1: 0.13 x 0.16 mm.

Supplied in cases of a hundred or two hundred units (depending on the code). Each case includes a bag of silica salt, and is supplied vacuum in tropical pack. Easy opening, no need for scissors.

code	dimensions mm	case quantity	case weight	case volume
D101818	18 x 18	5 x 200	0.52	0.0016
D102020	20 x 20	5 x 200	0.71	0.0016
D102222	22 x 22	5 x 200	0.65	0.0016
D102240	22 x 40	10 x 100	0.63	0.0016
D102424	24 x 24	5 x 200	0.72	0.0016
D102432	24 x 32	10 x 100	0.61	0.0016
D102440	24 x 40	10 x 100	0.65	0.0016
D102450	24 x 50	10 x 100	0.70	0.0016
D102460	24 x 60	10 x 100	0.61	0.0016





Microscope slide labels 23 x 23 mm

Write on these permanent adhesion labels and stick them onto microscope slides for better identification.

White colour.

Each sheet contains 40 labels.

Page size: 13.3 x 19 cm.

code	quantity	weight	volume
901600	1,000 labels	2.5	0.020





Mounting medium for microscopy

 $\mbox{VITROCLUD}{}^{\mbox{$^\circ$}}\mbox{ quick-hardening mounting medium for microscope preparations.}$ Neutral and colourless.

Short hardening time (20 min). Refractive index at room temperature: 1.5. Resistant to heat, cold (down to -17 °C), humidity and light including UV rays. It is recommended to close the container immediately after use.

If the medium becomes too viscous, add Xylene.

Supplied in aluminium bottles for a safer storage.

Store at maximum 20 °C.

code	description	case quantity	case weight	case volume
A20100	bottle 100 ml	1	0.12	0.001
A20250	bottle 250 ml	1	0.27	0.001
A20500	bottle 500 ml	1	0.54	0.002



Staining tray

Made of ABS.

It will accept 10 or 20 slides on plastic rails covered by a polymer strip to perfectly hold slides even if tray is held at an angle (see picture A).


When humidity is required, wells between rails may hold water, and raised rails will avoid water from touching the slides.

It features a draining plug to ease the emptying of the tray (see picture B). Four rubber feet ensure stability.

Clear lid made of PETG (also available in black ABS, for fluorescence work; ask for delivery time).

Cannot be used with acetone and chlorinated hydrocarbons.

	code	description	case quantity	case weight	case volume
MaW	989953	staining tray 10 slides	1	0.65	0.003
	989951	staining tray 20 slides	1	1.08	0.005

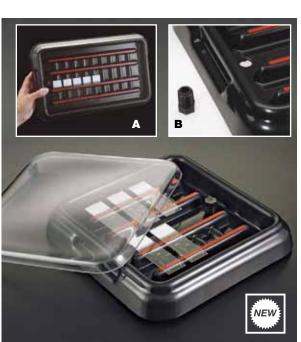





Stains

	code	description		case quantity	case weight	case volume
C€ (IVD)	808000	Eosin, solution, May Grunwald, 250 ml bottle Danger symbols: R11-23/24/25-39/23/24/25. S 7-16-36/37-45 Transport norms: UN: 1992. ADR: 3/19b.IMDG: 3.2/II. IATA: 3/II. PAX: 305. CAO: 307	&	48	14.40	0.045
C€ (IVD)	808001	Eosin, solution, May Grunwald, 1000 ml bottle Danger symbols: R11-23/24/25-39/23/24/25. S 7-16-36/37-45 Transport norms: UN: 1992. ADR: 3/19b.IMDG: 3.2/II. IATA: 3/II. PAX: 305. CAO: 307	48 12 48 48 48 48 48 48 48 48 48 4		13.00	0.045
C€ (IVD)	808100	Eosin, solution modified Giemsa, 250 ml bottle Danger symbols: R11-23/24/25-39/23/24/25. S 7-16-36/37-45 Transport norms: UN: 1992. ADR: 3/19b.IMDG: 3.2/II. IATA: 3/II. PAX: 305. CAO: 307	osin, solution modified Giemsa, 250 ml bottle symbols: R11-23/24/25-39/23/24/25. S 7-16-36/37-45 rt norms: UN: 1992. ADR: 3/19b.IMDG: 3.2/ll. IATA: 3/ll.			0.015
C€ (IVD)	808101	Eosin, solution modified Giernsa, 1000 ml bottle Danger symbols: R11-23/24/25-39/23/24/25. S 7-16-36/37-45 Transport norms: UN: 1992. ADR: 3/19b.IMDG: 3.2/II. IATA: 3/II. PAX: 305. CAO: 307		12	13	0.045
C€ (IVD)	808200	Eosin, solution Wright, 250 ml bottle Danger symbols: R11-23/24/25-39/23/24/25. S 7-16-36/37-45 Transport norms: UN: 1992. ADR: 3/19b.IMDG: 3.2/II. IATA: 3/II PAX: 305. CAO: 307		48	13.30	0.045
Gram	806030	Lugol bottle 250 ml		48	14.7	0.045
differential	805040	Violet crystal bottle 250 ml	<i>★</i>	48	14.4	0.045
staining	805001	Gram decolorizing bottle 250 ml		48	12.8	0.045
KIT: 805000*	805050	Safranin bottle 250 ml		48	14.4	0.045
	805001	Gram decolorzing, 250 ml bottle. Danger symbols: R-11-36-66-67. S 9-16-26 Transport norms: UN: 1993 ADR: 3/3b. IMDG: 3,2/II. IATA: 3/II PAX: 305. CAO: 307	*	16	4.12	0.014
	806030	Lugol lodine, 250 ml bottle. Danger symbols: R52		16	4.76	0.014
	805040	Violet crystal solution (Gentian violet), 250 ml bottle Danger symbols: R-21/22-36/38. S 28a Transport norms: UN: 2810 ADR: 6/25c. IMDG: 6.1/III. IATA 6.1/III. PAX: 611. CAO: 608	×	16	4.68	0.014
	805140	Violet crystal solution (Gentian violet), 1000 ml bottle Danger symbols: R-21/22-36/38. S 28a Transport norms: UN: 2810 ADR: 6/25c. IMDG: 6.1/III. IATA 6.1/III. PAX: 611. CAO: 618	×	12	13.00	0.045
Ziehl	805020	Fuchsin bottle 250 ml		16	14	0.045
differential staining	805030	Methylene blue bottle 250 ml	M X	48	14.40	0.045
KIT: 805010*	805011	Ziehl decolorizing 2 bottles 250 ml		48	13.40	0.045
	805011	Ziehl decolorizing, 250 ml bottle Danger symbols: R-11-36/37/38 S 7-16-26-45 Transport norms: UN: 2924 ADR: 3/26b. IMDG: 3,2/II. IATA: 3/II PAX: 305. CAO: 307	*	16	3.39	0.014
	805050	Safranin solution, 250 ml bottle		16	4.64	0.014
	805150	Safranin solution, 1000 ml bottle		12	13.00	0.045
	805120	Ziehl fuchsin solution, 1000 ml bottle Danger symbols: R10-21/22-36/38. S 36/37.	×	12	13.00	0.045
	805030	Methylene blue, 250 ml bottle	×	16	4.70	0.014
	805130	Methylene blue, 1,000 ml bottle	×	12	13.00	0.045
C€ (IVD)	805013	Fast stain for blood extensions, 250 ml containers, Kit with 2 containers A and 2 containers stain B		12 Kits	14.95	0.045

^{*}Please order the KITS 805000 and 805010 only for ground shipping. For maritime shipping orders, please order individual codes.

Minimum order quantity: 1 bottle or kit.



Staining jars with lid

Material: autoclavable TPX.

Two options available: vertical, "Hellendhall" type (Code **19335**) horizontal "Schifferdecker" type (Code **19351**).

code	description	case quantity	dimensions mm	case weight	case volume
19355	vertical, for 8 slides	4	58 x 53.5 x 86	0.27	0.001
19351	horizontal, for 10 slides	4	76 x 65 x 45	0.32	0.001





Staining jar

Made of **POM** (polyoxymethylene).
Colour: black. Ideal for slide staining.
Good resistance to alcohols and xylol (not to phenol).
Consists of a jar with a lid for slide staining.
Leakproof. Includes a white perforated slide partition for drying.
Not suitable for microwaves.

code	dimensions	case	case	case
	mm	quantity	weight	volume
191101	100 x 87 x 51	1	0.08	0.004



Staining rack

Made of **POM** (polyoxymethylene). Colour: black. Ideal for slide staining. Good resistance to alcohols and xylol (not to phenol). Accepts 25 slides. Includes a folding handle, and 25 numbered slots to facilitate classification. Slides insert easily in the rack. The rack ideally fits into the staining through presented above.

No suitable for microwaves.

code	dimensions	case	case	case
	mm	quantity	weight	volume
191100	91 x 79 x 38	1	0.03	0.001









Staining system

Slide staining system compound of a range of jars, and a 12 slides staining rack. Different jar colours allow the user to distinguish each stage of the staining process. Jar and rack made of **autoclavable** POM (polyoxymethylene), resistant to staining agents as alcohol, xylene, etc. (not phenol).

The jar embodies a hinged lid which allows storing the jar with the content. Reagent capacity: 80 ml.

Jars can be loosely joined to each other laterally (picture 1).

The rack features a handle in the hinged lid, and can be placed up when handled out of the jar. Since the slides are placed vertically in the rack, their writing area will not be stained, allowing their removal without the use of forceps (picture 2).

Rack and jar withstand from -170 °C and 121 °C.

Jar dimensions: 64 x 76 x 92 mm. Rack dimensions: 60 x 64 x 97 mm.

code	description	case quantity	case weight	case volume
191104	blue jar	6	0.72	0.0059
191105	green jar	6	0.72	0.0059
191106	white jar	6	0.72	0.0059
191107	yellow jar	6	0.72	0.0059
191108	191108 dark grey rack		0.32	0.0024



Staining jar and dish

Dish made of polypropylene and jar made of **autoclavable TPX**. The jar is supplied with two lids. The first lid provides evaporation free storage of staining liquid when the dish is not in use. The second includes a slot to allow staining dish insertion.

Code 19353: height (with lid) 70 mm, height (without lid) 65 mm, lenght 100 mm and width 75 mm.

Code **19354**: height 21 mm, length 83 mm, width 70 mm and lenght (handle) 160 mm.

code	description	case quantity	case weight	case volume
19353	jar with 2 lids	4	0.42	0.006
19354 staining dish for 20 slides		2	0.06	0.001



Staining racks

Stainless steel.

Models **S-004** and **S-002** can be leaning over to drain off the stain without falling the slide.

code	description	case quantity	case weight	case volume
S-004	for 12 slides length: 435 mm, width: 85 mm	1	0.24	0.003
S-002	for 24 slides lenght: 435 mm, width: 175 mm	1	0.34	0.003
S-003	without dividers lenght: 435 mm, width: 85 mm		0.19	0.003

Other stainless steel models can be made on request.



Staining jar and basket

Code 19360 is a stainless steel staining jar.

With a lid to avoid evaporation.

Made of stainless steel, this staining basket is designed for 30 slides of 26×76 mm for use in the above staining jar, code **19360**.

Includes an adjustable handle for more convenience.

code	description	dimensions mm	case quantity	case weight	case volume
19360	jar	115 x 88 x 77	1	0.42	0.001
19361	basket	110 x 85 x 55	1	0.13	0.001



Special baskets and trays for slides

Baskets are manufactured from stainless steel wire. Cuvettes are made of neutral glass and feature stainless steel lids.

Models CP-24, CP-30, CP-50 have folding handles.

code	description	dimensions mm	case quantity	case weight	case volume
CP-30 H	basket for 30 horizontal slides	172 x 86	1	0.111	0.0010
CP-45 H basket for 45 horizontal slides		235 x 85	1	0.159	0.0016
CP-30	basket for 30 vertical slides	70 x 90	1	0.10	0.0025



Coplin type staining jar

Jar made of opaque white, unbreakable and **autoclavable** polypropylene, designed for staining 5 standard microscope slides.

Rectangular container with rounded base and screw ribbed flat cap with an inner ring for preventing leakage.

Used as a storage system, it can hold up to 10 slides (2 per groove).

Dimensions: Base diameter: 60 mm.

Height of jar: 110 mm.

Height of jar with cap: 114 mm.

ı	code description		case quantity	case weight	case volume
	191087	Coplin jar for 5 slides	12	0.80	0.006









Slidefolder

Made of hight impact polystyrene. It will hold up to twenty standard slides. The case holding the slides horizontally offers numbered spaces for easy

Will resist temperatures between -80 °C and 100 °C.

Not autoclavable.

The two transparent doors cover the slides and offer a full view of each slide without having to remove them.

Easy opening.

Stackable.

External dimensions: 192 x 292 x 11 mm.

code slidefolder colour		case quantity	case weight	case volume
989945	blue	10	3.20	0.013





Slide tray

Made of high impact polystyrene. It will hold up to twenty standard slides. The slide tray is easily stackable and will take minimum space on any shelf or laboratory counter.

Not autoclavable.

External dimensions: 206 x 299 x 18 mm.

code slidefolder colour		case	case	case
		quantity	weight	volume
989956	white	10	0.91	0.013





See our slides on pages 96-97

Slide trays

Material: polystyrene. Stackable.

code	description	dimensions mm	case quantity	case weight	case volume
989910	slide tray for 10 slides	95x330x50	40	2.50	0.015
989920	slide tray for 20 slides	195x330x50	20	2.44	0.015







Slide mailer

Made of polypropylene.

Designed for slide transport and protection.

To send by mail we recommend to use special envelopes.

Compartment shape ensures tight fit: no risk of damaging slides.

Easy snap-lock. but not indicated for liquid samples.

mod.	code	slide capacity	dimensions mm	case quantity	case weight	case volume
1	19923	1	80 x 40 x 6	50	0.31	0.009
2	900025	2	87 x 47 x 16	100	1.30	0.009
3	979930	2	84 x 71 x 6	250	3.98	0.010
4	19924	3	84 x 99 x 6	10	0.29	0.0009
5	900028	5	82 x 17 x 29	100	1.02	0.012

Slide tube mailer

Tube in polypropilene for mailing, staining and storing of 4 standard slides.

Leakproof cap made of high density polyethylene in pink colour. Incorporates a **tamper evident locking mechanism** that can be

Incorporates a **tamper evident locking mechanism** that can be activated whenever user decides.

The cap has one space to insert a colour disc as identification (40911A, chapter 6. Tubes and microtubes).

Being transparent helps the user to notice if there are slides in it, an advantage in front of the popular white transparent tubes. Internal volume: 12 ml.

Dimensions: 35 x 87 mm.

code	description	case quantity	case weight	case volume
19926	tamper evident tube	100	1.90	0.017







Slide storage boxes

Manufactured in ABS, beige or blue colour. Economical boxes to store or transport 25, 50 or 100 slides. Safe closing and easy opening.

They include numbered sheets with a labeling space in the bottom of the case and in the inner side of the cover (excepting model 1, which does not include the one in the inner side of the cover).

Models 2 and 3 have a special modeling assuring a good stability when stacked.

mod.	code	slide capacity	minimum order quantity	dimensions mm	case quantity	case weight	case volume
1	19276.B	beige slide box for 25 slides	5	93 x 87 x 32	200	10.61	0.07
	19276.A	blue slide box for 25 slides	5	93 x 87 x 32	200	10.61	0.07
2	19277.B	beige slide box for 50 slides	2	200 x 89 x 32	100	11.28	0.07
	19277.A	blue slide box for 50 slides	2	200 x 89 x 32	100	11.28	0.07
3	19278.B	beige slide box for 100 slides	5	200 x 170 x 32	50	11.85	0.07
	19278.A	blue slide box for 100 slides	5	200 x 170 x 32	50	11.85	0.07



Cylindrical slide mailer

Polypropylene. Ideal for storing and transporting 5 to 10 slides in complete safety.

The deep threaded container provides a leakproof seal for transport.

code	description	case quantity	case weight	case volume
19922	for 10 slides	10	0.32	0.0028



Slide storage boxes

Manufactured in polystyrene for 100 slides.

Capacity can be doubled if two slides are placed into the same cell. The internal cork lining prevents damage to stored slides.

To ensure proper slide identification, each slot is numbered to correspond to the slide inventory sheet on the inside cover of the box. Stackable.

This model feature has a safety closure.

Dimensions: 208 x 162 x 32.

code	description	case quantity	case weight	case volume
19278.2	red, for 100 slides	1	0.30	0.002
19278.3	white, for 100 slides	1	0.30	0.002



Microscope slide boxes

Made of high impact polystyrene, those boxes are designed for slides. Each slot is numbered and can be cross-referenced with the corresponding index card supplied with each box.

Boxes are virtually unbreakable. Stackable.

Dimensions mm: 98 x 83 x 38 (25 slides), 230 x 97 x 35 (50 slides), 230 x 180 x 35 (100 slides).

code	description	case quantity	case weight	case volume
19276	for 25 slides	4	0.35	0.002
19277	for 50 slides	1	0.23	0.001
19278	for 100 slides	1	0.39	0.003



Slide boxes

Plain polyethylene slide box with clear polypropylene lid.

Not numbered.

Stackable.

Dimensions in mm: 105 x 88 x 32 (25 slides), 200 x 88 x 32 (50 slides)


cod	le	description	case quantity	case weight	case volume
1927	6.1	for 25 slides	48	3.56	0.019
1927	7.1	for 50 slides	24	3.40	0.018




Automatic slide dispenser

Made of ABS and PS, this dispenser can hold up to 50 standard slides. Slides are dispensed individually by rotating the knurled knob on either side of the moulded base.

code	description	case quantity	case weight	case volume
19540	for 50 slides 26 x 76 mm	1	0.20	0.002

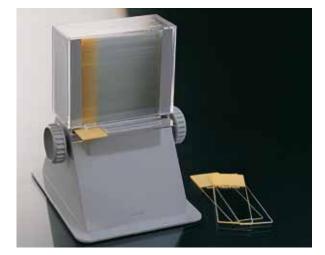



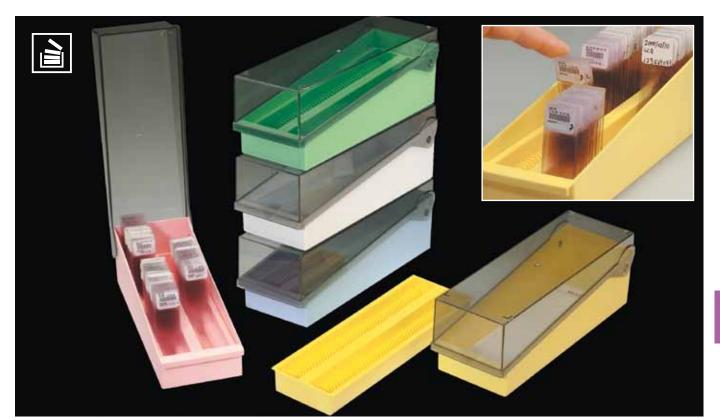

Slide storage system for 100 slides. Stackable

Made of high impact polystyrene. It can hold up to 200 standard slides.

Compound by a box with a transparent lid with an easy opening and an internal removable rack.

This removable rack is divided into two parts with 50 numbered slots in each one, so 100 slides can be separately placed and their strips can be read. For space saving purposes, you can double the amount of slides simply by storing 2 slides per slot, so each part can hold up to 100 slides resulting in a total capacity of 200 slides.

They can also be used without the internal rack, so the capacity reaches 400 slides.

Stackable. The internal rack allows discarding all the slides at the same time.

The frontal part of the box can be written.

External dimensions: 82 x 245 x 86 mm.

code	colour	case quantity	case weight	case volume
19279.5	blue	10	3.3	0.023
19279.7	green	10	3.3	0.023
19279.3	yellow	10	3.3	0.023

See minimum and delivery time for other colors.



Slide storage system for 50 slides. Stackable

Same characteristics as previous codes, it has half capacity. Each part of the inner rack can hold up to 25 separate slides with a total capacity of 50 slides, or for space saving purposes it has a total capacity of 100 slides by storing 2 slides in each slot.

External **dimensions**: 82 x 140 x 86 mm.

code	colour	case quantity	case weight	case volume
19274.5	blue	10	1.5	0.013
19274.7	green	10	1.5	0.013

Minimum order quantity: 1 storage system.







Whirl-Pak® bags with flat wires


Same features as those from the following page but with a flat thicker wire instead of rounded wire resulting in a safer closure, which is preferred in some markets.

With white write-on strip.

Steriles by ethylene oxide.

code	capacity	dimensions	thickness	case	case	case
	ml	cm	microns	quantity	weight	volume
200356	540	11.5 x 23	64	500	2.18	0.0074

Rounded wire
 Flat thicker wire

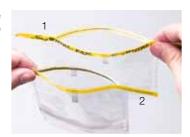





Whirl-Pak® stand-up bags

Same characteristics as Whirl-Pak® standard bags. However, they have a flat bottom which permits them to remain stand up when the sample is inside the bag.

This allows to have both hands free when working so you can add or to extract the content comfortably. No need of any rack.


Working as a beaker or bottle, when they are empty they take up much less space

They have a write-on white strip.

Steriles by ethylene oxide.

code	capacity ml	dimensions cm	thickness microns	case quantity	case weight	case volume
200361	540	11.5 x 23	64	500	2.50	0.017
200365	2,070	19 x 38	102	250	4.52	0.019



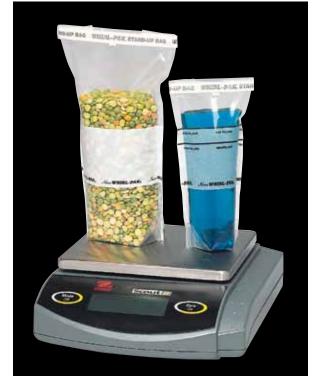













- 1. Tear off top of bag along perforation.
- 2. Use pull tabs to open bag.
- 3. Place sample in bag.
- 4a. Hold bag by wire ends and whirl three times to close.
- 4b. Or, fold the tab over tightly three times to close (used to close larger bags).
- 5. Bend wire ends onto bag.





Sterile Whirl-Pak® sample bags

Made of a low density polyethylene blend, resulting in a very resistant and transparent bags.

Suitable for both solid and liquid samples. The plastic thickness depends on the model (57 to 102 μ).

This special polyethylene blend makes unnecessary a higher thickness to achieve higher resistance.

Leakproof closure by several metallic rounded sticks when you turn closure bands three times (see scheme).

The Whirl-Pak® bags are the unique bags with a one-piece seam, avoiding the possible risks of the loss of the corners bags.

Steriles by ethylene oxide.

Made with materials suitables for alimentary use.

Suitable for liquid nitrogen.

Without strips

code	capacity ml	dimensions cm	thickness microns	case quantity	case weight	case volume
200325	60	7.5 x 12.5	57	500	0.84	0.0048
200340	120	7.5 x 18.5	57	500	1.06	0.0052
200341	210	9.5 x 18	76	500	1.50	0.0053
200342	390	13 x 19	76	500	2.50	0.0170
200329	540	11.5 x 23	64	500	2.04	0.0074
200332	720	15 x 23	76	500	2.94	0.017
200343	720	15 x 23	102	500	3.88	0.017
200345	1,080	12.5 x 38	76	76 500	3.90	0.017
200346	1,260	15 x 38	76	500	3.92	0.017
200347	2,070	19 x 38	76	500	5.62	0.017





Avec bande blanche pour identification

code	capacité ml	dimensions cm	épaisseur micron	quantité par carton	poids du carton	volume du carton
200326	60	7,5 x 12,5	57	500	0,78	0,0048
200349	120	7,5 x 18,5	57	500	1,02	0,0054
200364	390	13 x 19	64	500	2,02	0,017
200330	540	11,5 x 23	64	500	1,66	0,0071
200333	720	15 x 23	76	500	3,06	0,0170
200351	1.650	19 x 30	102	500	5,66	0,0182
200363	2.070	19 x 38	76	500	5,54	0,019
200357	2.700	25,4 x 38	102	250	5,32	0,022
200358	3.600	25,4 x 50,8	102	250	6,46	0,022
200359	5.400	38 x 50,8	102	100	4,68	0,022

La capacité des sachets a été calculée quand le sachet est fermé en pliant les tiges de fermeture métallique 3 fois. Les volumes sont approximatifs. Pas utiliser à une température supérieure à 82 °C.

Ils ne sont pas autoclavables.



See blender bags in Chapter 1. Microbiology



Zip-lock bags

Made of transparent polyethylene, very resistant to shocks and tearing. For all solid samples.

Suitable for laboratories (for tubes, swabs, etc.) as well as industries (to contain, spare pieces, small electronic components, buttons, sweets, etc.).

code	gauge	dimensions mm	case quantity	case weight	case volume
M 5555	200	55 x 55	2,000	1.13	0.0028
M 6080	200	60 x 80	2,000	1.35	0.0046
M 70100	200	70 x 100	2,000	1.57	0.0046
M 80120	200	80 x 120	2,000	2.54	0.0051
M 80160	200	80 x 160	2,000	2.80	0.0077
M 100150	200	100 x 150	2,000	3.90	0.0080
M 110110	200	110 x 110	2,000	2.90	0.0080
M 120180	200	120 x 180	2,000	4.56	0.0125
M 150220	200	150 x 220	2,000	6.87	0.0180
M 180250	200	180 x 250	2,000	9.16	0.0230
M 200300	200	200 x 300	2,000	13.94	0.0220
M 250330	200	250 x 330	2,000	16.66	0.0340
M 300400	200	300 x 400	2,000	23.10	0.0500





White strip zip-lock bags

Same features as the above products.

Feature white strip for writing for better sample identification.

code	gauge	dimensions mm	case quantity	case weight	case volume
M 5555B	200	55 x 55	2,000	1.13	0.0028
M 6080B	200	60 x 80	1,000	0.67	0.0026
M 80120B	200	80 x 120	1,000	1.27	0.0037
M 70110B	200	70 x 110	1,000	0.95	0.0037
M 100150B	200	100 x 150	1,000	1.59	0.0044
M 120180B	200	120 x 180	1,000	2.35	0.0065
M 160220B	200	160 x 220	1,000	4.32	0.0080
M 180250B	200	180 x 250	1,000	4.43	0.0125
M 200300B	200	200 x 300	1,000	6.82	0.0145



Specimen transport kangaroo bags

For tube transport. Double pouch design: one for specimen tubes and the other for documents. Includes a safety seal.

Reclosable (zip-lock design). Dimensions: 16x16.5 cm.


Made of polyethylene.

code	colour	case quantity	cases weight	case volume
V 160230	green	2 x 1,000	11.66	0.03
N 160230	orange	2 x 1,000	11.66	0.03





See blender bags in Chapter 1. Microbiology

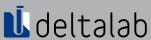



Whirl-Pak® specimen transport kangaroo bag

Non-sterile bag made of a low density polyethylene blend, resulting in a **very resistant and transparent bag.**

Designed for both solid and liquid samples.

Leakproof closure with several metallic rounded sticks.

Double pouch design: one for specimen tubes and the other one for documents.

It features the biohazard symbol.

Bag with a **one-piece seam**, avoiding the possible risks of the loss of the corners bags. **Suitable for liquid nitrogen**.

code	capacity ml	dimensions cm	thickness microns	case quantity	case weight	case volume
200372	720	15 x 23	64	500	3.38	0.0170





Mailing container

Mailing container with safety screw blue cap, both manufactured in polyethylene. Leak proof.

Dimensions uncapped 128×30 mm. Ideal for 10 ml tubes. The tube includes a piece of absorbent paper to prevent from any leakage.

Container and cap are sold separately.

code	description	case quantity	case weight	case volume
401301	01301 mailing container		7.70	0.09
401302	screw cap	500	1.90	0.009



Forceps

Made of **POM** (polyoxymethylene).


Colour: orange.

Self-sprung forceps with rounded ends.

Code **19503**, the inner part of the tip is ribbed to enhance sampling. Ribbed for a better handling.

Autoclavable. Flexible.

code	length mm	case quantity	case weight	case volume
19500	115	5	0.23	0.0006
19501	145	5	0.05	0.0005
19503	250	5	0.17	0.004

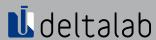

Slide mailers

Designed for slide transport and mailing.

Different models (tube, case, etc.) for one, two, three, four or even ten slides.



See more information in Chapter 4. **Histology, microscopy and staining**

Containers for sampling

Polystyrene or polypropylene containers. Volumes from 30 ml to 1,000 ml.



See more information in Chapter 7. **Sample Containers**



Exam gloves

(€ (EPI) **(**€ (MDD)

Wide range of single use gloves AQL 1.5.

Different options available:

- Powdered and non powdered latex.
- Powdered and non powdered vinyl.
- Non powdered nitrile.



100

EUROTUBO

See more information in Chapter 11. **Safety and general labware**

Swabs for microbiological sampling

C€ (MDD)

Sterile swabs with or without transport media. Tips made of viscose, cotton, dacron®, etc.

Supplied in tube, peel pack, flow-pack, etc. Transport media available:

- · Amies with and without charcoal
- Stuart
- Cary Blair
- Virus
- Chlamydia

See more information in Chapter 1. **Microbiology**



Autoclave resistant bags

Printed with the biohazard symbol and instructions of use. Two versions available:

- Transparent and translucent, suitable for routine autoclaving.
- Red color, suitable for autoclaving at 138 °C; they feature a sterilization indicator.






See more information in Chapter 11. **Safety and general labware**







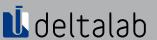





Cytology brush

Blue ABS shaft. For cervical specimen collection, designed for taking cell samples without damaging them.

Includes a soft tip to provide more comfort for patients.

Total length: 19.5 cm.

Sterile model is sterilised by ethylene oxide and includes instructions.

code	(description	case quantity	case weight	case volume
440150		non sterile	50 x 50	4.55	0.009
440151 STERILE EO individual		individual peel-pack	2 x 500	3.90	0.029



Vaginal speculum

(€ (MDD)

Manufactured in transparent polystyrene. Handle design that allows one-handed quick and convenient manner. Zip system for opening and closing. Aseptic. Rounded edges and ergonomically designed for single-handed use. Free of heavy materials as CE/94/62 directive.

Packaged individually with a peel-pack that facilitates the traceability and also indicates expiration date.

Available in 3 differents sizes.

code	shaft width mm	case quantity	case weight	case volume
440157	20	100	3.90	0.060
440158	30	100	4.22	0.060
440159	32	100	4.62	0.060



AYRE wood spatula

C€ (MDD)

Rounded edges.

Sterile model is sterilised by ethylene oxide and includes instructions.

code	C	description	case quantity	case weight	case volume
440142.0	2.0 non sterile		50 x 100	10.50	0.027
444150 STERILE EO individual peel-pac		individual peel-pack	2 x 500	5.10	0.029



AYRE plastic spatula

C€ (MDD)

Made of **autoclavable** polypropylene.

Rounded edges.

Sterile model is sterilised by ethylene oxide and include instructions.

code	C	description ca		case weight	case volume
440142	non sterile		5 x 500	7.0	0.025
444142	444142 STERILE EO individual peel-pack 444242 STERILE R individual flow-pack		2 x 500	4.6	0.029
444242			2 x 500	4.9	0.029

Minimum ordening quantity: 500





Isothermal bags



Isothermal bags for the transport of biological samples, manufactured from PVC, PE and polyurethane (foam). The isothermal bags have two functions: keeping the temperature of samples for 6 hours while being transported, whatever the climatic conditions; protecting the samples and their contents. Available in green.

Description: for added security, isothermal bags feature three different locks: by velcro (excepting 6 I model), by zip fastener and by a coded padlock.

Bags include the following elements:

- An inner pocket designed to hold the cold packs,
- A transparent outer pocket, at the front of the bag, for laboratory identification,
- A transparent pocket, over the top of the bag, to insert documents.
- Two cold packs (except **900073** which incorporates only 1) consisting of HDPE packs filled with a non toxic gel.

Before use, place cold packs in the freezer and keep at $-20\,^{\circ}\text{C}/-24\,^{\circ}\text{C}$ for 24 hours.

Quality features:

- Bags manufactured from sturdy and shockproof materials.
- Stand out from other isothermal bags existing on the market for their lightness and convenience.
- Offer a total security seal thanks to padlock.
- Space-saving and supplied folded in to reduce volume by 50 %.
- To clean, simply use standard detergents.
- Corners are rounded and reinforced by a double seam. Handles are sewn all around the bag for increased strength.
- Isothermal bags are traditionally manufactured (non-automated production), and each product undergoes rigid quality controls.

code	dimensions mm	capacity	case quantity	case weight	case volume
900073	900073 230x230x170		1	0.62	0.003
900074	220x220x260	10	1	1.10	0.010
900076	390x250x280	27	1	1.33	0.011
900078	450x300x300	40	1	1.46	0.014

Bottles for water sampling

Sterile bottles, designed according to ISO standards for sampling of consumption water.

Available with thiosulfate (microbiological analysis) and without thiosulfate (physico-chemical analysis).


Bottles are available in polyethylene or PET. Capacities vary from 500 to 1,500 ml.

Special configurations:

- Personalized labels,
- Individually bagged bottles,
- Special dosifications of thiosulfate,
- Other bottles or jars from the catalogue, etc.





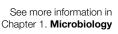





C€ (IVD)

High density polyethylene container and cap. Ideal for 24 hour urine collection and other liquid samples.

The ribbed cap features a pharmaceutical grade liner, which eliminates the need for an insert plug, allows a better handling and provides a watertight seal.

Containers are graduated in 50 ml increments up to 2.7 l.

Ergonomically designed for single-handed use. Stackable.

A ring on the neck of the containers avoids dripping when pouring.

Lenght x width: $130 \times 126 \text{ mm}$ ($\pm 1 \text{ mm}$).

Height, cap excluded: 249 mm.

Individually wrapped models are supplied capped.

Cases per palet: 12.

code	body colour	cap description	presentation	individually wrapped	sterile	internal mouth Ø mm	case quantity	case weight	case volume
407005	natural	white, with liner	unscrewed	no	no	64	34	5.70	0.152
407005/T	natural	white, with liner	screwed	no	no	64	34	5.70	0.152
407005.O	natural	white, with liner	screwed	yes	RADIATION	64	34	5.70	0.152
407006	brown	white, with liner	unscrewed	no	no	64	34	5.70	0.152
407006/T	brown	white, with liner	screwed	no	no	64	34	5.70	0.152



5 and 10 ml freestanding, conical bottom tubes

Freestanding tubes with screw caps. Conceived for sampling, transport and midterm storage. Feature molded graduations at 1mL increments and a large surface for labeling and identification. Round bottom tubes made of **autoclavable** polypropylene, leakproof caps made of polyethylene.

Withstand temperatures between 122 °C and –80 °C (–40 °C as the lowest working temperature).

Certified to the IATA 95kPA standard for safe shipping of samples.


code	description	case quantity	case weight	case volume
439905	5 ml PP thread tube	2 x 500	5.45	0.032
439910	10 ml PP thread tube	2 x 500	6.35	0.055


Dimensions (±0,5):

code	external cap Ø mm	internal tube Ø mm	length with cap mm	length without cap mm
439905	15.8	14.17	60.20	58.62
439910	15.8	14.17	94.74	92.96





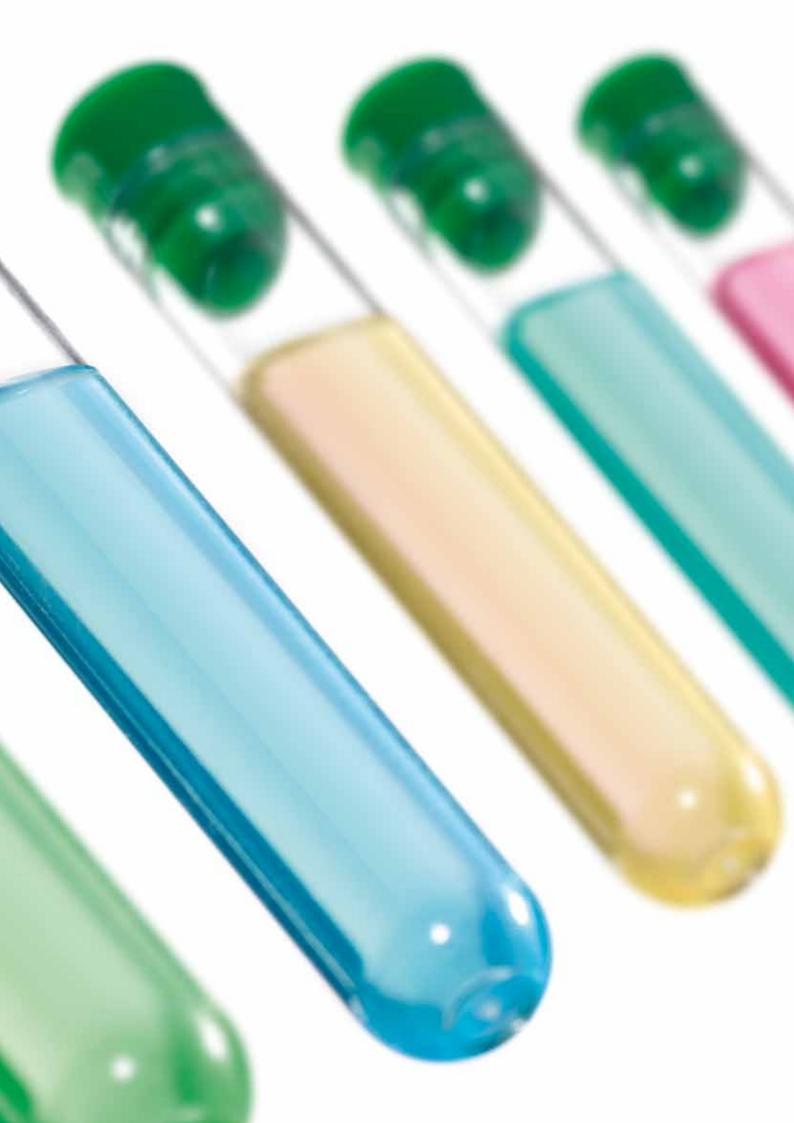

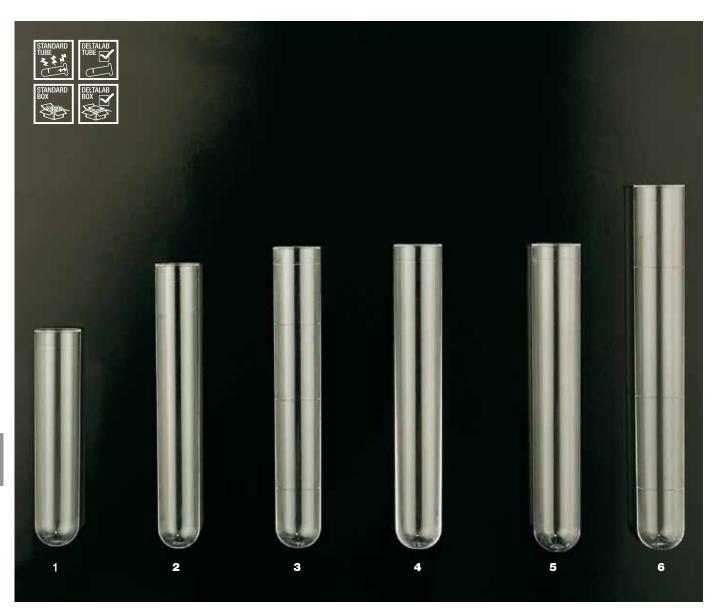

Round bottom polystyrene tubes

C€ (IVD)

Manufactured from high technology, new generation moulds to ensure the reproducibility of each tube. Tubes are made in one piece to ensure uniformity and dimensional accuracy from tube to tube.

mod.	code	dimensions mm	volume ml	case quantity	case weight	case volume	cases per pallet
1	300500	11 x 55	3	2 x 2,000	7.50	0.046	32
2	300300	11 x 70	4	3 x 1,000	6.20	0.046	32
3	300800	12 x 75	5	4 x 1,000	9.30	0.075	20
4	300800.2	12 x 75 thick wall	5	4 x 1,000	10.70	0.075	20
5	300800.1	12 x 75	5	4 x 1,000	9.60	0.075	20
6	300400	12 x 88	6	4 x 1,000	11.30	0.082	20

Codes **300800** and **300400** have graduation rings. Code **300800** (mod. 3) is graduated at 1, 2.5 and 4 ml; code **300400** (mod. 6) is graduated at 1, 2.5 and 5 ml.


Dimensions (±0,09)

code	internal diameter mm	external diameter mm	total length mm	maximum volume ml
300500	10 11.5		54.4	3.5
300300	9.1	10.8	69.4	4.0
300800	10.3	11.8	74.5	5.5
300800.1	10.2	11.8	74.5	5.5
300800.2	10	11.8	74.5	5.3
300400	10.4	12	87.6	5.9

Suitable for centrifugation up to 7,500 xg.

A complete range of test tube caps is detailed on page 125 of this catalogue.

All tubes are suitable for laboratory use.

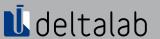



Round bottom polystyrene tubes

C€ (IVD)

Manufactured from high technology, new generation moulds to ensure the reproducibility of each tube. Tubes are made in one piece to ensure uniformity and dimensional accuracy from tube to tube.

mod.	code	dimensions mm	volume ml	case quantity	case weight	case volume	cases per pallet
1	300700	13 x 75	5	4 x 1,000	12.50	0.082	20
2	301700	13 x 100	7	6 x 500	11.60	0.082	20
3	300900	16 x 95	10	5 x 500	12.00	0.096	16
4	300705	16 x 100	10	5 x 500	14.80	0.110	16
5	300100	16 x 150	17	1,000	9,84	0.093	20

Code **300900** (mod. 3) has graduation rings at 2.5, 5 and 10 ml.

Dimensions (±0,09)

internal diameter mm	external diameter mm	total length mm	maximum volume ml	
10.9	12.9	74.7	6.7	
11.0	12.9	99.6	8.3	
13.8	15.8	94.6	12.2	
13.8	16.0	99.8	13.2	
14.6	16.3	151.0	21.1	
	10.9 11.0 13.8 13.8	diameter mm diameter mm 10.9 12.9 11.0 12.9 13.8 15.8 13.8 16.0	diameter mm diameter mm length mm 10.9 12.9 74.7 11.0 12.9 99.6 13.8 15.8 94.6 13.8 16.0 99.8	

Suitable for centrifugation up to 7,500 xg. (excepting code 300100, which does not withstand centrifugation). A complete range of test tube caps is detailed on page 125 of this catalogue.

All tubes are suitable for laboratory use.





Ultra clear polypropylene tubes

C€ (IVD)

Manufactured from high technology, new generation moulds to ensure the reproducibility of each tube. Tubes are made in one piece to ensure uniformity and dimensional accuracy from tube to tube. **Autoclavable** up to **121 °C**.

mod.	code	dimensions mm	volume ml	case quantity	case weight	case volume	cases per pallet
1	400500	11 x 55	3	2 x 2,000	6.50	0.046	32
2	400800	12 x 75	5	4 x 1,000	8.00	0.075	20
3	400800.1	12 x 75	5	4 x 1,000	8.00	0.075	20
4	400700	13 x 75	5	4 x 1,000	9.70	0.082	20
5	400400	12 x 88	6	4 x 1,000	9.60	0.082	20
6	401700	13 x 100	7	6 x 500	9.40	0.082	20
7	400900	16 x 95	10	5 x 500	10.10	0.096	16
8	400705	16 x 100	10	5 x 500	12.20	0.096	16
9	401100	15 x 50	5	5 x 1,000	10.40	0.082	20

Certain codes have graduation rings: **400800** (mod. 2): graduated at 1, 2.5 and 4.0 ml; code **400400** (mod. 5): graduated at 1, 2.5 and 5.0 ml; **400900** (mod. 7): graduated at 2.5, 5 and 10 ml; **401100** (mod. 9): graduated at 2.5 and 5 ml.

Dimensions (±0,09)

code	internal diameter mm	external diameter mm	total length mm	maximum volume ml
400500	9.9	11.4	53.9	3.4
400800	10.3	11.7	73.9	5.3
400800.1	10.2	11.7	73.9	5.3
400700	10.9	12.7	73.7	6.4
400400	10.3	11.9	86.9	5.8
401700	11.0	12.7	98.4	8.1
400900	13.7	15.5	93.7	11.9
400705	13.7	15.8	98.6	12.7
401100	13.7	15.5	48.9	6.5

Suitable for centrifugation up to 15,000 xg.

A complete range of test tube caps is detailed on page 125 of this catalogue.

All tubes are suitable for laboratory use.





Polystyrene tubes capped

Suitable for centrifugation up to 7,500 xg.

code	dimensions tube mm	height tube + cap mm	volume ml	case quantity	case weight	case volume	cases pallet
300702	13 x 75	83.4	5	2 x 1.000	7.5	0.046	32
300903	16 x 95	105.2	10	2 x 500	7.0	0.110	40
300907	16 x 100	109.6	11	4 x 500	15.0	0.110	16

Components:

code	tube		cap
300702	300700	+	307107
300903	310900	+	318107
300907	300705	+	308107



See Chapter 5. **Sampling and Transport.**

Polystyrene tubes, capped and labelled

(€ (IVD)

Suitable for centrifugation up to 3,000 xg.

Dimensions of the label: 35×25 mm.

code	dimensions tube mm	height tube + cap mm	volume ml	case quantity	case weight	case volume	cases pallet
300804	12 x 75	82.9	5	2 x 1,000	7.3	0.046	40
300704	13 x 75	83.4	5	2 x 1,000	7.5	0.046	40
300904	16 x 95	105.0	10	2 x 500	7.1	0.046	48
300908	16 x 100	110.4	11	4 x 500	14.0	0.096	16

Components:

code	tube		cap
300804	300800	+	305807
300704	300700	+	307107
300904	300900	+	308107
300908	300705	+	308107





Conical tubes (€ (IVD)

Manufactured in polypropylene (PP) and polystyrene (PS).

mod.		code	dimensions mm	material	volume ml	case quantity	case weight	case volume	cases per pallet
1	_	301200	16 x 102	PS	12	5 x 500	12.9	0.096	16
2	2 <u>₹</u>	401200	16 x 102	PP	12	5 x 500	10.6	0.096	16
3	_	301213	17 x 105	PS	12	1,500	8.80	0.088	20
4	_	301212	17 x 105	PS	12	1,500	8.80	0.088	20
5	_	301201	16 x 100	PS	12	5 x 500	12.5	0.096	16
6	2 <u>₹</u>	401201	16 x 100	PP	12	5 x 500	10.2	0.096	16
7	_	301202	16 x 102	PS	12	5 x 500	12.9	0.096	16
8	2 <u>₹</u> 20	401204	16 x 100	PP	12	5 x 500	10.2	0.096	16

Codes **301200** (mod. 1), **401200** (mod. 2) and **401204** (mod. 8) have internal graduations.

Codes 301213 (mod. 3) and 301212 (mod. 4) are flared at their top for easy pouring and have a special sedimentation base.

Codes **301200** (mod. 1), **401200** (mod. 2) and **301202** (mod. 7) have an external rim.

Dimensions (±0,09)

code	internal diameter mm	external diameter mm	total	maximum volume ml
004000			length mm	
• 301200	14.5	19.1/16.4	100.8	12.9
• 401200	14.4	18.8/16.2	99.8	12.5
• 301213	15.4/19.2	17.1/19.2	105.1	16.5
• 301212	15.4/19.2	17.1/22.0	105.0	16.4
301201	14.5	16.4	99.5	12.9
401201	14.4	16.2	98.4	12.5
• 301202	14.5	19.1/16.4	100.8	12.9
401204	14.4	16.2	98.4	12.5

• External diameter of the rim.

All models are suitable for centrifugation. PS tubes: 7,500 xg; PP tubes: 15,000 xg. Except code 301212 and 301213 to 3,000 xg.





Polyethylene caps

mad	códe	to fit tubes	bag	bag	bag	yellow	natural	malve	blue	green	red	black
mod.	code	to iii tubes	quantity	weight	volume			(to add at the end of the code)				
1	3053	300500, 400500, 300800.2	2 x 2,000	1.90	0.008		02	08		07	06	
2	3003	300300	5 x 1,000	3.40	0.016		02					
3	3058	300800, 400800, 300800.1, 300400, 400400, 400800.1	2,000	1.90	0.008		02	08	05	07	06	09
4	3071	300700, 301700, 400700, 401700, 300701, serial 300700.6	2,000	1.90	0.007	01	02	08	05	07	06	09
5	3072	300705, 300900, 400705, 400900, 401100, 300707	2,500	3.20	0.016		02		05	07	06	
6	3050	300705, 300900, 400705, 400900, 401100, 300907, 300904, 300908	2,500	4.10	0.020		02	08		07		
7	3081	300705, 300900, 400705, 400900, 401100, 300707, 300907, 300904, 300908	2,500	5.00	0.030		02			07	06	
8	3070	300100	1,000	1.20	0.010						06	
9	3055	301201, 401201, 301200, 401200, 301202, 401204	2,500	4.00	0.020		02			07		
10	3082	301201, 401201, 401204	2,500	5.10	0.028		02		05			
11	3066	301213, 301212	1,500	1.15	0.008	01						

Ask for other colours.



Re-caps

Ideal for re-capping blood collection tubes «vacuum type». Suitable for other glass and plastic tubes.


Manufactured in flexible polyethylene.

Easy to cap and remove.

code	description	case quantity	case weight	case volume
308502	Ø 12-13 mm natural colour	6 x 1,000	1.70	0.025
308506	Ø 12-13 mm red colour	6 x 1,000	2.18	0.025
308602	Ø 15,2-16 mm natural colour	6 x 1,000	3.31	0.025
308606	Ø 15,2-16 mm red colour	6 x 1,000	3.31	0.025

Other colours available on request, minimum order quantity is 6,000 units.



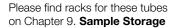



Internal diameter: 10.00 mm External diameter: 11.85 mm Total Length: 11.30 mm





Includes a small internal chamber. Light-weighted. Ergonomical. Ribbed for a better handling while using gloves. Cap includes two internal ribs to clamp firmly on the tube. Leakproof.

3003

Internal diameter: 9.25 mm External diameter: 11.35 mm Total length: 16.60 mm





Flat base. No chamber. Ribbed.

Cap includes three internal ribs to clamp firmly on the tube. Leakproof.





3058

Internal diameter: 10.60 mm External diameter: 15.00 mm Total length: 15.30 mm





Includes a small internal cavity designed to prevent possible drop spillage. Ergonomical.

Ribbed for a better handling while using gloves. Leakproof. Cap includes two ribs to clamp firmly on the tube. Top of the cap raised for easy-on, easy-off use and to avoid cap deformation from repeated removal.



Those three flange plug caps present significant advantages over other flange caps existing on the market, which are:

- a larger size, providing a larger surface for finger rest;
- flanges have a larger size than on other caps. More space between each flange increase the cap sealing qualities;
- a small internal cavity designed to prevent possible drop spillage;
- a raised top for easy-on, easy-off use and to avoid cap deformation from repeated removal;
- a ribbed surface;

Note on flange cap dimensions: using flange caps implies reducing tube capacity. If as much tube volume as possible is required, we recommend the use of those caps including a chamber, such as codes **3082**, **3050** and **3081**. To avoid possible liquid spillage on the benchtop, leave caps stand with flanges up.

code	internal Ø mm	external Ø mm	total length mm
3072	15.40	17.00	21.60
3055	14.50	17.05	21.77
3071	11.25	13.30	17.40











Internal diameter: 13.95 mm External diameter: 19.30 mm Total length: 13.15 mm





Includes an internal chamber for more liquid filling, thus compensating for the loss of space inside the tube caused by the insertion of a cap. Leakproof closing. Ergonomical. Ribbed surface. Cap shape is slightly conical for a better handling.

Also features an external flap so as to:

- avoid contact between liquid and fingers when opening the tube,
- · prevent dripping when cap is left on the benchtop,
- provide easy single-handed opening,
- offer a cap specially designed to fit our range of blood collection tubes, diameter 15 mm.



3081

Internal diameter: 13.80 mm External diameter: 19.88 mm

Total length: 19 mm






One of the caps offering the largest chamber on the market, thus compensating for the loss of space inside the tube caused by the insertion of a cap. Leakproof closing.

This cap allows an almost complete filling of the tube and eliminates the need for a fill-in line.

Ribbed for a better handling while using gloves.

The inner part of the cap, long-sized, consists of three rings to ensure a tight, safe leakproof closing. Top of the cap raised for easy-on, easy-off use and to avoid cap deformation from repeated removal. Ideal for use with urine.

To avoid possible liquid spillage on the benchtop, leave caps stand with flanges up.

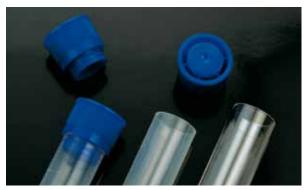

3082

Internal diameter: 14.55 mm External diameter: 20.45 mm Total length: 18.50 mm





It has the same characteristics than the previous cap.



3070

Internal diameter: 14.80 mm External diameter: 19.20 mm

Total length: 25.35 mm





Ribbed for a better handling while using gloves.

No chamber. Leakproof. Includes four internal rims to clamp firmly on tube. Reversible cap: use it as a watertight seal, or invert it for use in its "re-cap" position, as a simple closure, when no leakproof seal is required.



3066

External diameter: 23.30 mm

Total length: 5.90 mm

Light-weighed. Includes a lifting tab for easy opening.







Sterile culture tubes in polystyrene

Tubes supplied with either a two position ribbed polyethylene cap, which can be left loose for aerobic work or sealed for anaerobic cultures.

They are biologically inert, exempt from mold release agents, and withstand up to **1,400 xg** and **70 °C**.

Packaged in self-standing resealable zip-lock bags of 125 units.

code	dimensions mm	case quantity	case weight	case volume
300807	12 x 75	8 x 125	4.24	0.033
300808	17 x 100	8 x 125	7.14	0.060







PosSealed position for anaerobic cultures





EUROTUBO® 12 ml screw cap tubes, round bottom

C€ (IVD)

Made of autoclavable (121 °C) transparent polypropylene or polystyrene. Green cap made of high density polyethylene. Dimensions: 15 x 102 mm. The external skirt allows the tubes to remain free-standing. Supplied screwed. Recommended volume: 12 ml. The sterile model (ethylene oxide) is supplied individually in flow-pack bag, with indication of batch, expiration date, etc. According to the guidelines for sterile products. Resistance to centrifugation: PS: 7,500 xg. PP: 15,000 xg.

Attention: For autoclaving, the cap should be loose on the thread and not tightly fitted.

	mod.	code	description	sterile	maximum volume ml	volume recom. ml	case quantity	case weight	case volume	cases pallet
	1	301403	polystyrene	no	14.4	12	6 x 250	12.66	0.082	20
	2	301402	polystyrene individually wrapped	STERILE EO	14.4	12	6 x 250	13.00	0.096	16
21.0 121.0	3	401402	polypropylene individually wrapped	STERILE EO	14.2	12	6 x 250	11.80	0.096	16
121°C	4	401403	high transparency polypropylene	no	14.2	12	6 x 250	11.00	0.082	20

Dimensions (±0,09):

code	external cap Ø mm A	external tube Ø mm B	internal tube Ø mm C	length with cap mm D	length without cap mm E	D — D	E
301402, 301403	20.9	16.3	14.4	103.9	102.5	A	\\\\ ↑c в
401402, 401403	20.9	16.2	14.3	102.9	101.5		<u> </u>





15 ml and 50 ml centrifugal tubes

C€ (IVD)

 $\label{eq:made_of_polypropylene} \mbox{Made of polypropylene, suitable for both clinical and research applications.}$

DNAsa, RNAsa and pyrogen free. They are also free from natural rubber and heavy metals.

High transparency of the material for a clear visualization during experiments specially for molecular biology and animal tissue culture.

Tube and cap designed with the system of flat threads for a complete leakproof.

Highly smooth hydrophobic surface for minimum disturbance during centrifugation.

Silk-screen blue graduation in the tube and large white frosted portion for easy writing.

Autoclavable at 121 °C.

Centrifugation resistance: 14.000 xg.

Available models: 15 ml non-skirted and 50 ml skirted and non-skirted.

code	description	presentation	sterile	case weight	case volume	cases per pallet
15 ml tubes						
429940	non-skirted tube	500 tubes in bulk	NO	500	4.5	0.0264
429945	non-skirted tube	20 bag to 25 tubes.	NO	500	4.6	0.0264
429942	non-skirted tube	20 bag to 25 tubes.	STERILE R	500	4.6	0.0264
50 ml tubes						
429930	non-skirted tube	20 bag to 25 tubes.	NO	500	8	0.0663
429931	non-skirted tube	20 bag to 25 tubes.	STERILE R	500	8	0.0663
429950	skirted tube	20 bag to 25 tubes.	NO	500	9	0.0663
429951	skirted tube	20 bag to 25 tubes.	STERILE R	500	9	0.0663

5 and 10 ml freestanding, conical bottom tubes

Freestanding tubes with screw caps. Conceived for sampling, transport and midterm storage. Feature molded graduations at 1 ml increments and a large surface for labeling and identification. Round bottom tubes made of **autoclavable** polypropylene, leakproof caps made of polyethylene.

Withstand temperatures between 121 °C and -80 °C (-40 °C as the lowest working temperature).

Certified to the IATA 95kPA standard for safe shipping of samples.

code	description	case quantity	case weight	case volume
439905	5 ml PP thread tube	2 x 500	5.45	0.032
439910	10 ml PP thread tube	2 x 500	6.35	0.055

Dimensions (±0,5):

code	external cap Ø mm	internal tube Ø mm	length with cap mm	length without cap mm
439905	15.8	14.17	60.20	58.62
439910	15.8	14.17	94.74	92.96





15 ml and 50 ml conical high resistance tubes. Sterile

C€ (IVD)

Tubes made of transparent, copolymer polypropylene.

Green caps made of polyethylene with an internal liner which ensures leakproofness.

Tubes feature a solvent resistant white panel and black graduations for use both for clear or dark samples.

DNAse, RNAse, endotoxins and metal free.

They are sterile by radiation and withstand temperatures down to **-90 °C** (15 ml) and **-80 °C** (50 ml)...

Resistance to centrifugation: 17,000 xg (15 ml) and 20.000 xg (50 ml).

Autoclavable tube (the liner of the cap is not autoclavable)

Supplied in bags of 50 units.

Bags are printed with instructions for use and feature a double closure: a first tamper evident seal that helps guaranteeing sterility, and a secondary zip-lock, resealable closure.

	code	presentation	capacity ml	dimensions mm	case quantity	case weight	case volume
1	409920	bag 50 units	15	17 x 118	10 x 50	4.12	0.035
2	409922	bag 50 units	50	29.6 x 114.6	10 x 50	8.16	0.076





EUROTUBO® 15 ml and 50 ml conical tubes

C€ (IVD)

Tubes made of transparent polypropylene, conical bottom. Suitable for centrifugation tests in immunology laboratories, microbiology, etc. With continuous thread and molded external graduation each 0.5 ml in 15 ml tubes, and each 5 ml in 50 ml tubes.

Blue polyethylene cap with internal elastic obturation to assure an hermetical closure.

Resistance to centrifugation:

15 ml tubes, non sterile models until 7.000 xg, sterile models 5.000 xg.

50 ml tubes, non sterile models until 12.000 xg, sterile models 7.000 xg.

Sterile models are not suitable for autoclaving.


Non sterile models are autoclavable at 121 °C, with the cap positioned on the thread but uncapped.

50 ml tubes have skirted and non skirted models.

	code	description	sterile	external cap Ø mm	internal tube Ø mm	external tube Ø mm	length with cap mm	case quantity	case weight	case volume	cases per pallet
1	5 ml conica	l tubes									
	429910	non sterile	no	20.9	14.35	17	120	500	3.8	0.029	54
	429920	100 units in bag	STERILE R	20.9	14.35	17	120	5 x 100	3.9	0.029	54
	429946	individualy bagged	STERILE R	20.9	14.35	17	120	500	3.9	0.038	40
5	i0 ml conica	l tubes									
	429900	non skirted	no	34.4	27.2	29.5	117.5	500	7.9	0.075	20
	429900SP	non skirted, unscrewed cap	no	34.4	27.2	29.5	117.5	500	7.9	0.082	20
	429901	skirted	no	34.4	27.2	29.5	117.5	500	8.5	0.082	20
	429926	non skirted, individually wrapped	STERILE R	34.4	27.2	29.5	117.5	500	8.2	0.082	20
	429926.25	non skirted, bag 25 units	STERILE R	34.4	27.2	29.5	117.5	20x25	8.2	0.082	20
	429926.10	non skirted, bag 100 units	STERILE R	34.4	27.2	29.5	117.5	5x100	7.9	0.082	20
	429927	skirted, individually wrapped	STERILE R	34.4	27.2	29.5	117.5	500	8.8	0.082	20
(NaW)	429927.25	skirted, bag 25 units	STERILE R	34.4	27.2	29.5	117.5	20x25	8.8	0.082	20
	429927.10	skirted, bag 100 units	STERILE R	34.4	27.2	29.5	117.5	5x100	8.5	0.082	20





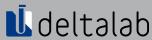

GLASS TUBE RANGE

The high quality of our range of disposable glass culture tubes and screw threaded tubes is reflected in the uniformity of their wallthickness, shape of their bottoms, close tolerances, well formed bottoms and well moulded screw-neck.

All of our tubes are fully annealed at temperatures above **500 °C** to eliminate contamination, reduce breakage and increase strength.

The majority of our tubes are packaged in plain outer cases.

Each tube is polished individually.

Culture tubes:

They are available in either borosilicate or soda lime glass. Both glasses differ primarily in their chemical composition and their chemical resistance.

Borosilicate glass is a type of glass with a higher resistance to thermal shock. For this reason borosilicate glass is more suitable when tubes are intended to be heated or cooled rapidly.

Soda lime glass is a more economical alternative and suits for most culture medium, bacteriological and laboratory applications where the tube will not be subjected to extreme temperature changes.


Screw-neck tubes are ideal for tissue culture and general bacteriological work. Tubes are supplied with standard GPI finishes.

The screw – neck vials are made of borosilicate glass and are basically used for diagnosis, pharmaceutics, chemicals, laboratory analysis, storage of liquids and powders. The caps presented are manufactured in autoclavable polypropylene with internal silicone joint.

Different packagings are available, depending on the model of tube:

Boxes: sturdy cardboard boxes. Each box contains small shrink-wrapped cases. Tubes are presented horizontally.

Trays: sturdy shrink-wrapped cardboard trays. Tubes are presented vertically.

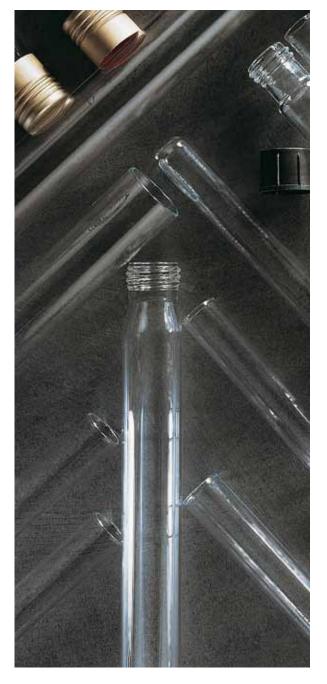


Culture tube cap

Made of polypropylene, **autoclavable** up to **121 °C**. Suitable with codes **916100, 816100, 916150, 816150** and **816160**. Vented caps (with inner flats). Recommended for glass tubes.

code	recommended for tubes Ø mm	colour	case quantity	case weight	case volume
617402	16	natural	4 x 1,000	5.3	0.08



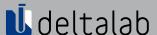




Round bottom glass tubes

Made of **borosilicate** or **soda** glass.

The high quality of those tubes is reflected in the uniformity of their wall thickness and of their diameter and height dimensions. Supplied in small quantities per case for a more convenient use in laboratory.

code boro	code soda	total capacity ml	Ø int. mm tube	Ø ext. mm tube	height mm	thickness mm	case quantity	case weight	case volume	
Supplied in boxes (1)										
901075	801075	4	8.20	9.75	75	0.60	4 x 250	3.60	0.010	
901275	801275	6	10.20	11.60	75	0.60	4 x 250	4.50	0.014	
913100	813100	10	11.10	12.70	100	0.60	4 x 250	5.90	0.022	
916100	816100	15	13.95	15.75	100	0.60	4 x 250	9.10	0.034	
916150	816150	22	13.55	16.00	150	0.70	4 x 250	13.60	0.049	
	816160	27	14.40	16.00	160	0.55	500	5.50	0.018	
918150	818150	28	15.00	18.00	150	0.85	2 x 250	7.30	0.030	
	820150	34	17.20	20.00	150	0.85	100	1.92	0.006	
	820200	47	17.15	19.25	200	0.85	250	6.30	0.020	
Supplied in tra	ays (2)									
	801175T	6	10.10	11.60	75	0.50	500	1.89	0.005	



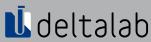



Screw threaded glass tubes

Made of **borosilicate** or **soda glass**. The high quality of those screw threaded tubes is reflected in the uniformity of their wall thickness and of their diameter and height dimensions.

Supplied in small quantities per case for a more convenient use in laboratory.

Tubes are supplied without caps. To order caps see below.

Polypropylene caps are designed for the round bottom tubes.

Aluminium caps are designed for flat bottom tubes.

	'											
code	total capacity ml	screw	Ø int. mm tube	Ø ext. mm tube	height mm	thickness mm	case quantity	case weight	case volume	suitable cap		
Round bottom borosilicate tubes, supplied in boxes.												
713100	5	13 - 415	10.65	12.75	100	1.00	4 x 250	9.56	0.025	617100		
716100	11	15 - 415	14.00	16.00	100	1.05	4 x 250	13.60	0.045	617200		
716125	14	15 - 415	14.00	16.00	125	1.05	4 x 250	14.64	0.045	617200		
716150	18	15 - 415	14.00	16.00	150	1.05	4 x 250	17.42	0.047	617200		
720150	34	18 - 415	17.70	20.00	150	1.15	4 x 125	11.80	0.036	617300		
Flat bottom soda glass tubes supplied in trays.												
617145	24	15 - 415	14.20	17.75	144.5	1.05	252	5.54	0.010	617000		

Screw caps for the above tubes

For round bottom borosilicate glass tubes, choose **polypropylene caps** (black colour).


For flat bottom soda glass tubes, choose **aluminium** caps, with an internal liner made of caoutchouc.

Pharmaceutical quality.

Both models of caps are fully **autoclavable**.

code	suitable tube	type of cap	bag quantity	bag weight	bag volume
617000	617145	aluminium	500	0.4	0.005
617100	713100	polypropylene	1,000	0.9	0.001
617200	716100 / 716125 / 716150	polypropylene	1,000	1.2	0.002
617300	720150	polypropylene	500	1.1	0.003



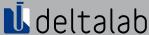



Flat bottom screw threaded vials

Made of ultra clear borosilicate glass. Fully **autoclavable**.

Supplied in trays individually shrink-wrapped. Thickness: 1 mm. Caps are supplied separately. To order caps see below.

mod.	code	dimensions mm	volume ml	screw	tray quantity	tray weight	tray volume	suitable cap
1	900204	15 x 45	4	13-425	100	0.4	0.001	917000
2	900205	18 x 50	5	15-425	250	1.1	0.003	918000
3	900211	22 x 48	10	18-400	100	0.9	0.001	922000
4	900212	19 x 65	12	18-400	100	1.1	0.002	922000
5	900220	28 x 65	20	22-400	2 x 150	5.4	0.010	928000
6	900225	28 x 70	24	22-400	100	1.9	0.005	928000

Caps for the above vials

Made of autoclavable polypropylene. Black caps feature an internal silicone joint.

letter	code	screw	bag quantity	bag weight	bag volume
Α	917000	13-425	100	0.10	0.001
В	918000	15-425	250	0.18	0.001
С	C 922000	18-400	250	0.22	0.001
D	928000	22-400	300	0.28	0.001



Racks for containers up to 23 mm Ø

Made of stainless steel. 25 mm square holes. Features a strong woven mesh base.

265 x 265 x 85	100 (10 x 10)	1	0.425	0.0059
219 x 166 x 85	48 (8 x 6)	1	0.238	0.0030
65 x 110 x 85	24 (6 x 4)	1	0.145	0.0015
111 x 84 x 85	12 (4 x 3)	1	0.095	0.0008
	219 x 166 x 85 65 x 110 x 85	219 x 166 x 85 48 (8 x 6) 65 x 110 x 85 24 (6 x 4)	219 x 166 x 85 48 (8 x 6) 1 65 x 110 x 85 24 (6 x 4) 1	219 x 166 x 85 48 (8 x 6) 1 0.238 65 x 110 x 85 24 (6 x 4) 1 0.145



See other racks for tubes in Chapter 9. **Sample Storage**




0.4 ml microtubes, Beckman® type

Manufactured in polyethylene.

Dimensions: 5.50 x 47.50 mm (Ø x h), capped.

Withstand centrifugation up to 11.000 xg and can be used at temperatures from $-10 \,^{\circ}\text{C}$ to $130 \,^{\circ}\text{C}$.

	code	description	case quantity	case weight	case volume
	900007	Beckman® type microtube 0.4 ml	18 x 1,000	11.32	0.070

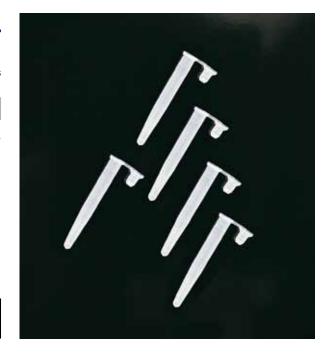

See our wide range of racks for microtubes in Chapter 9. **Sample Storage**



0.5 ml microtubes, Vitatron® type

Manufactured in polypropylene.

Autoclavable. Graduated.

Easy handling caps: they can be opened and closed with one hand.

Will withstand centrifugation up to: 14,000 xg.

Dimensions: 7.87 x 31.8 mm (diameter x height) capped.

code	description	case quantity	case weight	case volume
900118	Vitatron® type microtube 0.5 ml	12 x 1,000	7.0	0.052



0.5 ml microtubes, Vitatron® type

Made of **autoclavable** polypropylene. Ideal for little sample volumes, microtechniques, microsedimentation, etc.

Suitable for transport, storage and freezing of biological samples.

Watertight closing system.

Easy handling caps: they can be opened and closed with one hand. Easily pierceable.

Will withstand centrifugation up to 11,000 xg.

The rim in the inner part of the cap assures an optimum watertightness.

Resist temperatures from -10 °C to 130 °C.

Dimensions: 7,50 x 31,25 mm (Ø x h) capped.

code	description	case quantity	case weight	case volume
900117	Vitatron® type microtube 0.5 ml	14 x 1,000	7.30	0.069







0.5 ml microtubes

Made of ultra clear polypropylene. **Autoclavable**.

Tubes have a frosted writing surface on side and on closure for better sample identification.

Caps are easy to handle and can be opened and closed with one hand. Lids are easily pierceable. Tubes have moulded-in graduations in 100 ul increments from 0.1 to 0.6 ml. Withstand centrifugation up to **17,000 xg**.

Dimensions: 30 x 8 mm. Can be used from -80 °C to 121 °C.

RNAse, DNAse and PCR inhibitors free.

Model **4092.1NS**, with a low adhesion surface, is specially designed for research procedures such as protein work and nucleic acid amplifications. Made from special resins to minimise liquid retention and ensure optimum sample yield, thus eliminating the use of lubricants that may be harmful to samples.

The proprietary formulation is completely **non-reactive.**

С	ode	description	case quantity	case weight	case volume
409	92.1N	natural graduated	10 x 1.000	7.28	0.053
409	2.1NS	low adherence natural graduated	10 x 500	3.42	0.028



Standard 1.5 ml microtubes

Made of ultra clear polypropylene. **Autoclavable**. Attached flat caps.

Tubes can be easily opened and closed with one hand. Ideal for the storage and freezing of any kind of biological material (cells, sperm, bacteria, etc.). Can be used at temperatures down to $-100~^{\circ}\text{C}$.

Withstand centrifugation up to 11,000 xg.

code	dimensions Ø x h mm	colour	case quantity	case weight	case volume
200400	10.86 x 39.06	natural	5 x 1,000	5.40	0.051
200401	10.86 x 39.06	yellow	5 x 1,000	5.40	0.051
200405	10.86 x 39.06	blue	5 x 1,000	5.40	0.051
200407	10.86 x 39.06	green	5 x 1,000	5.40	0.051
200410	10.86 x 39.06	pink	5 x 1,000	5.40	0.051



1.5 ml microtubes

 $\label{thm:manufactured} \mbox{Manufactured in polypropylene. } \textbf{Autoclavable}.$

Easy opening and closing.

Graduated tube.

It features a frosted area for writing. Flat cap.

Will withstand centrifugation up to: 12,500 xg.

code	dimensions Ø x h mm	case quantity	case weight	case volume
200400P	10.2 x 39.9	12 x 500	7.42	0.056

1.5 ml safety lock microtubes

Made of ultra clear polypropylene.

Feature a special "snap-seal" secure lock.

A special locking clasp with a distinctive audible snapseal design ensures tubes will not open during centrifugation. Tubes have a frosted writing surface on side and on attached closure for better sample identification.

Caps are easy to handle and can be opened and closed with one hand. Lids are easily pierceable.

Tubes have moulded-in graduations in 0.5 ml increments.

Autoclavable.

Can be used from -90 °C to 121 °C.

Withstand centrifugation up to 17,000 xg.

code	dimensions Ø x h mm	case quantity	case weight	case volume
4092.5N	11 x 39	10 x 500	5.76	0.055





1.5 microtubes with secure lock pick-up tab

Same characteristics as code **4092.5N** but these microtubes include a pickup tab for easier handling without actually touching the tube. Pierceable cap.

Withstand centrifugation up to 17,000 xg.


Autoclavable.

Withstand temperatures from -175 °C to 121 °C.

code	colour	case quantity	case weight	case volume
4091.1N	natural	10 x 500	4.94	0.045





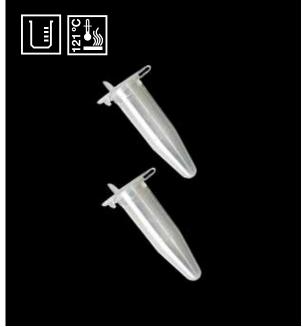

1.5 ml graduated microtubes

Made of ultra clear polypropylene.

Autoclavable.

Feature a frosted writing surface on side and on closure for better sample identification. Caps are easy to handle and can be opened and closed with one hand. Lids are easily pierceable.

Can be used from -80 °C to 121 °C.

Dimensions: 39 x 11 mm.

Caps specially designed to reduce popping when centrifugation, freezing, boiling, storing or shipping is required.


- when the objective is to avoid liquid remaining on the wall of the tube, just push the cap slightly and a little centrifugation time is required.
- when centrifugation is required, push the cap firmly to achieve complete closure.

Tubes withstand centrifugation up to 13,000 xg.

code	colour	case quantity	case weight	case volume
4092.2N	natural	10 x 500	6.82	0.048

Other colours under request.



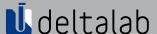



1.5 ml microtubes. DNASE, RNASE free

Made of ultra clear polypropylene. **Autoclavable**. Attached flat caps.

Graduated from 100 μl.

Tubes can be easily opened and closed with one hand. Ideal for the storage and freezing of any kind of biological material (cells, sperm, bacteria, etc.). Can be used at temperatures from **-90 °C** to **121 °C**.

Withstand centrifugation up to 20,000 xg.

DNAse, RNAse, DNA and PCR inhibitors free.

Dimensions (Ø x h): 10.7 x 39.1 mm.

code	description	case quantity	case weight	case volume
4092.3N	natural colour	10 x 500	7.30	0.053
4092.3NS	natural colour, clear siliconised	10 x 250	3.70	0.026



2 ml microtubes with attached caps

Made of ultra clear polypropylene.

Autoclavable.

Certified RNAse and DNAse free.

Feature a frosted writing surface on side and on closure for better sample identification.

Easy handling: tubes can be opened and closed with one hand.

Lids are easily pierceable.

Tubes have moulded-in graduations in 100 μ l increments up to 2 ml.

Dimensions: 11 x 40 mm.

Can be used at temperatures from **-90 °C** up to **121 °C**.

Withstand centrifugation up to 17,000 xg.

Model **4092.6NS** has a low adhesion surface for special applications such as protein work.

code	description	case quantity	case weight	case volume
4092.6N	natural colour	10 x 500	7.27	0.053
4092.6NS	natural colour, clear siliconised	10 x 250	4.09	0.036



2 ml graduated microtube

Microtubes made of ultra clear polypropylene. **Autoclavable**.

2 ml graduated tubes with frosted area for writing.

Flat cap that features a secure lock with safety shield.

Easy handling with just one hand.

They can be used at temperatures from -20 °C up to 121 °C.

Withstand centrifugation up to 12,000 xg.

code	dimensions Ø x h mm	case quantity	case weight	case volume
4092.7N	10.3 x 40.5	10 x 500	6.50	0.055



Screw thread microtubes

Made of **autoclavable** polypropylene, they can be used at extreme temperatures from **-190 °C** to **121 °C**. Two versions available: in transparent polypropylene, or opaque brown **(UV resistant**, designed to be used with light sensitive samples).

Certified RNAse, DNAse and pyrogen free. Tubes and caps can be centrifuged at 17,000 kg. Caps are supplied separately, see below. Dimensions: $11 \times 44 \text{ mm.}$

mod.	code	volume ml	skirt	case quantity	case weight	case volume		
Transparent polypropylene								
1	409110.1	0.5	yes	1,000	1.46	0.010		
2	409110.2	1.5	yes	1,000	1.45	0.010		
3	409110.3	1.5	no	1,000	1.13	0.010		
4	409110.4	2.0	yes	1,000	1.34	0.010		
Brown po	Brown polypropylene							
5	409113.1	0.5	yes	1,000	1.53	0.010		
6	409113.2	1.5	yes	1,000	1.34	0.010		
7	409113.3	1.5	no	1,000	1.34	0.009		
8	409113.4	2	yes	1,000	1.34	0.009		



Caps for screw thread microtubes

Made of polypropylene.

Caps are available either with an attachment loop or without loop.

Both models have a sealing O-ring (red) of silicone to ensure a positive leakproof seal.

For sample identification, colour coding inserts can be placed upon caps (made of polypropylene).

Cap dimensions: 13 x 8 mm.

mod.	code	description	case quantity	case weight	case volume				
Transpare	Transparent polypropylene								
1	409007.N	cap with loop	1,000	0.40	0.003				
2	409008.N	cap without loop	1,000	0.40	0.003				
Brown po	Brown polypropylene								
3	409007.M	brown cap with loop	1,000	0.40	0.003				
4	409008.M	brown cap without loop	1,000	0.40	0.003				
Inserts	Inserts								
5	409111A	blue insert	500	0.58	0.003				
6	409111AM	yellow insert	500	0.58	0.003				
7	409111B	white insert	500	0.58	0.003				
8	409111R	red insert	500	0.58	0.003				



Screw cap microtubes, with cap. Sterile

Tubes and caps in medical grade, transparent polypropylene.

The cap embodies a non-reactive ethylene-propylene o-ring.

Suitable for autoclave, liquid nitrogen (gaz) and boiling processes.

Perfect for long term storage.

Withstand temperatures from -190 °C. Centrifugation Resistance: 20,000 xg.

DNAse, RNAse, DNA, and PCR inhibitors free.


Graduated models feature a frosted area for writing.


The non-graduated model incorporates an external grip for an easy handling. Microtubes are supplied capped, in bags of 50 units.

code	description		graduation	case quantity	case weight	case volume
409115/4	0.5 ml skirted	STERILE	×	50 x 50	5.90	0.030
409115/2	1.5 ml non skirted	STERILE	\checkmark	50 x 50	5.90	0.030
409115/6	2 ml skirted	STERILE	\checkmark	50 x 50	5.90	0.030
409115/3	2 ml non skirted	STERILE	\checkmark	50 x 50	5.90	0.030





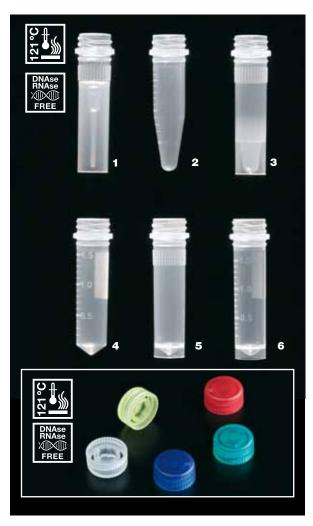

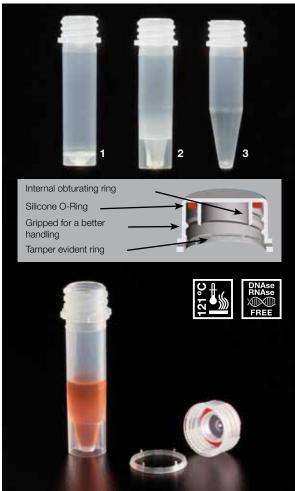

Screw thread microtubes

Made of transparent polypropylene.

Suitable for use in liquid nitrogen, autoclave and for boiling applications, and can be used at temperatures down to -190 °C.

Certified RNAse, DNAse and PCR inhibitors free.

Products ideal for long term sample storage.

Withstand centrifugation at 20,000 xg. The codes 409111/4, 409111/5 and 409111/6, with external moulded slots for better handling with gloves.

Dimensions: 10.3 x 44.5 mm (except code **409111/2**: 10.3 x 43.6 mm).

Caps are supplied separately, see below.

mod.	code	description	graduation	case quantity	case weight	case volume
1	409111/4	0.5 ml skirted	×	500	0.70	0.002
2	409111/2	1.5 ml	\checkmark	500	0.60	0.002
3	409111/5	1.5 ml skirted	×	500	0.70	0.002
4	409111/3	2.0 ml	\checkmark	500	0.70	0.002
5	409111/6	2.0 ml skirted	×	500	0.70	0.002
6	409111/7	2.0 ml skirted	\checkmark	500	0.70	0.002

Caps for screw thread microtubes

Made of medical grade polypropylene.

Feature an internal O-ring to ensure leakproof seal.

Dimensions 13.0 x 6.0 mm.

Certified RNAse, DNAse and PCR inhibitors free.

code	colour	case quantity	case weight	case volume
409112/0	natural	500	0.024	0.0001
409112/1	blue	500	0.024	0.0001
409112/2	green	500	0.024	0.0001
409112/4	red	500	0.024	0.0001
409112/6	yellow	500	0.024	0.0001





Screw cap until locking ring clicks over serrated tube neck.

Contents are now protected until cap is removed and tamper-evident ring is

Screw thread tamper evident microtubes

Microtubes and caps are made of autoclavable ultra clear polypropylene. Ribbed cap with internal silicone O-ring for a positive leakproof seal. Super fast 1/4 turn thread design. Tamper evident seal which notices if microtube has been opened.

Rnase, Dnase and Pyrogen free.

They can be used at extreme temperatures from -196 °C to 121 °C.

Resistance to centrifugation: 17,000 xg.

Microtube dimensions: 11 x 44 mm.

Cap dimensions (with tamper-evident ring): 15 x 9 mm.

- Used on: Test of fertility and DNA testing
- Packaging of diagnostic kits and reagents
- Forensic laboratories

	mod.	code	volume ml	case quantity	case weight	case volume
	1	409110.4T	2	1,000	2.06	0.013
	2	409110.2T	1.5	1,000	2.06	0.013
EW	3	409110.5T	1.5	1,000	0.65	0.013





CRYOINSTANT: Porous beads for microbiological culture preservation

CryoInstant is an excellent preservation system for storing microbiological culture (for example fungus on sporulation phase), consisting on **2 ml cryovials** containing 25 porous beads and cryopreservative-added broth.

With these vials we can:

- Get a perfect freezing storage system.
- Obtain up to 25 replicates from the original microorganism, to use them during many years
- Ease the recovery of the stored strain streaking an inoculated bead directly onto solid media
- Prevent changes and repetitive thawing / freezing of the rest of beads
- **Avoid** the **growth** of **ice** when recovering,
- **Improve recovery rates** of frozen cultures by adding cryopreservative to broth
- Minimise the risk of cross contamination,
- Save space in the freezers when storing.

The external thread cryovials are made of polypropilene (page 145). Their smooth inner surface prevents from contamination. Their caps embody a silicone washer to ensure a positive leakproof seal. Tubes and caps withstand up to **-190 °C**.

In order to ease the classification of the samples, we offer the possibility to get cryovials with caps and beads in five different colours: every code displays caps and beads from a single and same colour (excepting code **409113/6**, which is an assortment), allowing a fast and easy identification of the samples, and providing a system whereby users are able to code different bacterial species, different samples, different laboratories, etc. And even more: with colour coded inserts (described on page 144) a further classification can be made.

Cryovials are presented in carton boxes (resistant to **-100 °C**; dimensions: 150x150x55 mm) with 100 units; boxes are described in page 206. Each box is supplied plastic-wrapped and labelled with code reference, batch, expiration date and specifying the colour of the caps and beads.


code	beads and cap colour	case quantity	case weight	case volume
409113/1	white	100	0.53	0.002
409113/2	blue	100	0.53	0.002
409113/3	yellow	100	0.53	0.002
409113/4	red	100	0.53	0.002
409113/5	green	100	0.53	0.002
409113/6	assorted: 5 colours x 20 cryovials		0.53	0.002

Way of use:

- 1 Take the sample from the strain with a loop (please see our loops on pages 32-33)
- 2 Inoculate the cryovial by putting the loop in contact with the broth
- 3 Close the cryovial and shake it gently in order to permeate the strain into the beads
- **4** Extract the cryopreservative-added broth with a Pasteur pipette (see our plastic pipettes on pages 189-193)
- 5 Close the cryovial and take it to the freezer
- 6 Every time we need to reproduce the strain, we will take a bead with a loop or a forceps
- 7 We will put it down on a Petri dish with culture medium, attempting to put in contact all the surface of the bead with the medium





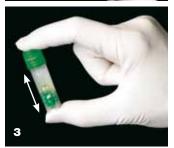

















Solutions for cryopreservation

Wilmut offers a **wide range of biological simple storage tubes** that fit in the 96-well plate SBS standard format, and single tubes which are compatible with any freezing box format.

With volumes in between 0.5 ml and 1.2 ml and different options of codification that allows Wilmut system **adapt to the needs of any laboratory or biobank**. They are products designed and oriented exclusively to the user.





Main features:

- Tubes and racks for microbiological sample storage
- It is the tube with more capacity in the minimum space: maximum optimization of the freezers (from 15 % to 75 % of space saving)
- Maximum traceability: 2D datamatrix code, legible code at the tube, and bar code in the rack
- Traceability files available at www.wilmut.es or using software management
- Big quantity of accessories: caps, lockers, pisckers, cappers, decappers, labels, scanners, etc.



Sterile cryovials

Tubes made of **autoclavable** polypropylene.

Designed for the storage of biological material at temperatures as low as -196 °C. (Liquid Nitrogen but in gas state).

Cryovials have printed graduations and feature a white marking area for sample identification.

Caps and tubes have both the same coefficient of expansion, which further enhances their leakproof qualities at changing temperatures.

Caps are manufactured in polyethylene* and embody a silicone O-ring that guarantees leakproof.

Colour code inserts can be fitted into the caps for identification.

All cryovials are certified RNAse, DNAse and pyrogen free.

Cryovials are sterilized by radiation and are packaged in safety-lock bags of 100. Supplied capped.

External threads:

minimize liquid retention. Non skirted versions with stand centrifugation up to 17,000 kg.

* 10 ml cryovial (code **401410**) features a polyethylene cap which does not accept inserts.

Internal threads:

Non skirted versions withstand centrifugation up to 14,000 xg.

Cryovials with external threads:

code	volume ml	skirt	dimensions mm*	case quantity	case weight	case volume
409105.1	1.2	yes	12.5 x 42	10 x 100	2.55	0.014
409103.1	2.0	no	12.5 x 47	10 x 100	2.70	0.016
409106.1	2.0	yes	12.5 x 49	10 x 100	2.68	0.021
409107	3.0	yes	12.5 x 71	10 x 100	3.88	0.023
409108	4.0	yes	12.5 x 77	10 x 100	3.90	0.028
409109	5.0	yes	12.5 x 92	10 x 100	4.60	0.023
401410	10.0	yes	17.0 x 84	10 x 50	2.80	0.020

^{*} Capped.

Cryovials with internal threads:

code	volume ml	skirt	dimensions mm*	case quantity	case weight	case volume
409001	1.2	yes	12.5 x 41	10 x 100	1.94	0.015
409002	2.0	no	12.5 x 48	10 x 100	2.22	0.016
409002.1	2.0	yes	12.5 x 49	10 x 100	2.24	0.014
409003	4.0	no	12.5 x 70	10 x 100	3.79	0.028
409003.1	4.0	yes	12.5 x 72	10 x 100	3.90	0.028
409003.2	5.0	no	12.5 x 90	10 x 100	4.60	0.024

^{*} Capped.



See racks at page 209



Colour coded inserts


Made of polypropylene. Very convenient for colour identification in freezers, on benchtops etc.

Fit precisely into the caps of cryovials. They feature a little hole in order to help removing them from the cap. Not suitable for code **401410.**

code	colour	case quantity	case weight	case volume
409003A	blue	500	0.08	0.0004
409003AM	yellow	500	0.08	0.0004
409003R	red	500	0.08	0.0004
409003V	green	500	0.08	0.0004





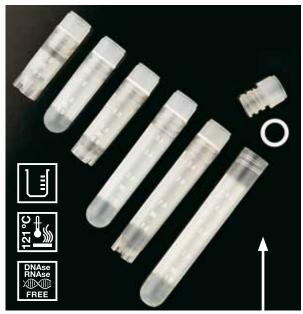













Non sterile cryovials with external threads

Made of polypropylene. Designed for the storage and transportation of biological material. The external thread design provides a smooth and uniform inner surface, thus reducing the risk of contamination.

Tubes can be **autoclaved** to **121 °C** and withstand temperatures down to **-190 °C**. The closures and tubes are both manufactured of polypropylene having the same coefficient of expansion, which further enhances their leakproof qualities at changing temperatures.

Cryovials withstand centrifugation up to **17,000 xg**. Skirted tubes are not recommended for centrifugation.

Available with or without graduation and white band.

Caps suitable for these cryovials are supplied separately, see code **409110** and series.

Cryovials printed without graduation:

code	volume ml	skirt	dimensions mm*	case quantity	case weight	case volume
409102.1	1.2	yes	12.5 x 43	1,000	1.38	0.007
409104.1	2.0	no	12.5 x 48	1,000	1.58	0.007
409106.2	2.0	yes	12.5 x 49	1,000	1.68	0.006
409107.1	3.0	yes	12.5 x 72	1,000	2.48	0.020
409108.1	4.0	yes	12.5 x 76	1,000	2.75	0.017
409109.1	5.0	yes	12.5 x 93	1,000	3.66	0.023

^{*} Capped

Cryovials with graduation and band:

code	volume ml	skirt	dimensions mm*	case quantity	case weight	case volume
409132	2.0	yes	12.5 x 49	1,000	2.68	0.010
409133	3.0	yes	12.5 x 72	1,000	3.77	0.010
409135	4.0	yes	12.5 x 76	1,000	3.90	0.020

^{*} Capped

Caps for the above tubes

Made of polypropylene. Feature a long skirt and a special thread design allowing them to be removed or sealed with a single turn and assuring a better protection of the content. A silicone washer ensures a positive leakproof seal at all temperatures.

Caps are made of polypropylene having the same coefficient of expansion as the above cryovials. See the previous page for colour inserts to fit into these caps.

code	colour	case quantity	case weight	case volume
409110	natural	1,000	0.631	0.003
409110/1	blue	1,000	0.082	0.003
409110/4	red	1,000	0.802	0.003

Canes for cryovials

Made of aluminium.

Suitable for the storage of five or six 1.2 or 2.0 ml tubes in liquid nitrogen freezers such as Dewar flasks.

1	code	dimensions	case quantity	case weight	case volume
	409120	290 mm	12	0.20	0.0001











30 ml universal polystyrene containers (25 x 90 mm)


C€ (IVD)

Polypropylene caps. 30 ml containers. These containers are manufactured with an innovative process resulting in one of the most resistant containers available on the market.

The standard cap colour is white.

Other colours can be supplied (ask for quantities).

mod.	code	description	sterile	case quantity	case weight	case volume	cases pallet
1	309202	polystyrene	no	400	5.50	0.045	32
2	309222	polystyrene individually wrapped	no	400	5.70	0.045	32
2	309222.O	polystyrene individually wrapped	ETHYLENE-OXIDE	400	5.70	0.045	32

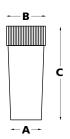

Dimensions:

A: 25.1 mm; B: 30.8 mm; C: 92.5 mm. Maximum volume: 31.7 ml. Recommended volume: 25 ml. **Withstand centrifugation up to: 7,000 xg.**

In an **ASEPTIC PRODUCTION**, manufacturing and packaging methods and environment are maintained to avoid or significantly reduce bacterial contamination (bioburden).

Our **EUROTUBO**® containers are injected at a temperature of **250** °C and are automatically connected to the assembling and packaging processes, without manual intervention, so as to guarantee an aseptic production.





30 ml universal polypropylene containers (25 x 90 mm)


C€ (IVD)

Made in Polypropylene. 30 ml containers. These containers are manufactured with an innovative process resulting in one of the most resistant containers available on the market. The screw closure can assure an almost leakproof closing system. Polypropylene cap in white colour, except **409602** (blue). Other colours and labeled models can be supplied (ask for quantities).

Autoclavable (except the code 409602, the polystyrene spatula must be removed before autoclaving).

mod.	code	description	sterile	case quantity	case weight	case volume	cases pallet
1	409202	polypropylene	no	400	4.90	0.046	32
2	409222	polypropylene individually wrapped	no	400	5.10	0.046	32
2	409222.O	polypropylene individually wrapped	ETHYLENE-OXIDE	400	5.10	0.046	32
3	409602	polypropylene with PS spoon	no	400	5.30	0.046	32

Dimensions: A: 24.8 mm; B: 30.8 mm; C: 91.6 mm. Maximum volume: 31.1 ml. Recommended volume: 25 ml. Withstand centrifugation up to: 10.000 xg.

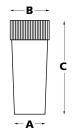



30 ml container with boric acid (25 x 90 mm)

This container includes a boric acid preservative which allows the storage and transport of the urine sample to the laboratory.

The boric acid supports the viability of the bacteria while preventing further multiplication. This preservative does not alter the characteristics of the urine and will not affect analysis. The bacteria remain viable for 24-48 hours or more, depending on the type of microorganisms.

Boric acid quantity: 0.4 gr.

Sterile by ethylene oxide.

code	description	case quantity	case weight	case volume	cases pallet
409202.B	polypropylene with boric acid	400	5.70	0.046	32








50 ml containers (38 x 70 mm)

(€ (IVD)

With polyethylene caps. Volume: 50 ml.

code	description	case quantity	case weight	case volume	cases pallet
409802	ultra clear polypropylene	700	9.30	0.096	16
409852	ultra clear polypropylene with PS spoon	700	9.30	0.096	16
409826	ultra clear polypropylene individually wrapped	600	8.40	0.096	16

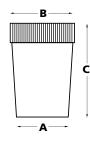

Dimensions:

A: 35.1 mm; B: 41.5 mm; C: 72.3 mm. Maximum volume: 56 ml. Recommended volume: 45 ml.

40 ml containers (30 x 70 mm)

C€ (IVD)

With polyethylene caps. Graduated. Volume 40 ml.

code	description	case quantity	case weight	case volume	cases pallet
409402	ultra clear polypropylene	8 x 100	8.35	0.096	16
409426	ultra clear polypropylene, individually wrapped	800	10.50	0.110	16
309402	polystyrene	8 x 100	11.00	0.096	16

Dimensions:

A: 29.1 mm; B: 35.7 mm, C: 69.8 mm. Maximum volume: 39 mm.

Recommended volume: 39 mm.







50 ml tamper evident containers (38 x70 mm)


(€ (IVD)

Ultra clear polypropylene containers with a specially designed tamper evident ring. Polyethylene caps. Standard colour of caps: red. Volume: 50 ml.

code	description	case quantity	case weight	case volume	cases pallet
419802	ultra clear polypropylene	700	9.30	0.096	16
419826	ultra clear polypropylene individually wrapped	600	8.40	0.096	16

Containers are supplied capped but not fully locked so that the container is not fully sealed.

Once samples (urine, milk, water, etc.) are collected, screw the cap in a normal way and then screw again to fully position the tamper evident ring.

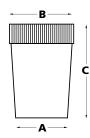

Dimensions:

A: 35.1 mm; B: 41.5 mm; C: 72.3 mm. Maximum volume: 56 ml. Recommended volume: 45 ml.










60 ml polystyrene containers (38 x 65 mm)

C€ (IVD)

New generation of 60 ml containers with high safety closure biflex caps. Caps made of polyethylene.

Standard colour of the cap: red.

mod.	code	description	sterile	case quantity	case weight	case volume	cases pallet
1	309502	polystyrene	aseptic	600	9.10	0.082	20
2	309552	polystyrene with spoon	no	600	8.20	0.096	16
3	309526	polystyrene individually wrapped	STERILE A	750	11.55	0.140	12

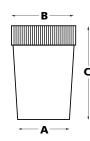

Dimensions:

A: 39.3 mm; B: 44.2 mm; C: 66.0 mm.

Maximum volume: 63 ml. Recommended volume: 50 ml.

We can supply these containers labeled. Please ask quantities. We can supply these containers sterilised by radiation.

Racks for containers up to 47 mm Ø

Stainless steel. 50 mm square holes. With a strong woven mesh base.

code	dimensions mm	distribution	quantity	weight	volume
I-240	315 x 210 x 45	24 (6 x 4)	1	0.320	0.0042
I-250	210 x 160 x 45	12 (4 x 3)	1	0.228	0.0022

Wide range of racks and baskets, see chapter 9.











60 ml polypropylene containers (38 x 65 mm)


C€ (IVD)

New generation of 60 ml sterile containers with high safety closure biflex caps. Specially designed for pneumatic transport of specimens. Caps made of polyethylene. Standard colour of the cap: red.

mod.	code	description	sterile	case quantity	case weight	case volume	cases pallet
1	409501	polypropylene labelled	aseptic	600	8.20	0.096	16
2	409502	polypropylene	aseptic	600	8.00	0.096	16
3	409552	polypropylene with spoon	aseptic	600	8.00	0.096	16
4	409526	polypropylene individually wrapped	STERILE A	750	10.30	0.140	12
-	409526.O	polypropylene individually wrapped	ETHYLENE-OXIDE	750	10.30	0.140	12
5	409526.1	polypropylene individually wrapped, labelled	STERILE A	750	10.45	0.140	12

Dimensions:

A: 38.8 mm; B: 44.2 mm; C: 66.1 mm. Maximum volume: 62 ml. Recommended volume: 50 ml.

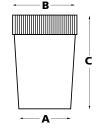



Graduated 60 ml containers 38 x 65 mm

(€ (IVD)

Made of transparent popypropylene.

Moulded graduation in 5 ml increments.

Red polyethylene cap, biflex.

For other cap colours ask the minimum order quantities.

Ideal for pneumatic transport because it has a watertight closing system. Supplied capped.

code	description	sterile	case quantity	case weight	case volume	cases pallet
409502G	polypropylene	aseptic	600	8.00	0.096	16
409526G	polypropylene, individually wrapped	STERILE A	750	10.30	0.140	12

Dimensions:

A: 38.8 mm; B: 44.2 mm; C: 66.1 mm.

Maximum volume: 62 ml. Recommended volume: 50 ml.








150 ml containers (57 x 73 mm) with spoon

(€ (IVD)

Ultra clear polypropylene containers with polyethylene leak proof screw caps. Spoon made of Polystyrene. Graduated in 100 ml increments. Frosted panel on side for easy writing. Standard colour of the cap: red.

mod.	code	description	sterile	case quantity	case weight	case volume	cases pallet
1	409752	polypropylene with PS spoon	aseptic	450	7.95	0.140	16
2	409756	polypropylene individually wrapped with PS spoon	aseptic	350	6.50	0.140	16
2	409756.B	polypropylene individually wrapped with PS spoon	RADIATION	350	6.50	0.140	16

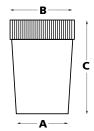

Dimensions:

A: 56.3 mm; B: 62.1 mm; C: 74.3 mm.

Maximum volume: 150 ml.
Recommended volume: 125 ml.

These containers can be supplied sterilized by radiation.

150 ml screw cap containers (58 x 69 mm)


Manufactured in polystyrene or polypropylene depending on the model. Cap manufactured in polyethylene.

Total volume 150 ml.

Code 309752 is supplied with a polystyrene spoon inside.

Code **409700** is sterile by radiation.


mod.	code	description	sterile	case quantity	case weight	case volume
1	309752	polystyrene white cap and spoon	no	360	13.45	0.106
2	409700	polypropylene red cap	RADIATION	360	11.58	0.106

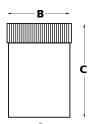

Dimensions:

A: 58 mm; B: 62 mm; C: 69 mm.



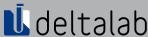





150 ml containers (57 x 73 mm)

(€ (IVD)


Ultra clear polypropylene containers with polyethylene leak proof screw caps. Moulded graduation in 20 ml increments up to 100 ml. Frosted panel on side for easy writing. Standard colour of the cap: red.

mod.	code	description	sterile	case quantity	case weight	case volume	cases pallet
1	409701	polypropylene labelled	aseptic	450	7.95	0.140	16
2	409702	polypropylene	aseptic	450	7.80	0.140	16
3	409726	polypropylene individually wrapped	STERILE A	350	6.50	0.140	16
-	409726.G	polypropylene individually wrapped	RADIATION	400	7.60	0.140	16



Dimensions:

A: 56,3 mm; B: 62,1 mm; C: 74,3 mm. Maximum volume: 150 ml. Recommended volume: 125 ml.

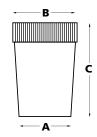

150 ml containers uncapped (57 x 73 mm) (6 (IVD)

With the same characteristics as the previous model, but these containers are supplied uncapped.

Containers and caps are supplied in the same box.

code	description	case quantity	case weight	case volume	cases pallet	
409703	polypropylene uncapped	620	10.50	0.082	20	



C€ (IVD)

With the same characteristics as codes **409702** (p. 153) and **409502** (p. 151) but in blue colour (containers and cap). Ideal for photosensitive samples or when it is required not to show the content. Ask for other labeling or for not labeled containers to our commercial department.

Minimum order quantity: 30,000 units.

Minimum order quantity: 1 pallet = 12,000 units (code **409511**) and 7,600 units (code **409711**).


code	description	case quantity	case weight	case volume	cases pallet
409511	60 ml container white label	600	8.20	0.096	16
409711	150 ml container white label	450	7.95	0.140	16



200 ml graduated container (57 x 116 mm) (6 (IVD)

Containers made of ultra clear polypropylene and caps made of polyethylene. Those containers have double moulded-in graduations, (in 20 ml increments up to 200 ml and in 1 oz increments up to 7 oz), and feature a frosted panel on their side for easy writing. Standard colour of caps: red.

code	max. capacity	case quantity	case weight	case volume	cases pallet
409902	240 ml	275	7.10	0.140	16

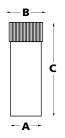

Dimensions:


A: 56.5 mm; B: 62.2 mm; C: 115.5 mm; Maximum volume: 240 ml; Recommended volume: 210 ml.



200 ml screw cap container (58 x 95 mm)


Manufactured in polystyrene, cap made of white polyethylene. Total volume: 220 ml.

code	max. capacity	case quantity	case weight	case volume	cases pallet
309903	220 ml	300	9.6	0.090	16

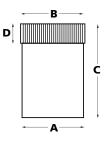

Dimensions:

A: 58 mm; B: 62 mm; C: 95 mm; D: 22 mm.

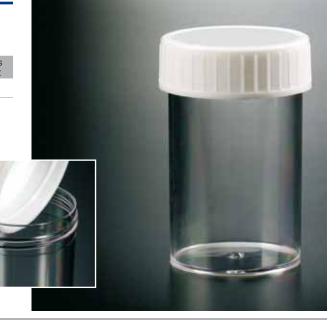

Urine collection system designed to save space in the lab, and to be used when urine samples need to be collected into tubes for a subsequent centrifugation. Designed to be user friendly.

The patient urinates in the 150 ml plastic cup and then pours the urine into the tube. Once the tube is filled up, the patient caps it and delivers the sample to the lab, where the health assistant will fill in the label with the required data. The patient just has to follow the hospital or laboratory instructions. Manufactured under aseptic conditions. There are several versions and custom-designed versions can be manufactured. Bad printed with bacth number and expiry date.





Easy pouring



We can manufacture sets with boric acid.



New bag with instructions.



mod.	code	bag	cup 150 ml	tube	cap	capped	cap colour	label 25 x 35 mm	labelled	case quantity	case weight	case volume
1	309324.E3	1	/	400705	308102	X	natural	/	X	300	4.50	0.140
2	309316.E	1	1	400705	308107	X	green	/	X	300	4.50	0.140
3	309316	1	/	400705	308107	1	green	✓	1	300	4.50	0.140
4	309324.E	1	1	401204	308202	X	natural	/	X	300	4.70	0.140
5	309324	/	1	401204	308207	1	green	/	1	300	4.70	0.140
6	309318	1	1	301212	306601	1	yellow	/	1	300	4.70	0.140




Disposable cups

Disposable sample cups, for general sample collection. Cups are supplied stacked.

code	volume ml	material	mouth Ø x height	case quantity	case weight	case volume
202248	100	polypropylene	57.5 x 62	50 x 100	7.40	0.080
202317	150	polystyrene	70 x 70	30 x 100	8.32	0.080





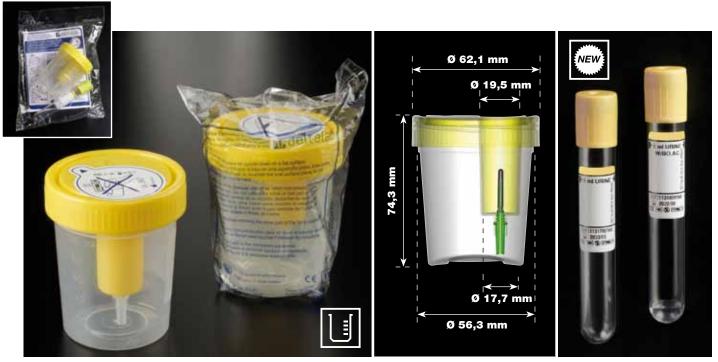

Containers and tubes for vacuum system

C€ (IVD)


120 ml containers designed to use with vacuum tubes. It allows a clean transfer of the urine from the container to the vacuum tube. Manufactured in transparent Polypropylene. Leakproof cap made of yellow polyethylene, with an internal obturating ring. It includes a cannula with a needle inside protected by a rubber cap and a plastic suction nozzle. The cap features a label that seals the cannula in order to prevent contact with the collection needle. The container features a wide marking area for sample identification and a graduation that helps to quantify the sample when uncapped. Once capped, the graduation is inaccurate as the volume is altered by the cannula. Maximum volume with cap: 120 ml, recommended volume: 60 ml.

Tubes in PET produced with the inner cap in caoutchouc and outer cap in yellow polyethylene. Incorporate label for sample identification, with white background and a transverse line in yellow to easily see that it is a tube for urine. Printed with the lot, expiration date, CE mark, volume of suction, manufacturer, and only one use symbol.

Label dimensions: 36 x 20 mm. Supplied in foam racks of 100 units, shrink-wrapped and labeled.

code	presentation	sterile	case quantity	case weight	case volume	cases pallet
Container (dimensions: 57	x73 mm)					
408702	in bulk	no	475	9.60	0.14	16
408726	individually wrapped	no	350	7.50	0.14	16
408702.G	in bulk	RADIATION	250	5.05	0.08	24
408726.G	individually wrapped	RADIATION	350	7.50	0.14	16
Vacuum tube (dimensions	: 16x100 mm)					
408900	vacuum tube, 9 ml	RADIATION	10 x 100	12.60	0.062	10
408910	vacuum tube, 9 ml with boric acid	RADIATION	10 x 100	13.70	0.062	10
KIT container and vacuu	m tube					
408720	KIT 408702.G + 408900		250	3.90	0.10	20
408721	KIT 408702.G + 408910		250	4.10	0.11	20

Process of sample transfer:

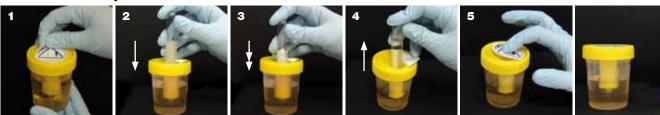

- 1. Remove the cap label.
- 2. Place evacuated tube into cavity on cap, stopper down.
- 3. Advance the tube over puncture point to pierce stopper. The urine tube will be filled. Repeat if another tube is to be collected.
- 4. Hold tube in position until filled and remove tube from device.
- 5. Replace label over hole and reseal.

This is a completely hygienic and safe system.

It prevents sample loss and avoid the contamination of the sample, the user or the work environment.







Rectangular graduated until 3 litres container with handle

Product designed for collection of 24 hours urine samples. Container and cap made of polyethylene.

Brown opaque body and cap depending on the model:

- 408600 Yellow cap with vacuum system
- 408610 Standard white cap

Both are completely leakproof thanks to its sealing ring and ensure a total sealing of the sample.

The vacuum system ensures a clean transfer of the urine from the container to the vacuum tube and prevents the risk of cross contamination in the handling of the sample. It also avoids sample loss by splashing or outpouring.

The container has a handle with ergonomic empty space to facilitate a much more comfortable dump and the manual transport of the sample. It incorporates a translucent vertical stripe on one side where the graduation allows the control of the level of sample containing. Graduated each 100 ml, to 3,000 ml.


Total capacity: 3.3 l.

Unit weight (without cap): 115 g (\pm 3 g). Cap weight (vacuum system): 8.2 g (\pm 0.2 g). Cap weight (standard system): 5.5 g (\pm 0.2 g).

	code	presentation	sterile	case quantity	case weight	case volume	cases pallet
WEW	408600*	3 I rectangular container with vacuum system	no	30	4.10	0.15	12
WEW	408610	standard 3 I rectangular container	no	30	3.90	0.15	12

^{*}Suitable with vacuum tubes 408900 and 408910.



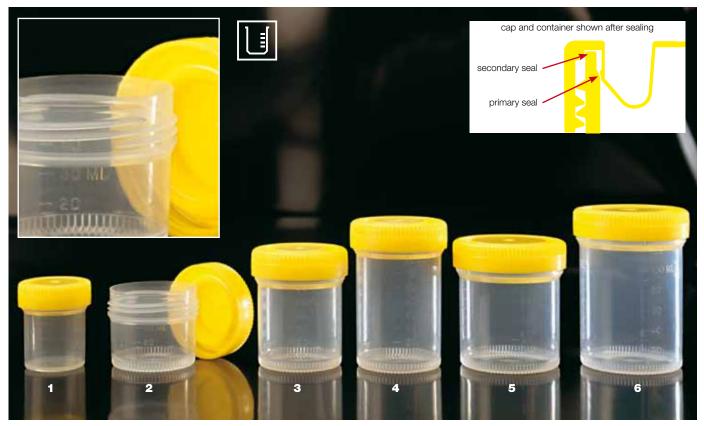

Security screw cap containers

C€ (IVD)

Ultra clear polypropylene containers with yellow polyethylene cap.

They have moulded graduations in ml.

Supplied uncapped with the caps and containers in the same box.

Ideal for the pneumatic transport of liquid samples, including those that have phormaline, because they have a watertight closure. The cap, with an innovative inner shape, has a double closure. (see the figure).

The base and the cap are ribbed so the handling with gloves is comfortable and reliable.

There is a ring in the upper part of the body, in order to avoid the dripping in case of liquid decanting.

For other cap colours (blue, orange, magenta), ask for the minimum order quantities.

Certified as suitable for Formalin biopsies

mod.	code	grad. volume ml	total volume ml	mouth Ø mm	height mm	case quantity	case weight	case volume
1	202840	20	26	34	43	1,000	8.84	0.090
2	202841	40	50	45	43	600	8.75	0.094
3	202842	60	76	45	60	500	8.20	0.110
4	202843	80	92	45	75	400	7.46	0.104
5	202848	90	108	53	68	300	6.35	0.102
6	202844	120	146	53	78	300	6.66	0.109

EUROTUBO® Leak proof screw cap containers


C€ (IVD)

Specially designed for histology, they are provided with double internal security closure. Manufactured in translucent polypropylene with polyethylene ribbed yellow cap. They have moulded graduations.

Their leakproof closure makes them ideal for liquid or solid sample storage.

They meet the standard UNE-EN 14401:2005 Rigid plastic containers. Methods to test the efectiveness of closures.

Supplied uncapped.

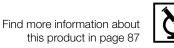









Cylindrical container graduated up to 2 litres with handle

C€ (IVD)

High density polyethylene container and cap. Ideal for 24 hour urine collection and any other liquid samples.

We recommend references with white and ribbed cap, that features a pharmaceutical grade liner, which eliminates the need for an insert plug, allows a better handling and provides a watertight seal.

Containers are graduated in 50 ml increments up to 2 l. Include a convenient handle for an easy single-handed use.

Base diameter: 108 mm. Height, cap included: 290 mm.

Cases per palet: 16.

Codes **408001.1** and **408001.0** are individually wrapped by retractile system. Ask for other versions cap colours, capped containers and sterile containers.

code	body colour	cap description	presentation	individually wrapped	sterile	internal mouth Ø mm	case quantity	case weight	case volume
408001	natural	white, with liner	unscrewed	no	no	65	34	7.00	0.140
408001.O	natural	white, with liner	screwed	yes	RADIATION	65	34	7.00	0.140
408001-86	natural	black + plug	unscrewed	no	no	86	33	7.00	0.140
408001OBT	natural	black + plug	unscrewed	no	no	66.5	33	4.25	0.140
408001O/T	natural	black + plug	screwed	no	no	66.5	33	4.59	0.140
408001.1	natural	white, with liner	screwed	yes	no	65	34	7.00	0.140







Rectangular containers, graduated up to 2 litres

C€ (IVD)

Containers and caps made of high density polyethylene.

Ideal for 24 hour urine collection and other liquid samples.

Containers are in transluscent colours and moulded graduation in 50 ml increments up to 2 l.

Ergonomically designed for single-handed use with a lateral tap of 66 mm length Height, cap included: 253 mm.

Base: 84 x 125 mm.

Cases per pallet: 16.

Models with white cap, have pharmaceutical grade liner, which **eliminates the need of a plug**, allows a better handling and provides a leakproof closure.

Individually wrapped models are supplied capped.

Ask for other versions: colour (both, cap or body), and sterile containers.

code	body colour	cap description	presentation	individually wrapped	sterile	internal mouth Ø mm	case quantity	case weight	case volume
407001	natural	white, with liner	unscrewed	no	no	67.8	43	6.00	0.140
407001/T	natural	white, with liner	screwed	no	no	67.8	43	6.00	0.140
407008	natural	white, with liner	screwed	yes	no	67.8	43	6.40	0.140
407008.O	natural	white, with liner	screwed	yes	RADIATION	67.8	43	6.40	0.140
407001OBT	natural	black + plug	unscrewed	no	no	67.8	41	6.50	0.140
407008OBT	natural	black + plug	screwed	yes	no	67.8	41	6.40	0.140
407003	brown	white, with liner	unscrewed	no	no	67.8	43	6.00	0.140





Squared containers, graduated up to 2.7 I

C€ (IVD)

High density polyethylene container and cap. Ideal for 24 hour urine collection and other liquid samples.

The ribbed cap features a pharmaceutical grade liner, which **eliminates the need for an insert plug**, allows a better handling and provides a watertight seal.

Containers are graduated in 50 ml increments up to 2.7 l.

Ergonomically designed for single-handed use. Stackable.

A ring on the neck of the containers avoids dripping when pouring.

Lenght x width: 130 x 126 mm (±1 mm).

Height, cap excluded: 249 mm (±1 mm).

Individually wrapped models are supplied capped.

Cases per palet: 12.

code	body colour	cap description	presentation	individually wrapped	sterile	internal mouth Ø mm	case quantity	case weight	case volume
407005	natural	white, with liner	unscrewed	no	no	64	34	5.70	0.152
407005/T	natural	white, with liner	screwed	no	no	64	34	5.50	0.146
407005.O	natural	white, with liner	screwed	yes	RADIATION	64	34	5.70	0.152
407006	brown	white, with liner	unscrewed	no	no	64	34	5.70	0.152
407006/T	brown	white, with liner	screwed	no	no	64	34	5.70	0.152





Security screw cap containers

C€ (IVD)

Natural colour high density polyethylene containers. Capacities: 30, 60 and 125 ml. Ideal for liquids. 60 and 125 ml containers: Bottles with screw cap, leakproof closure with child-proof system and split seal.

The cap is compound by two pieces that turn around themselves avoiding the cap to be opened, working as a security closure avoiding the accidental manipulation. To open the cap push it down following the figure on the cap. The sealed caps need to be unscrewed to break the seal, showing that it has not been manipulated. 30 ml container **444603202** has a security seal, not the child-proof system closure. Pallet cases quantity: 16.

code	cap colour	body colour	height mm	Ø base mm	body weight g	case quantity	case weight	case volume			
1. Capacity: 30 r	nl. Ø Pilfer mout	th 25, security cl	osure tamper ev	ident							
444603202	white	natural	69.0	31.0	6.2	1,100	11.32	0.103			
1. Capacity: 30 r	nl. Ø Pilfer mout	th 25, security cl	osure childproof	system							
444603204	white	white	60.5	35.6	7.0	600	10.00	0.078			
2. Capacity: 60 r	nl. Ø Pilfer mout	th 28, security cl	osure childproof	system							
444603300	white	natural	93.7	37.6	9.5	600	8.00	0.110			
444602903	white	white	93.7	37.6	9.5	600	8.00	0.110			
444602901	blue	natural	87.5	37.6	9.5	600	7.50	0.110			
444602801	blue	white	87.5	37.6	9.5	600	7.50	0.110			
3. Capacity: 125	3. Capacity: 125 ml. Ø Pilfer mouth 28, security closure childproof system										
444603402	white	natural	114.8	48.0	13.6	300	5.15	0.110			

We can supply other colours if required.

25 ml cups

Ultra clear polypropylene cup with cover in natural polyethylene, graduated up to $25\ \mathrm{ml}.$

Dimensions: Ø mouth: 35 mm, Ø base: 27 mm, height 40 mm.

The cups are supplied bagged in stacks of 100 units.

code	description	case quantity	case weight	case volume
430130	standard graduated cup	60 x 100	8.60	0.065
430130.1	reinforced graduated cup	60 x 100	8.60	0.065
430100.1	push cover	8 x 600	4.90	0.034











Hinged lid container

Specially designed for **milk, wine and agroalimentary oil,** although it is also suitable for other liquid or solid samples.

Recommended for its use with FOSS® and DELTA® analytical systems, among others.

Made of **autoclavable** polypropylene. Tamper evident Lid, attached to the body by means of a hinge.

The lid is conceived to maintain 90° after its opening, enabling this pot to be used with just one hand. It has a frosted squared surface for identification. Excellent mechanical and chemical resistance.

For **coloured versions**, please ask our commercial department

Hinged lid container, 50 ml


C€ (IVD)


Graduated up to 50 ml.

Code **410046** is supplied uncapped.


Code 410056 is supplied uncapped inside an individual bag.

code	description	case quantity	case weight	case volume	cases pallet
410046	50 ml container	800	7.5	0.14	12
410056	50 ml container individually wrapped	800	7.5	0.14	12

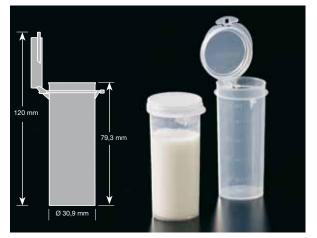

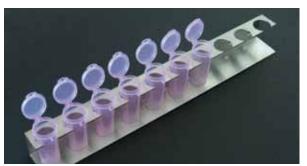

We can manufacture special racks in aluminium or stainless steel under request. Please ask the commercial department.

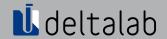






The perfect aid for a safe, easy, quick and accurate liquid handling.

- Low retention
- Autoclave resistant
- Good subjection
- High transparency
- Uniformity
- Free of metals
- Traceability

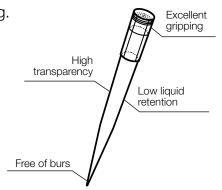


Tips with filter

The filter tip avoids aerosols and prevents pipette and / or sample contamination.

STANDARD


- 1. It purifies the air that goes through the filter, removing aerosol and contaminants
- 2. It partially absorbs the sample in case of overaspiration (10-20%)
- 3. It protects the sample from an accidentally non sterile pipette

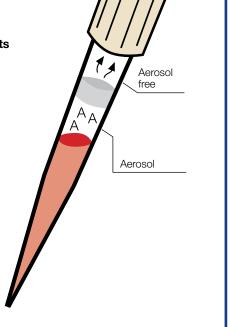

HYDROPHOBIC POLYPROPYLENE

EXTRA

- 1. It purifies the air that goes through the filter, removing aerosol and contaminants
- **2.** In case of overaspiration, **the filter compacts when touching the liquid,** and stops the pipetting
- 3. It protects the sample from an accidentally non sterile pipette

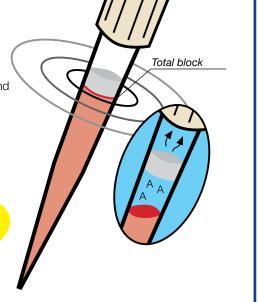

HYDROPHOBIC POLYPROPYLENE

Allows sample recovering!!

Our range of tips meet almost every liquid handling need, as they are available in a wide scope of options:

- With or without filter in bulk or rack (with hinged, sliding, or free lid).
- Sterile, DNAse and RNAse free, with Standard or Extra filter, in rack.
- Sterile, without filter, packed in individual bag.




• KIT: cardboard box with a plastic rack with hinged lid and 10 refills of 96 tips.





Our racks are made of autoclavable polypropylene, and have been developed according to an identifying colour code: racks with tips without filter are blue or translucent colour, depending on the code; racks with filter tips are red colour and racks with Extra filter tips are white.

*Excepting codes 327-26 (page 175) and 200076RF (page 178), in blue and translucent colour respectively.

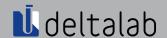



0.1 - 10 µl tip

Natural colour universal crystal tip, graduated. Gilson® type. El Rack made of ${\bf autoclavable}$ polypropylene up to ${\bf 121~^{\circ}C}$.

Racks with transparent lid and base in different colours:

Rack without filter: blue base.

STANDARD filter: red base / transparent lid. EXTRA FILTER: white base / transparent lid.

Suitable (not exclusive)

Gilson®, Eppendorf®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Nichiryo®, HTL®.

C	code	description	sterile	case quantity	case weight	case volume
•	200024	tips in bulk	no	20 x 1,000	3.80	0.026
•	199084RC	96 tips rack	no	10 x 96	1.92	0.009
•	326-19	filter tips in bulk	no	20 x 1,000	3.70	0.026
•••	327-19	96 tips rack filter	RADIATION	10 x 96	1.82	0.009
	327-19X	96 tips rack EXTRA filter	RADIATION	10 x 96	1.84	0.007

Code 326-19 minimum unit sale: 1 bag.



0.1 - 10 µl tip

Natural colour universal crystal tip, graduated. Gilson® type. Long shape, perfect to pipette samples in deep vessels.

Suitable (not exclusive)

Gilson®, Socorex®, Eppendorf®, Biohit®, Nichiryo®, Finnpipette®, Dsg®, Elkay®, Genex®, Jencons®, HTL®.

Rack made of autoclavable polypropylene up to 121 °C.


Racks with transparent lid and base in different colours:

Rack without filter: blue base. STANDARD filter: red base. EXTRA filter: white base.


code	description	sterile	case quantity	case weight	case volume
901-01	tips in bulk	no	15 x 1,000	3.00	0.026
901-01R	96 tips rack	no	10 x 96	1.80	0.009
901-01F	filter tips in bulk	no	15 x 1,000	2.92	0.016
999 301-01FR	96 tips rack filter	RADIATION	10 x 96	1.74	0.009
●●● 301-01FRX	96 tips rack EXTRA filter	RADIATION	10 x 96	1.80	0.009

Tips in bulk: minimum unit sale: 1 bag.

Code **301-01F** suitable for Eppendorf® Research Plus.

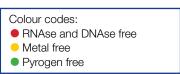



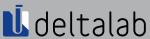



0.5 - 10 µl tip

Natural colour universal crystal tip, bevelled. Eppendorf® type.

Suitable (not exclusive)

Eppendorf®, Gilson®, Biohit®, Brand®, Dsg®, Elkay®, Finnpipette®, Jencons®, Nichiryo®, Socorex®, HTL®.

Rack made of **autoclavable** polypropylene up to **121 °C**.

Racks with transparent lid and base in different colours:

Rack without filter: blue base; STANDARD filter: red base; EXTRA filter: white base;

	code	description	sterile	case quantity	case weight	case volume
•	301-11	tips in bulk	no	25 x 1,000	4.80	0.050
•	301-10	96 tips rack	no	10 x 96	1.82	0.009
•••	301-09	96 tips rack	RADIATION	10 x 96	1.80	0.009
•	326-37	filter tips in bulk	no	25 x 1,000	5.40	0.035
•••	327-37R	96 tips rack filter	RADIATION	10 x 96	1.82	0.009
•••	327-37RX	96 tips rack EXTRA filter	RADIATION	10 x 96	1.80	0.009

Tips in bulk: minimum unit sale: 1 bag.



2 - 20 µl tip

Natural colour universal tip. Eppendorf® type. Available without filter (2-200µl).

See our codes 200072.

Suitable (not exclusive)

Eppendorf®, Biohit®, Finnpipette®, Gilson®, Socorex®, Brand®, HTL®.

Rack made of **autoclavable** polypropylene up to **121 °C**.


Racks with transparent lid and base in different colours:

Rack without filter: blue base.

STANDARD filter: red base.

code	description	sterile	case quantity	case weight	case volume
••• 200010FR	96 tips rack filter	RADIATION	10 x 96	2.0	0.009

Tips in bulk: minimum unit sale: 1 bag.

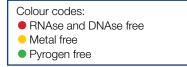

2 - 30 µl tip

Natural colour universal filter tip, bevelled. Gilson® type.

Suitable (not exclusive)

Gilson® P20 P100 P200, Eppendorf®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Socorex®, Nichiryo®, HTL®, etc.

Rack made of **autoclavable** polypropylene up to **121 °C**.

Racks with transparent lid and base in different colours:

Rack with STANDARD filter: red base.

	code	description	sterile	case quantity	case weight	case volume
•	327-36	filter tips in bulk	no	25 x 1,000	9.5	0.071
•••	327-36R	96 tips rack filter	RADIATION	10 x 96	2.0	0.009

Tips in bulk: minimum unit sale: 1 bag.



2 - 100 µl tip

Natural colour universal filter tip, bevelled. Gilson® type.

Suitable (not exclusive)

Gilson® P100 P200, Eppendorf®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Nichiryo®, HTL®, etc.

Rack made of autoclavable polypropylene up to 121 °C.

Racks with transparent lid and base in different colours:

Rack with STANDARD filter: red base.

EXTRA filter: white base.

	code	description	sterile	case quantity	case weight	case volume
•	327-23	filter tips in bulk	no	25 x 1,000	9.0	0.071
•••	327-24	96 tips rack filter	RADIATION	10 x 96	2.0	0.009
•••	327-24X	96 tips rack EXTRA filter	RADIATION	10 x 96	2.1	0.009

Tips in bulk: minimum unit sale: 1 bag.



Colour codes:

RNAse and DNAse free

Metal free





2 - 160 µl tip

Natural colour universal filter tip, bevelled. Gilson® type. Identical tip as code **200078**, but having the filter placed up to 160 μ l.


Suitable (not exclusive)

Gilson® P20 P100 P200, Eppendorf®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Socorex®, Nichiryo®, HTL®, etc.

Rack made of **autoclavable** polypropylene up to **121 °C**. Racks with transparent lid and base in different colours:

Rack with STANDARD filter: red base.

d	code	description	sterile	case quantity	case weight	case volume
	327-15	96 tips rack filter	RADIATION	10 x 96	2	0.009

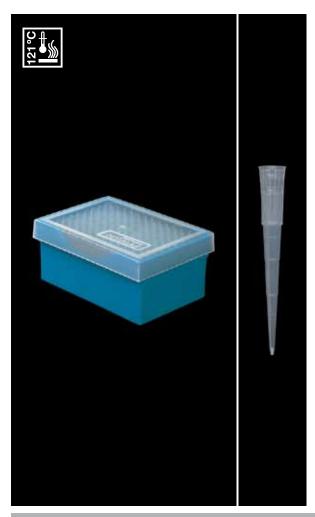

2 - 200 µl tip

Natural colour universal tip, bevelled. Gilson type[®].

Suitable (not exclusive)

Gilson® P20 P100 P200, Eppendorf®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Socorex®, Nichiryo®, HTL®, etc.

Rack made of ${\bf autoclavable}$ polypropylene up to ${\bf 121~^{\circ}C};$ with transparent lid and blue base.

	code	description	sterile	case quantity	case weight	case volume
•	200078	tips in bulk	no	25 x 1,000	9.2	0.071
•	199078RC	96 tips rack	no	10 x 96	1,96	0.009
•••	199078N	96 tips rack	RADIATION	10 x 96	1,96	0.009

Tips in bulk: minimum unit sale: 1 bag.

Colour codes:

RNAse and

RNAse and DNAse free

Metal free

2 - 200 (100) µl tip

Natural colour universal tip. Eppendorf® type.

Suitable (not exclusive)

Eppendorf®, Gilson®, Brand®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Kartell®, Nichiryo®, Socorex®, HTL®.

Rack made of **autoclavable** (up to 121 °C) polypropylene.

Racks with transparent lid and base in different colours:

Rack without filter: blue base.

STANDARD filter: red base.

code	description	sterile	case quantity	case weight	case volume			
2 - 200 µl tip	2 - 200 μl tip							
200072	tips in bulk	no	25 x 1,000	9.1	0.071			
99014RC	96 tips rack	no	10 x 96	1.9	0.009			
99014N	96 tips rack	RADIATION	10 x 96	2	0.009			
2 - 100 µl filter tip								
927-33	tips in bulk	no	25 x 1,000	10	0.071			
999 327-33R	96 tips rack filter	RADIATION	10 x 96	1.97	0.009			

Tips in bulk: minimum unit sale: 1 bag.



2 - 200 µl tip

Universal tip, bevelled. Gilson® type. Tips with filter: natural colour. Tips without filter: yellow colour.

Suitable (not exclusive)

Gilson®, Eppendorf®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Nichiryo®, HTL®, Socorex® serie Akura.

Rack made of autoclavable (up to 121 °C) polypropylene.

Racks with transparent lid and base in different colours:

Rack without filter: blue base. STANDARD filter: red base. EXTRA filter: white base.

	code	description	sterile	case quantity	case weight	case volume			
Natura	Natural colour:								
•	327-34	filter tips in bulk	no	25 x 1,000	9.2	0.083			
•••	327-34R	96 tips rack filter	RADIATION	10 x 96	1.98	0.009			
•••	327-34RX	96 tips rack EXTRA filter	RADIATION	10 x 96	1.98	0.009			
Yellow	colour:								
•	200016	tips in bulk	no	25 x 1,000	8.5	0,071			
•	200016B	single tips individually wrapped	RADIATION	200	0,19	0.002			
	200016R	96 tips rack	no	10 x 96	1.94	0.009			

Code **327-34** minimum unit sale: 1 bag.




RNAse and DNAse free

Metal free





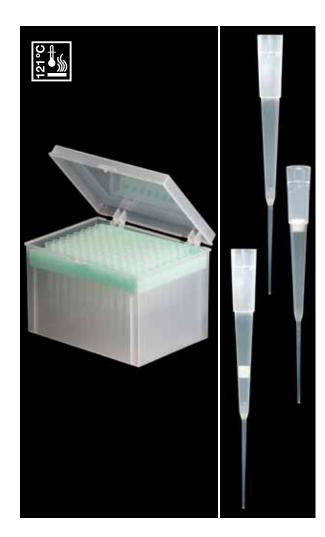

Universal natural colour tip with collar for gel.

Suitable (not exclusive)

Pipetman®, Eppendorf®, Titertek®, Biohit®, Socorex®, Nichiryo®, Costar Corning®, Oxford®, Rainin®, Finnpipette®.

Translucent rack made of **autoclavable** polypropylene with hinged lid.

code		description	sterile	case quantity	case weight	case volume
•	200079	tips in bulk	no	10 x 1,000	9.2	0.050
•••	200079R	96 tips rack	RADIATION	10 x 96	1.2	0.010
•••	200079F	96 tips rack filter 1-20 μl	RADIATION	10 x 96	1.8	0.009
•••	200079FS	96 tips rack filter 1-100 µl	RADIATION	10 x 96	1.8	0.009

Tips in bulk: minimum unit sale: 1 bag.



1 - 200 µl tip

Universal natural colour tip with 91 mm length.


Suitable (not exclusive)

Pipetman®, Eppendorf®, Nichiryo®, Costar Corning®, Oxford®, Volac®, Socorex®, Titertek®, SMI®, Rainin®, Finnpipette®, etc.

Rack made of polypropylene **autoclavable** up to **121 °C**. Transparent lid and blue base, with 204 tips in 12×17 distribution.

code	description	sterile	case quantity	case weight	case volume
0 200026	tips in bulk	no	5 x 1,000	3.50	0.050
327-26	204 tips rack filter	RADIATION	8 x 204	2.03	0.013

Tips in bulk: minimum unit sale: 1 bag.

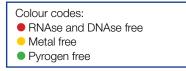

Natural colour universal tip. Designed to work with viscous fluids like glycerol or thick tissue homogenates. Suitable also for fragile cellular structures, minimizing the risk of tearing or rupture.

Suitable (not exclusive)

Eppendorf®, Socorex®, Nichiryo®, Finnpipette®, Biohit®, Oxford®, Pipetman®, Rainin®, Costar Corning®, Volac®.

Rack made of polypropylene **autoclavable** up to **121 °C**. Transparent lid and blue base.

code	description	sterile	case quantity	case weight	case volume
927-40	tips in bulk	no	10 x 1,000	4.20	0.020

Tips in bulk: minimum unit sale: 1 bag.



5 - 200 µl tip

Eppendorf® universal tip.

Suitable (not exclusive)

Eppendorf®, Gilson®, Brand®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Nichiryo®, Socorex®, HTL®.

Rack made of ${\bf autoclavable}$ polypropylene up to ${\bf 121~^{\circ}C}.$

Transparent lid and blue base.


	code	description	sterile	case quantity	case weight	case volume
•	200009	tips in bulk	no	25 x 1,000	8.82	0.071
•	200009B	single tips individually wrapped	RADIATION	200	0.19	0.002
•	200009R	96 tips rack	no	10 x 96	1.98	0.009



Colour codes:

RNAse and DNAse free

Metal free

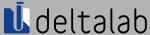



Yellow universal tip, Gilson® type.

Suitable (not exclusive)


Gilson®, Eppendorf®, Biohit®, Dsg®, Elkay®, Finnpipette®, Genex®, Jencons®, Nichiryo®, HTL®.

Rack made of polypropylene **autoclavable** up to **121 °C**. Transparent lid and blue base.

code		description	sterile	case quantity	case weight	case volume
•	199083	tips in bulk	no	25 x 1,000	9.0	0.071
•	199620	96 tips rack	no	10 x 96	1.96	0.010
NEW	KIT199083R	kit: rack + 10 refills of 96 tip	no	1	1.00	0.0089

Tips in bulk: minimum unit sale: 1 bag.



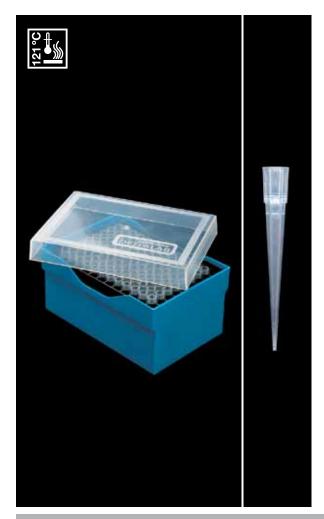

2 - 300 µl tip

Natural colour universal tip. $Biohit^{\otimes}$ type.

Suitable (not exclusive)

Gilson® P200, Finnpipette® (old model), Labsystems®, etc.

Rack made of ${\bf autoclavable}$ polypropylene up to ${\bf 121~^{\circ}C}.$ Transparent lid and blue base.

	code	description	sterile	case quantity	case weight	case volume
•	200074	tips in bulk	no	24 x 1,000	8.3	0.067
•	199076RC	96 tips rack	no	10 x 96	1.8	0.009

Tips in bulk: minimum unit sale: 1 bag.



RNAse and DNAse free

Metal free

Yellow universal tip with crown and with a 50 µl margin to prevent any accidental overaspiration. Filter tip is neutral instead of yellow.

Translucent rack made of autoclavable polypropylene up to 121 °C; with hinged lid.

Suitable (not exclusive)

Eppendorf®, Socorex®, Gilson® except NEO model, Nichiryo®, Biohit®, Brand®, Jencons®, Thermo®, DSG®, Elkay®, HTL®, Finnpipette®, Genex®,

code		description	sterile	case quantity	case weight	case volume
	200076	tips in bulk	no	24 x 1,000	8.82	0.071
	200076R	tips in rack	no	10 x 96	1.98	0.009
NEW	KIT200076R	kit: rack + 10 refills of 96 tip	no	1	1.00	0.0089
	200076RF	filter tips in rack	RADIATION	10 x 96	2.04	0.009

Tips in bulk: minimum unit sale: 1 bag.





5 - 200 µl tip


Natural colour universal tip, bevelled. MLA® type.

Translucent rack made of autoclavable polypropylene up to 121 °C; with hinged lid.

	code	description	sterile	case quantity	case weight	case volume
•	200017	tips in bulk	no	25 x 1,000	10.5	0.071
•	200017R	96 tips rack	no	10 x 96	1.8	0.010

Tips in bulk: minimum unit sale: 1 bag.



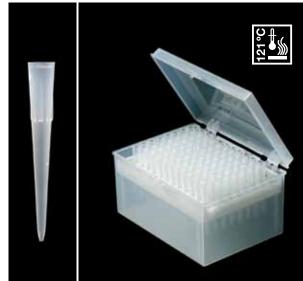

5 - 200 µl tip


Natural colour universal tip, bevelled. Oxford® type.

Suitable for the DiaMed-Micro typing System machine.

Translucent rack made of autoclavable polypropylene up to 121 °C; with hinged lid.


	code	description	sterile	case quantity	case weight	case volume
•	200019	tips in bulk	no	12 x 1,000	6.9	0.071
•	200019R	96 tips rack	no	10 x 96	1.8	0.010

Tips in bulk: minimum unit sale: 1 bag.

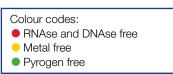





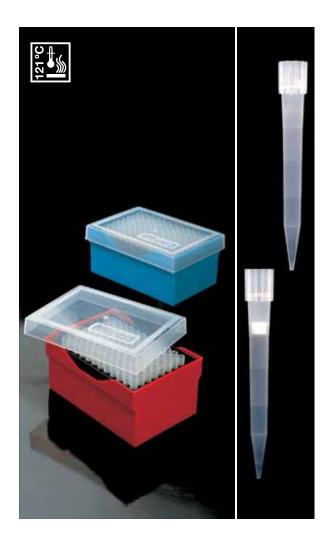

20 - 300 µl tip

Natural colour universal tip. Eppendorf® type.

Suitable (not exclusive)

Eppendorf®, Gilson®, Brand®, Biohit®, Elkay®, Jencons®, Finnpipette®, Nichiryo®, Socorex®.

Rack made of **autoclavable** (up to **121 °C**) polypropylene.

Racks with transparent lid and base in different colours:

Rack without filter: blue base. STANDARD filter: red base.

code		description	sterile	case quantity	case weight	case volume
•	200077	tips in bulk	no	20 x 1,000	9.7	0.071
•	200077R	96 tips rack	no	10 x 96	1.94	0.009
•	200077F	filter tips in bulk	no	20 x 1,000	9.9	0.071
	200077FR	96 tips rack filter	RADIATION	10 x 96	1.92	0.009

Tips in bulk: minimum unit sale: 1 bag.



100 - 1,000 µl tip

Universal tip with collar. Blue or natural colour.

Suitable (not exclusive)

Eppendorf®, Gilson®, Biohit®, DSG®, Elkay®, Finnpipette®, Genex®, Jencons®, Nichiryo®. Suitable for Menarini®.

Rack made of polypropylene, **autoclavable** up to **121 °C**.

Racks with transparent lid and base in different colours depending on the model.

Tips without filter: blue base. STANDARD filter: red base. EXTRA filter: white base.

	code	description	sterile	case quantity	case weight	case volume		
Blue c	olour:							
	200012	tips in bulk	no	10 x 1,000	8.4	0.085		
•	200012B	single tips individually wrapped	RADIATION	175	0.27	0.0035		
	200012R	96 tips rack	no	10 x 96	2.36	0.009		
Natura	Natural colour:							
•	200070	tips in bulk	no	10 x 1,000	8.0	0.083		
•	200070R	96 tips rack	no	10 x 96	2.4	0.009		
•••	200070NR	96 tips rack	RADIATION	10 x 96	2.38	0.010		
•	327-30	filter tips in bulk	no	10 x 1,000	3.5	0.083		
•••	327-30R	96 tips rack filter	RADIATION	10 x 96	2.38	0.013		
•••	327-30RX	96 tips rack EXTRA filter	RADIATION	10 x 96	2.4	0.015		

Codes 200070 and 327-30 minimum unit sale: 1 bag.

Colour codes:

- RNAse and DNAse free
- Metal free
- Pyrogen free

100 - 1,000 µl tip

Gilson® Universal tip.

Suitable (not exclusive)

Gilson®, Eppendorf®, Biohit®, DSG®, Elkay®, Finnpipette®, Genex®, Jencons®, Nichiryo®, Socorex®, HTL®.


Rack made of polypropylene, autoclavable up to 121 °C.

Racks with transparent lid and base in different colours depending on the model.

Tips without filter: blue base. STANDARD filter: red base.

code	description	sterile	case quantity	case weight	case volume				
Blue colour:									
9 200080.1	tips in bulk	no	10 x 1,000	9.86	0.082				
9 200080.1R	96 tips rack	no	10 x 96	2.62	0.015				
Natural colour:	Natural colour:								
9 200082	tips in bulk	no	10 x 1,000	10.8	0.083				
900082R	96 tips rack	no	10 x 96	2.66	0.015				
• 200082NR	96 tips rack	RADIATION	10 x 96	2.62	0.013				
900082F	filter tips in bulk	no	10 x 1,000	11.2	0.083				
327-16	96 tips rack filter	RADIATION	10 x 96	2.64	0.013				
••• 327-16X	96 tips rack EXTRA filter	RADIATION	10 x 96	2.7	0.013				

Tips in bulk: minimum unit sale: 1 bag.

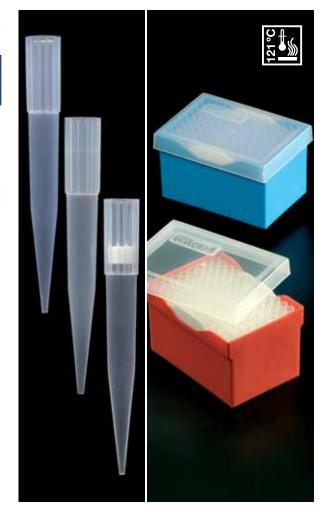

100 - 1,000 µl tip

Universal tip with collar. Natural colour, Eppendorf® type.

Suitable (not exclusive)

Eppendorf®, Nichiryo®, Biohit®, Brand®, Dsg®, Gilson®, Elkay®, Finnpipette®, Genex®, Jencons®, HTL®. Suitable for Menarini®.

Translucent rack made of ${\bf autoclavable}$ polypropylene up to ${\bf 121}~{}^{\bf c}{\bf C};$ with hinged lid.

	code	description	sterile	case quantity	case weight	case volume
	200029	tips in bulk	no	10 x 1,000	8.1	0.083
•	200029R	96 tips rack	no	10 x 96	2.4	0.010
•	200029F	filter tips in bulk	no	10 x 1,000	8.1	0.083
	200029RF	96 tips rack filter	RADIATION	10 x 96	2.4	0.013

Tips in bulk: minimum unit sale: 1 bag.




Colour codes:

RNAse and DNAse free

Metal free

Pyrogen free

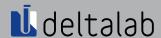



5 - 500 µl tip

Universal natural colour tip. Beckmann® type.

For use (but not limited to)

Sherwood®, Lancer®.

code	description	sterile	case quantity	case weight	case volume
000021	tips in bulk	no	10 x 1,000	6.9	0.071

Tips in bulk: minimum unit sale: 1 bag.



100 - 1,000 μl tip

Universal tip. Oxford® type.

For use (but not limited to)

Finnpipette®, Labsystem®.

	code	description	sterile	case quantity	case weight	case volume
•	200020	tips in bulk	no	8 x 1,000	10.40	0.071

Tips in bulk: minimum unit sale: 1 bag.



100 - 1,000 µl tip

Natural colour universal tip special for MLA®. Bevelled.

Suitable (not exclusive)

Finnpipette® colour.

code	description	sterile	case quantity	case weight	case volume
200018	tips in bulk	no	8 x 1,000	2.75	0.071

Tips in bulk: minimum unit sale: 1 bag.





100 - 1,000 µl tip

Natural colour universal tip, bevelled and graduated.

For use (but not limited to)

Gilson®, Rainin®, Finnpipette®, Nichiryo®, Eppendorf®, Biometra®.

code	description	sterile	case quantity	case weight	case volume
200022	tips in bulk	no	10 x 1,000	10.9	0.080

Tips in bulk: minimum unit sale: 1 bag.



Colour codes:

- RNAse and DNAse free
- Metal free
- Pyrogen free

Pipetting reservoir

Reservoir specially designed for use with multichannel pipettes. Made of polystyrene. Code **M-210** is supplied sterile and individually wrapped.

code	description	sterile	case quantity	case weight	case volume
M-200	150 ml	no	50	0.73	0.0095
M-210	150 ml	RADIATION	50	0.73	0.0095
M-203	25 ml	no	5 x 10	12.58	0.098



See our plates for liquid handling on chapter 2. **Molecular Biology**

Pipetting reservoir

8 and 12 channel reservoir, 40 ml, PS, sterile. Individually wrapped.

Dimensions:

330089: 83 x 126 x 13 mm **330129**: 126 x 83 x 13 mm

code	description	sterile	case quantity	case weight	case volume
330089	8 channel. 40 ml	RADIATION	10	0.23	0.0035
330129	12 chanel. 40 ml	RADIATION	10	0.23	0.0035





1 - 5 ml macrotip

Natural colour universal macro tip, with collar. Eppendorf $\ensuremath{^{\$}}$ type.

For use (but not limited to)

Biohit®, Dsg®, Eppendorf®, Genex®, Jencons®, Socorex® Acura 835.

Rack made of polypropylene **autoclavable** up to **121 °C**. Polycarbonate transparent lid and marble base. Sliding lid.

code		description	sterile	case quantity	case weight	case volume
•	200075	tips in bulk	no	8 x 250	5.9	0.084
•	199075	50 tips rack	no	4 x 50	1.98	0.015
	200075RF	50 tips rack filter	RADIATION	4 x 50	2.5	0.015



1 - 5 ml macrotip

Natural colour universal macro tip with collar. Socorex® type.

For use (but not limited to)

Nichiryo®, Oxford®.

code	description	sterile	case quantity	case weight	case volume
831-312	tips in bulk	no	12 x 250	9.76	0.096
831-312B	single tips individually wrapped	RADIATION	75	0.3	0.003

Colour codes:

- RNAse and DNAse free
- Metal free
- Pyrogen free



1 - 5 ml macrotip

Blue colour universal macro tip with collar. Socorex $\ensuremath{^{^{\$}}}$ type.

For use (but not limited to)

Gilson®, Nichiryo®.

Rack made of polypropylene ${\bf autoclavable}$ up to ${\bf 121~^{\circ}C}.$

Transparent lid and marble base.

Sliding lid.

code 200083		description	sterile	case quantity	case weight	case volume
•	200083	tips in bulk	no	8 x 250	5.78	0.085
•	200083R	50 tips rack	no	4 x 50	2.1	0.015



Natural colour universal macro tip with collar. Finnipipette® type.

For use (but not limited to)

Brand®, Labsystem®.

Rack made of polypropylene autoclavable up to 121 °C.

Transparent lid and marble base.

Sliding lid.

code	description	sterile	case quantity	case weight	case volume
0 200028	tips in bulk	no	8 x 250	6.76	0.085
900028R	50 tips rack	no	4 x 50	2.2	0.015
•0• 200028RF	50 tips rack filter	RADIATION	4 x 50	2.50	0.015



RNAse and DNAse free

Metal free

Pyrogen free





2 - 10 ml macrotip

Natural colour universal macro tip with collar. Gilson® type.

For use (but not limited to)

Labsystem®, Socorex®.

Rack made of polypropylene **autoclavable** up to **121 °C**. Transparent lid and marble base. Sliding lid.

code		description	sterile	case quantity	case weight	case volume
200015		tips in bulk	no	10 x 100	7.24	0.085
o 200015R		25 tips rack	no	4 x 25	2.25	0.015
	200015RF	25 tips rack filter	RADIATION	4 X 25	2.2	0.015

- Colour codes:
- RNAse and DNAse free
- Metal free
- Pyrogen free
- Please contact our commercial department to get more information about our wide range of Socorex® products:
- * Pipettes: adjustable and fixed volume; digital, analogical and electronical; multichannel, etc.
- * Microdispensers, table top dispensers and compact dispensers
- * Self-refilling syringes



FILTER TIPS DS-BARRIER

Filter tips specially designed for PCR and other amplification techniques or for preventing aerosol contamination from infectious materials.

Tips are made of low retention virgin polypropylene with a filter made of ultrahigh molecular weight polyethylene with a pore smaller than 20 microns.

DS-Barrier filter tips racked on hinge racks, sterilized by ionizing radiation, **DNase, Rnase, Metals and pyrogens free**, or in a thermosealed bag non sterile.

DS-Barrier racks are **autoclavable at 121 °C** for 15 minutes and are identified by color code (every tip size has a different rack color for easy identification).









DS-Barrier Filter tip for volumes between 10 and 20 μl for micro tip cone pipettes. Universal compatibility.

Designed to fit with in most market pipettes.

code	description	volume µl	presentation	sterile	case quantity	case weight	case volume
249010	fiter tip	10	rack 96 u.	RADIATION	10 x 96	1.23	0.009
169010	fiter tip	10	bag 960 u.	no	72 x 960	13.4	0.144
249021	fiter tip	20	rack 96 u.	RADIATION	10 x 96	1.25	0.009
169021	fiter tip	20	bag 960 u.	no	72 x 960	13.6	0.144



DS-Barrier filter tip for volumes between 10 and 100 μ l for standard tip cone pipettes. Universal compatibility.

Designed to fit with in most market pipettes.

code	description	volume µl	presentation	sterile	case quantity	case weight	case volume
249011	fiter tip	10	rack 96 u.	RADIATION	10 x 96	1.31	0.009
169011	fiter tip	10	bag 960 u.	no	40 x 960	13.2	0.144
249020	fiter tip	20	rack 96 u.	RADIATION	10 x 96	1.31	0.009
169020	fiter tip	20	bag 960 u.	no	40 x 960	13.2	0.144
249100	fiter tip	100	rack 96 u.	RADIATION	10 x 96	1.31	0.009
169100	fiter tip	100	bag 960 u.	no	40 x 960	13.2	0.144

DS-Barrier 200 - 300 µl

DS-Barrier filter tip for volumes between 200 and 300 μl for standard tip cone pipettes. Universal compatibility.

Designed to fit with in most market pipettes.

code	description	volume µl	presentation	sterile	case quantity	case weight	case volume
249200	fiter tip	200	rack 96 u.	RADIATION	10 x 96	1.35	0.009
169200	fiter tip	200	bag 960 u.	no	32 x 960	11.25	0.144
249300	fiter tip	300	rack 96 u.	RADIATION	10 x 96	1.35	0.009
169300	fiter tip	300	bag 960 u.	no	32 x 960	11.25	0.144

DS-Barrier 1.000 µl


DS-Barrier filter tip for volumes to 1000 μ l for standard tip cone pipettes. Universal compatibility.

Designed to fit with in most market pipettes.

code	description	volume µl	presentation	sterile	case quantity	case weight	case volume
249001	fiter tip	1,000	rack 60 u.	RADIATION	8 x 60	1.42	0.009
169001	fiter tip	1,000	bag 960 u.	no	12 x 960	13.52	0.144







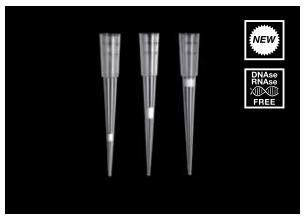





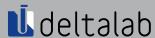

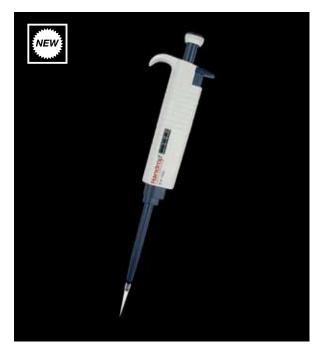

HANDROP DIGITAL - SC

Single channel variable volume digital pipettes with digital display. Precise volume selector.

Easy calibración and maintenance.

Lightweight and ergonomic design with tip ejector.

code	volume range µl	increments µl	accuracy %	precision %	case weight	case volume
770005	0.1 - 2.5	0.05	2.5 - 12	2 - 5	0.21	0.0009
770000	0.5 - 10	0.1	1.0 - 2.5	0.8 - 1.5	0.21	0.0009
770080	2 - 20	0.5	0.9 - 3.0	0.4 - 2.0	0.21	0.0009
770050	10 - 100	1	0.8 - 3.0	0.15 - 1.5	0.21	0.0009
770070	20 - 200	1	0.6 - 3.0	0.15 - 1.0	0.21	0.0009
770060	100 - 1,000	5	0.6 - 2.0	0.2 - 0.7	0.21	0.0009
770040	200 - 1,000	5	0.6 - 0.9	0.2 - 0.3	0.21	0.0009
770110	1 - 5 ml	0.05 ml	0.5 - 0.7	0.15 - 0.3	0.21	0.0009



HANDROP DIGITAL - MC

8 and 12 channel variable volume digital pipettes with digital display. Precise volume selector.

Easy calibración and maintenance.

Lightweight and ergonomic design with tip ejector.

code	volume range µl	increments µI	accuracy %	precision %	case weight	case volume
8 chanels						
770210	0,5 - 10	0,1	1,0 - 4,0	1,5 - 4,0	0,34	0.0017
770220	5 - 50	0,5	1,0 - 3,0	0,5 - 2,0	0,34	0.0017
770240	50 - 300	5	0,7 - 1,5	0,25 - 0,8	0,34	0.0017
12 chanels	;					
770310	0,5 - 10	0,1	1,0 - 4,0	1,5 - 4,0	0,34	0.0017
770320	5 - 50	0,5	1,0 - 3,0	0,5 - 2,0	0,34	0.0017
770340	50 - 300	5	0,7 - 1,5	0,25 - 0,8	0,34	0.0017



HANDROP STEPPER

Single channel stepper.

Allows to dispense up to 40 different volumes and up to 48 doses.

Compatible with tips on the next page.

Volumes from 1 µl to 5 ml can be dispensed.

code	description	case weight	case volume
721008	Handrop Stepper	0.22	0.0013



Ripette Pro® mechanical stepper

Very practical to use mechanical stepper.

Compatible with 12 tip sizes and 10 volume adjustments that offers up to 120 program steps.

It is a very easy-to-use product because volume adjustments and dosing volumes for different tips are shown in a vision panel.

Moreover, the dosing button is very easy-moving to dispense the liquid conveniently.

It has an ergonomic and light design.

Compatibility with most market tips.

code	description	case weight	case volume
316N	mechanical stepper	0.30	0.002



Ritips® professional

Different capacities. Black graduations.


Each box of 25 ml and 50 ml contains a free of charge adapter.

For use (but not limited to)

Biohit® eLine® Lite and Pro, Brand® HandyStep electronic, Eppendorf® Multipette Plus, Eppendorf® Repeater Plus, Multipette Pro, and Multipette 4780, Nichimate® Stepper (Nichiryo®), Gilson® Distriman, StepMate®, Ripette®, etc.

When used with Multipette plus and Multipette Pro the display dispenser is inoperative.

	code	description	case quantity	case weight	case volume
	317-314.1	Combitip plus 1.0 ml	100	0.365	0.006
ERILE	317-315.1	Combitip plus 2.5 ml	100	0.478	0.006
띮	317-316.1	Combitip plus 5 ml	100	0.511	0.006
N STI	317-317.1	Combitip plus 10 ml	100	0.617	0.006
8	317-318.1	Combitip plus 25 ml	25	0.413	0.006
	317-319.1	Combitip plus 50 ml	25	0.511	0.006

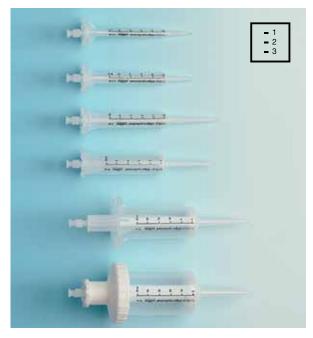

Ritips® standard


Different capacities, standard and sterile models. Red graduations.

For use (but not limited to)

Eppendorf® Multipette 4780, Brand® Handystep and Handystep Electronic, Nichimate® Stepper, Ripette® and Ripette® GenX, Minilab® 100/101, Stepmate®, Easystep®, Gilson® Distriman®, etc.

Sterile models are **Rnase**, **Dnase** and pyrogen free. Individually wrapped.

	code		description	case quantity	case weight	case volume
빌	316-314.1		Combitip 1,25 ml	100	0.447	0.006
STERILE	316-315.1		Combitip 2,5 ml	100	0.532	0.006
S	316-316.1		Combitip 5 ml	100	0.690	0.006
NON	316-317.1		Combitip 12,5 ml	100	0.874	0.006
	316-314.2	••	Combitip 1,25 ml	100	0.610	0.009
RILE	316-315.2	••	Combitip 2,5 ml	100	0.660	0.009
STE	316-316.2	••	Combitip 5 ml	100	0.820	0.009
	316-317.2	••	Combitip 12,5 ml	100	1.020	0.009

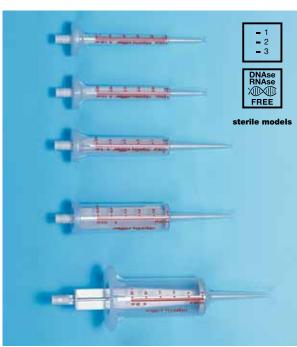








POLYETHYLENE Transfer Pipettes

- Sterile and non-sterile
- Low-affinity surface
- Attached bulb avoids cross-contamination
- One-step biohazard disposal
- Shatterproof, low density polyethylene which is non-toxic and inert
- Can be frozen in liquid nitrogen

You can ask for our pipettes in three kinds of presentation:

- Bag made (steriles or not): of polyethylene (picture 1).
- Peel-pack (steriles): packed individually with its code, lot number and expiration date printed on the paper (picture 2).
- Flow-pack (steriles): packed individually with its code, lot number and expiration date printed on polyethylene bag (picture 3).

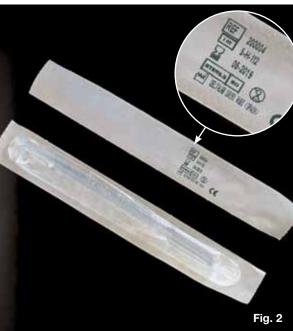

Used in Lab Procedures with

- Beakers
- Volumetric flasks
- Vials
- Centrifuge tubes
- Cryogenic vials
- Dishes
- Microtiter plates
- Microscope slides

Life Sciences Applications


- Cells in flasks and tubes: large volume, tong pipettes
- Isolating Cells: large aperture pipettes
- Control techniques: narrow stem and fine tip

Industrial Research and Testing

- Chromatography
- Viscous liquid samples and transfers
- Field Testing
- Spectrophotometry Cuvettes
- Reagent Handling

Clinical Applications

- Blood Banking
- Small volume dispensing
- Sedimentation and Storage
- Bottle Sampling
- Slides

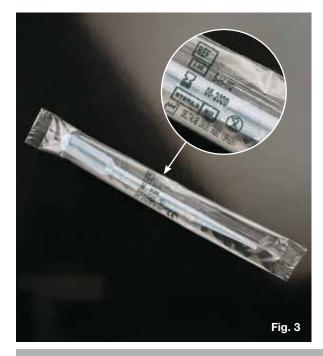





Micropipettes

Made in transparent, inert and low density polyethylene. They are used as volume measurers, to incorporate in kits of reagents, pregnancy tests, etc. They can be frozen in liquid nitrogen.

Before use, press the upper bulb, filling all the stem, and the rest remains in the inferior bulb, which it is not due to give. Volumetrical tolerance: ± 10 %.

mod.	code	presentation	lenght mm	dispense volume µl	case weight	case volume
Α	201110	500	83,0	100	0.30	0.002
В	201250	500	102,0	250	0.39	0.005



Special purpose pipettes

Made in transparent, inert and low density polyethylene. They can be frozen in liquid nitrogen. Volumetrical tolerance: \pm 10 %.

mod.	code	presentation	lenght mm	Ø stem mm	total capacity ml	bulb draw µl	case weight	case volume
Α	210030	8 x 500	130	3.8	0.5	0.2	2.9	0.023
В	210005	500	50	1.19	1.3	1.0	1.9	0.022






Graduated pipettes

Made in transparent, inert and low density polyethylene.

Sterile models are ethylene oxide sterilized (except codes **200037** and **200034** that are radiation sterilized). Volumetrical tolerance: ±10%.

mod.	code	presentation	sterile	lenght mm	Ø stem mm	graduated until ml	total capacity ml	bulb draw ml	graduations ml	case weight	case volume	pallet quantity
Α	200002	10 x 500	no	150	7.0	2	5.0	1.8	0.50	7.86	0.099	100,000
В	200006	10 x 500	no	154	7.8	3	7.0	3.2	0.50	9.50	0.098	100,000
В	200062.0	10 x 500 bag 20 u.	OXIDE	160	7.8	3	7.0	3.2	0.50	12.60	0.130	35,000
В	200007	6 x 500 individually wrapped peel-pack	OXIDE	154	7.8	3	7.0	3.2	0.50	7.20	0.165	27,000
В	200037	6 x 500 individually wrapped flow-pack	RADIATION	154	7.8	3	7.0	3.2	0.50	7.20	0.162	27,000
С	200003	6 x 500	no	154	5.0	1	5.0	3.1	0.25	4.88	0.058	90,000
С	200032.0	6 x 500 bag 20 u.	OXIDE	158	5.0	1	5.0	3.1	0.25	6.48	0.067	48,000
С	200004	6 x 500 individually wrapped peel-pack	OXIDE	154	5.0	1	5.0	3.1	0.25	7.20	0.165	27,000
С	200034	6 x 500 individually wrapped flow-pack	RADIATION	154	5.0	1	5.0	3.1	0.25	7.40	0.162	27,000



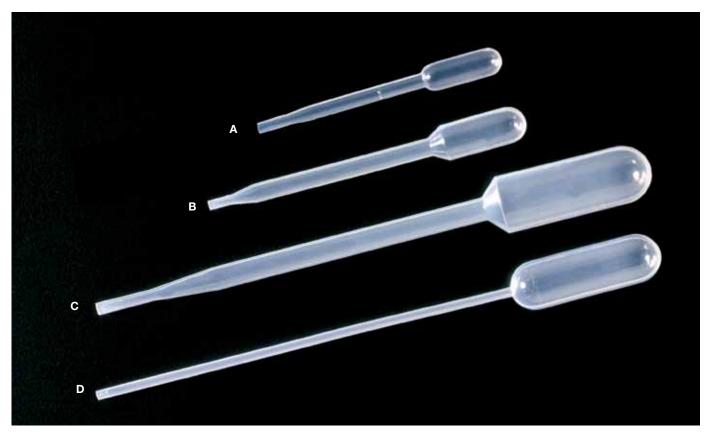

General purpose pipettes

Made in transparent polyethylene.


Model 200038 is sterile by radiation and code 200008 is sterile by ethylene oxide.

Volumetrical tolerance: ±10 %.

mod.	code	presentation	sterile	lenght mm	Ø stem mm	total capacity ml	bulb draw ml	case weight	case volume
Α	200121	10 x 1,000	no	67	3.9	1.1	0.5	4.2	0.051
В	200001	8 x 500	no	87	4.7	1.7	0.9	3.5	0.038
С	200000	10 x 400	no	150	6.3	8.0	4.6	8.7	0.096
D	200005	6 x 500	no	150	2.5	4.0	3.1	4.2	0.057
D	200008	6 x 500 individually wrapped peel-pack	OXIDE	150	2.5	4.0	3.1	13.0	0.162
D	200038	6 x 500 individually wrapped flow-pack	RADIATION	150	2.5	4.0	3.1	13.0	0.162





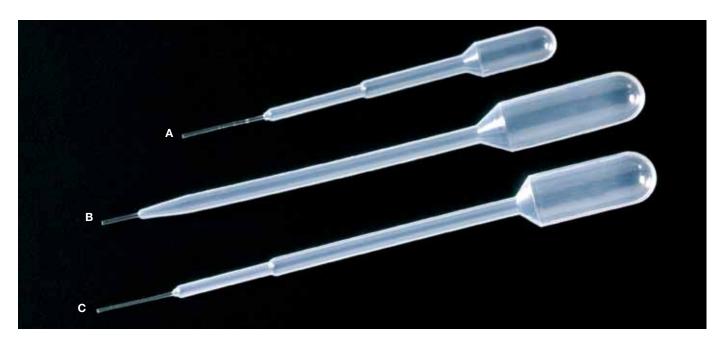

General purpose pipettes

Made in transparent an inert polyethylene.

Sterile models are sterilized by ethylene oxide. Volumetrical tolerance: ± 10 %.


mod.	code	presentation	sterile	lenght mm	Ø stem mm	total capacity ml	bulb draw ml	case weight	case volume
Α	210003	8 x 400	no	104	3.0	1.5	1.0	2.4	0.022
А	210023	500 individually wrapped peel-pack	OXIDE	104	3.0	1.5	1.0	1.38	0.021
В	210002	10 x 500	no	144	5.0	5.0	3.3	7.2	0.099
В	210022	500 individually wrapped peel-pack	OXIDE	144	5.0	5.0	3.3	1.65	0.03
С	210004	8 x 250	no	153	3.0	5.5	3.0	3.30	0.042
С	210024	500 individually wrapped peel-pack	OXIDE	153	3.0	5.5	3.0	1.62	0.03



Extra long pipettes

Made in transparent polyethylene. Model 200050 is graduated each ml up to 5 ml. Volumetrical tolerance: ±10 %.

mod.	code	presentation	lenght mm	Ø stem mm	total capacity ml	bulb draw ml	case weight	case volume
Α	200050	10 x 250	216	-	5.0	3.5	8.90	0.105
В	200060	10 x 100	300	9.0	23.0	7.1	7.3	0.105
С	210006	10 x 400	225	5.0	6.0	2.3	9.04	0.098

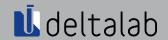



Sterile serological pipettes

Single use, made of crystal polystyrene. Accuracy > 98 %.

Sterilised by radiation. Made in one piece for an optimum precision, except code 90043 (50 ml).

No liquid retention: as tip is not sealed, liquid is not held at the joint. Total dispensing guaranteed.

Non pyrogenic (ISO 10993-11), non cytotoxic (ISO 10993-5) and non haemolytic (ISO 10993-4). **DNase and RNase free.** Volumes range from 1 to 50 ml. Pipettes include a coloured cotton plug for better identification. The peel-pack have a **fiber-free** packaging including a strip for easy opening. Black graduations are bright, black and **permanent**. 5, 10 and 25 ml models have negative graduations for extra capacity, and a second reverse scale (ascending and descending).

	code	capacity ml	presentation	cotton colour	case quantity	tip	graduation ml	negative graduation ml	total capacity ml	case weight	case volume	cases per pallet
	900030	1	1 peel pack	•	4 x 500	Α	0/0.9	until –0.3	1.3	10.00	0.068	24
	900031	1	bag of 25	•	4 x 1,000	Α	0/0.9	until –0.3	1.3	14.52	0.071	30
	900130	1	1 peel pack	•	4 x 500	С	0/0.9	until –0.3	1.3	10.10	0.071	30
	900131	1	bag of 25	•	4 x 1,000	С	0/0.9	until –0.3	1.3	14.00	0.071	25
	900032	2	1 peel pack	•	4 x 500	А	0/1.8	until -0.4	2.4	11.60	0.068	30
	900033	2	bag of 25	•	4 x 1,000	Α	0/1.8	until –0.4	2.4	17.20	0.071	30
	900034	5	1 peel pack	•	4 x 250	А	0/4	until –3	8	12.20	0.071	24
	900038	5	bag of 25	•	4 x 375	Α	0/4	until –3	8	15.00	0.071	30
NEW	900144N	5	1 peel pack	•	4 x 125	В	0/4	until –3	8	7.00	0.066	30
NEW	900134N	5	1 peel pack	•	4 x 250	С	0/4	until –3	8	7.20	0.071	30
	900036	10	1 peel pack	•	4 x 200	Α	0/9	until –3	13	11.07	0.071	24
	900037	10	bag of 25	•	4 x 250	Α	0/9	until –3	13	13.08	0.071	30
NEW	900146N	10	1 peel pack	•	4 x 200	В	0/9	until –3	13	7.86	0.067	30
NEW	900136N	10	1 peel pack	•	200	С	0/9	until –3	13	8.20	0.071	30
	900041	25	1 peel pack	•	4 x 100	А	0/23	until –10	35	11.20	0.071	24
	900043	50	1 peel pack	•	4 x 25	А	1/2	until –10	60	2.50	0.018	60







Pipetting bulb

Made of rubber. Used to avoid mouth pipetting and contamination risk. Can be opened, cleaned and **autoclaved**. Ideal for Wintrobe and Westergren pipettes.

code	description	bag quantity	bag weight	bag volume
19201	red pippeting bulb	1	0.03	0.001



Pipette pumps

Several models for various pipette volumes.

Designed for fast and efficient pipetting with simple, one handed operation. Pipettes fit smoothly into collar.

Rotate the knurled thumb wheel on the side for precision filling or dispensing, and press the fast release lever for quick emptying.

Easy to use and easily disassembled for cleaning.

Sizes are colour coded.

Pipette pumps resist acids and alkalies.

Attention: New models coming soon.

code	description	colour	bag quantity	bag weight	bag volume
W-100	up to 2 ml	blue	1	0.05	0.0001
W-110	up to 10 ml	green	1	0.05	0.0001
W-120	up to 25 ml	red	1	0.05	0.0001



See our Chapter 9. Sample storage



See our Chapter 11. Hygiene and Safety







Non graduated serological pipettes. Sterile

Single use, made in polystyrene, without graduation. In sterile packing by radiation.

DNAse and RNAse free, non cytotoxic, non haemolitic, apyrogen. Without cotton.

Standard tip, see picture (A) at the page 194.

Suitable for Integra VacuSafe Comfort and others.

Dimensions:

External tip diameter, between 5.45 - 5.7 mm, internal diameter tip end: between 1.8 - 2.2 mm.

code	capacity ml	presentacion	case quantity	case weight	case volume
900032.1	2	peel pack de 1	4 x 500	11.30	0.071



Sterile single use Straw-like pipettes

Single use straw-like non graduated pipettes, ideal for food industry use. In polyethylene.

Suitable for automatic pipettes. **Sterilized** by radiation.

Packaged in flow packs of 25 units.

Outer diameter: 4.35 mm. Inner diameter: 3.65 mm.

code	capacity	lenght	case	case	case
	ml	mm	quantity	weight	volume
900050	1.0	190	40 x 25	1.19	0.007

Pipette trays


Made of white PVC.

Resistant to temperatures from **-20 °C** to **80 °C**.


Model 1, code **19252**, is perfect to fit inside drawers. It includes four compartments that can hold up to thirty pipettes of 1, 2, 5, or 10 ml.

Model 2, code **19996**. Sixteen pipettes up to 10 mm diameter can be accommodated longitudinally; laterally, seven pipettes up to 20 mm diameter. Tray edges are ergonomically designed for a better handling. It is also useful with other instruments.

mod.	code	dimensions mm	case quantity	case weight	case volume
1	19252	426 x 300 x 30	1	0.53	0.006
2	19996	283 x 216 x 40	6	1.38	0.017

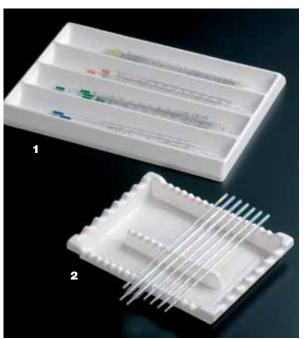



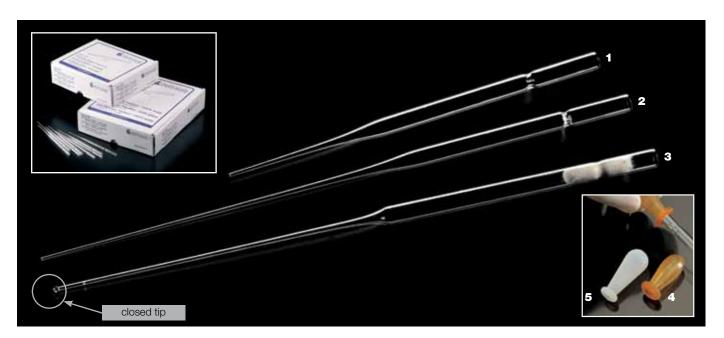

Glass Pasteur pipettes, open mouth

Soda glass. Dimensions: Ø mouth: 6.95 ± 0.15 mm; Ø tip: 1.2 ± 0.15 mm; 0.53 ± 0.03 thick.

mod	code	total lenght mm	stem lenght mm	case quantity	case weight	case volume
1	701	150	59	4 x 250	3.33	0.012
2	702	230	130	4 x 250	3.50	0.016

Glass Pasteur pipettes, closed tip


Made of neutral soda glass, with the tip closed and cotton in the mouth. Easy opening of the closed part by pushing it slightly. Dimensions: mouth: 7.00 mm; tip: 1.4 to 1.6 mm; thickness: 0.50 mm.

mod	code	total lenght mm	stem lenght mm	case quantity	case weight	case volume
3	712	230	130	4 x 250	4.00	0.015

Pipetting bulbs

Designed for pipetting with plastic and glass pipettes.

mod	code	material	total length mm	Ø max. mm	capacity pipetting ml	bag quantity	bag weight	bag volume
4	19194	latex	35	14	2	100	0.15	0.0007
5	19195	silicone	40	15	1,5 - 1,8	50	0.26	0.0014

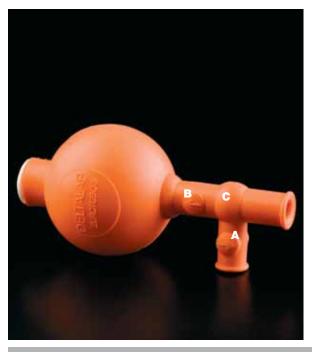

EUROTUBO® Pipetting bulb

This one-handed design is the simplest safety pipette filler to use available. Manufactured in natural orange rubber.

Approx. drawn capacity: 25 ml.

Single hand use, only two operating points.

Evacuate via the automatic valve.

Standard model, accommodates all pipettes. Ability to clean inside of bulb by removing patented valve and rinsing out.

code	description	bag quantity	bag weight	bag volume
19200	pippeting bulb	1	0.10	0.002



1. Evacuate the air by 2. Intake by pressing pressing the bulb, as the drawing indicates.



on point B




3. Drain the liquid by 4. pressing on point . A 🕠




Blow out the pipette by pressing as the drawing indicates (point C).





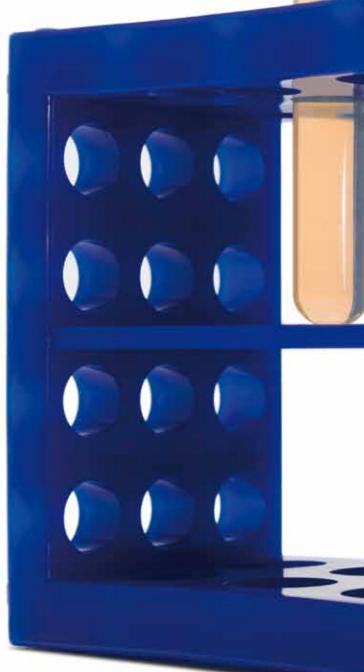



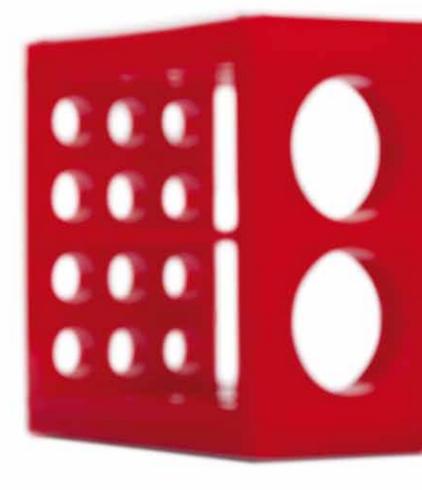



Polycarbonate racks for cryovials

These racks provide ultra low temperature storage (from -196 °C to 121 °C), holding 2 or 5 ml cryovals (depending on the model). Compatible with other storage systems, they can be sterilised using **autoclave** protocols.

They can also be immersed in liquid nitrogen, where its design means minimal issues with excess liquid on their removal.

code	capacity*	colour		dimensions mm	max. tubes capacity	max. tubes height	case quantity	case weight	case volume
M-061NA	25	orange	D	75 x 75 x 52	2 ml	48	5	0.33	0.003
M-061AM	- 25	yellow	Е	75 X 75 X 52	2 1111		5	0.33	0.003
M-062B	81	blue		132 x 132 x 53	2 ml	48	5	1.02	0.006
M-062R	01	red	С	132 X 132 X 53		40	5	1.02	0.006
M-063B	81	blue	В	132 x 132 x 95	5 ml	93	5	1.36	0.010
M-063R	01	red		132 X 132 X 93	51111	93	5	1.36	0.010
M-064V	100	green	А	132 x 132 x 53	2 ml	48	5	0.90	0.006
M-064AM	100	yellow		102 X 102 X 00	2 mi	40	5	0.90	0.006

^{*} The capacity of codes M-064V and M-064AM is halved when used with capped cryovials with external threads.

Polystyrene racks for cryovials

Made of polystyrene, those boxes are designed to hold 1 and 2 ml cryovials at temperatures from -80 °C to 80 °C. Stackable.

Boxes with white base and numbered inner rack. Boxes are stackable and have air vents.

Covers are keyed to the base to prevent misalignment.


Each box is supplied with an innovate vial picker for easy vial removal.

Dimensions: 133 x 133 x 52 mm.

Capacity: 100.

Ø wells: 12.8 mm.

code	colour	capacity	case quantity	case weight	case volume
M-036B	blue	81	4	0.62	0.0045
M-036P	pink	81	4	0.62	0.0045
M-036AM	yellow	81	4	0.62	0.0045

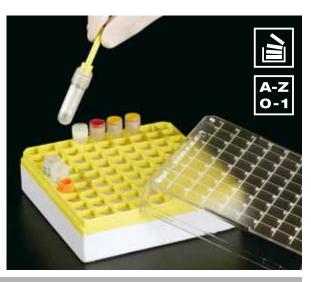







Polycarbonate racks for cryovials

Made of extra strong polycarbonate, these boxes are designed to be used at temperatures between **-190 °C** and **121 °C** and are **autoclavables**. Racks can hold either 81 or 100 tubes from 1 ml to 5 ml.

They can be immersed in liquid nitrogen. Transparent covers are printed with an alphanumerical grid.

Boxes made to accept 81 tubes have a white base and a coloured grid (in which each well is numbered) (A type).

Boxes made to accept 100 tubes* have a coloured base, with each place numbered in the lid, instead of a grid (B type).

Each box is supplied with an innovative vial picker for easy vial removal (see C).

Boxes are stackable and have air vents.

Covers are keyed to the base to prevent misalignment.

colour	case quantity	case weight	case volume				
For 1.2 and 2 ml cry	ovials. 81 holes. Dimensions	: 133 x 133 x 52 mm					
yellow	4	0,73	0,0043				
red	4	0,73	0,0043				
green	4	0,73	0,0043				
For 3.4 and 5 ml cryovials. 81 holes. Dimensions: 133 x 133 x 95 mm							
yellow	10	2,28	0,020				
red	10	2,28	0,020				
green	10	2,28	0,020				
For 1.2 and 2 ml cryo	ovials. 100 holes. Dimensions	s: 133 x 133 x 52 mm					
yellow	4	0,72	0,0043				
red	4	0,72	0,0043				
green	4	0,72	0,0043				
	For 1.2 and 2 ml cry yellow red green For 3.4 and 5 ml cry yellow red green For 1.2 and 2 ml cryo yellow red	For 1.2 and 2 ml cryovials. 81 holes. Dimensions yellow 4 red 4 green 4 For 3.4 and 5 ml cryovials. 81 holes. Dimensions yellow 10 red 10 green 10 For 1.2 and 2 ml cryovials. 100 holes. Dimensions yellow 4 red 4	For 1.2 and 2 ml cryovials. 81 holes. Dimensions: 133 x 133 x 52 mm yellow				

^{*} The capacity of codes M-035AM, M-035R and M-035V, is halved when used with capped cryovials with external threads.



Polycarbonate rack for cryovials

Made of polycarbonate, this box is designed to store up to 42 cryovials (10 ml, code 401410) at temperatures between **-190 °C** and **121 °C**.

Autoclavable and stackable.

Transparent cover is printed with the indication of the well number from 1 to 42. White base and blue grid, with numbered wells.

They can be immersed in liquid nitrogen.

Upon removal of the bath, excess liquid is discharged through vents incorporated.

Cover is keyed to the base to prevent misalignment.

Dimensions: 133 x 133 x 95 mm

code	description	case quantity	case weight	case volume
M-036A	PC rack for 10 ml cryovials	5	1.22	0.010



Rack 9x9 height 45

200401

Polypropylene rack. Dimensions with lid: 133 x 133 x 45 mm.

Available in fluorescent colours with transparent lid to help identification.

It resists -196 °C and is autoclavable up to 121 °C.

Sections will hold 81 microtubes (0.5, 1.5 or 2 ml) or cryotubes (1.2 ml) up to 12 mm diameter and 42 mm height.


The capacity is reduced when storing external screw cryotubes.

It embodies an alpha numeric grid printed on the base and over the upper rims for a clearer identification of samples.

Suitable with tubes:

200405 200407 200410 4092.2N 4092.5N

4092.6N	4092.6NS	409001			
code	e	colour	quantity	weight	volume
M-020AM		yellow	1	0.17	0.0011
M-020	DN	natural	1	0.17	0.0011
M-020R		red	1	0.17	0.0011
M-020V		green	1	0.17	0.0011

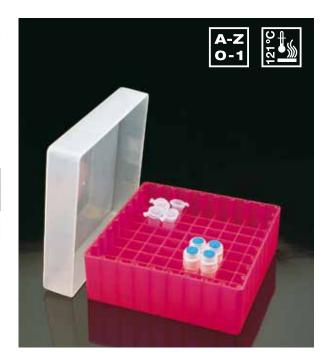

Rack 9x9 height 50

Polypropylene rack.

200400

Dimensions with lid: 133 x 133 x 53 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196° C and is autoclavable up to 121 °C.

Sections will hold 81 microtubes (0.5, 1.5 or 2 ml) or cryotubes (1.2 or 2 ml) up to 12.5 mm diameter and 50 mm height.

The capacity is reduced when storing external screw cryotubes.

It embodies an alpha numeric grid printed on the base and over the upper rims for a clearer identification of samples.

Suitable with tubes:							
200400	200400 4092.2N 4092.5N 4092.6N 4092.6NS 409001 409002						
409002.1 409110.1 409111/2							

code	colour	quantity	weight	volume
U-9950CAM	yellow	1	0.19	0.0011
U-9950CN	natural	1	0.11	0.0011
U-9950CNA	orange	1	0.11	0.0011
U-9950CR	red	1	0.11	0.0011
U-9950CV	areen	1	0.11	0.0011



Rack 9x9 height 95

Polypropylene rack. It can be freezed down to -196 °C, and autoclaved up to 121 °C.

Dimensions with lid: 133 x 133 x 95 mm

Available in fluorescent colors with translucent-transparent lid.

natural


Ideal for 4 or 5 ml cryotubes.

81 sections will hold microcontainers up to 12 mm diameter and 93 mm high: screw cap microtubes, cryovials (even those of 4 and 5 ml), snap cap microtubes, and short-medium sized tubes.


Suitable with tubes:								
300300 300500 300700 300800 300800.2 400500 4007								
400800								
code	•	colour		bag quantity	bag weight	bag volume		

0.0026

0.23

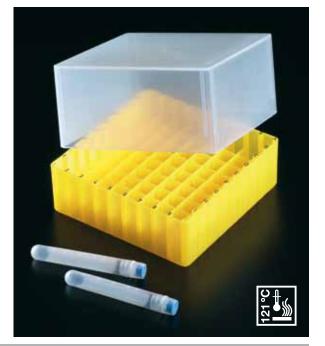



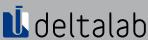



Rack 8x8 round hole height 45

Polypropylene rack, **autoclavable** up to **121 °C**.

Dimensions with lid: 133 x 133 x 53 mm.


Available in fluorescent colours with transparent lid to help identification. **Resists -196° C**.

Sections will hold 64 microtubes (0.5, 1.5 or 2 ml) or cryotubes (1.2 ml) up to 11 mm diameter and 42 mm height.

Alphanumeric identification.

Suitable with tubes:							
4092.1N 4092.1NS 4092.2N 4092.5N 4092.6N 4092.6NS 200400							
409111/1							

code	colour	quantity	weight	volume
U-8845A	blue	1	0.15	0.0009
U-8845N	natural	1	0.15	0.0009

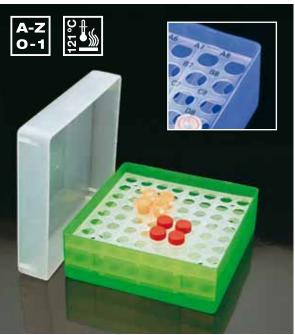

Rack 8x8 round hole height 50

Polypropylene rack, autoclavable up to 121 °C.

Dimensions with lid: 133 x 133 x 53 mm.

Available in fluorescent colours with transparent lid to help identification. $\textbf{Resists-196}^{\circ}$ C.

Sections will hold 64 microtubes (0.5, 1.5 or 2 ml) or cryotubes (1.2 or 2 ml) up to 11 mm diameter and 50 mm height.

Alphanumeric identification.

Suitable with tubes:							
200400 408003 409110.1 409110.2 409110.3 409110.4 4091							
409111/1	409111/2	409111/3	409111/4	409111/5	409111/6	4092.2N	
4092.5N 4092.6N 4092.6NS							

code	colour	quantity	weight	volume
U-8850A	blue	1	0.15	0.0011
U-8850N	natural	1	0.15	0.0011
U-8850NA	orange	1	0.15	0.0011






Polypropylene rack, ${\bf autoclavable}$ up to ${\bf 121~^{\circ}C}.$

Dimensions with lid: 133 x 133 x 53 mm.

Available in fluorescent colours with transparent lid to help identification. Resists -196° C.

Sections will hold 49 microtubes (0.5, 1.5 or 2 ml), cryotubes (1.2 or 2 ml), o low volume tubes, up to 16 mm diameter and 50 mm height.

	Suitable with tubes:					
409002 409002.1 409104.1 409106.1 409106.2						
cod	е	colour		quantity	weight	volume
U-775	60A	blue		1	0.15	0.0011
U-7750AM U-7750N		yellow natural		1	0.15	0.0011
				1	0.15	0.0011

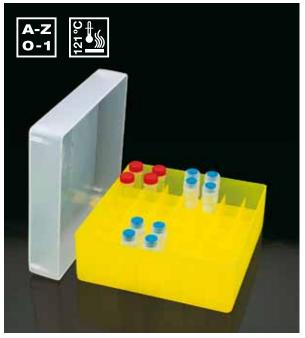



Freezing cardboard storage box, pre-assembled compartments

Designed for storage and preservation of frozen samples.

Made of white treated cardboard and covered by special white paper handstuck and hand-assembled for a better quality.

Pre-assembled compartments.

Resist temperatures up to -100\,^{\circ}\text{C}. Ideal for microtubes and cryovials.

There are alphanumeric identifications printed on each box.

Suitable for the racks of freezer chests.

The models **M-600** and **M-610** are silk screen printed on the base and lid. The other models have no printings.

These box can be custom-designed made, ask to our commercial department.

code	external dimensions mm	grid num.	Ø holes mm	case quantity	case weight	case volume
M-600	150 x 150 x 55	10 x 10	12	10	1.91	0.021
M-610	150 x 150 x 100	10 x 10	12	10	2.66	0.043
M-607	134 x 134 x 47	10 x 10	10	10	1.50	0.013
M-606	134 x 134 x 57	10 x 10	10	10	1.50	0.013
M-603	134 x 134 x 75	9 x 9	12	10	2.30	0.020
M-604	134 x 134 x 79	9 x 9	12	10	2.35	0.020



Rack 10x10 height 37

Polypropylene rack.

Dimensions with lid: 133 x 133 x 37 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196 °C and **autoclavable** up to **121 °C**. Sections will hold 100 microtubes up to 10 mm diameter and 34 mm height.

	Suitable with tubes:							
4092.1N	4092.1N 4092.1NS 900117							
code	е	colour		quantity		weight	volume	
U-053	7A	blue		1		0.12	0.00062	
U-0537N		natural		1		0.12	0.00062	
U-0537NA		orange		1		0.12	0.00062	





Rack 4x4 height 45

Polypropylene rack.

Dimensions with lid: 67 x 67 x 45 mm.

Available in fluorescent colors with transparent lid. It can be freezed up to -196 °C, and autoclaved up to 121 °C.

16 sections (15x15 mm) will hold screw cap microtubes (range **409110.1** for example), cryovials (as code **409001**), snap cap microtubes (see codes **200400** or **4092.2N**), and short tubes (codes **300500**, **400500**, **401100**).

Black color rack includes a black lid and is perfect to store light sensitive samples; it can be combined with our opaque microtubes (range **409113.1**) It may be used with taller tubes as the lid lies on their caps (up to 60 mm high).

code	colour	bag quantity	bag weight	bag volume
U-4445N	natural	5	1.17	0.0012
U-4445A	blue	5	1.17	0.0012
U-4445AM	yellow	5	1.17	0.0012
U-4445NE	black	5	1.17	0.0012







Racks with lid

Polypropylene rack, **autoclavable** and designed for storage of samples down to **-90 °C**. The lid is jointed to the rack by hinges. Available in fluorescent colours.

The racks are alphanumerically matrixed and are designed for 1.5 and 2 ml microtubes or 2.0 ml cryotubes; and we can supply two different models which can accommodate 50 and 100 tubes respectively.

Dimensions: 71 x 153 x 53 mm (50 wells).

140 x 153 x 53 mm (100 wells).

code	colour	holes	case quantity	case weight	case volume
M-500	natural	50	20	2.48	0.020
M-501	neon green	50	20	2.48	0.020
M-502	neon pink	50	20	2.48	0.020
M-503	neon yellow	50	20	2.48	0.020
M-510	natural	100	20	4.14	0.028
M-511	neon green	100	20	4.14	0.028
M-512	neon pink	100	20	4.14	0.028
M-513	neon yellow	100	20	4.14	0.028

Minimum order quantity: 5 and multiples of 5 up to 20. From 20, multiples of 20.



Microtube racks with lift-off telescopic lid

Made of **autoclavable** polypropylene, those racks are designed for long-term storage of samples at temperatures down to -90 °C.

Model **M-450** can accommodate 50 microtubes of 0.5 ml (hole diameter 10 mm).

Model **M-410** can accommodate 50 microtubes of 1.5/2.0 ml (hole diameter 12 mm).

Lids have a moulded-in alphanumeric location identification system with a clear centre panel for easy tube viewing.

Four frosted side panels and a frosted rim on the edges of the lid provide additional labeling options.

code	dimensions mm	case quantity	case weight	case volume
M-450	120 x 60, height 37,5	5	0.47	0.0021
M-410	140 x 73, height 37,5	5	0.56	0.0035





Made of white polypropylene. Designed to hold 24 or 48 microtubes or cryovials of $1.5/2.0 \, \text{ml}$ sizes.

Tubes are presented on four levels for better sample viewing.

Each well is numbered.

Model M-052 consists of one 24 hole-rack supplied in a white hinged box with transparent cover; model **M-053** consists of two 24-hole racks, also housed in a white hinged box with transparent cover.

Racks are removable from boxes and are stackable.

They float in water bath and can incorporate dry ice for use with frozen samples on the bench.

Fully autoclavable (121 °C) and freezable down to -90 °C.

Dimensions: **M-052**: 120 x 85 x 60 mm.

M-053: 235 x 85 x 60 mm.

code	description	case quantity	case weight	case volume
M-052	rack for 24 microtubes	10	1.34	0.0089
M-053	rack for 48 microtubes	5	0.74	0.0087






Reversible racks with removable handles


Made of polypropylene, autoclavable up to 121 °C.

The racks have two sides for different tubes: one side holds up to 80 tubes of 12 mm diameter (for example 10×75 mm or 12×75 mm tubes, microtubes, etc.) and the other side holds up to 60 microtubes or PCR (8 mm \emptyset , 0.2 or 0.5 ml).

By removing the handles, units can be firmly anchored laterally to one another and more units can be added on top of each other, thanks to small pins supplied with the racks.

code	colour	case quantity	case weight	case volume
M-563A	blue	10	1.60	0.008
M-563NA	orange	10	1.60	0.008
M-563P	pink	10	1.60	0.008

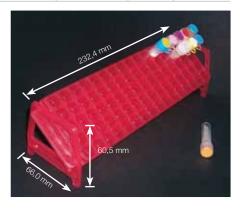

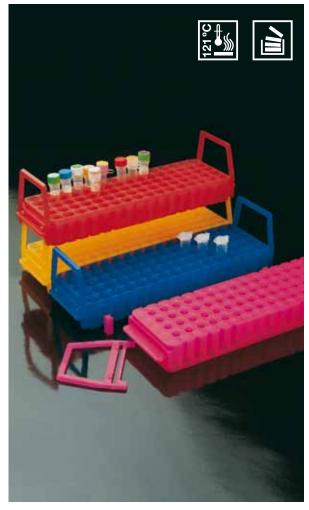

80 well racks

Made of autoclavable polypropylene, those racks are designed for freezer storage and are rated to $\ -80\ ^{\circ}C.$

Racks are designed to accommodate microtubes of diameter up to 11 mm. They are space-saving, as they feature 80 wells in a 5×16 configuration. An alpha-numeric matrix allows accurate sample identification.


Lateral handles allow easy handling of the racks. Racks can be used alone or placed into the box **M-011.1** (made of translucent plastic).

Racks are compatible with the Gilson FC80 dispenser, for example.


code	description	dimensions mm	case quantity	case weight	case volume
M-010.1	yellow rack	225 x 65 x 27	5	0.85	0.0022
M-010.2	natural rack	225 x 65 x 27	5	0.85	0.0022
M-010.5	blue rack	225 x 65 x 27	5	0.85	0.0022
M-010.6	red rack	225 x 65 x 27	5	0.85	0.0022
M-010.7	green rack	225 x 65 x 27	5	0.85	0.0022
M-011.1	box with lid	230 x 85 x 60	5	2.00	0.029





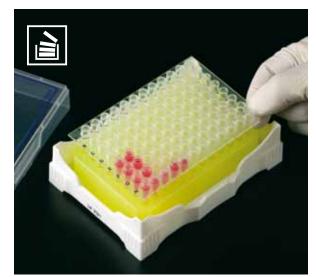

Rack at room temperature

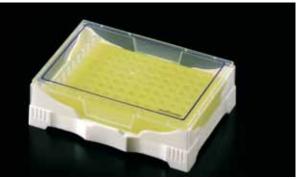

Rack at -20 °C





Isofreeze Rack with temperature indicator

Manufactured from polypropylene, this rack is filled with a non-toxic gel that maintains **4 °C** temperature for almost three hours and a half when the rack is removed from the freezer (or in case of freezer failure), protecting PCR reagents, cells, enzymes, etc, from fast temperature changes.


Moreover, the plastic changes its colour from green to yellow as long as the temperature of the rack changes (see figure). It can hold up to 96 0,2 ml single microtubes, strips, and plates.

Dimensions with lid: 141 x 99 x 44 mm.

Dimensions without lid: 141 x 99 x 38 mm.

ATTENTION: Freeze for a minimum of eight hours at **-20 °C** before use. Keep it covered while being used.

Store at the freezer while not being used

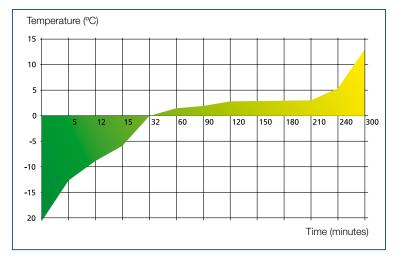

code	description	case quantity	case weight	case volume
M-572	yellow-green isofreeze rack	2	0.64	0.0005



See our PCR tubes and plates in Chapter 2. **Molecular biology**

96 well PCR rack

Base and cover made of autoclavable polypropylene.

Designed for storing and working 0,2 ml PCR tubes, separate or in strips. 96-well plates can also be accommodated.

Both base and cover have an alphanumeric identification.

Once covered they are stackable.

Racks can be horizontally attached to each other in order to build-up any configuration.

Dimensions with cover: $132 \times 91 \times 32$ mm It withstands from **-90 °C** to **121 °C**.

code	description	case quantity	case weight	case volume
M-585	blue PCR rack	20	3.0	0.02
M-581	yellow PCR rack	20	3.0	0.02

Minimum order quantity: 5 and multiples of 5 up to 20.

From 20, multiples of 20



24 well incubation rack

Made of autoclavable polypropylene.

This floating rack will incubate up to 24 microtubes in a water bath and can accommodate either 1.5 ml or 0.5 ml tubes.

Tubes are held snugly and then released when the rack is pressed down on a benchtop.

The rack is alphanumerically labelled for better sample identification. Dimensions: $119 \times 78 \text{ mm}$.

code	description	case quantity	case weight	case volume
M-022	incubation rack	4	0.29	0.0024



Round incubation racks

Made of **autoclavable** polypropylene, those racks are designed to hold either 8 or 20 microtubes (1.5 ml) in 10.8 mm holes.

Even when the rack is fully loaded, tubes are completely immersed while tops rest above the rack, floating in a beaker of water.

Perfect for use in benchtop incubation.

The rack is easily unloaded by lightly pressing it down against the bench until the legs touch, releasing the tubes.

Racks are 6.4 mm thick with 19 mm long legs and a central handle.

code	capacity	Ø rack	case quantity	case weight	case volume
M-024	8 tubes	66	4	0.14	0.0015
M-026	20 tubes	105	4	0.20	0.0075





IsoFreeze reversible racks with lid

Manufactured from virtually unbreakable polycarbonate and filled with a non-toxic gel, those Isofreeze temperature maintenance racks provide consistent thermal protection for PCR reagents, cells, enzymes or other temperature sensitive agents when removed from the freezer or refrigerator for extended periods of time, and provide safe storage on the event of freezer failure. Each rack can accommodate up to 20 1.5/2.0 ml microtubes on one side, or 20 0.5 ml microtubes on the reverse side.

Racks feature a transparent lid for better sample viewing. Each well in the rack is numbered.

Available in two temperature ranges.

Model M-570 (white) will maintain a temperature of $0~^{\circ}C$ on the bench top for 5 hours.


Model **M-571** (blue) will maintain a temperature between **-10 °C** and **-20 °C** on the bench top up to 3 hours.

Dimensions with lid: 108 x 172 x 104.

code	description	case quantity	case weight	case volume
M-570	white isofreeze rack	1	1.30	0.005
M-571	blue isofreeze rack	1	1.30	0.005





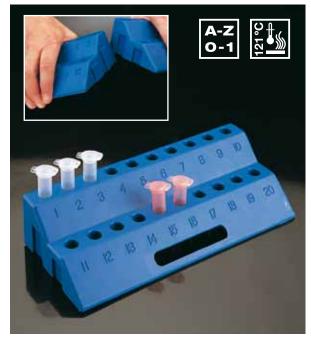

20 well interlocking tube racks


Made of blue autoclavable polypropylene.

These racks slide together on both sides to allow the user to attach several for larger experiments.

They are designed to accommodate 1.5/2.0 ml Eppendorf type tubes on two levels of 10 wells each, with easy-to-read numbers for better tube viewing and identification.

The innovative front finger slot makes them easy to handle when working with gloves. Dimensions: $210 \times 90 \times 48$ mm.

code	description	case quantity	case weight	case volume
200420	blue 20 hole rack	6	0.80	0.0040

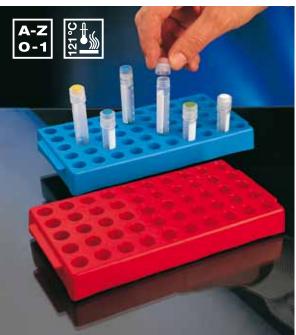


50 well microtube racks

Made of **autoclavable** polypropylene, these racks are designed for use with our tubes presented on page 144 and successive. They can hold 50 tubes in a 10 x 5 configuration. These racks allow easy one handed opening and closing of tubes thanks to an innovative locking system which ensures that tubes will securely lock in each well and will not turn. Each place is identified with an alphanumeric index for better sample location. Racks have strong handles for better handling and are stackable.

Dimensions: 100 x 250 x 25 mm.

code	description	case quantity	case weight	case volume
M-028.5	blue rack	4	1.0	0.004
M-028.6	red rack	4	1.0	0.004

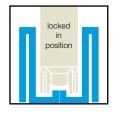



Styrofoam racks with lid

Made of white heavy Styrofoam, these racks are designed for long-term storage of microtubes and cryovials in the freezer.

Both racks accommodate tubes maximum 44 mm high.


Diameter of the well: 13mm.

Model **G-050** (50 wells) is suitable for 1.5/2.0 ml microtubes and/or 1.2 ml cryovials. Dimensions: $21 \times 10.7 \times 6.5$ cm.

Model **G-100** (100 wells) is suitable for 1.5/2.0 ml microtubes. Dimensions: $33.5 \times 9.5 \times 6.8$ cm.

Racks are alpha-numerically matrixed for better sample location. Stackable, light-weighted and economical.

code	description	case quantity	case weight	case volume
G-050	50 positions	10	0.50	0.028
G-100	100 positions	10	0.57	0.029

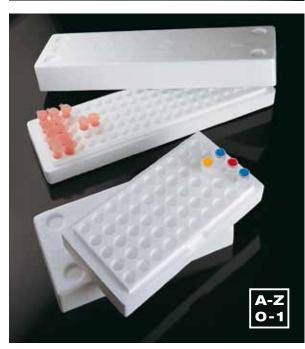

Reversible 32+8 microtube racks with lid

This rack made of **autoclavable** at **121 °C** polypropylene and resistant to **-90 °C.** Is supplied with a clear polypropylene lid for better sample viewing. Side one has 32 wells for 0.5 ml tubes plus an additional 8 wells across for 1.5/2.0 ml tubes. Side two has the opposite configuration of 32 wells for 1.5/2.0 tubes and 8 wells for 0.5 ml tubes. Both sides are alphanumerically labeled for better sample location.

Racks have frosted writing patches on all four sides for smear resistant labeling. Small side handles.

Dimensions: 170 x 80 x 27 mm.

code	description	case quantity	case weight	case volume
M-019AM	yellow rack	5	0.80	0.007
M-019R	pink rack	5	0.80	0.007
M-019V	green rack	5	0.80	0.007



PCR workstations

Autoclavable polypropylene. This storage/working rack, with lid, is designed to work with different PCR volumes at the same time. It can accommodate up to 32 individual 0.2 ml PCR tubes, 24 1.5/2.0 ml and 16 0.5 ml microtubes. Tube wells are alpha-numerically matrixed with raised bold type for accurate sample identification. The rack also has a transparent lid for easy viewing of samples. Frosted writing panels on four sides allow labelling or writing. The rack is supplied with a removable rack accommodating 0.2 ml PCR tubes. The removable rack can be purchased separately: see product code below. Dimensions: main rack: 215 x 118 x 50 mm.

Removable rack: 125 x 85 x 15 mm.

Both racks withstands from -90 °C to 121 °C.

code	description	case quantity	case weight	case volume
M-540	green rack	5	1.50	0.007
M-541	blue rack	5	1.50	0.007
M-542	yellow rack	5	1.50	0.007

Minimum unit sale: 5



96 well PCR racks

Autoclavable polypropylene.

Our 96 well plate, individual tubes or strips of 8 or 12 tubes can be used with this rack. An alphanumeric grid helps sample identification.

The included lid is transparent and has a non-slip surface to allow safe storage of multiple racks. This rack can be used alone or inserted in the above workstation, codes M-540 to M-543.

Dimensions: 125 x 85 x 15 mm. It withstands from -90 °C to 121 °C.

code	description	case quantity	case weight	case volume
M-550	blue rack	5	0.4	0.002
M-551	yellow rack	5	0.4	0.002
M-552	pink rack	5	0.4	0.002
M-553	green rack	5	0.4	0.002







Stacking PCR work up rack

Autoclavable polypropylene. The stacking PCR work-up rack is a robust and convenient workstation or storage rack for plates, strips or tubes.

It withstands from -90 °C to 121 °C. The conveniently small footprint ensures it is an equally useful tool whether it is used on the bench, in the freezer or some other laboratory situation.

The positive click ensures the lid will not become easily dislodged. Yet only a gentle squeeze on each side will allow you to lift off the clear lid. It accommodates every different type of 0,2 ml PCR plate; including fully skirted types, strips or tubes.

The imaginative shaping of the rack's feet make this system easy to use, providing a stable platform as well as enabling stacking with or without the **lid** and whether the vessels are open or sealed.

The racks are manufactured in distinctive colours and are packaged as an assortment.

code description		case	case	case
		quantity	weight	volume
M-554	stacking PCR rack	5	0.45	0.0035



COMBI Rack for PCR

The box is made of PVC and the rack is made of **autoclavable** polypropylene. Rack and box can be used both together and separately.

The box: can be used not only as a storage rack but also as a work station. The white base accepts both 96 and 384 well plates and a transparent cover. It is stackable, and its low height helps to save space. Lid dimensions: 108.0 x 148.6 x 35.6 mm. Box dimensions: 108.0 x 148.6 x 24.9 mm.

The rack: can hold up to 96 PCR tubes or plates or 12 strips of 8 tubes with caps.

Each hole is identified with an alpha numeric numbering system.

It features a writing area.

The grid stands on 4 legs and can be placed on a lab counter or in a refrigerator. It is stackable.

Can withstand temperatures from **-80 °C** to **121 °C**.

Dimensions: 84.6 x 125.5 x 28.9 mm

code	description	case quantity	case weight	case volume
M-640 white box with lid		5	0.20	0.0004
M-641	96 hole blue rack	5	0.20	0.0004






Made of durable autoclavable polypropylene. This reversible rack with a clear polypropylene lid is designed with wide spacing to store up to 96 microtubes on either side. One side accommodates 1.5/2.0 ml microtubes; the reverse side accommodates 0.5 ml microtubes. The rack is alphanumerically matrixed for better sample location and has frosted writing patches on all four sides for smear resistant labelling. Operating temperature range: -90 °C to 121 °C. With ergonomical lateral handles.

Dimensions: 213x114x26 mm.

code	description	case quantity	case weight	case volume
M-021AM	yellow rack	5	2.05	0.008
M-021R	pink rack	5	2.05	0.008
M-021V	green rack	5	2.05	0.008

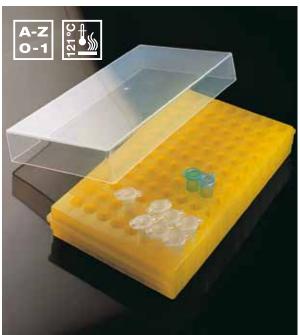




Microtube racks

Those easy to assemble racks are made of white **autoclavable** polypropylene. These racks have open sides that make them adaptable for use in a tray of crushed ice. Designed for 1.5 ml 42 (6 x 7) microtubes. Racks are supplied dismantled.

Dimensions: 128 x 105 x 43 mm.

code	description	case quantity	case weight	case volume
M-110	microtube racks	8	0.97	0.006

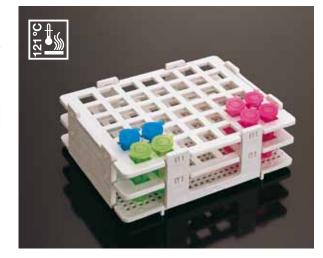


42 well cuvette rack

Made of white **autoclavable** polypropylene, this easy to assemble rack provides gentle no-scratch support for fragile spectrophotometer cells. Forty-two 10 mm light path cells fit in this 6 x 7 row layout (See page 81). Also suitable for tubes.

Racks are supplied dismantled.

Dimensions: 128 x 105 x 43 mm.

code description		case	case	case
		quantity	weight	volume
M-100	well cuvette rack	8	1.04	0.01

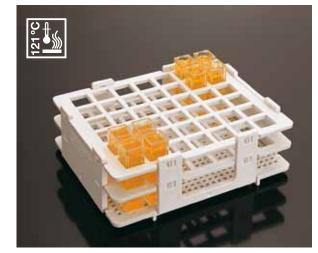

Test tube baskets

For rinsing, transport and storage.

Made of **autoclavable** polypropylene, these baskets with sliding lids are ideal for rinsing or carrying glassware.

code	description	case quantity	case weight	case volume
M-040	basket 127 x 102 x 102 mm	6	0.76	0.018
M-050	-050 basket 152 x 152 x 152 mm		1.39	0.003
M-060	basket 229 x 229 x 229 mm	6	3.12	0.080



Pipette rack

Made of **autoclavable** polypropylene.

For orderly, space-saving pipette storage. This rack holds 50 pipettes in a vertical 5 x 10 arrangement. This system ensures that pipettes will not contact table surface.

The end plates have handles for easy carrying.

Hole diameter: 16 mm.

Dimensions: 213 x 114 x 222 mm.

code	description	case quantity	case weight	case volume
19973	pipette rack	1	0,54	0,006







Universal racks

Made of durable polypropylene, this versatile rack is specially **designed to attach to another by pressing the two sides together**. Two racks or more can then be attached to present tubes in varied configurations.

The rack can accomodate the following:

- 4 x 50 ml conical tubes
- 12 x 15 ml conical tubes
- 32 x Ø12 mm tubes or 1.5/2.0 ml microtubes
- 32 x 0.5 ml microtubes

Alpha-numeric matrix on the four sides for better sample location. Can be autoclaved and can be placed in freezers and water baths.

Dimensions: 172 x 50 x 90 mm.

code	colour	case quantity	case weight	case volume
M-520	natural	5	0.81	0.005
M-521	fluorescent blue	5	0.81	0.005
M-522	fluorescent pink	5	0.81	0.005
M-523	fluorescent yellow	5	0.81	0.005



Universal racks

Made of **autoclavable** polypropylene, it withstands from **-90 °C** to **121 °C**. This tube rack system is designed to support a variety of tube sizes and styles:

- Ø 30 mm for 8 x 50 ml tubes

- Ø 20 mm for 20 tubes
- Ø 16/17 mm for 24 x 15 ml tubes
- Ø 12/13 mm for 40 tubes

Dimensions: 198 x 118 x 96 mm.

code	colour	case quantity	case weight	case volume
M-530	fluorescent blue	5	2.10	0.014
M-531	M-531 fluorescent green		2.10	0.014
M-532	-532 fluorescent pink		2.10	0.014
M-533	fluorescent yellow	5	2.10	0.014





Made of **autoclavable** polypropylene.

It with stands from $-90~^{\circ}\text{C}$ to $121~^{\circ}\text{C}.$

This design consists of a unit of four independent cubes packaged in a clear support tray.

Each rack may be independently rotated to allow the simultaneous storage and use of multiple sizes of tubes.

Each side of the racks is alpha-numeric matrixed for better sample location. The following tubes can be accommodated in each rack:

- 2 x 50 ml conical tubes
- 6 x 15 ml conical tubes
- 9 x 12 mm diameter tubes or 1.5/2.0 ml microtubes
- 12 x 0.5 ml microtubes

Dimensions:


Tray (handles included): 242 x 120 x 40 mm.

Rack: 48 x 92 x 92 mm.

code	description		case weight	case volume
330	combi rack	1	0.60	0.0030





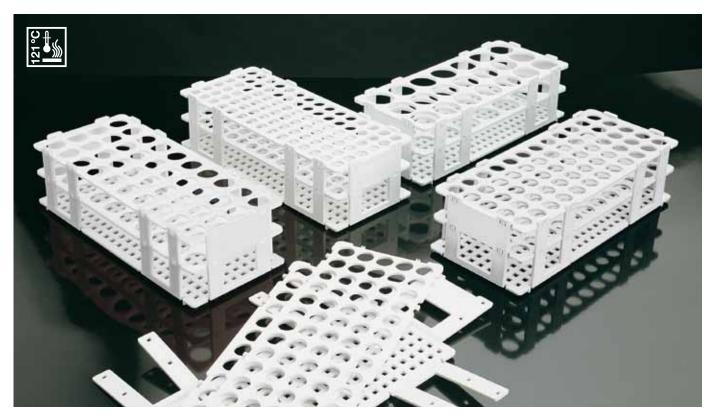

Universal test tube racks

Round hole test tube racks made of white **autoclavable** polypropylene.

A side panel provides a labelling option. Very easy to assemble.

Racks are supplied dismantled.

Dimensions: 240 x 103 x 65 mm.

code	max. diameter tubes	tube capacity	bag quantity	bag weight	bag volume
W-014	13	6 x 15 (90)	1	0.19	0.012
W-015	16	5 x 12 (60)	1	0.19	0.012
W-016	20	4 x 10 (40)	1	0.19	0.012
W-017	25	3 x 8 (24)	1	0.018	0.010
W-018	30	3 x 7 (21)	1	0.016	0.008



Square hole test tube racks with grip


These high quality polypropylene racks are **autoclavable**, stackable when unloaded and can be dismantled for space saving.

A unique grip feature ensures that tubes are firmly held in racks. An alphanumeric grid reference ensures an easy sample location. Dimensions: $250 \times 105 \times 68$ mm.

Very easy to assemble.

Racks are supplied dismantled.

code	max. diameter tube	tube capacity	bag quantity	bag weight	bag volume
W-011	13	6 x 15 (90)	1	0.17	0.0012
W-012	16	5 x 12 (60)	1	0.17	0.0012
W-013	20	4 x 10 (40)	1	0.09	0.0009

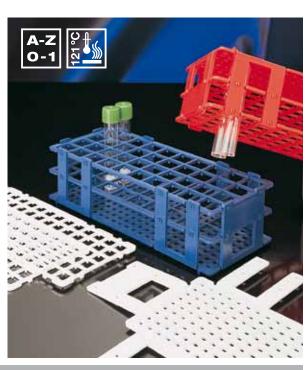








Universal test tube racks

Made of white **autoclavable** polypropylene. Specially designed to ensure that all test tubes will be incubated at the same temperature. An alpha-numeric grid reference on top tier of the rack allows an easy sample location. Racks are easily assembled, and are stackable when unloaded. Each rack is supplied dismantled. Easy to assemble.

Also available in colour: Blue: add .04 to the code.

Yellow: add .06 to the code. Red: add .10 to the code.

code	places	Ø test tube mm	dimensions mm	bag quantity	bag weight	bag volume
19563	100 (4 x 25)	micro	262 x 108 x 45	5	0.80	0.007
19564	90 (6 x 15)	13	246.5 x 104 x 60	5	0.80	0.007
19565	60 (5 x 12)	16	246.5 x 104 x 70	5	0.90	0.006
19566	40 (4 x 10)	20	246.5 x 104 x 70	5	0.90	0.006
19567	40 (4 x 10)	25	296.5 x 124 x 85	5	1.00	0.007
19568	24 (3 x 8)	30	300 x 111 x 83	5	1.00	0.007

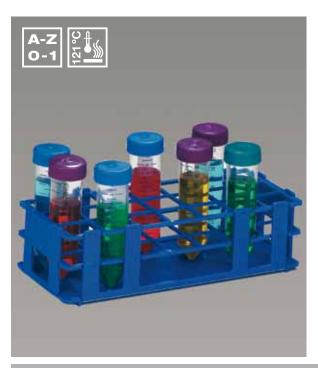

Assembling rack

This rack is moulded of blue **autoclavable (121°C)** polypropylene.

Suitable for daily usage and storage.

It holds up to 21 tubes of 50 ml (30 mm diameter) each, and features end plates for labelling, as well as numbered and lettered rows for tube identification. Interlocking feet allows stacking of racks for storage.

This rack is supplied flat, individually bagged.

Fast and easy assembling (instructions included).

code	dimensions mm	rows	case quantity	case weight	case volume
19570	246 x 104 x 64	3 x 7	8	1.50	0.020

Acetal polymer racks

Made of acetal polymer. Autoclavable up 121 $^{\rm o}{\rm C}$ for 20 minutes. Stackable.

Compact and lightweight.

They remains stable if submerged.

Color blue except model M-562L, which is lilac.

With handles on each side for better handling.

External dimensions for all codes: 293 mm x 115 mm x 65 mm H.

Recommended for tubes: 12×75 , 15×95 , 16×100 , centrifugation tubes, urine collection tubes, etc.

code	for tubes Ø	capacity	case quantity	case weight	case volume
M-560A	up to 13 mm	(6 x 14) 84	10	3.27	0.032
M-561A	up to 16 mm	(5 x 12) 60	10	3.00	0.029



Single use racks

Low-price, flexible racks made of polystyrene.

Dimensions:

901407: 200 x 240 x 30 mm **901410**: 160 x 170 x 25 mm

code	Ø hole	n° holes	case quantity	case weight	case volume
901407	16 mm Ø	10 x 12	25	1.90	0.033
901410	13 mm Ø	10 x 10	25	0.66	0.012





Combi rack height 128

Polypropylene rack.

Dimensions with lid: 133 x 133 x 128 mm.

Available in fluorescent colours with transparent lid to help identification.

Resist -196 °C and autoclavable up to 121 °C.

Sections will hold 10 tubes up to 30 mm diameter, 2 tubes up to 16 mm Ø and 1 tube up to 11 mm Ø.

All of them up to 125 mm height.

	Suitable with tubes:									
10100T	300400	300705	300707	300900	300907	300908				
301200	301201	301202	301212	301213	301403	301700				
400400	400705	400900	401200	401201	401204	401403				
401700	409920	429900	429901	429910	429926	429927				
429930	429931	713100	716100	716125	720125	901075				
913100	916100									

code	colour	bag quantity	bag weight	bag volume
U-50128A	blue	1	0.20	0.0035
U-50128N	natural	1	0.20	0.0035
U-50128R	red	1	0.20	0.0035
U-50128V	green	1	0.20	0.0035







Rack 1x1 height 45

Polypropylene rack.

Dimensions with lid: 133 x 133 x 45 mm.

Available in fluorescent colours with transparent lid to help identification. **Resists -196 °C** and **autoclavable** up to **121 °C**.

Only 1 cavity.

Suitable for all kind of tubes and other objects.

code	colour	quantity	weight	volume
U-1145A	blue	1	0.12	0.00011
U-1145AM	yellow	1	0.12	0.00076
U-1145N	natural	1	0.12	0.00076
U-1145NA	orange	1	0.12	0.00076



Rack 1x1 height 70

Polypropylene rack.

Dimensions with lid: 133 x 133 x 70 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196° C and autoclavable up to 121 °C.

Only one cavity. Suitable for all kind of tubes and other objects.

code	colour	quantity	weight	volume
U-1170A	blue	1	0.12	0.0014
U-1170AM	yellow	1	0.12	0.0014
U-1170N	natural	1	0.12	0.0014





Polypropylene rack.


Dimensions with lid: 133 x 133 x 70 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196 °C and autoclavable up to 121 °C.

Sections will hold 81 tubes up to 12 mm diameter and 67 mm height. Ideal for our references **300500** and **400500**, etc.

code	colour	quantity	weight	volume
U-9970A	blue	1	0.18	0.0014
U-9970N	natural	1	0.18	0.0014
U-9970R	red	1	0.18	0.0014

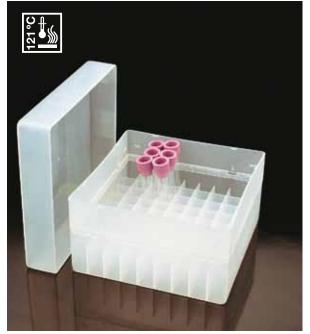

Rack 7x7 height 80

Polypropylene rack.

Dimensions with lid: 133 x 133 x 80 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196 °C and autoclavable up to 121 °C.

Sections will hold 49 tubes up to 16 mm diameter and 77 mm height.

Suitable with tubes:									
401100	4090	03.1	409107.1	4091	08.1				
code	•		colour		qu	antity	weight	volume	
U-7780)A		blue			1	0.18	0.0014	
U-7780	ON		natural			1	0.18	0.0014	
U-7780	NA		orange			1	0.18	0.0014	

Cuitable with tubes.



Rack 9x9 height 80

Polypropylene rack.

Dimensions with lid: 133 x 133 x 80 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196 °C and autoclavable up to 121 °C.

Sections will hold 81 tubes up to 12 mm diameter and 77 mm height.

400500	409003.2	409108	409108.1	801175T	901075
409003.1	409107	409107.1	801075	801275	
code		colour		weight	volume
U-9980A		blue	1	0.18	0.0021
U-9980N		natural	1	0.18	0.0021

Suitable with tubes:





Rack 7x7 height 96

Polypropylene rack.

Dimensions with lid: 133 x 133 x 96 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196° C and autoclavable up to 121 °C.

Sections will hold 49 tubes up to 16mm diameter and 93 mm height. Ideal for tubes 16×90 -92 mm or similar tubes.

code	colour	quantity	weight	volume
U-7796A	blue	1	0.17	0.0021
U-7796V	green	1	0.17	0.0021







Rack 9x9 height 96

Polypropylene rack.

Dimensions with lid: 133 x 133 x 96 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196° C and autoclavable up to 121 C.

Sections will hold 81 tubes up to 12 mm diameter and 93 mm height.

Suitable with tubes:								
300700	300701	300702	300700.4	300704	300700.5	300700.6		
300700.7	300800	300800.1	300801	300800.2	300804	301500		
400700	400800	400800.1	409003.1	409003.2	409107	409108		
409107.1	409109	409108.1	409109.1	801175T	801275	901275		

code	colour	quantity	weight	volume
U-9996A	blue	1	0.20	0.0021
U-9996N	natural	1	0.20	0.0021



Rack 5x5 round hole height 105

Polypropylene rack.

Dimensions with lid: 133 x 133 x 105 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196° C and autoclavable up to 121 °C.

Sections will hold 25 tubes up to 17 mm diameter and 102 mm height.

	Siutable with tubes:							
300705	301200	301700	401200	401700	810100T			
300707	301201	400705	401201	713100	913100			
300900	301202	400900	401204	716100	916100			

code	colour	quantity	weight	volume
U-105105A	blue	1	0.18	0.0023
U-105105Y	yellow	1	0.18	0.0023



Rack 5x5 round hole height 128

Polypropylene rack.

Dimensions with lid: 133 x 133 x 128 mm.

Available in fluorescent colours with transparent lid to help identification.

Resists -196° C and autoclavable up to 121 °C.

Sections will hold 25 tubes up to 17 mm diameter and 125 mm height.


	Suitable with tubes:								
300705	300707	300900	300904	300907	300908	301200			
301201	301202	301212	301213	301403	400705	400900			
401200	401201	401204	401403	401700	409920	429910			
429940	429942	429945	713100	716100	810100T	913100			
916100									


code	colour	quantity	weight	volume
U-15128A	blue	1	0.20	0.0032
U-15128AM	yellow	1	0.20	0.0032
U-15128N	natural	1	0.20	0.0032







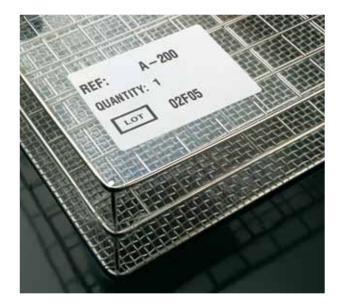

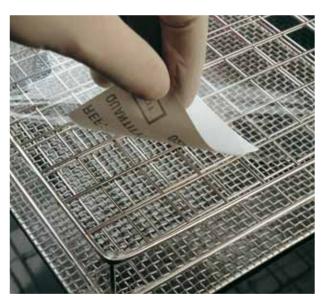

Stainless steel racks and baskets

For over 25 years Deltalab has been producing a range of fine handmade stainless steel Racks and Baskets.

Materials of the highest quality are used to ensure final products have top finishings and are sturdy and durable.

The base of all racks and baskets are manufactured from a very strong stainless steel mesh. Handmade production enables Deltalab to manufacture **Racks and Baskets customised** to your own specifications.

Just indicate the product size you need and Deltalab will provide a product made to your own requirements.

Most racks can be manufactured with a stainless steel side identification plate upon request.

Racks also available with band.

Also can be made with safety lid for transport.









Aluminium racks 10 x 5 holes

Stackable. Alpha-numerical index for better sample location. Supplied encased and interlocked, wrapped in a protecting film.

Code Z-300 suitable for 12 ml screw tubes, series 301402, and for 15 ml centrifuge tubes, series 429900 (See chapter of tubes)


code	dimensions mm	Ø hole mm	case quantity	case weight	case volume
Z-100	165 x 85 x 45	12	10	1.26	0.010
Z-250	195 x 100 x 45	15	10	1.52	0.010
Z-300	220 x 110 x 45	17	10	1.86	0.010

Aluminium racks 5 x 5 holes

Stackable. Alpha-numerical index for better sample location. Supplied encased and interlocked, wrapped in a protecting film. Code Z-075 suitable for 12 ml screw tubes, series 301402, and for 15 ml centrifuge tubes, series **429900** (See chapter of tubes)

code	dimensions mm	Ø hole mm	case quantity	case weight	case volume
Z-025	35 x 85 x 45	12	10	0.78	0.010
Z-050	100 x 100 x 45	15	10	0.89	0.010
Z-075	110 x 110 x 45	17	10	1.60	0.010

Minimum unit sale: 1

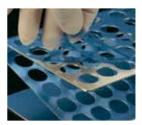

Stainless steel racks 10 x 5 holes

Stackable. Alpha-numerical index for better sample location. Supplied encased and interlocked, wrapped in a protecting film. Code **Z-600** suitable for 12 ml screw tubes, series 301402, and for 15 ml centrifuge tubes, series 429900 (See chapter of tubes)

code	dimensions mm	Ø hole mm	case quantity	case weight	case volume
Z-500	165 x 85 x 50	12	10	3.04	0.012
Z-550	195 x 100 x 50	15	10	3.94	0.012
Z-600	220 x 110 x 50	17	10	4.48	0.012




Stainless steel racks 5 x 5 holes

Stackable. Alpha-numerical index for better sample location. Supplied encased and interlocked, wrapped in a protecting film. Code **Z-800** suitable for 12 ml screw tubes, series 301402, and for 15 ml centrifuge tubes, series 429900 (See chapter of tubes)

code	dimensions mm	Ø hole mm	case quantity	case weight	case volume
Z-700	85 x 85 x 50	12	10	1.90	0.010
Z-750	100 x 100 x 50	15	10	2.06	0.010
Z-800	110 x 110 x 50	17	10	2.45	0.010

Minimum unit sale: 1



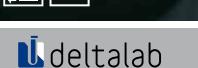

Stainless steel racks

For tubes of 12 mm \emptyset , hole 13 mm \emptyset .


code	lenght x width x height mm	tubes capacity	quantity	weight	volume
F-200	155 x 155 x 35	100 (10 x 10)	1	0.150	0.0008
F-210	155 x 80 x 35	50 (10 x 5)	1	0.100	0.0004
F-240	95 x 63 x 35	24 (6 x 4)	1	0.070	0.0002
F-211	155 x 80 x 70	50 (10 x 5)	1	0.140	0.0009



Stainless steel racks

For tubes of 14 mm Ø, hole 15 mm Ø.

code	lenght x width x height mm	tubes capacity	quantity	weight	volume
B-200	175 x 175 x 60	100 (10 x 10)	1	0.225	0.0018
B-210	175 x 85 x 60	50 (10 x 5)	1	0.145	0.0009
B-220	347 x 87 x 60	100 (20 x 5)	1	0.270	0.0018
B-240	110 x 75 x 60	24 (6 x 4)	1	0.090	0.0005

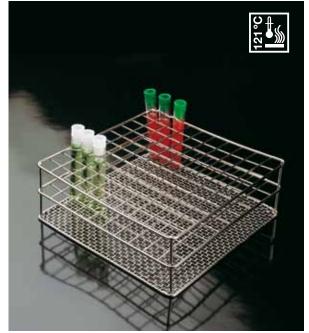


Stainless steel racks


For tubes of 18 mm \emptyset , hole 19 mm \emptyset .

Codes A-202, A-212, A-222 and A-242 are suitable for 12 ml screw tubes, series **301402**, and for 15 ml centrifuge tubes, series **429900** (See chapter of tubes)

code	lenght x width x height mm	tubes capacity	quantity	weight	volume
A-201	210 x 210 x 40	100 (10 x 10)	1	0.290	0.0017
A-211	210 x 110 x 40	50 (10 x 5)	1	0.185	0.0009
A-231	163 x 125 x 40	48 (8 x 6)	1	0.170	0.0008
A-241	125 x 85 x 40	24 (6 x 4)	1	0.160	0.0004
A-251	87 x 65 x 40	12 (4 x 3)	1	0.085	0.0002
A-200	210 x 210 x 60	100 (10 x 10)	1	0.306	0.0026
A-210	210 x 110 x 60	50 (10 x 5)	1	0.200	0.0013
A-230	165 x 125 x 60	48 (8 x 6)	1	0.185	0.0012
A-240	125 x 85 x 60	24 (6 x 4)	1	0.177	0.0006
A-250	87 x 65 x 60	12 (4 x 3)	1	0.100	0.0003
A-260	310 x 210 x 59	150 (15 x 10)	1	0.445	0.0040
A-202	210 x 210 x 85	100 (10 x 10)	1	0.320	0.0037
A-212	210 x 110 x 85	50 (10 x 5)	1	0.200	0.0019
A-222	415 x 110 x 85	100 (20 x 5)	1	0.335	0.0038
A-242	125 x 85 x 85	24 (6 x 4)	1	0.190	0.0009

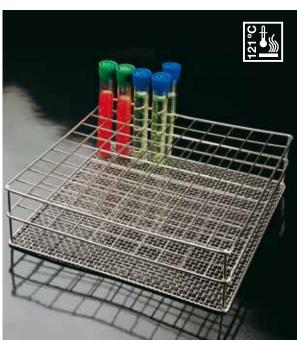



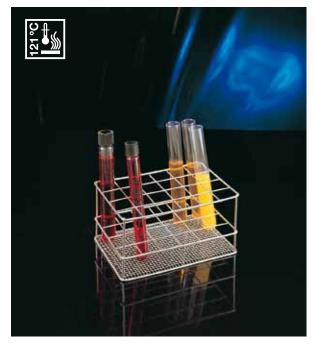

Racks for containers up to 23 mm Ø

Made of stainless steel. 25 mm square holes. Features a strong woven mesh base.

code	lenght x width x height mm	tubes capacity	quantity	weight	volume
D-200	265 x 265 x 85	100 (10 x 10)	1	0.425	0.0059
D-230	219 x 166 x 85	48 (8 x 6)	1	0.238	0.0030
D-240	165 x 110 x 85	24 (6 x 4)	1	0.145	0.0015
D-250	111 x 84 x 85	12 (4 x 3)	1	0.095	0.0008



Suitable for glass tubes shown on pages 133-134

Racks for containers up to 28 mm Ø

Made of stainless steel. 30 mm square holes. Features a strong woven mesh base.

code	lenght x width x height mm	tubes capacity	quantity	weight	volume
C-211	315 x 160 x 60	50 (10 x 5)	1	0.280	0.0030
C-256	200 x 70 x 60	12 (6 x 2)	1	0.115	0.0008





Suitable for containers shown on pages 148-149

Racks for containers up to 28 mm Ø

Made of stainless steel. 30 mm square holes. Features a strong woven mesh base.


code	lenght x width x height mm	tubes capacity	quantity	weight	volume
C-200	315 x 315 x 40	100 (10 x 10)	1	0.530	0.0040
C-254	195 x 65 x 40	12 (6 x 2)	1	0.100	0.0005





Suitable for containers shown on pages 148-149



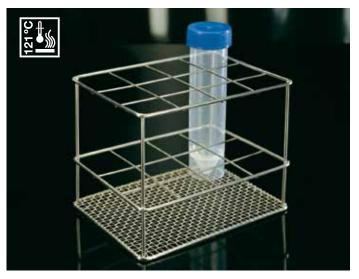



Racks for containers up to 30 mm Ø

Made of stainless steel. 33 mm square holes. Features a strong woven mesh base.

code	lenght x width x height mm	containers capacity	quantity	weight	volume
R-292	137 x 104 x 85	12 (3 x 4)	1	0,10	0,001
R-293	171 x 171 x 85	25 (5 x 5)	1	0,18	0,002
R-281	138 x 70 x 92	8 (4 x 2)	1	0,50	0,001
R-282	138 x 138 x 92	16 (4 x 4)	1	0,50	0,002
R-283	206 x 206 x 89	36 (6 x 6)	1	0,60	0,005



Suitable for containers shown on pages 35 and 37

50 ml conical tubes made of transparent polypropylene. Suitable for centrifugation. Sterile and non steril version skirted and not skirted.



You can find our 50 ml centrifuge tubes at chapters 1 and 6.





Racks for containers up to 38 mm Ø

Stainless steel. 40 mm square holes. With a strong woven mesh base.

code	lenght x width x height mm	containers capacity	quantity	weight	volume
L-250	168 x 128 x 43	12 (4 x 3)	1	0.140	0.0009
L-240	255 x 170 x 43	24 (6 x 4)	1	0.246	0.0019



Suitable for containers shown on pages 150 and 151

Racks for containers up to 47 mm Ø

Stainless steel. 50 mm square holes. With a strong woven mesh base.

code	lenght x width x height mm	containers capacity	quantity	weight	volume
I-240	315 x 210 x 45	24 (6 x 4)	1	0.325	0.0029
I-250	210 x 160 x 45	12 (4 x 3)	1	0.192	0.0015






Suitable for containers shown on pages 152 and 153

Racks for containers up to 58 mm Ø

Stainless steel. 63 mm square holes.

code	lenght x width x height mm	containers capacity	quantity	weight	volume
N-240	393 x 262 x 45	24 (6 x 4)	1	0.450	0.0046
N-250	393 x 132 x 45	12 (6 x 2)	1	0.282	0.0023



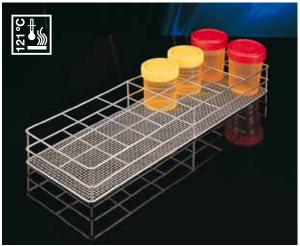



Suitable for containers shown on pages 154, 155, 156 and 158

Racks for containers up to 58 ml Ø, with lid

Made of stainless steel with a lid. Hole of 63 mm.

Ideal for carrying containers, as they stand firmly in the rack thanks to the lid.

code	lenght x width x height mm	containers capacity	quantity	weight	volume
NT-250	392 x 131 x 80	12 (6 x 2)	1	0.40	0.004





Suitable for containers shown on pages 154, 155, 156 and 158



Bottle carrier with folding handlers

Manufactured in stainless steel and designed for the transport of six 1 liter bottles such as chemical reagent bottles.

code	lenght x width x height mm	quantity	weight	volume
19700	380 x 240 x130	1	0.35	0.018



See our range of bottles and containers in chapter 12.



Baskets with folding handles

Made of stainless steel. Features a strong woven mesh base. For bottles and containers.

code	lenght x width x height mm	quantity	weight	volume					
With dividers every 12 cm. 4 departments									
H-261	480 x 300 x 150	1	1.84	0.021					
Without dividers	Without dividers								
H-260	480 x 298 x 150	1	1.18	0.021					






Instrument trays

Made of stainless steel thick woven mesh. Reinforced construction to withstand weight. With not folding handles.

code	lenght x width x height mm	quantity	weight	volume
T-200	367 x 238 x 77	1	0.98	0.0070
T-210	360 x 260 x 80	1	1.12	0.0075



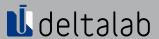



Baskets for autoclave

Stackable baskets made of stainless steel.

code	lenght x width x height mm	quantity	weight	volume
H-290	H-290 470 x 300 x 200		1.095	0.0285
H-295	570 x 300 x 200	1	1.228	0.0345



Square baskets with optional lid

Made of stainless steel. Feature a strong woven mesh base. Baskets have a folding handle for more convenience and are available with or without lid. To order baskets with lids, simply add ${\bf T}$ to the end of product code (ex: **H-200T** to order the first model with a lid). Minimum order quantity: 5 units. Other models with internal divisions are available on request. Please ask our commercial department for more details.

code	lenght x width x height mm	quantity	weight	volume
H-200	118 x 118 x 144	1	0.252	0.0020
H-210	160 x 160 x 140	1	0.324	0.0035
H-215	139 x 139 x 139	1	0.270	0.0027
H-220	210 x 210 x 180	1	0.524	0.0080
H-240	300 x 300 x 200	1	0.778	0.0180





Made of stainless steel, the base of these baskets is reinforced by a thick stainless steel mesh for added strength. Baskets feature a folding handle for more convenience. Other models with inner crossed dividers are available on request. Please ask our commercial department for more details.

code	height mm	diameter mm	quantity	weight	volume
R-200	125	150	1	0.265	0.0020
R-215	150	200	1	0.350	0.0035
R-220	200	200	1	0.528	0.0040
R-230	180	208	1	0.550	0.0040
R-240	180	230	1	0.560	0.0040
R-250	180	240	1	0.570	0.0045



Stainless steel basket for blood collection

Manufactured in stainless steel. Ideal to collect blood samples at hospital chambers, and carry them to and inside the laboratory. The basket includes:


Picture A: A-295, rack for 64, 10 ml tubes (95 x 15 mm).

Picture B: B-240, rack for 24, 5 ml tubes (75 x 12 mm).


Picture C: A-261, rack for 28, 5 ml tubes (50 x 15 mm flat bottom).

Picture D: CP-24, small-sized basket for 24 slides vertically placed.

code	description	lenght x width x height mm	quantity	weight	volume
H-401	basket with racks	340 x 290 x 85	1	1.340	0.0083

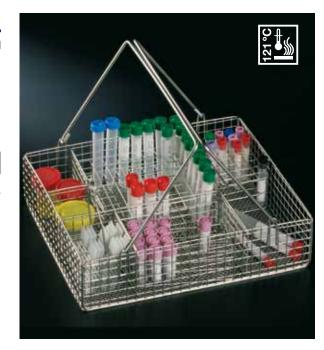

Code **H-401** is compound of:

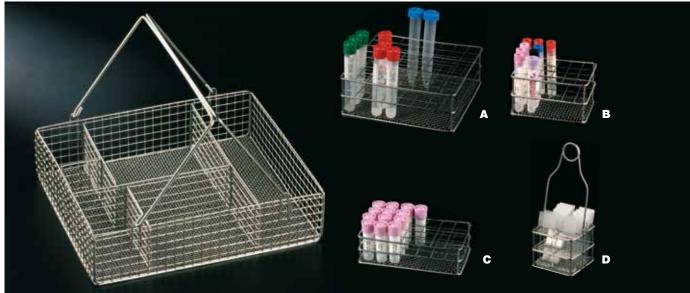




Plastic tray for blood collection tubes

Made of sturdy autoclavable polypropylene. Colour: white Ideal to carry material from lab to patient. Designed for flexibility and versatility, this tray can be arranged to meet your specific needs. It includes two 40-place test tube supports (for tubes up to \varnothing 16 mm) and one 18-place slide holder (for 76 x 26 mm slides), a large requisition section, and adjustable partitions to hold varied materials. The handle is formed by an end mounted, detachable loop that is 178 mm high.

code	dimensions mm	quantity	weight	volume
M-300	356 x 292 x 51	1	0.26	0.030





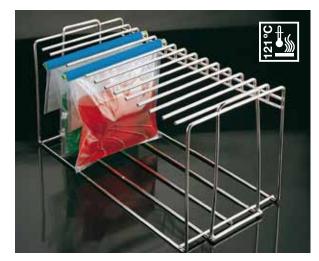


Racks for bags

Made of stainless steel, specially designed for Stomacher type homogeneizer bags. Dimensions $390 \times 200 \times 240$ mm. See bags on chapter 15.

code	description	quantity	weight	volume
15006	rack for 14 bags	1	1.1	0.140
15007	clips to adjust bags	10	0.2	0.001

Customer-designed models can be produced. Please ask the commercial department.

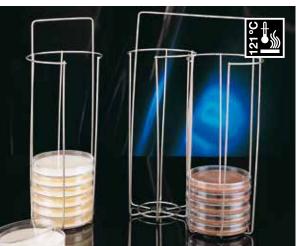


Petri dish baskets

These stainless steel baskets are suitable for dishes of diameter up to 10 cm. Two models available to hold 10 or 20 Petri dishes

code	description	dimensions* mm (height without handle)	height handle mm	case quantity	case weight	case volume
H-600	1 section, holds up to 10 dishes	105x110x250	105	1	0.13	0.005
H-601	2 sections, holds up to 20 dishes	105x215x250	100	1	0.28	0.009

*Lenght x width x height. Heigh with handle: 360 mm.

We can supply customer-designed models.

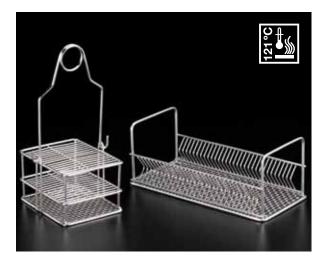

Special baskets and trays for slides

Baskets are manufactured from stainless steel wire. Cuvettes are made of neutral glass and feature stainless steel lids.

Models CP-24, CP-30, CP-50 have folding handles.

code	description	dimensions mm	case quantity	case weight	case volume
CP-30 H	basket for 30 horizontal slides	172 x 86	1	0.150	0.0011
CP-45 H	basket for 45 horizontal slides	235 x 85	1	0.186	0.0015
CP-30	basket for 30 vertical slides	70 x 90	1	0.110	0.0004







See more staining systems in chapter 4. **Histology, microscopy and staining**

Stainings racks for microscope slides

Made of stainless steel.


code	description	quantity	weight	volume
S-004	for 12 slides lenght: 435 width: 85	1	0.210	0.0020
S-002	for 24 slides lenght: 435 width: 175	1	0.410	0.0030
S-003	without divisions lenght: 435 width: 85	1	0.170	0.0020

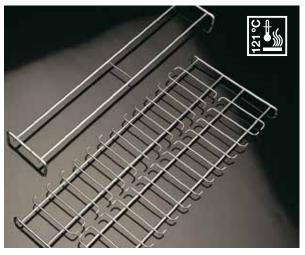

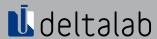





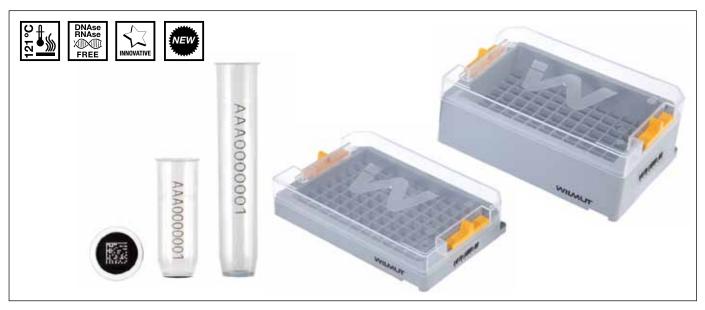

TUBES WITH 2d DATA MATRIX CODIFICATION

PP high quality tubes with U bottom made in clean room **ISO Class 7** designed to give the best traceability and maximum security for the samples. They accomplish the following:

2D datamatrix code on the bottom of the tube and lateral human readable alphanumeric codification on the side. Laser engraved Laser engraved code.

Solid one-piece construction.

Maximum capacity in the minimum space. Optimization of te capacity of the freezer: from 15% plus of samples than with other brands in the same freezer.

Low adherence PP,low cytotoxicity and extractable levels, and **DNAse, RNAse, human DNA and Endotoxines/Pyrogens free**. Allows storage up to **-86 °C**.

Autoclavables at 121 °C.


Availability of traceability files in the web site: ww.wilmut.es.

Wilmut 2D standard tubes (W-2DST) *

The tubes have a unique and randomized code.

code	gross volume ml	net volume ml	presen- tation	min. order quantity	case quantity	case weight	case volume
W053100	0.65	0.5	bag	960 u.	40x960	12.9	0.096
W123100	1.4	1.2	bag	960 u.	20x960	11.9	0.096
W053960	0.65	0.5	coded rack	20x96	6x20x96	12	0.055
W123960	1.4	1.2	coded rack	20x96	6x12x96	13	0.055

*Availability in sterile presentation and usp grade pp. For more information contact our sales dept.

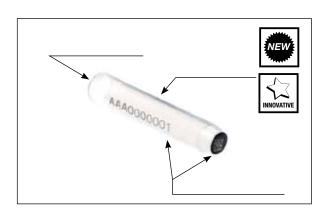




Wilmut 2D Right Position (W-2DRP) *

The tubes are visually linked to the rack and each tube has its position inside the rack.

ı	code	gross volume ml	net volume ml	presen- tation	min. order quantity	case quantity	case weight	case volume
	W052960	0.65	0.5	coded rack	20x96	6x20x96	12	0.055
	W122960	1.4	1.2	coded rack	12x96	6x12x96	13	0.055

*Availability in sterile presentation and usp grade pp. For more information contact our sales dept.

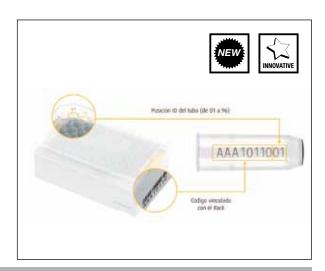




Made with high quality PP, with volume indicator in the tube, autoclavables at 121 °C and allows storage up to -86 °C.



Wilmut Seroteca (W-SER)

Low adherence PP, low cytotoxicity and extractable levels, and **DNAse, RNAse, human DNA and Endotoxines/Pyrogens free.**

code	gross volume ml	net volume ml	presen- tation	min. order quantity	case quantity	case weight	case volume
W051100	0.65	0.5	bag	960 u	40x960	12.9	0,096
W121100	1.4	1.2	bag	960 u	20 x 960	11.9	0,096
W051960	0.65	0.5	rack	20x96	6x20x96	12	0,055
W121960	1.4	1.2	rack	12x96	6x12x96	13	0,055
W059960	0.65	0.5	coded rack	20x96	6x20x96	12	0,055
W129960	1.4	1.2	coded rack	12x96	6x12x96	13	0,055

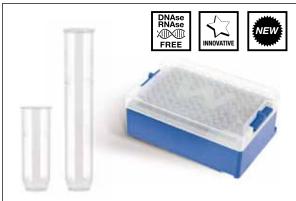

Wilmut Eco (W-ECO)

These racks are made of recycled materials.

code	gross volume ml	net volume ml	presen- tation	min. order quantity	case quantity	case weight	case volume
W058100	0.65	0.5	bag	960 u	40x960	12.9	0.096
W128100	1.4	1.2	bag	960 u	40x960	11.9	0.096
W058960	0.65	0.5	rack	18x96	18x96	1.56	0.069
W128960	1.4	1.2	rack	10x96	10x96	1.55	0.069

WILMUT ACCESSORIES

Wilmut offers a wide range of accessories to complement all the range of wilmut tubes.



Plug Style caps W-SEAL*

Made with high quality TPE that assures a perfect adjustment with the tube and an hermetic sealing at temp. up to -86 °C.

Available in differnet colors for visual sample identification.

code	presentation	colour	case quantity	case weight	case volume
W999600	sheet 10x96	natural	35x10x96	6.08	0.026
W999601	sheet 10x96	blue	35x10x96	6.08	0.026
W999602	sheet 10x96	black	35x10x96	6.08	0.026
W999603	sheet 10x96	red	35x10x96	6.08	0.026
W999604	sheet 10x96	yellow	35x10x96	6.08	0.026
W990000	bag 960 u	natural	70x960	9.22	0.026
W990001	bag 960 u	blue	70x960	9.22	0.026
W990002	bag 960 u	black	70x960	9.22	0.026
W990003	bag 960 u	red	70x960	9.22	0.026
W990004	bag 960 u	yellow	70x960	9.22	0.026
W990082	bag 960 u - ECO	black	70x960	9.22	0.026

PET 8 strip caps polyethylene

High quality polyethylene 8 strip caps.

code	presentation	case quantity	case weight	case volume
W992960	bag 120x8	35	6.2	0.026



PP WILMUT CRYOBOXES

PP Cryoboxes with standard SBS 96-well format with security closers easy to open even when are frozen.

W-RACK 96 Cryoboxes*

Made with high quality PP, resistant up to -86 °C and autoclavable at 121 °C.

code	description	case quantity	case weight	case volume
W000050	cryobox for 0.65 ml tube	18	1.14	0.0069
W000120	cryobox for 1.2 ml tube	10	1.16	0.0069
W000059	coded cryobox for 0.65 ml tube	18	1.14	0.0069
W000129	coded cryobox for 1.2 ml tube	10	1.16	0.0069

^{*}Available with cryolabel compatible with liquid handling robots. Please ask the commercial department.



Cryoboxes W-ECO RACK 96

Made with recicled PP.

code	description	case quantity	case weight	case volume
W000051	cryobox for 0.65 ml tube	18	1.03	0.0069
W000121	cryobox for 1.2 ml tube	10	1.03	0.0069




Clousers W-LOCK

Clousers for Wilmutw-RACK and W-ECO.

Available ujin different colors for visual identification of the sample. Made with POM with high resistence to impact and deformation.

code	description	case quantity	case weight	case volume
W700000	bolsa 50 - natural	35x50	2.64	0.026
W700001	bolsa 50 - blue	35x50	2.64	0.026
W700002	bolsa 50 - black	35x50	2.64	0.026
W700003	bolsa 50 - red	35x50	2.64	0.026
W700004	bolsa 50 - yellow	35x50	2.64	0.026

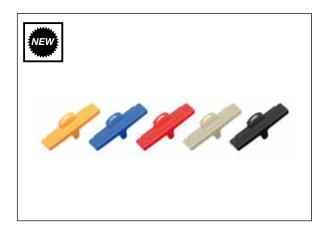




PICKERS W-PICKER AND W-MULTIPICKER

Tools designed for easily pick tubes from W-rack, one at at time or in strips of 8. Compatible with the whole range of Wilmut tubes, capped or non capped.

	code	tubes capacity	presentation	case quantity	case weight	case volume
1	W000040	1	1	40	1	0.005
2	W000048	8	1	6	6	0.004

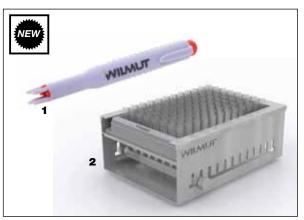





HAND OPERATED CAPPER W-CAPP10

Enables to capp up to 96 tubes at a time of any existing brand in SBS format. Compatible with the whole range of Wilmut tubes, 1,2 ml and 0,65 ml with an adapter.

Made of top quality materials.

code	description	case quantity	case weight	case volume
W000010	capper	1	3.5	0.013



AUTOMATIC CAPPER W-CAPP96

Enables to cap up to 96 tubes at a time of any existing brand in SBS format. Compatible with the whole range of Wilmut tubes, 1,2 ml and 0,65 ml with an adapter.

Compatible with robotic stations.

code	description	case quantity	case weight	case volume
W000020	capper	1	15	0,056



DECCAPERS W-DECAP

Made of aluminium.

Autoclavables and resistant to UV decontamination. Enable to easily remove caps from the tubes, one at a time or in rows of 8 caps.


code	description	presentation	case quantity	case weight	case volume
W000011	decapper for 1 tube	5	30x65	4,8	0,005
W000018	decapper for a row of 8 tubes	5	30x65	4,8	0,005



CRYOLABELS W-LABELS*

Resistant to low temperature and immersion in liquid Nitrogen. Available for all kinds of tubes, cryoboxes and cryoracks.

code	sizes (mm)	presentation (rolls of units)
W092170	10x25	4,500
W092171	19x25	2,500
W092221	11x55 + cap	3,000
W037100	11x37	2,000
W025100	32x25	2,000
W092006	6x56	3x3,000

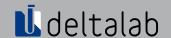

WILMUT TRACEABILITY MANAGEMENT

Wilmut offers all the tools for an easy and complete traceability of the samples. From our web site: www.wilmut.es you can download all the files of the complet 96 tubes cryobox that a bed sancanner would generate.

This allows to work with 2D code tubes without an scanner.



2D Bed Scanner W-DATAPAQ 96

Allows to scann 96 tubes rack in just 10 seconds.

Integrates fully into robotic platforms. Surface treated specially for avoiding errors while reading freezed tubes.

code	description	case quantity	case weight	case volume
W009680	bed scanner for 96 tubes	1	2.8	0.014



2D and 1D manual code reader W-RADER

Scan of tubes, boxes, or cryoracks one by one.

They are supplied with a support that allows their use in an automatic way. Allows the reading of codes on shiny surfaces.

code	description	case quantity	case weight	case volume
W001400	USB code reader	1	0,8	0,0031
W003500	wireless code reader	1	0,8	0,0041





SAMPLE MANAGEMENT SOFTWARE*

Powerful tool for sample management that allows to have a complete traceability of samples that can include: information of the initial sample, codes, processing, storage, search tools, etc.

*For more inforamtion contact our sales dept.







Space optimization and organization of biological samples.

Made with high quality materials compatible and adaptable to any storage format: boxes, racks, tubes, moulds, etc., and all types of cryostorage equipment: refrigerators, freezers, deep freezers, nitrogen tanks, walk in cool

We produce all sort of metalic racks: vertical, horizontal, cell system racks, tray racks, etc., including also customized solutions to adapt to the customer

Three cryorack ranges:

- Aluminum Range: Made of aluminum
- Standard Range: Made of Stainless Steel type 304
- Premium Range: Made of Stainless Steel type 304 with extra advantages to improve the traceability and the cool transmition to the samples.

For more inforamtion contact our sales dept

CRYOGENIC STORAGE BOXES



W-Coat

Vegetal fiber storage boxes for biological samples to temp down to -196 °C. Made with a spetial moisture repelent coating to ensure greater durability. They are supplied assembled and include divider.

code	colour	case size mm	dividers	cells Ø mm	case quantity	case weight	case volume
W350900	white	50x133x133	9x9	13	36	3.4	0.034
W370900	white	75x133x133	9x9	13	36	4.0	0.053
W310900	white	100x133x133	9x9	13	36	4.9	0.073

Available in diferent box colours and diferent dividers for all sort of tubes. Contact our commercial department.






W-Coat AIY

Assemble it yourself vegetal fiber storage boxes for biological samples to temp down to -196 °C.

Made with a spetial moisture repelent coating to ensure greater durability.

code	description	case quantity	case weight	case volume
Boxes				
W350990	50x133x133 mm	36	1.8	0.037
W370990	75x133x133 mm	36	3.2	0.055
W310990	100x133x133 mm	36	3.7	0.073
Dividers				
W030930	9x9 de Ø 13 mm, height 30 mm	36	0.8	0.004
W031030	10x10 de Ø 12 mm, height 30 mm	36	0.8	0.004
W030965	9x9 de Ø 13 mm, height 65 mm	36	1.5	0.007

Dividers availables in other sizes. Contact our commercial department.

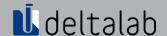







Safety containers

Containers manufactured with virgin raw materials enabling incineration; four and six liters containers are made in autoclavable polypropylene, the rest in polyethylene. Suitable for solid waste type II and III (see page 243).

They are compound by a body bearing the hazard warning label, and two lids:

- The first offers an entire opening, taking advantage of the total diameter of the mouth of the container; very useful to throw away bigger volume waste. May be closed hermetically and definitively
- The secondary one presents a partial opening, which may be closed temporarily or definitively. It embodies two devices in order to easy the extraction/disconnection of dental and analytical needles, and a 70 x 42 mm rectangular orifice. All containers, excepting code **241500**, embody a handle

In accordance with major European and international standards (ONU ADR, ISO 6706, BS 5404, etc.)

	mod.	code	capacity 	material (body and lid)	Ø mouth mm	Ø base mm	height mm	case quantity	case weight	case volume	cases per pallet
	1	241500	1.5	PE	148	122	134	60	10.26	0.120	12
	2	242000	2.0	PE	148	122	162	60	10.10	0.116	12
121°S	3	244000	4.0	PP	195	175	172	40	11.68	0.147	12
121°C 121°C	4	246000	6.0	PP	220	193	201	20	7.82	0.146	12
	5	240007	7.0	PE	210	178	240	20	8.92	0.144	12

Needle collection container

Made of high density polyethylene. Instructions for use (in English and Spanish) are printed on the containers.

Suitable for waste type III (see page 243).

Capacity: 1 I.


It features two caps.

The first one has a mouth which embodies three different devices to remove needles.

The other one is designed to close the container temporarily or definitively. Can be incinerated without release of toxic fumes, like all our containers.

code	colour	dimensions mm	case quantity	case weight	case volume
241400	yellow	Ø mouth 70; Ø base 100; height with cap 173	76	7.9	0.15



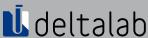



Safety containers

These containers are manufactured in virgin raw materials enabling incineration. Suitable for solid waste type II and III.

They are compound by a body bearing the hazard warning label, a handle, and two lids, allowing two using options:

- Total opening, taking advantage of the total diameter of the mouth of the container; very useful to throw away bigger volume waste
- Partial opening, 120 mm diameter, usable for syringes, contact plates, and midsize waste

In accordance with ONU ADR, ISO 6706, BS 5404 standards and major European standards.

mod.	code	capacity I	material (body and lid)	Ø mouth mm	Ø base mm	height mm	sales unit	sel.ut. weight	sel.ut. volume	cases per pallet
1	240011	11	PE PP	260	226	255	10	4.69	0.050	40
2	240021	20	PE PP	323	277	328	10	9.39	0.090	19
3	240028	25	PE PE	323	283	387	10	9.71	0.100	13



Safety labelling

Permanent adhesion labelling to alert from risk or danger.

Labels are shiny and fluorescent.


In compliance with international safety standards.

Recommendable temperature of use:

- Codes 901531, 901533: from -50 °C to 105 °C
- Code 901538: from 4 °C to 105 °C

mod.	code	description	bag quantity	bag weight	bag volume
1	901531	roll of 500 labels (labels size 25 x 25 mm)	1	0.06	0.0002
2	901533	bag of 10 units 200 x 250 mm	10	0.06	0.001
3	901538	roll of 13 m x 25 mm	12	0.57	0.001



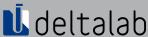



Safety containers

Manufactured in polypropylene, they are also shockproof, perforation and solvents resistant, and can be incinerated. Two ways of locking: provisional; and permanent. **Autoclavables**. Portable models **260000** and **232806** are designed to be easy opened and closed with just one hand.

- 1. **232810**, 0.20 I capacity. Approximate capacity for 100 needles. Suitable for needles and lancets. Device **type A**. Easy opening and closure with just one hand.
- 2. **260000**, 0.45 I capacity. Designed for those people who need regular injections and have to carry a needle container. With a device **type B** and a handle. Black colour.
- 3. 232811, 0.45 I capacity. Similar to the previous model, but in yellow colour, and featuring device C.
- 4. **232806**, 0.6 I capacity. Usable for insulin needles, it is also right for complete syringes. Features a device **type A** and a handle on one side.
- 232809, 1 I capacity. Suitable in places where waste quantity to be disposed of is small. With a handle and devices A and C.
- 6. 232808, 2 I capacity. Features devices type D.
- 7. **232801**, 4 I capacity. The most used container for discarding residues in general. Embodies devices **type D** and a handle.
- 8. 232802, 7 I capacity. Its height makes it ideal for contaminated pipettes. With devices type D.

All models, excepting **232810**, **232811**, **260000**, and **232806**, embody a socket to place needle caps in vertical position. All of them are in accordance with ADR (UN). Excepting portable models (**260000**, and **232806**) the others are also in accordance with BS 7320, TRBA250 and AFNOR NF X 30 500 standards. These models are constantly updated; design can be slightly modified. Suitable for solid waste type II and III



mod.	code	capacity	container dimensions mm	opening devices	useful capacity l	container weight g	case quantity	case weight	case volume
1	232810	0.20	79 x 37 x 146	А	0.20	64	100	6.2	0.040
2	260000	0.45	105 x 50 x 167	В	0.35	88	100	10.0	0.080
3	232811	0.45	105 x 50 x 167	С	0.35	88	100	9.8	0.080
4	232806	0.60	110 x 55 x 220	А	0.47	112	100	12.3	0.112
5	232809	1.00	107 x 90 x 190	A+C	0.81	125	100	13.7	0.093
6	232808	2.00	195 x 120 x 170	D	1.60	216	50	12.5	0.092
7	232801	4.00	175 x 175 x 248	D	3.00	310	50	17.4	0.130
8	232802	7.00	175 x 175 x 382	D	5.60	434	50	23.0	0.190

DEVICES:

- A. Oval aperture for syringes and lancets.
- B. Removal of needles with non return system.
- C. Removal of luerlock, luerslip and pen needles.
- D. Five kinds of aperture: a wide one to discard elements of medium size; another one designed to remove scalpel blades; one to remove luerlock and blood collection systems; to remove butterfly and luerslip needles and the last one for insulin pen.




mod. 1, 4




mod. 2



mod. 3

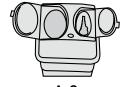

A+C mod. 5

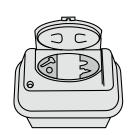

D mod. 6, 7, 8





Cardboard waste containers

Cardboard containers with lid, designed for storage and disposal of group II solid and semi-solid waste. With an integrated low density polyethylene bag (stuck in their interior). They can be incinerated without release of toxic fumes. The lid features two seals. One temporary that will avoid unpleasant smells and eliminate the risk of contamination while using the container. The second one is a definitive, positive seal: closing the bag with a belt already included, sealing the lids, and reinforcing them with an adhesive tape.

These containers embody lateral handles.

Model 3, code **250050**, 50 I capacity, is made of cardboard 2.7 mm minimal thick. PEBD 60 μ bag, adhered to the base and the walls. **Its height allows users not to need to lean forward as they throw away waste**. Easy to assemble thanks to its auto-mountable bottom. A maximum filling line is printed on its body.

Model 4, code **270045**, 50 I capacity, is manufactured in **reinforced cardboard, double thickness** (4 mm minimum). PEBD bag 54 μ adhered to the base. Safe closure with the help of raised edges. Both are printed with the biohazard anagram and text; assembling and closing instructions and drawings; identification formulary, and standards. In accordance with UN ADR (road transport), and manufactured pursuant to ISO 9001 standards. Containers are supplied folded up for space-saving. Those containers are not suitable for needle disposal. For this purpose see the special containers on previous pages.

mod.	code	quantity	capacity I	container dimensions mm	maximum load kg	weight ut. g	pallet dimensions	pallet quantitty
3	250050	10	50	263 x 263 x 756	12.5	642	100 x 120 x 180	40 x 10
4	270045	10	50	290 x 377 x 490	18.0	930	120 x 100 x 220	300

Clasificación of waste by type:

Type I waste: medical waste treated as municipal waste. Easy management (cardboard, paper, office equipment, kitchens, garages, gardening and non-infectious patients waste).

Type II waste: non specific medical waste. Preventive measures in the handling, collection, storage and transport in the field of health center (dressing materials, plasters, clothing and single-use material contaminated with blood, secretions and/or excretions).

Type III waste: specific health waste or risk waste. Preventive measures in the handling, collection, storage, transportation, treatment and disposal, as they may pose a health risk to workers and the public. They can be divided into 5 groups: medical waste or infectious anatomical waste, blood and blood products in liquid, needles, sharps and live attenuated vaccines.

Type IV waste: waste specified in singular regulations. Subject to especial requirements from the point of view of hygiene and environtment. Include: cytostatic waste,

chemical residues, expired pharmaceuticals, minerals and synthetic oils, waste metals, radioactive waste and human anatomical remains entity.

Source: NTP 372 – Spanish National Institute of Safety and Health at Work.



Rectangular waste containers (PP)

Made of yellow, perforation resistant, virgin polypropylene. Suitable for solid or semi-solid waste, II and III groups. Lids feature a central handle for a better handling in the daily use and also while carrying containers.

Secure closure thanks to their fourteen closing points. Lateral handles for further help carrying. Can be incinerated with no release of toxic furnes. Stackable when empty, or full and closed. Minimum wall thickness: 2.5 mm.

In accordance with UN ADR and AFNOR NF X30 500. Manufactured following ISO 9002 and ISO 14001.

mod	code	quantity	capacity I	container dimensions mm	maximum load (ADR) kg	weight ut. g	pallet dimensions	pallet quantity
1	240035	10	30	415 x 314 x 373	14	1,400	80 x 120 x 200	130
2	240065	10	60	415 x 314 x 575	25	2,000	80 x 120 x 200	115



Bags for Autoclave

Printed with the biohazard graphic symbol and instructions of use in five languages. Do not close the bags hermetically when placing them in the **autoclave**.

code	dimensions cm	aprox. capacity I	tickness microns	bag quantity	bag weight	bag volume
200100	50 x 75	24	30	100	1.61	0.006
200310	30 x 60	27	35	500	6.00	0.011
200311	40 x 75	45	40	300	6.35	0.022
200312	60 x 75	68	40	200	6.35	0.020
200318	75 x 90	100	40	100	4.40	0.011

Minimum order quantity: 50.000 pieces.



Autoclave resistant red bags

Heavy duty, high impact biohazard bags made of a 0.055 mm thick high molecular weight polypropylene blend.

They are autoclavable to 138 °C, saving autoclave time.

Bags printed with the standard biohazard warning symbol and precautionary procedures in four languages: English, Spanish, French and German.

They feature a **sterilization indicator** patch that darkens when exposed to steam sterilization. Do not close the bags hermetically when placing them in the autoclave.

code	dimensions cm	volume I	case quantity	case weight	case volume
200320	48 x 58	24.4	200	5.10	0.011
200321	65 x 91	75.1	200	10.14	0.020
200322	78 x 96	121.0	200	13.25	0.029



Bottles for sterilization

Made of soda glass 3.3. Autoclavable up to 140 °C.

Graduations printed in white colour. Screw caps made of blue polypropylene. Ribbed screw caps for a better handling.

Drip ring on bottle neck to avoid dripping. Ideal for chemical reagents and culture media.

They can be used fo freeze substances up to -40 °C.

It is recommended to frezze them in an inclined position (ca 45°) and filled up to max. 3/4.

code	capac. ml	screw	Ø ext. mm	height mm	case quantity	case weight	case volume
402001	100	GL45	56	100	10	1.74	0.005
402002	250	GL45	70	138	10	2.60	0.009
402005	500	GL45	86	176	10	3.86	0.016
402010	1,000	GL45	101	225	10	6.12	0.029
402020	2,000	GL45	136	260	10	10.54	0.059









Self-adhesive sterilisation indicator tapes made of semicreped paper. It sticks on different surfaces such as plastic, metal, fabrics, paper, cardboard, etc. Thickness range from 0.15 mm to 0.17 mm depending on the model. Tape will remain affixed for 18 months.

191050: indicator for sterilisation by ethylene oxide gas.

The word GAS printed in brown changes to green once it has been sterilized. Following the specifications:

- exposition to 500 mg/l EO and 60% RH for minimum 30 min at 30 °C
- exposition to 12% EO and 60% RH for minimum 48 min at 54 °C.

191051: indicator for sterilisation with poupinel method (dry heat). The green strips turn to brown by the effect of time and temperature as per the following specification:

• exposition a minimum 20 min at 160 °C or more in dry air.

191052: indicator for sterilisation with steam.

The white strips turn to dark brown by the effect of time, temperature and saturated steam as per the following specification:

- exposition at 121 °C for 10 min
- exposition at 134 °C for 2 min.

mod.	code	type of sterilisation	lenght. mm	width mm	case quantity	case weight	case volume
1	191050	ethylene oxyde gas	50	19	8	0,88	0,0026
2	191051	dry heat (poupinel)	50	19	8	0,75	0,0032
3	191052	autoclave	50	19	8	0,75	0,0032



Sterilization indicator tapes

For steam autoclave.

Resistance to cold: -40 °C.

Model 1: The tape has printed on it **CONTAMINATED**.

After sterilisation to 121 °C the word STERILIZED appears.

Model 2: The tape has printed on it **-Batch**, **-Con**. (Control) and **-Date**. After sterilisation to **121 °C** the word **STERILIZED** appears.

mod.	code	descrtiption	case quantity	case weight	case volume
1	191222	roll of 19 mm x 13 m	8	0.32	0.0007
2	191223	roll of 19 mm x 13 m	8	0.32	0.0007



Alcohol pads

Desinfecting pads with 70% isopropyl alcohol, made from a smooth, fiber-free fabric.

Ideal for use as an antiseptic for cleaning small areas.

Supplied in aluminium bags.

	code	dimensions cm	case quantity	case weight	case volume
iew	192928	6.3 x 3	24 x 100	2.86	0.010





Coloured adhesive tapes for sample identification

Ideal for use in freezers, autoclaves, etc. Peel off without leaving a sticky residue. Withstand temperatures from $-50\,^{\circ}\text{C}$ to $130\,^{\circ}\text{C}$.

High resistance to oils, acids, water, etc.

Suitable for glass, plastic, paper, metal and rubber surfaces.

Can be marked with pencil, ball-point pen, etc.

Tape roll inside diameter: 13 m rolls: 26 mm, 55 m rolls: 76 mm.

All tape models are 19 mm wide.

The code 191220 is a 12 roll tape metal dispenser that accepts 13 m rolls.

code	colour	tape length m	case quantity	case weight	case volume
191221/1	white	13	16	0.64	0.001
191221/2	red	13	16	0.64	0.001
191221/3	yellow	13	16	0.64	0.001
191221/4	green	13	16	0.64	0.001
191221/5	orange	13	16	0.64	0.001
191221/6	blue	13	16	0.64	0.001
191221.1	white	55	12	2.00	0.005
191221.2	red	55	12	2.00	0.005
191221.3	yellow	55	12	2.00	0.005
191221.4	green	55	12	2.00	0.005
191221.5	orange	55	12	2.00	0.005
191221.6	blue	55	12	2.00	0.005
191220	dispenser	-	1	2.75	0.005



Adhesive labels

White colour.

code	description	bag quantity	bag weight	bag volume
901031	54 x 36 mm printed	2 x 5,000	2.34	0.010
901031.2	54 x 36 mm plain	2 x 5,000	1.20	0.003
901041	50 x 30 mm printed	1 x 10,000	2.70	0.002
901016	50 x 30 mm plain	1 x 10,000	2.34	0.005
901030	35 x 25 mm printed	2 x 7,000	2.25	0.006








Pipetting bulb

Made of rubber. Used to avoid mouth pipetting and contamination risk. Can be opened, cleaned and **autoclaved**. Ideal for Wintrobe and Westergren pipettes.

code	description	bag quantity	bag weight	bag volume
19201	red pippeting bulb	1	0.03	0.0002





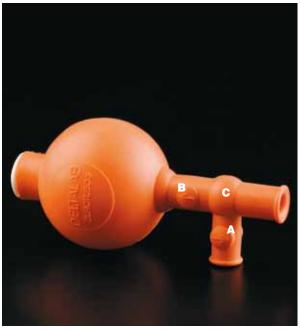








3. Drain the liquid by 4. Blow out the pressing on point pipette by pressing as the drawing indicates (point C).

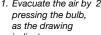

on point B 🕥 A 🔱

EUROTUBO® Pipetting bulb

This one-handed design is the simplest safety pipette filler to use available. Manufactured in natural orange rubber.

Approx. drawn capacity: 25 ml.

Single hand use, only two operating points.

Evacuate via the automatic valve.

Standard model, accommodates all pipettes. Ability to clean inside of bulb by removing patented valve and rinsing out.

code	description	bag quantity	bag weight	bag volume
19200	pippeting bulb	1	0.05	0.0004



Pipette pumps

Several models for various pipette volumes.

Designed for fast and efficient pipetting with simple, one handed operation. Pipettes fit smoothly into collar.

Rotate the knurled thumb wheel on the side for precision filling or dispensing, and press the fast release lever for quick emptying.

Easy to use and easily disassembled for cleaning.

Sizes are colour coded.

Pipette pumps resist acids and alkalies.


Attention: New models coming soon.

code	description	colour	bag quantity	bag weight	bag volume
W-100	up to 2 ml	blue	1	0.05	0.0001
W-110	up to 10 ml	green	1	0.05	0.0001
W-120	up to 25 ml	red	1	0.05	0.0001

Pipette pump support rack

Made of acrylic resin and designed to keep in position 4 pipette pumps. Suitable for the above pipette pumps.

code	dimensions mm	bag quantity	bag weight	bag volume
W-150	109 x 208 x 144	1	0.38	0.005

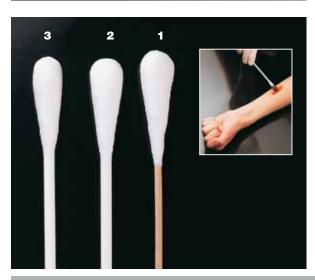


Hospital swabs

(€ (MDD)

Cotton tip Ø 11 mm.

Swabs designed for cleaning and desinfection when tip is soaked in a desinfectant, thus reducing the risk of a direct contact between user's hands and the wound or surface to clean.

mod.	code	description	case quantity	case weight	case volume
1	300234	wood + cotton 150 x 2.5 mm	20 x 50	2.20	0.012
2	300236	polypropylene + cotton 150 x 4 mm	20 x 50	2.30	0.012
3	300237	polypropylene + cotton 200 x 4 mm	10 x 50	3.30	0.012

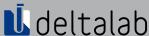





Safety glasses

(€ (PPE)

Made of polycarbonate (1 mm thickness). Scratch resistant, antifogging, antistatic lenses.

Protect against UV rays.

With lateral vents. Conform to EN 166 standard.

code	description	quantity	weight	volume
19385	glasses	5	0.47	0.002



Safety glasses

(€ (PPE)

Made of polycarbonate (2 mm thickness).

Can be worn over conventional corrective glasses.

Conform to EN 166, EN 167, EN 168, EN $\overset{1}{1}$ 69, EN 170, EN 171 and EN 172. Standard.

code	code description		weight	volume
19386	glasses	12	0.69	0.004



Emergency eyewash bottle


Made of low density polyethylene, the bottle is flexible and easy to squeeze. The bottle has an eyecup with a snap-on cap. The eyecup has a perforated plate to break up the stream providing a gentle spray to the eye. A capped air vent valve prevents drawback of contaminated eye wash into bottle. Contaminated wash exits through a drain tube, which rotates in any direction.

Instructions for use printed in various languages.

Bottle capacity: 1 l.

Supplied empty.

code	capacity ml	case quantity	case weight	case volume
19384	1,000	1	0.25	0.008

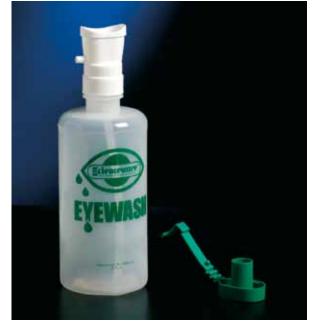



Emergency eyewash bottle with holder

Made of low density polyethylene, the bottle is flexible and easy to squeeze. The bottle features an eyecup with a cap. The eyecup has a perforated plate to break up the stream providing a gentle spray to the eye. Contaminated eyewash drips into the container attached to the cap, thus preventing any liquid from entering the bottle. Includes a green 440 x 370 mm polystyrene holder.

Supplied with screws. Bottle capacity: 1 I.

Instructions for use enclosed.

Supplied empty.

code	capacity ml	case quantity	case weight	case volume
19387	1,000	1	0.87	0.022







Wash bottles

Made a combination of high and low density polyethylene. Transluscent bottle, with HDPE cap (blue), delivery tube and seal cap (natural).

Bottle caps are ribbed for a better handling.

High resistance to all liquids, acids and bases.

Include a large neck for safe and easy filling.

code	capac. ml	neck mm	base mm	height mm	case quantity	case weight	case volume
44330000	250	30	65	155	90	5,30	0,078
44300000	500	30	75	185	60	4,85	0,080
44310000	1,000	30	90	222	45	5,12	0,106



Wide mouth wash bottles

Transluscent low density polyethylene bottles, with screw cap and delivery tube available in 3 different colours (blue, yellow and red) for easy identification. Feature a one piece cap and delivery tube assembly, thus eliminating air or liquid absorption through the cap. Include a wide neck for safe and easy filling. Ribbed caps for a better handling while opening and closing.

Moulded-in graduations.

 $250~\mathrm{ml}$ model is graduated in $25~\mathrm{ml}$ increments, $500~\mathrm{ml}$ model, in $100~\mathrm{ml}$ increments.

code	capac. ml	cap colour	int. neck Ø mm	base Ø mm	height mm	case quantity	case weight	case volume
191637.04	250	blue	38	60	140	10	0.47	0.009
191637.06	250	yellow	38	60	140	10	0.47	0.009
191637.10	250	red	38	60	140	10	0.47	0.009
191638.04	500	blue	38	75	165	10	0.97	0.019
191638.06	500	yellow	38	75	165	10	0.97	0.019
191638.10	500	red	38	75	165	10	0.97	0.019



Integral wash bottles

Made of transluscent polyethylene. Feature a one piece cap and delivery tube assembly, which eliminates air or liquid absorption through the cap and ensures a positive leaktight seal. To increase liquid flow, cut the tip of the delivery tube; close the tip by means of the cap attached to the delivery tube. Caps are ribbed for a better handling while opening and closing.

Tagging points on neck and cap for a tamper proof closure. Bottles include a fill in line.

code	capac. ml	height mm	neck d/n std.	case quantity	case weight	case volume
191633	250	143	GL32	10	0.68	0.008
191634	500	181	GL32	10	0.84	0.015





Wood tongue depressors

(€ (MDD)

Tongue depressor with smooth rounded edges. Made from high-quality wood to provide a smooth, splinter-free finish. Odour and taste-free.

Maintains a high tensile strength.

In addition to their traditional use, depressors can be used for spreading balms and ointments, for stirring liquids or mixing medicine, etc.

The sterile models are supplied with "instructions for use".

Size: 15 x 1.8 cm.

code	description		case quantity	case weight	case volume	cases per pallet
441142	non sterile		50 x 100	13.78	0.028	40
441242.2	STERILE R individual flow-pack		4 x 500	12.30	0.070	24

Sterile models are supplied with "instructions for use".



Plastic tongue depressors

(€ (MDD)

Made of polystyrene.

Colour: cream. Tongue depressor with smooth rounded edges.

Odour and taste-free. Maintains a high tensile strength.

More flexible than wood depressors.

Features a ribbed finger rest for easy handling. In addition to their traditional use, depressors can be used for spreading balms and ointments, for stirring liquids or mixing medicine, etc.

Size: 16 x 2 cm.

code	description		case quantity	case weight	case volume	cases per pallet
441142.1	non sterile		6 x 250	8.30	0.014	108
441142.2	STERILE R	individual flow-pack	4 x 500	12.66	0.070	24

Sterile models are supplied with "instuctions for use".



Absorbent paper

This 120 gsm absorbent paper provides a total protection thanks to its two layers. The top side consists of a layer of cellulose of great liquid absorption power; the reverse side consists of a layer of polyethylene ensuring that it is completely waterproof.

For use in the following applications:

- laboratory trays and tables protection
- pathological anatomy laboratories
- seed germination (between-paper method)
- spillage recovering using its leaktight surface
- particularly suitable for working with materials that are valuable (precious metals) or hazardous (toxic or corrosive materials)

Paper thickness: 0.244 mm.

Absorption: 170 g/m². Klemm absorption: 86 mm.

code	description	case quantity	case weight	case volume
120003	50 m x 50 cm reel	1 reel	3.50	0.012
120002	50 x 50 cm sheets	1x100 sheets	2.90	0.008

Filter paper

Our range of filter paper reams is made from high quality cellulose fibres to ensure a good resistance to humidity and a high absorption capacity.


Pore diameter: \pm 30-40 microns. For use in the following applications:

- laboratory worktop protection
- simple filtering operations for varied products
- labware sterilisation
- pharmaceutical product preparation

code	sheet size cm	grammage	thickness mm	Klemm absorp. (mm)	case quantity	case weight	case volume
120010	50 x 50 cm	73 g/m²	0,170 ±0,02	89	1x100 sheets	2.04	0.006





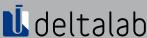



Concentrated deodorant Nilodor

Scientifically elaborated concentrated deodorant.

Ideal for hospitals (in microbiology, autopsy and oncology applications), veterinary clinics, private homes and farms to avoid unpleasant smells from animals.

This product does not cover odours.

On the contrary it quickly neutralises and eliminates them.

Non toxic.

Does not contaminate.

Does not leave stains.

2 versions available: spray or concentrated liquid.

Spray comes in a plastic bottle; liquid is supplied in a glass bottle. Expiration date: 5 years.

code	description	case quantity	case weight	case volume
910000	liquid, 7.5 ml bottle (200 drops)	12	0.54	0.0014
911000	spray, 100 ml bottle	15	1.82	0.009



Urine collection bottle

(€ (MDD)

Made of polyethylene.

Colour: natural. Collection bottle for male urine, including a 110 mm handle. The bottle has moulded-in graduations in 100 ml increments up to 1,000 ml. Total capacity: 1.1 l.

Rounded rim on the neck of the bottle. Inner diameter of the neck is 50 mm, outer diameter is 69 mm. Supplied **individually wrapped**.

Can be sterilised by ethylene oxide.

code	total capacity	case quantity	case weight	case volume	cases per pallet
19591	11	60	6.25	0.160	12

Urine collection bottle with lid

(€ (MDD)

With the same characteristics than the previous model this one features a plug in lid.

Individually wrapped.

code	total capacity	case quantity	case weight	case volume	cases per pallet
19591/T	11	60	6.54	0.160	12



Bed pan

(€ (MDD)

Made of white polypropylene.

Especially designed for disabled people, for the collection of urine and faeces. Features **rounded edges** for more comfort.

L-shaped handle for a better handling.

Allows a secure one-handed grip and avoids accidental spillage.

Useful for hanging on handles, beds, etc.

Total capacity: 2.5 l.

Autoclavable.


Can be washed in automatic machines.

Supplied individually wrapped.

code	dimensions	case	case	case	cases per
	mm	quantity	weight	volume	pallet
19586	440 x 270 x 95	20	10.00	0.140	12





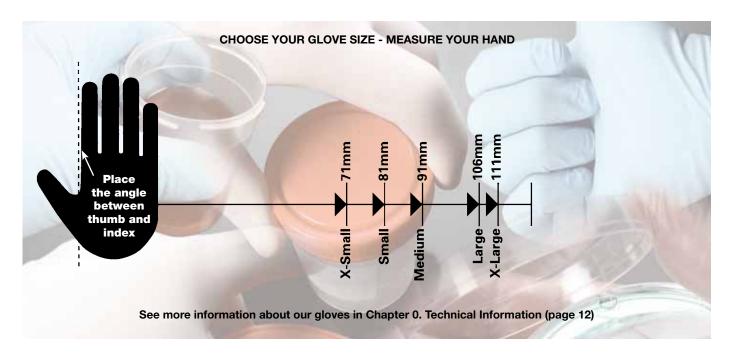

Nitrile gloves

(€ (PPE) (€ (MDD)

Single use gloves made of Butadiene Nitrate (Nitrile). Medical and laboratory use, non sterile gloves. Ambidextrous. Hypoallergenic. Microtexture.

Gloves are available **powder-free**. 100 units per case.

Minimum length: 280 mm.

AQL 1.5. In compliance with the EN 455, EN 374 and EN 420 European standard; and with the ASTM D 6319 American standard.

code	size	description	case quantity	case weight	case volume
402NP	small (6-7)	powder-free	10 x 100	4.23	0.017
402NM	medium (7-8)	powder-free	10 x 100	4.50	0.020
402NG	large (8-9)	powder-free	10 x 100	5.05	0.020

Cases per pallet: 63




Vinyl gloves

(€ (PPE) (€ (MDD)


Single use gloves made of top quality polyvinyl. Non sterile gloves. Designed for medicale use. Ambidextrous. Hypoallergenic (powder-free gloves). Gloves are available **powder-free** or **pre-powdered**. Both models feature beaded cuffs. 100 units per case.

Minimum length: 240 mm.

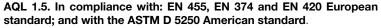

code	size	description	case quantity	case weight	case volume
202VP	small (6-7)	powder-free	10x100	4.30	0.020
202VM	medium (7-8)	powder-free	10x100	5.06	0.014
202VG	large (8-9)	powder-free	10x100	5.11	0.014
200VP	small (6-7)	pre-powered	10x100	4.14	0.013
200VM	medium (7-8)	pre-powered	10x100	5.06	0.014
200VG	large (8-9)	pre-powered	10x100	4.62	0.014

Cases per pallet: 77

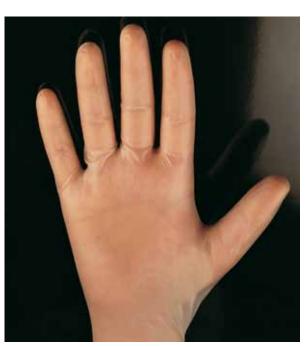



Latex gloves




Single use gloves made of top quality latex. Non sterile gloves. Designed for medical and laboratory use. Ambidextrous. Gloves are available **powder-free** or **pre-powdered**. Both models feature beaded cuffs. 100 units per case. Minimum length: 240 mm, average thickness 0.08 mm.

AQL 1.5. In compliance with the EEN 455, EN 374 and EN 420 European standard and with the ASTM D 3578 American standard.

code	size	description	case quantity	case weight	case volume
102LP	small (6-7)	powder-free	10x100	5.90	0.017
102LM	medium (7-8)	powder-free	10x100	6.28	0.020
102LG	large (8-9)	powder-free	10x100	6.76	0.017
100LSP	extra-small (5-6)	pre-powered	10x100	5.62	0.016
100LP	small (6-7)	pre-powered	10x100	5.24	0.017
100LM	medium (7-8)	pre-powered	10x100	6.40	0.017
100LG	large (8-9)	pre-powered	10x100	6.90	0.017

Cases per pallet: 63

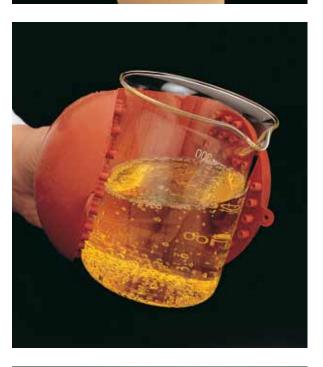

Safety gloves

Provide superior protection when handling hot, cold or damp objects. Made of strong non toxic silicone rubber which resists tearing and withstands temperatures from **-55** °C up to **260** °C. Thumb and fingers fit into end pockets. The gripping surfaces have multiple concave tipped studs for non-slip grip (48+48 studs).

code	dimensions	case	case	case
	mm	quantity	weight	volume
19270	190 x 100	1	0.19	0.0010



Thermic safety gloves

(€ (EPI)

All cotton, asbestos free gloves. They are ideal for taking hot objects out not only from the autoclave but also from the oven, etc., and handling them. Gauntlet with black sewn reinforcements at the edges for long lasting use. Heat resistant up to **230 °C**.

They may be both used with the washing and the drying machine.

code	total lenght mm	sleeve lenght mm	bag quantity	bag weight	bag volume
400S	470	280	1 pair	0.36	0.002






Measuring cylinders with pentagonal base

Material: polypropylene. **Autoclavable** up to **121 °C**. In compliance with ISO 6706 and BS 5404 standards.

code	capac. ml	grad. ml	subdiv. ml	O.D. Ø mm	height mm	case quantity	case weight	case volume
408025	25	5	0.5	19	195	30	0.78	0.010
408050	50	10	1.0	26	200	30	1.10	0.010
408100	100	10	1.0	31	250	30	1.47	0.018
408250	250	20	2.0	41	315	12	1.40	0.015
408500	500	50	5.0	56	360	12	1.94	0.029
481000	1,000	100	10.0	66	440	6	1.58	0.029
482000	2,000	200	20.0	80	530	6	3.75	0.048

Minimum order quantity: 1.

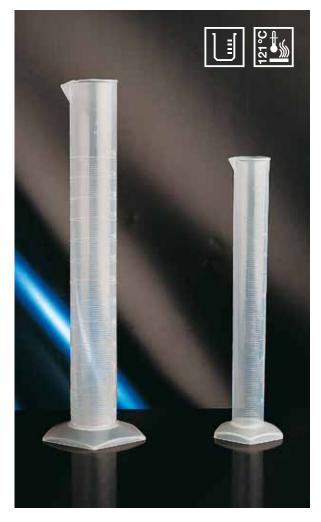

Graduated measuring cylinders pentagonal base

Made of polypropylene. Autoclavable up to 121 $^{\circ}$ C. ISO 6706 and BS 5404 standards.

						POL	YPROF	YLENE
code	capac. ml	grad. ml	subdiv. ml	O.D. Ø mm	height mm	case quantity	case weight	case volume
192562*	50	10	1.0	26	200	30	1.15	0.011
192563*	100	10	1.0	31	250	30	1.31	0.021
192564	250	20	2.0	41	315	12	1.50	0.014
192565	500	50	5.0	56	360	12	2.50	0.028
192566*	1,000	100	10.0	66	440	6	2.15	0.029
192567*	2,000	200	20.0	80	530	6	3.25	0.048

^{*}Minimum order quantity: 1.

Made of **TPX**. **Autoclavable** up to **170 °C**.


150 6706 and B5 5404 standards.							TPX	
code	capac. ml	grad. ml	subdiv. ml	O.D. Ø mm	height mm	case quantity	case weight	case volume
192571	25	5	0.5	19	195	30	0.85	0.008
192572	50	10	1.0	26	200	30	1.15	0.011
192573	100	10	1.0	31	250	30	1.40	0.018
192574	250	20	2.0	41	315	12	1.45	0.014
192575	500	50	0.5	56	360	12	2.45	0.028
192576	1,000	100	10.0	66	440	6	2.15	0.030
192577	2,000	200	20.0	80	530	6	3.25	0.048

Minimum order quantity: 1.







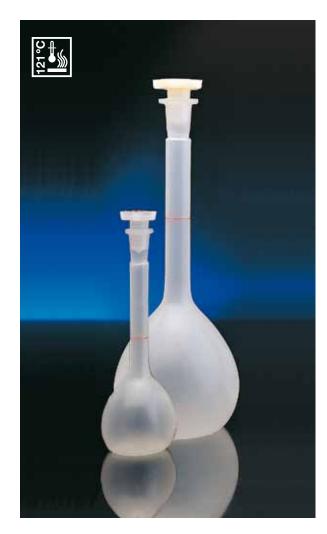

Volumetric flasks with stopper

Material: polypropylene. Autoclavable up to 121 °C.

code	capac. ml	neck	case quantity	case weight	case volume
191484	25	12 / 21	20	0.24	0.002
191485	50	14.5 / 23	20	0.32	0.005
191486	100	14.5 / 23	20	0.39	0.012
191487	250	19 / 26	10	0.52	0.009
191488	500	19 / 26	10	0.60	0.018
191489	1,000	19 / 26	5	0.55	0.019

Minimum order quantity: 1.



NOTE: NS Neck:

The first measure regards to the inner \varnothing of the lower neck. The second one refers to the inner \varnothing of the upper neck.



Material: polypropylene. **Autoclavable** up to **121 °C**. Graduated.

code	capac. ml	grad. ml	neck	quantity	weight	volume
191460	50	10	14.5 / 23	1	0.02	0.001
191461	125	25	14.5 / 23	1	0.03	0.001
191462	250	25	19 / 26	1	0.03	0.001
191463	500	50	24 / 29	1	0.10	0.0002
191464	1,000	50	29 / 32	1	0.20	0.0002
191465	2,000	100	34.5 / 35	1	0.30	0.0003





NOTE: NS Neck:

The first measure regards to the inner \varnothing of the lower neck. The second one refers to the inner \varnothing of the upper neck.



Graduated beakers

Made of TPX. **Autoclavable** to **170 °C**. ISO 6706 and BS 5404 standards.

code	capac. ml	subdiv. ml	Ø mm	height mm	case quantity	case weight	case volume
191720	25	1	35	49	20	0.17	0.001
191721	50	2	42	60	20	0.28	0.003
191722	100	5	52	72	12	0.70	0.004
191723	250	10	71	96	16	1.45	0.011
191724	500	10	90	120	12	2.10	0.015
191725	1,000	20	110	149	4	1.20	0.010
191726	2,000	50	135	184	4	2.00	0.018
191727	3,000	500	160	200	2	1.80	0.014

Minimum order quantity: 1.



Graduated beakers

Made of transparent polypropylene. Autoclavable a 121 $^{\circ}$ C. ISO 6706 and BS 5404 standards. Blue graduation.

•								
code	capac. ml	subdiv. ml	Ø mm	height mm	case quantity	case weight	case volume	
191822	50	2	42	60	20	0.31	0.004	
191823	100	5	52	72	12	0.35	0.004	
191824	250	10	71	96	16	0.75	0.010	
191825	500	10	90	120	12	1.05	0.015	
191826	1,000	20	110	149	4	0.49	0.010	
191827	2,000	50	135	184	4	0.60	0.018	
191828	3,000	500	160	200	2	0.90	0.014	
191829	5,000	500	191	230	2	1.40	0.022	




Graduated beakers

Material: transparent polypropylene. Autoclavable up to 121 $^{\circ}$ C. In compliance with ISO 6706 and BS 5404 standards.

code	capac. ml	Ø mm	height mm	case quantity	case weight	case volume
409025	25	35	50	20	0.22	0.002
409050	50	42	60	20	0.29	0.004
409100	100	53	72	12	0.33	0.004
409250	250	71	96	16	0.72	0.011
409500	500	90	120	12	0.88	0.015
491000	1,000	110	149	4	0.49	0.009
492000	2,000	135	184	4	0.86	0.017
493000	3,000	160	200	2	0.88	0.022
495000	5,000	191	230	2	2.00	0.022





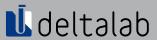



Beakers with three spouts

Graduated beakers made of ultra clear polypropylene Autoclavable up to 121 °C.

Disposable.

Moulded-in graduations.

code	volume ml	minimum order quantity	case quantity	case weight	case volume
490050	50	10	100	1.02	0.011
490100	100	25	100	1.38	0.011
490250	250	10	100	2.50	0.028
490400	400	10	100	3.30	0.044
490900	900	10	100	5.89	0.060



Conical measuring flasks

In polypropylene. Autoclavable up to 121 $^{\circ}\text{C}.$

code	volume ml	height mm	case quantity	case weight	case volume
19423	250	160	18	1.14	0.019
19424	500	180	8	1.00	0.028
19425	1,000	270	3	1.00	0.016





Analytical funnels

Material: polypropylene. **Autoclavable**. 60° angle.


code	Ø mm	stem Ø mm	lenght stem mm	total lenght mm	case quantity	case weight	case volume
19145	25	4	37	56.3	20	0.03	0.0008
19146	35	5	37	63.4	20	0.05	0.0019
19147	45	5	43	79.0	20	0.07	0.0035
19148	65	9	62	112.0	20	0.20	0.0110
19150	80	10	70	132.0	20	0.34	0.0135
19152	100	11	82	161.5	1	0.03	0.0007
19153	120	11	86	184.0	1	0.03	0.0037



High speed funnels

 $\label{eq:material:TPX} \textit{Material: TPX}. \textit{With interior helicoid channels for a quicker filtration. Prevents the filter paper from sticking. Transparent. \textbf{Autoclavable}.$


code	Ø mm	stem Ø mm	lenght mm	case quantity	case weight	case volume
19175	51	3	150	6	0.71	0.013
19172	70	3	150	6	0.13	0.004
19173	100	7	108	4	0.18	0.006
19174	140	10	132	6	0.95	0.025

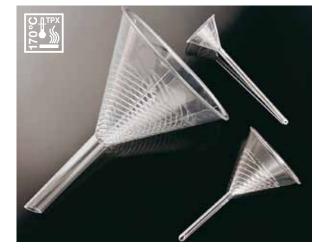

Long-stem funnels

Material: polypropylene. 60°C angle. **Autoclavable**.

code	Ø mm	stem Ø mm	lenght stem mm	total lenght mm	case quantity	case weight	case volume
19162	60	8	150	194	10	0.16	0.010
19164	80	8	143	206	10	0.19	0.016
19165	100	8	145	244	10	0.33	0.016

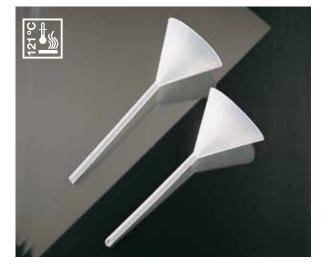



Funnels

Material: polypropylene.

Autoclavable.

code	stem internal Ø mm	stem external Ø mm	bouche mm	quantity	weight	volume
194171	5.7	8.1	80	1	0.027	0.0001
194172	9.1	11.7	100	1	0.022	0.0005
194173	9.6	13.0	120	1	0.032	0.0005
194175	13.9	18.2	180	1	0.080	0.0012
194176	18.5	21.8	220	1	0.160	0.0028
194177	21.7	26.4	260	1	0.20	0.0081
194178	24.7	29.7	300	1	0.34	0.0110







Powder funnels

Material: polypropylene. **Autoclavable**. 60° angle.

code	Ø ext. mm	stem Ø mm	case quantity	case weight	case volume
19171	60	15	20	0.17	0.002
19167	80	15	20	0.28	0.013
19168	100	25	20	0.56	0.010
19169	120	30	1	0.05	0.001
19170	150	35	1	0.05	0.003
19166	180	40	1	0.09	0.008



«IMHOFF» sedimentation cone

Material: SAN.

Moulded graduations according to DIN standards 12672.

Resistant to temperatures from -20 °C to 85 °C.

mod.	code	capac. mm	Ø mm	height mm	case quantity	case weight	case volume
1	191055	1,000	125	480	4	1.75	0.019



Cone

Rack for sedimentation cones

Material: PMMA.

For 2 cones. Resistant to temperatures from **–40 °C** to **85 °C**.

1	mod.	code	dimensions mm	case quantity	case weight	case volume
	2	191056	150 x 300 x 290	1	1.20	0.001




Büchner funnels

Material: polypropylene.

Unbreakable and resistant to corrosion and heat.

Easily dismantled for cleaning. **Autoclavable**.

code	capac. ml	filter Ø mm	upper Ø mm	bottom Ø mm	holes mm	case quantity	case weight	case volume
19437	40	42.5	13.8	9.6	1	1	0.05	0.006
19438	70	55	15.9	10.7	1	1	0.05	0.007
19440	285	80	17.7	11.4	1.5	1	0.13	0.002
19441	390	90	17.7	11.2	2	1	0.20	0.002
19442	810	110	30.1	22.6	2	1	0.30	0.004
19443	2,100	160	30.4	19.7	2.5	1	0.60	0.012
19445	6,000	240	37.7	23.7	3	1	1.00	0.028

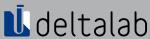



Burette clamps

Made of autoclavable polypropylene.

Clamps can fit onto rods \emptyset 8 to 14 mm. Stainless steel centre arm includes a spring and firmly holds in place the burette by means of rubber pads. Burette clamps do not obscure graduations.

Suitable for all burettes. Base is supplied separately, see details below.

mod.	code	description	case quantity	case weight	case volume
1	19139	for 1 burette	5	0.21	0.006
2	19140	for 2 burette	5	0.31	0.013

Rectangular base for burettes

Base made of **autoclavable** polypropylene, with chromium plated steel rod. Rods can be screwed one onto the other to obtain different heights.

All models come with 2 rods of Ø 250 x 12 mm.

2 models available: base with central hole, and base with off centre hole.

mod.	code	description	case quantity	case weight	case volume
3	19266	base with central hole and + 2 rods	1	0.75	0.001
4	19267	base with off-centre hole and + 2 rods	1	0.80	0.004

Rods available separately.

Please contact our Commercial department.



Separation funnel holder

Made of autoclavable polypropylene.

Supplied with two adaptors to fit support rods diameter 8 to 14 mm. Designed for funnels Ø 75 mm minimum.

code	description	case quantity	case weight	case volume
19970	for 1 funnel	5	0.05	0.0004





Funnel holders

Made of **autoclavable** polypropylene.

Double holder, for funnels from 25 mm (adaptor supplied) up to 170 mm. Without adaptor, holders can be used with powder funnels with diameter not exceeding 35 mm.

Suitable for use on \emptyset 8 to 14 mm rods (which include codes **19266** and **19267** for example) by means of adaptors (included).

code	description	case quantity	case weight	case volume
19269	for 2 funnels	5	0.65	0.001







Pipette stand

Made of polypropylene. Consists of a rounded base with a rod, and a circular pipette stand to hold up to 28 pipettes.

Base and pipette stand are supplied separately.

- Code **19262**, rounded base with rod, is specially designed to hold the pipette stand.
- \bullet The stand code **19265** accomodates 18 pipettes Ø up to 10 mm, and 10 pipettes Ø up to 15 mm.

The lower disc has little holes for liquid draining. The ring included with the stand allows its raising and lowering on the rod.

code	description	case quantity	case weight	case volume
19262	rounded base	1	0.18	0.002
19265	18 small and 10 large holes stand	1	0.35	0.007



Rotary pipette stand

This stand which can hold 94 pipettes consists of 2 discs rotating on a central vertical axis for convenient pipette selection. The lower rotating disc has tapered sides and open bottoms to protect tips and allow water run off. Rotating discs are 18 cm diameter.

The pipette stand comes unassembled with instructions for quick and easy assembly.

Made of autoclavable polypropylene.


Can accomodate all our serology pipettes presented on page 39-194, except the 25 ml model.

code	dimensions mm (Ø x h)	case quantity	case weight	case volume	
19254	230 x 480	1	0.39	0.004	





See Chapter 1. Microbiology

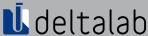

Plastic dessicator

Base: polypropylene. Cover: transparent polycarbonate. Will hold a vacuum of 740 mm Hg for 24 hours, with loss not exceeding 20 mm Hg. The cover seals onto a neoprene O-ring set into the flange of the base to ensure a positive seal. Non-return valve (polycarbonate) and vacuum release tap (polyethylene) assembly included; a specially grooved stopper in this assembly allows air to be re-admitted very slowly to avoid disturbing the contents. A little grease spread on the O-ring will assist in establishing initial vacuum. The dessicator includes a removable pan for holding calcium chloride or other drying agents.

code	aprox. capac. I	Ø mm	h. ext. mm	h. int. mm	case quantity	case weight	case volume
19230	4.5	200	230	175	1	2.00	0.004
19554	9.2	250	300	225	1	1.38	0.017



Dessicator discs

(€ (IVD)

Material: white polypropylene. 25 holes. Autoclavable.

code	Ø	for	case	case	case
	mm	dessicator	quantity	weight	volume
19553	239	19554	10	1.57	0.010

Minimum order quantity: 1



Glass dessicator

Base and cover made of 3.3 soda glass. High resistance to chemicals. Can be used with drying agents such as calcium chloride, silica gel, phosphor pentoxide and even sulfuric acid.

Before use, clean the sealing surfaces and lubricate them with grease for vacuum.

Will hold a vacuum of 500 mm Hg for 24 hours.

Dimensions: 200 mm Ø, 270 mm high. Capacity approx. 2.5 l.

code	description	case quantity	case weight	case volume
19232	dessicator with a stainless steel plate	1	5.10	0.026
19235	dessicator with a stainless plate and PTFE valve	1	5.10	0.026









Stainless steel water still

Automatic stainless steel water distiller.

It supplies an output of 6 - 7 l/h of high purity water (conductivity: $2.5~\mu S \cdot cm^{-1}$ a $20~^{\circ}C$) free of metallic ions and pyrogen. The water quality is considered ASTM III/IV (laboratory grade). In order to always obtain the maximum degree of water purity, we recommend to control periodically the quality of the water supplied to the still and to maintain the inner part of the still in good cleaning conditions.

It embodies two security systems in order to protect still from overheating and save energy: the first system controls the temperature level and turns off the still in case the temperature rises more than required. The second one interrupts the process in case water supply fails. Easy maintenance.

Supplied assembled.

Electrical consumption: 6,000 W

Heaters: 2 x 3 KW Voltage: 230 V Intensity: 30 A

Water consumption: 60 - 70 l/h

code	code dimensions mm		case weight	case volume
19580	260 x 260 x 627	1	16	0.12



Flask weight ring

Heavy vinyl covered ring ideal to improve stability of Erlenmeyer flasks in water bath.

High chemical resistance.

code	erlenmeyer capacity	internal Ø mm	quantity	weight	volume
192470	125 - 500 ml	48	1	0.32	0.0003
192471	250 - 1,000 ml	51	1	0.59	0.001
192472	500 - 2,000 ml	57	1	0.62	0.0003




Vacuum pump

Material: polypropylene. This vacuum pump combines excellent suction capacity with limited water consumption, under low pressure (0.5-1 kg/cm²) as well as high presssure (10 kg/cm²). The built-in, non return valve protects the pump from possible water back-flow. **Autoclavable**.

Can be dismounted for easy cleaning.

code description		quantity	weight	volume
191395	water jet vacuum pump	1	0.025	0.00001

Minimum order quantity: 1.

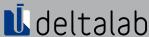



Test tube tongs


Made of brown plastic.

High resistance to heat up to $70\ ^{\circ}C$. Ideal for handling hot tubes or to heat test tubes.

With a longer arm for a better handling and to avoid close contact between fingers and the flames of the burner.

Total lenght: 17.5 cm


code	code description		bag weight	bag volume
19509	plastic tongs	10	0.278	0.001

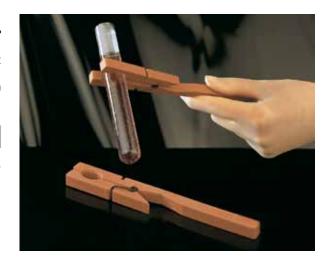

Laboratory spatulas

Material: polystyrene.

code	description	length mm	bag quantity	bag weight	bag volume
19593	spatula-spoon	210	10	0.067	0.0004
19596	spatula-spoon	180	10	0.012	0.0670

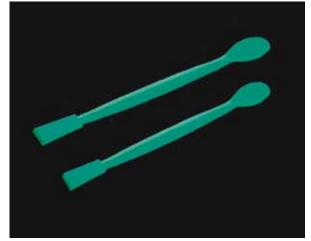

Measuring scoops

Material: white polypropylene. **Autoclavable**.

code	capacity ml	length mm	bag quantity	bag weight	bag volume
191057	10	100	12	0.055	0.004
191058	25	135	12	0.095	0.001
191059	50	160	12	0.162	0.004
191060	100	200	12	0.290	0.005
191062	250	260	6	0.271	0.007
191063	500	315	6	0.46	0.013
191064	1,000	400	6	0.72	0.010






Measuring scoops

Material: high density polyethylene.

code	capacity ml	length mm	bag quantity	bag weight	bag volume
191065	25	135	12	0.076	0.0014
191066	65	185	12	0.20	0.0011
191067	110	215	12	0.35	0.0036
191068	150	250	12	0.42	0.0040
191070	350	310	6	0.39	0.005
191071	750	350	6	0.69	0.019
191069	1,250	400	6	0.99	0.0115





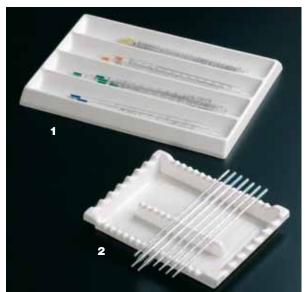







PVC Pipette tray

Made of white PVC, resistant to temperatures from **-20 °C** to **80 °C**. **Model 1**, code **19252**, is perfect to fit inside drawers. It includes four compartments that can hold up to thirty pipettes of 1, 2, 5, or 10 ml.

Model 2, code **19996**. Sixteen pipettes up to 10 mm diameter can be accommodated longitudinally; laterally, seven pipettes up to 20 mm diameter. Tray edges are ergonomically designed for a better handling. It is also useful with other instruments.

mod.	code	dimensions mm	case quantity	case weight	case volume
1	19252	426 x 300 x 30	1	0.53	0.006
2	19996	283 x 216 x 40	6	1.38	0.016



See our range of pipettes on pages 39

PVC Antiacid trays

Resistant to temperatures from **-20 °C** to **80 °C**. Ideal for photographic purposes thanks to their ribbed interior surface.

code	dimensions mm	quantity	weight	volume
19280	200 x 150 x 50	1	0.13	0.0007
19281	250 x 200 x 60	1	0.23	0.0012
19282	320 x 260 x 70	1	0.47	0.0020
19283	350 x 300 x 85	1	0.51	0.0026
19284	430 x 330 x 95	1	0.78	0.004
19285	520 x 420 x 95	1	1.14	0.008

ABS trays

They withstand temperatures from **-40 °C** to **85 °C**.

code	dimensions mm	body weigth g	quantity	weight	volume
195701	200 x 150 x 20	78	20	1.80	0.007
195702	300 x 150 x 20	115	20	2.50	0.008
195703	350 x 250 x 20	217	20	5.72	0.014
195705	200 x 150 x 40	91	20	1.88	0.007
195706	300 x 150 x 40	134	20	2.76	0.008
195707	350 x 250 x 40	254	20	6.00	0.019
195709	200 x 150 x 80	110	20	2.36	0.008
195711	350 x 250 x 80	299	20	6.64	0.022
195713	400 x 300 x 20	345	10	3.50	0.014

Drying rack

PVC coated steel wire drying rack for general labware.

It embodies 55 pegs, and 5 brackets specially designed to hold erlenmeyer flasks and beakers. Pegs are plastic coated for better protection.

Two rings on the rack allow the rack to be mounted on the wall.

It includes a tray in order to avoid drops on the surfaces where the rack is placed.

code	dimensions	case	case	case
	mm	quantity	weight	volume
19212	640 x 360 x 140	1	2.12	0.130





Pipette and burette rinsing set

Made of polyethylene and polypropylene.

Can hold pipettes and burettes of length up to 60 cm.

Ensures a fast, safe washing.

Ideal for use in labs where caustic and toxic materials are used.

code	description	height mm	Ø mm	case quantity	case weight	case volume
19215	burette washer	990	165	1	3.50	0.110
19217	pipette washer	734	165	1	3.35	0.080
19218	pipette jar	650	165	1	1.80	0.042
19219	pipette basket	650	130	4	2.00	0.026
19221	pippete jar	500	165	1	1.50	0.030
191219	extension for 19219	650	130	2	0.56	0.00001

Suggested sets:

• For pipettes: 1 x **19217**

2 x **19219** 1 x **19218**

• For shorter pipettes: 1 x 19217

1 x **19221**

• For burettes: 1 x **19215**

2 x 19219 2 x 191219 1 x 19218



Drying rack

Wall mountable high impact polystyrene drying rack for general labware.

This rack consists of a single mould 4 mm plate with 72 peg sockets.

Push-in pegs can be placed in any configuration on the rack, and are easily removable to accomodate odd shaped items.


Peg sockets have closed ends to avoid dripping and eliminate the risk of biological contamination.

Rack suitable for items with diameter of the mouth larger than 15 mm.

A **drip channel** collects dripping that drains via a drain tube provided. Includes 11 additional peg sockets (Ø 6 mm) for smaller parts.

Comes complete with a mounting kit and fixing template.

code	dimensions mm	case quantity	case weight	case volume
19213	450 x 63 x 110	1	5.00	0.050

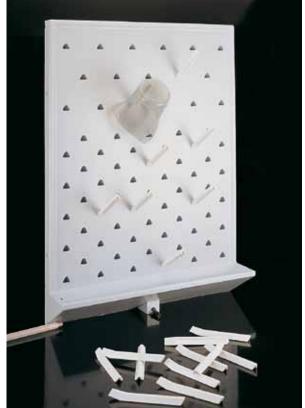







Hexagonal weighing dishes

Made of white translucent polystyrene. Easily bent into pouring spouts, the dishes enable non-spill transfer. Biologically inert, contaminant-free, they are resistant to diluted acids, aqueous solutions, alcohols and bases.

Antistatic. Suitable for microwave.

Will withstand temperatures from 80 °C to -10 °C.

code	capacity ml	upper Ø (a) x lower Ø (b)	case quantity	case weight	case volume
1911106	20	35 x 25	500	0.25	0.001
1911107	50	64 x 47	500	1.00	0.003
1911108	200	110 x 85	500	2.30	0.007
1911109	350	130 x 89	500	2.80	0.013



Squared weighing dishes

Made of white translucent polystyrene.

Easily bent into pouring spouts, the dishes enable non-spill transfer. Biologically inert, contaminant-free, they are resistant to diluted acids, aqueous solutions, alcohols and bases. **Not antistatic**. Suitable for microwave.

Will withstand temperatures from 80 °C to -10 °C.

code	dimensions	weight (g)	bag quantity	bag weight	bag volume
1911103	41 x 41 x 8	0.34	500	0.31	0.001
1911104	89 x 89 x 25	1.07	500	1.32	0.005
1911105	140 x 140 x 22	3.27	500	3.31	0.013



Weighing bottles

Light-weight, re-usable and **autoclavable**, made of translucent polypropylene. Push on lids ensure an airtight seal.

High resistance to chemicals. Tops can be squeezed for more convenient pouring. It is recommended to remove lids before autoclaving.

code	tare (g)	capacity ml	Ø mm	height mm	bag quantity	bag weight	bag volume
19340	8	30	30	50	10	0.07	0.002
19345	10	50	50	30	10	0.11	0.004



Watch glasses

Made of **autoclavable** polypropylene. Moulded base for more stability. Easy **cleaning**.

	-				
code	weight (g)	Ø mm	case quantity	case weight	case volume
19452	3.7	60	10	0.05	0.0002
19454	5.5	80	10	0.10	0.0002
19455	10	100	10	0.09	0.0018





Reinforced silicone tubing. Non toxic

Reinforced silicone tubing with a tough polyester inner-braid imbedded in the wall, it is much more resistant to pressure than standard silicone tubing. Very flexible. Its clear wall permits positive flow monitoring and inspection. It has a smooth inner wall that will not trap bacteria-harboring particles, making it ideal for food and beverage use.

Indicated for elevated working pressure.

It is not toxic.

Resist temperatures from -50 °C to 200 °C (beaks of 250 °C).

Minimum order quantity: one roll of 15 meters.

code	dimensions mm Ø int. x Ø ext.	bag weight	bag volume
350613	6 x 13	2.10	0.0043
350815	8 x 15	2.33	0.0130
359517	9.5 x 17	3.06	0.015



Silicone sheets

 50×50 cm silicone sheets ideal as a laboratory benchtop protection or for use in autoclaves to protect any material.

For cleaning simply use a wet cloth.

Resistant up to 260 °C.

code	thickness mm	bag quantity	bag weight	bag volume
350002	2	1	0.60	0.030
350004	4	1	1.20	0.0012






Silicone mat


Made of red silicone, this mat will resist temperatures up to **260 °C**. Features a high adherence and anti-skid surface.

For use over the benchtop or over laboratory stirrers to prevent beakers from slipping.

Dimensions: 25 x 25 x 0.4 cm.

code	description	bag quantity	bag weight	bag volume
900000	silicone mat	1	0.132	0.004

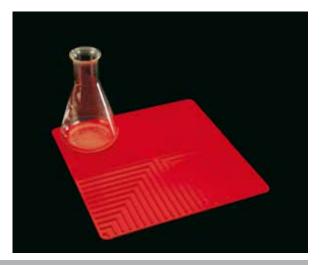

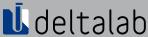



Medical grade silicone tubing. Non toxic

Medical grade silicone. For use in a variety of medical, pharmaceutical, industrial and clinical applications. Translucent tubing. Very high mechanical and chemical resistance. Silicone tubing has a chemical composition similar to quartz and glass and share their excellent properties.

Autoclavable.

Technical features:

Specific gravity: 1.2 + 0.1 Toughness: 40 - 80 Sh

Resistance to traction: 77-95 kg/cm²

Extension: 275-780%

Resistance to tearing: 20-50 kg/cm²

Temperature range: -50 °C to 200 °C (beaks of 250 °C)

		·		
code	dimensions mm Ø int. x Ø ext.	bag quantity (meters)	bag weight	bag volume
350051	0.5 x 1	1 x 100	0.08	0.0011
350115	1 x 1.5	1 x 100	0.21	0.004
350013	1 x 3	1 x 15	0.12	0.001
350024	2 x 4	1 x 15	0.21	0.021
350034	3 x 4	1 x 15	0.21	0.004
350035	3 x 5	1 x 15	0.25	0.002
350037	3 x 7	1 x 15	0.49	0.003
350046	4 x 6	1 x 15	0.50	0.002
350048	4 x 8	1 x 15	0.77	0.004
350057	5 x 7	1 x 15	0.34	0.003
350059	5 x 9	1 x 15	0.66	0.004
350069	6 x 9	1 x 15	0.61	0.01
350610	6 x 10	1 x 15	0.82	0.009
350612	6 x 12	1 x 15	0.21	0.004
350079	7 x 9	1 x 15	0.57	0.009
350710	7 x 10	1 x 15	0.69	0.007
350810	8 x 10	1 x 15	0.62	0.009
350812	8 x 12	1 x 15	1.09	0.005
350814	8 x 14	1 x 15	1.84	0.0007
350912	9 x 12	1 x 15	0.81	0.012
350915	9 x 15	1 x 15	0.21	0.004
351014	10 x 14	1 x 15	1.39	0.008
351420	14 x 20	1 x 15	2.72	0.024
351622	16 x 22	30	3.99	0.064
352030	20 x 30	30	14.00	0.098

 $\textbf{NOTE} : \text{Tube thickness} = \underbrace{\text{external } \varnothing \text{ - internal } \varnothing}_{}$




Tubing clamps

Material: autoclavable acetal.

These clamps will not distort, rust or corrode. Different positions to control flow from full open to complete closure.

code	description	bag quantity	bag weight	bag volume
19505	for tubes of 3.2 to 11 mm (15 positions)	12	0.06	0.0005
19504 for tubes of 11 to 19 mm (12 positions)		6	0.09	0.0007

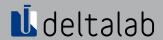

«T» and «L» connectors

Material: polypropylene. **Autoclavable**.

code	fits tubing Ø mm	Ø mm* valley/crest	Ø mm interior mouth	shape	bag quantity	bag weight	bag volume
19460	6	4.5 / 5.4	2.8	Т	20	0.01	0.0001
19461	8	6.9 / 7.6	4.4	Т	20	0.23	0.0001
19462	10	8.7 / 9.5	6.3	Т	20	0.54	0.001

^{*} See the scheme






«Y» and «+» connectors

Material: polypropylene. Autoclavable.

code	fits tubing Ø mm	Ø mm* valley/crest	Ø mm interior mouth	shape	bag quantity	bag weight	bag volume
191411	6	4.8 / 5.6	2.7	+	20	0.20	0.0001
191412	8	6.9 / 7.4	4.6	+	20	0.30	0.0001
191414	12	10.4 / 11.5	8.0	+	20	0.08	0.0006
19465	6	4.7 / 5.7	2.7	Υ	20	0.11	0.0001
19467	10	8.9 / 9.9	5.5	Υ	20	0.03	0.00002
19468	12	10.5 / 11.6	7.3	Υ	20	0.08	0.00003

^{*} See the scheme



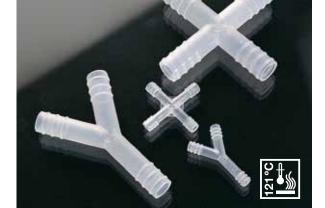

«Straight» and 3 way connectors

 $\label{eq:material:$

code	fits tubing Ø mm	Ø mm* valley/crest	Ø mm interior mouth	shape	bag quantity	bag weight	bag volume
19510	4-5-6	3/4.5/5.5	2.0	straight line	20	0.16	0.00002
19512	9-8-10	6.3/9/11	4.5	straight line	20	0.06	0.0003
19513	10-11-12	8.3/10.5/13	6.5	straight line	20	0.73	0.001
19530	10-11-12	9.9/11.5	5.5	three branches	20	0.13	0.001

^{*} See the scheme








Laboratory holders

Ideal for holding tubing or thermometers onto the edge of beakers or other vessels. **Autoclavable**.


code	description	bag quantity	bag weight	bag volume
19971	for tubing of up to 9.5 mm Ø Polycarbonate	6 x 6	0.13	0.001
19972	for thermometers from 3.2 mmto 11 mm Ø Polypropylene	12 x 3	0.16	0.001

Minimum order quantity: 19971 - 6 units

19972 - 3 units

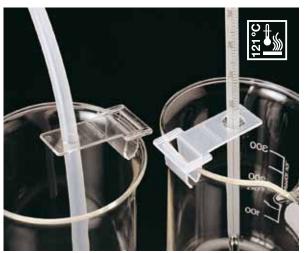



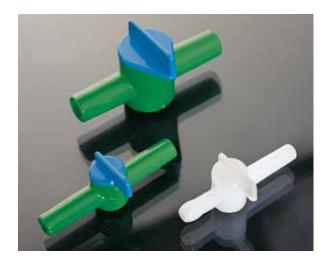

Stopcock, straight bored

Ideal for liquids and gases under low pressure or low vacuum. Material: polypropylene and high density polyethylene that gives a smooth turning action requiring no lubricants.

All are straight bored.

code	fits tubing Ø mm	Ø mm interior mouth	bag quantity	bag weight	bag volume
19374	8	4.9	10	0.04	0.0005
19371	10	7.0	10	0.06	0.0006
19379	16	12.7	10	0.13	0.001




Composed Disconnectors / connectors

Material: polyethylene. Very useful for connecting rubber or plastic tubing of different diameter. Composed by two parts.

Connectors fit tightly together.

code	fits tubing Ø mm	Ø mm crest max/min	Ø mm interior mouth	bag quantity	bag weight	bag volume
19434	04-05-06	3.5 / 6.0	1.6	20	0.49	0.001
19435	06-07-08	5.3 / 7.6	2.7	20	0.61	0.001
19436	08-09-10	7.2 / 9.7	3.6	20	0.83	0.001
19535	10-11-12	9.3 / 12.2	5.5	20	1.06	0.001



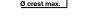



Material: polypropylene. **Autoclavable**.

Maximum work pressure: 2 BAR (2 kg/cm²). Minimum work pressure: 0,07 BAR, at 20 °C.

code	Ø ext	Ø mm valley	bag	bag	bag
	valves	max/min	quantity	weight	volume
19418	8-9-10	7.8 / 10.1	10	0.04	0.0002





Flow indicator (1)

Material: SAN. Perfectly transparent.

The rotation of the red ball allows an immediate visual indication of glass or liquid flow. Very useful for connecting tubes of different diameter from 6.5 to 10 mm.

code	lenght	Ø	bag	bag	bag
	mm	mm	quantity	weight	volume
19538	102	57	5	0.07	0.002



Flow indicator-paddle wheel (2)

Material: SAN.

Not for use with corrosive fluids or temperatures over 85 °C.

code	dimensions	tubing	bag	bag	bag
	mm	mm	quantity	weight	volume
19539	8 x 17 x 40	6-11	5	0.08	0.0004





Wash bottles

Made a combination of high and low density polyethylene.

Transluscent bottle, with HDPE cap (blue), delivery tube and seal cap (natural). Bottle caps are ribbed for a better handling.

High resistance to all liquids, acids and bases.

Include a large neck for safe and easy filling.

code	capac. ml	neck mm	base mm	height mm	case quantity	case weight	case volume
44330000	250	30	65	155	90	5.30	0.078
44300000	500	30	75	185	60	4.85	0.080
44310000	1,000	30	90	222	45	5.12	0.106



Wide mouth wash bottles

Transluscent low density polyethylene bottles, with screw cap and delivery tube available in 3 different colours (blue, yellow and red) for easy identification. Feature a one piece cap and delivery tube assembly, thus eliminating air or liquid absorption through the cap. Include a wide neck for safe and easy filling. Ribbed caps for a better handling while opening and closing.

Moulded-in graduations.

250 ml model is graduated in 25 ml increments, 500 ml model, in 100 ml increments.

code	capac. ml	cap colour	int. neck Ø mm	base Ø mm	height mm	case quantity	case weight	case volume
191637.04	250	blue	38	60	140	10	0.47	0.009
191637.06	250	yellow	38	60	140	10	0.47	0.009
191637.10	250	red	38	60	140	10	0.47	0.009
191638.04	500	blue	38	75	165	10	0.97	0.019
191638.06	500	yellow	38	75	165	10	0.97	0.019
191638.10	500	red	38	75	165	10	0.97	0.019



Integral wash bottles

Made of transluscent polyethylene. Feature a one piece cap and delivery tube assembly, which eliminates air or liquid absorption through the cap and ensures a positive leaktight seal. To increase liquid flow, cut the tip of the delivery tube; close the tip by means of the cap attached to the delivery tube. Caps are ribbed for a better handling while opening and closing.

Tagging points on neck and cap for a tamper proof closure. Bottles include a fill in line.

code	capac. ml	height mm	neck d/n std.	case quantity	case weight	case volume
191633	250	143	GL32	10	0.60	0.008
191634	500	181	GL32	10	0.83	0.018







Special materials: PTFE

This section details special products that are used in critical situations where glass is not suitable. PTFE is a registered brand which refers to a type of polyethylene where hydrogen atoms have been replaced by fluorine ones, commonly called «fluocarbons».

The main characteristic of these chemical compounds is their high resistance to nearly all chemical components including nitric acid. Some variants of PTFE are: FEP, PFA, etc.

PTFE main features:

- Unbreakable, resistant to almost all temperatures, resistant to almost all chemical products, except fluorine and other alkaline metals during fusion
- Non flammable
- These PTFE products are highly recommended for those cases where extreme conditions of transport and storage are required



Low form graduated beakers

Material: PFA.

High chemical and temperature resistance (-250 °C / 270 °C).

code	capacity ml	div. ml	height ml	Ø mm	quantity	weight	volume
191537	250	50 / 1	96	67	1	0.04	0.001



Wash bottles high resistance

Material: FEP.

High chemical and temperature resistance (-250 °C / 270 °C).

code	capacity ml	height mm	Ø mm	quantity	weight	volume
191693	250	157	61	1	0.06	0.001
191694	500	189	76	1	0.09	0.00005




Beakers

Material: PTFE. No graduation.

High temperature and chemical resistance (–200 °C / 270 °C).

code	capac. ml	Ø body mm	height mm	quantity	weight	volume
19718	100	54	68	1	0.03	0.001



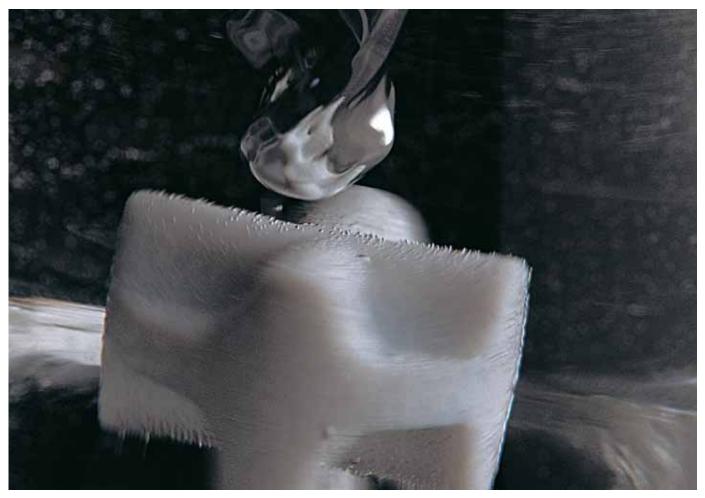

Magnetic stirring bars

This section details our vast range of PTFE magnetic stirring bars. A variety of different sizes and colours is available.

Micro stirring bars: for use in the smallest containers.

Cylindrical stirring bars: these are the most common and lend themselves to many general applications. The diameter of the bar chosen will depend on the diameter of the flask used although it is recommended to select the largest stirring bar possible to achieve optimum results.

Octagonal stirring bars with pivot ring: ideal where the base of the container may be slightly curved or irregular.

Triangular stirring bars: these are used primarily to help dissolve solids or sediments. Also used for higher viscous liquids.

Ovoid stirring bars: used in rounded bottom or Erlenmeyer flasks.

Special shape stirring bars: cruciform stirring bars are used in round bottom tubes.



Stirring bar retriever

Polypropylene coated.

For the removal of stirrer magnets from vessels of all kinds. Includes a magnetic tip, so it is easy to remove the magnets.

code	diameter mm	lenght mm	bag quantity	bag weight	bag volume
19919	10	300	1	0.03	0.00001








Micro stirring bars

code	diameter mm	lenght mm	quantity	weight	volume
19851	3	3	1	0.04	0.0001
19852	2	5	1	0.04	0.0001
19853	3	6	1	0.04	0.0001
19749	2	7	1	0.04	0.0001
19854	1.5	8	1	0.04	0.0001
19855	3	8	1	0.04	0.0001
19856	3	10	1	0.04	0.0001
19857	3	13	1	0.04	0.0001

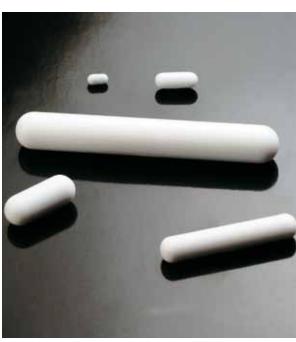

Cylindrical stirring bars


code	Ø mm	lenght mm	quantity	weight	volume
19753	3	6	1	0.04	0.0001
19748	3	8	1	0.04	0.0001
19750	4.5	12	1	0.04	0.0001
19751	6	15	1	0.04	0.0001
19756	8	20	1	0.06	0.0001
19757	8	25	1	0.06	0.0001
19758	8	30	1	0.06	0.0001
19763	6	35	1	0.06	0.0001
19764	8	40	1	0.06	0.0001
19698	8	50	1	0.06	0.0001
19705	10	70	1	0.08	0.0001



Triangular stirring bars

code	dimenions mm	quantity	weight	volume
19782	9 x 35	1	0.08	0.0001
19783	12 x 50	1	0.08	0.0001



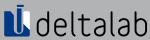


Octahedral stirring bars with pivot ring

code	diameter mm	lenght mm	quantity	weight	volume
19809	8	13	1	0.04	0.0001
19786	8	15	1	0.04	0.0001
19810	8	22	1	0.06	0.0001
19787	8	25	1	0.06	0.0001
19793	8	28	1	0.06	0.0001
19788	8	38	1	0.06	0.0001
19794	8	51	1	0.06	0.0001
19797	10	51	1	0.08	0.0001
19811	13	75	1	0.08	0.0001



Colored octahedral stirring bars

code	colour	diameter mm	lenght mm	quantity	weight	volume
19787.4	blue	8	25	1	0.06	0.0001
19788.4	blue	8	38	1	0.06	0.0001
19788.10	red	8	38	1	0.06	0.0001

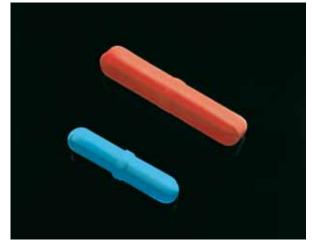

Ovoid stirring bars

code	dimensions mm	quantity	weight	volume
19806	50 x 20	1	0.004	0.0002






Star-shapped stirring bars

code	diameter mm	lenght mm	quantity	weight	volume
19784	10	8	1	0.08	0.0001



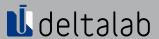



Cylindrical or octahedral stirring bar assortment pack

Assorted cylindrical or octahedral stirring bars supplied in a convenient transparent plastic box including dividers for magnet classification. Box includes 2 units of the following stirring bars:

Code 19400: (diameter x lenght in mm): 8x13, 8x15, 10x25, 10x38, 10x51 and 10x64

Code 19400.1: (diameter x lenght in mm) 6x10, 4,5x15, 6x20, 6x25, 6x30, 8x40, 8x50, 10x60, 10x70.

code	description	quantity	weight	volume
19400	box with 12 octahedrals stirring bars	1	0.25	0.0005
19400.1	box with 18 cylindrical stirrings bars	1	0.25	0.0005



One-hand timer

60 Minute times model "Onehand" non-slip surface. Can be rewound with one hand only.

	code	description	Ø x h mm	quantity	weight	volume
NEW	19901	one-hand timer	72.5 x 29.5	1	0.22	0.0002



Digital timers

Digital timer with hours, days and months, 24 h. programme facility, 60 seconds, alarm, includes clip and small bench stand.

These digital timer count forward and down. Both include a LR44 battery. Code **900400** includes instructions in English and German an code **900600** includes instructions in Spanish, English and French.

mod.	code	description	dimensions mm	quantity	weight	volume
1	900600	digital timer with 2 memories	64 x 61 x 14	1	0.06	0.0003
2	900400	digital timer with 1 memory	52 x 67 x 12	1	0.04	0.0002











Polyethylene wide mouth jars with lid and insert plug

High density polyethylene jars with screw cap and insert plug. Jars and plugs are natural translucent colour; lids are black. Supplied uncapped (except sterile models, supplied capped and bagged).

Upon request, jars can be supplied individually wrapped for a minimum order of 3 cases by model, and sterile for a minimum order of 6 cases by model.

Ask for minimum quantity and delivery time for caps and jars in other colours.

	code	capacity ml	sterile, individual wrap	no sterile, individual wrap	Ø internal mouth mm	Ø body mm	unitary weight g	height with lid mm	case quantity	case weight	case volume	pallet quantity
	202811	30	_	yes	25	32	4.2	53	100	1.00	0.011	14.700
	202810	60	_	yes	38	52	8.2	48	200	3.30	0.046	8.000
	202808	90	_	yes	51	60	8	51	150	3.80	0.046	4.800
	202809	125	_	yes	51	60	12	62	200	6.00	0.082	4.000
	202801	170	_	yes	51	60	14	80	160	6.45	0.082	3.200
	202802	250	_	yes	54	67	20	100	125	5.25	0.078	2.500
	202802.1	250	yes	_	54	67	20	100	125	5.50	0.078	2.500
NEW	202802B	250	_	yes	54	67	20	100	125	5.50	0.078	2.500
	202814	400	_	yes	60	74	27.3	124	130	8.00	0.140	1.560
	202814.1	400	yes	_	60	74	27.3	124	130	8.00	0.140	1.560
	202803	500	_	yes	67	80	35	131	120	8.40	0.140	1.920
	202823.1	500	yes	_	67	80	35	131	120	8.10	0.140	1.920
NEW	202821.1	500	_	yes	86	103	44.9	93	95	9.20	0.140	1.520









	code	capacity ml	sterile,	no sterile, individual wrap	Ø internal mouth mm	Ø body mm	unitary weight g	height with lid	case	case weight	case volume	pallet quantity
	202821	500	—	ves	86	103	44.9	93	quantity 95	9.20	0.140	1.520
	202821.O	500	yes	_	86	103	44.9	93	95	9.20	0.140	1.520
Nav	202823	500	_	yes	67	80	35	131	120	8.10	0.140	1.920
***************************************	202813	750	_	yes	86	103	48	127	75	7.70	0.140	1.200
	202813.10	750	yes	_	86	103	48	127	75	7.40	0.140	1.200
	202818	1,000	_	yes	86	103	56.7	157	50	7.00	0.140	800
	202828G	1,000	yes	_	86	103	56.7	157	50	7.00	0.140	800
NEW	202828	1,000	_	yes	86	103	56.7	157	50	7.00	0.140	800
	202804	1,250	_	yes	86	103	56.7	173	50	6.10	0.140	800
	202815	1,550	_	yes	86	113	64.6	181	45	6.00	0.140	720
	202815B.O	1,550	yes	_	86	113	64.6	181	45	6.00	0.140	720
	202820	2,000	_	yes	86	103	83.1	283	35	5.80	0.140	560
	202820.O	2,000	yes	_	86	103	83.1	283	35	6.10	0.140	560
	202839*	2,000	_	yes	105	128	73	205	32	4.00	0.140	512





Polyethylene beakers

Volumes available range from 60 ml to 1,000 ml.

All models come without lid, except the 1,000 ml beaker, which is available with or without lid. For use with liquid and solid samples.

Manufactured in a flexible material for easy pouring.

code	volume ml	body Ø mm	mouth Ø mm	height mm	weight body g	case quantity	case weight	case volume
242801	60	40	50	48	5.0	1,100	5.95	0.034
242802	100	40	55	58	6.5	1,000	7.20	0.034
242803	230	55	70	78	15.0	700	10.82	0.042
242804	275	53	72	98	22.7	500	11.66	0.042
242805	500	70	90	108	38.7	200	8.50	0.045
242806	1,000	90	120	132	63.2	120	8.22	0.042
242811*	1,000	90	120	132	65.0	65	5.93	0.042

^{*} With lid. Beaker and lid are supplied in the same bag.

Disposable cups

Disposable sample cups, for general sample collection. Cups are supplied stacked.

code	volume ml	material	case quantity	case weight	case volume
202248	100	PP	50 x 100	7.40	0.080
202317	150	PS	30 x 100	8.32	0.080











Tamper evident buckets, large capacity

 $\label{eq:decomposition} \mbox{Disposable containers made of polypropylene } \textbf{autoclavable}. \mbox{ Include tight fitting and leakproof lids.}$

All models include a sturdy and convenient white plastic handle for easier carrying.

Buckets are supplied with lid apart. Made with materials suitables for alimentary use.

code	capacity 	Ø ext. lid mm	Ø int. mouth mm	height mm	*body weight g	case quantity	case weight	case volume
222802	3	200	184.00	137.50	105.00	80	9.00	0.140
222803	4.5	225	210.00	156.00	134.00	48	9.50	0.140
222804	5.6	225	211.00	194.00	152.00	45	10.00	0.140
222805	10.6	268	251.00	263.00	314.00	20	10.50	0.140

^{*} Weight with no lid and handle.



Jars with lid

Body and lid in white ${\bf autoclavable}\ {\bf polypropylene}.$

Self sealing security lid.

Made with materials suitables for alimentary use.

Code 241014 embodies a plastic handle.

Jars are supplied uncapped.

Cases per pallet: 24, excepting code **241010**, palletized in 40 cases.

code	nominal Ø int. mouth x capacity ml Ø base x height mm		weight g	case quantity	case weight	case volume
241010	250	84.5 x 79.5 x 62	18	150	3,00	0,042
241011	550	107 x 98 x 80	32	250	5,00	0,100
241013	1,000	112.5 x 112 x 137.5	46	150	7,68	0,078
241014	1,560	133.5 x 134 x 140	58	150	4,40	0,130





Buckets with spout

Natural low density polyethylene buckets with a stainless metal handle. Include a spout for easier pouring. Moulded-in raduations. Very sturdy construction.

These buckets may be printed with your name and logo.

code	capacity I	dimensions mm	quantity	weight	volume
19905	9	270 x 200 x 250	1	0.50	0.009
19906	12	300 x 220 x 280	1	0.70	0.013
19907	15	340 x 250 x 310	1	1.00	0.019

The specified data may vary when changing to new models.

Industrial jugs


Made of thick rigid polyethylene. Natural colour. Recommended for industrial use.

code	capacity I	mouth Ø mm	height mm	weight ut. g	case quantity	case weight	case volume
470006	3	150	180	215	5	1.43	0.026
470007	5	190	240	400	5	2.21	0.028



Round jugs with handle

Polypropylene. Autoclavable to 121 °C.

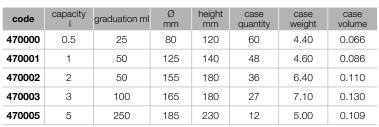

Minimum order quantity: 1.







Self sealing security jars

Excellently finished, body in autoclavable polypropylene.

Screw cap in high density polyethylene with strips. Security sealing with tag for easy opening. Tamper evident system. External moulded ring. Smooth internal side.

Each jar specifies on its bottom its capacity, and the material it is made of. Supplied uncapped.

	code	volume ml	body Ø mm	lid external Ø mm	height mm	weight ut. g	case quantity	case weight	case volume
NEW	202911N	100	54.1	61.5	53	11	150	3.40	0.042
NEW	202912N	250	71.2	76.6	74	25	115	4.80	0.082
NEW	202913N	350	70.5	76.6	101.5	30	80	3,64	0.078
NEW	202914N	500	90.5	100.0	92.1	45	110	6,60	0.140
NEW	202915N	1,000	105.3	116.2	130.2	80	60	7.13	0.140
NEW	202916N	1,250	105.0	116.2	162.5	90	48	5.86	0.130



Square containers with star-shaped caps

Colour: translucent white.

Containers and insert plugs made of high density polyethylene. Star-shaped screw caps made of polypropylene.

Tagging points on neck and cap for a tamper proof closure or for labeling.

All models have moulded-in graduations.

Containers are supplied plugged and capped. Not recommended for liquids transport.

Containors are supplied plugged and supped that recommended for inquide transport.										
code	capacity ml	graduation each ml	dimensions mm	neck Ø mm	body weight g	case quantity	case weight	case volume	pallet quantity	
444608	25	12.5	32 x 34 x 52	18	6	500	6.60	0.090	12,000	
444609	50	10	38 x 38 x 70	23	15	400	7.50	0.074	9,600	
444610	100	20	42 x 48 x 90	33	22	600	16.94	0.170	9,600	
444611	250	50	57 x 60 x 110	33	32	300	13.80	0.170	4,800	
444612	500	100	70 x 80 x 138	44	55	150	11.50	0.170	2,400	
444613	1.000	100	80 x 103 x 176	58	97	80	11.16	0.170	1,280	
444614	2.000	100	100 x 134 x 210	58	177	40	9.16	0.170	640	





Wide neck containers with star-shaped cap

Made of high density polyethylene. Natural colour.

Cylindrical container featuring a star-shaped screw cap.

Tagging points on neck and cap for a tamper proof closure, or for labeling.

Containers are supplied non capped.

code	volume ml	internal mouth Ø mm	body Ø mm	height without lid mm	body weight g	case quantity	case weight	case volume	pallet quantity
44420000	60	23.2	42	80.5	13	200	3.90	0.043	8,000
44430000	125	29.5	52	103.5	18	100	3.18	0.045	4,000
44440000	250	37	65	128	24	120	5.18	0.087	2,880
44450000	500	47.4	82	152	40	60	4.60	0.080	1,440
44460000	1,000	59.4	102	183.2	90	55	7.38	0.140	880
44470000	2,000	59.4	123.5	244.5	128	30	5.00	0.140	480



Wide neck graduated containers with star-shaped cap

Colour: translucent white.

Bottles and insert plugs made of high density polyethylene.

Star-shaped screw caps made of polypropylene.

Tagging points on neck and cap for a tamper proof closure.


Sloping shoulders for easy pouring.

Bottles are supplied plugged and capped.



code	volume ml	graduation ml	internal mouth Ø mm	body Ø mm	height mm	body weight g	case quantity	case weight	case volume	pallet quantity
19408	50	10	24	40	75	13.20	100	1.65	0.0185	9,600
19409	100	20	24	50	90	18.30	100	2.37	0.027	8,400
19410	250	25	34	60	125	35.80	50	2.16	0.037	3,600
19411	500	50	34	75	160	50.70	20	1.20	0.024	1,120
19412	1,000	100	45	95	200	92.20	20	1.95	0.046	960



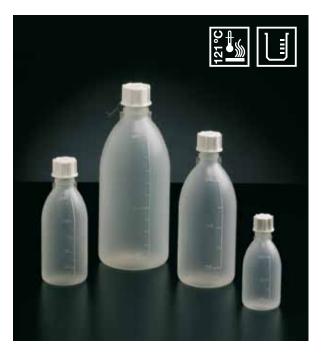

Bottles with standard, narrow neck

Bottles made of **autoclavable**, natural colour polypropylene. Hermetical closure.

Easy opening and closing caps thanks to its size. Gripped, in order to ease its usage with gloves. Cap and body with holes for sealing or labeling. Neck and cap with vents.

Graduation: 100 ml model, every 20 ml; 250 ml, every 25 ml; 500 ml and 1,000 ml, every 100 ml.

code	volume ml	neck DIN	mouth Ø mm	body Ø mm	height mm	body weight g	case quantity	case weight	case volume
191597	100	GL18	13	48	106	18.5	100	1.90	0.040
191599	250	GL25	19	60	148	33.1	5 x 50	10.00	0.180
191601	500	GL25	19	75	180	49.3	150	9.00	0.180
191603	1,000	GL32	23	95	222	100.9	80	8.80	0.180



Bottles with standard, wide neck

Wide neck bottles made of **autoclavable** polypropylene, or polyethylene (not autoclavable).

Body and caps in natural colour. Hermetical closure.

Easy opening and closing caps thanks to its size. Gripped, in order to ease its usage with gloves.

Neck and cap with vents. Graduation: 100 ml model, every 20 ml; 250 ml, every 25 ml; 500 ml, 1,000 ml, and 2,000 ml, every 100 ml.

Made of polyethylene, not autoclavable:

code	volume ml	neck DIN	mouth Ø mm	body Ø mm	height mm	body weight g	case quantity	case weight	case volume
191610	100	GL32	24	48	103	25.7	100	2.78	0.036
191612	250	GL45	38	60	140	39.2	5 x 50	2.06	0.180
191614	500	GL45	38	75	168	55.9	5 x 25	9.30	0.180
191616	1,000	GL63	55	95	206	109.3	68	9.30	0.180
191618	2,000	GL63	55	120	249	145.2	30	6.30	0.180



Made of polypropylene, autoclavable:

code	volume ml	neck DIN	mouth Ø mm	body Ø mm	height mm	body weight g	case quantity	case weight	case volume
191623	100	GL32	24	48	103	26.6	100	2.98	0.040
191625	250	GL45	38	60	140	40.8	250	12.40	0.180
191627	500	GL45	38	75	168	56.5	125	8.50	0.180
191629	1,000	GL63	55	95	206	100.2	68	8.84	0.180
191631	2,000	GL63	55	120	249	175.5	30	6.30	0.180





Graduated narrow neck bottles

Made of transluscent polyethylene.

Designed for the storage of acids and bases. Include a ribbed screw cap for a better handling, and an insert plug.

Feature sloping shoulders for easy pouring.

Tagging points on neck and cap for a tamper proof closure and for labeling. Bottles have the following moulded-in graduations: code 19317, every 10 ml; 19319, every 20 ml; 19323, every 25 ml, 19324, every 50 ml, 19325 and 19326, every 100 ml.

Supplied capped.

code	volume ml	neck Ø mm	bottom Ø mm	height mm	case quantity	case weight	case volume
19317	50	18	40	75	4 x 100	6.70	0.073
19319	125	18	50	110	6 x 100	15.50	0.151
19323	250	22	60	145	6 x 50	12.80	0.151
19324	500	22	75	175	6 x 25	9.00	0.195
19325	1,000	33	95	220	4 x 20	9.00	0.184
19326	2,000	33	120	270	2 x 20	8.40	0.151



Graduated narrow neck opaque bottles

Opaque cylindrical bottles made of dark grey polyethylene, designed for the storage of light sensitive liquids or solid substances. Include a black screw cap and an insert plug for a positive leakproof seal.

Ribbed caps for a better handling.

Heavy duty bottles, with sloping shoulders for easy pouring.

Tagging points on neck and cap for a tamper proof closure. Body and cap with holes for sealing or labeling.

Supplied capped.

code	volume ml	neck Ø mm	bottom Ø mm	height mm	case quantity	case weight	case volume
193241	500	22	75	175	6 x 25	9.00	0.195
193251	1,000	33	95	220	4 x 20	5.00	0.192

Minimum order quantity: 1 bag.





Narrow neck, high resistance bottles

Bottles made of translucent, white, high density polyethylene, designed for the storage of acids and bases.

Blue screw cap on a combination of high and low density polyethylene. Screw caps are ribbed for a better handling.

Sloping shoulders for easy pouring. Heavy duty bottles.

Supplied uncapped.

code	volume ml	neck Ø mm	bottom Ø mm	height mm	case quantity	case weight	case volume
44260000	250	30	65	155	90	4.90	0.076
44270000	500	30	75	185	70	4.86	0.076
44280000	1,000	30	90	222	50	6.70	0.120















Tamper evident, polyethylene bottles

Translucent cylindrical bottles made of high density polyethylene with PP screw cap and LDPE insert plug, designed for the storage of alcohol. Caps are ribbed for a better handling while working with gloves. Insert plugs are dispensed pierced. Tamper evident cap.

No graduation. Supplied uncapped.

code	capac. ml	int. neck Ø mm	bottom Ø mm	height mm	body weight g	case quantity	case weight	case volume	pallet quantity
292300	75	19	39.4	97.6	8	220	2.92	0.042	6,160
292303	125	19	45	106	10	140	2.44	0.043	5,600
292304	250	19	58	128	18	80	2.10	0.044	3,200
292305	500	19	73	160	31	100	4.30	0.110	1,200
292306	1,000	19	82	230	46	60	5.00	0.110	720



Tamper evident bottles, tall models

Translucent cylindrical bottles made of natural colour low density polyethylene, with blue tamper evident closures.

Caps are ribbed for a better handling while working with gloves.

Bottles specially designed with tamper evident rings to ensure a positive seal. Bottles are supplied uncapped.

To fully seal them, screw cap in a normal way and then screw again to fully position the tamper evident ring. To open, unscrew cap.

The ring, automatically broken, will remain on the neck of the bottle.

code	capac. ml	body Ø mm	neck Ø mm	height mm	body weight g	case quantity	case weight	case volume	pallet quantity
282303	125	48	37	105	16.2	125	3.10	0.045	5,000
282304	250	59	37	129	21.9	75	2.48	0.045	3,000
282305	500	66	37	200	35.0	100	5.00	0.110	1,600
282306	1,000	87	37	230	52.0	60	4.10	0.140	960



Screw cap containers with seal and dropper

Made of white polypropylene.

White watertight screw closing system with dropper-plug and blue security seal. Hermetical closure. Mouth diameter: 10.2 mm. Supplied uncapped.

code	volume ml	Ø base mm	height mm	body weight g	case quantity	case weight	pallet quantity
53020	20	24.6	65	3.50	2,000	12.70	32,000
53125	125	42.2	134.7	15.50	275	5.88	4,400
53250	250	50.8	167.3	21.50	200	5.95	3,200
53500	500	68.3	196.5	41.00	100	5.50	1,600

Case volume: 0.114 m³





Rectangular bottles

Body and cap made of polyethylene.

Bottles with blue screw caps, including seal ring and internal joint.


Tamper evident system.

Body in natural colour, graduated (moulded):

250 ml model, graduated each 50 ml.

500 and 1,000 ml models, graduated each 100 ml.

Supplied uncapped.

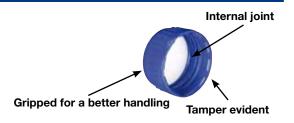

code	volume ml	internal mouth Ø mm	dimensions (A x B x C) mm	weight g	case quantity	case weight	case volume	pallet quantity
292819G	250	30.3	58 x 58 x 115	21	210	6.0	0.096	3,360
292820G	500	30.3	77 x 77 x 130	35	153	9.7	0.140	2,448
292821G	1,000	30.3	96 x 75 x 184	62	85	6.4	0.140	1,360

Nozzle bottles

Cylindrical bottles with screw caps, made of high density polyethylene, with a narrow nozzle design. Natural colour.

Cut the nozzle tip with scissors before liquid dispensing.

No graduation.

Bottles and caps are supplied in separated bags.

code	volume ml	neck Ø mm	bottom Ø mm	height mm	case quantity	case weight	case volume
272301	50	20	35	75	160	2.10	0.026
272303	125	20	48	100	135	2.52	0.045
272304	250	20	59	143	75	2.24	0.047









Squared bottles

Bottles with screw cap and insert plug, made of polyethylene. Bottles and insert plugs are white (transluscent), caps are red. Ribbed caps for a better handling when using gloves. Supplied capped.

code	volume ml	mouth Ø mm	dimensions mm	weight g	case quantity	case weight	case volume	pallet quantity
292818	150	37	50.5 x 50.5 x 88	15	300	8.40	0.096	4,800
292819	250	37	58 x 58 x 103	21	210	6.00	0.096	3,360
292820	500	50.5	77 x 77 x 116	35	165	8.18	0.140	2,640
292821	1,000	50.5	96 x 75 x 172	62	90	7.00	0.140	1,440



Polyethylene bottles


Rectangular body and cap made of polyethylene. The smaller sides are gripped for a better handling. The bigger sides are flat for labelling. Self-sealing blue cap with internal security ring. Hermetical closure. Tamper evident system. Supplied uncapped.

code	volume ml	dimensions mm	int. mouth Ø mm	body weight g	case quantity	case weight	case volume	pallet quantity
282310	500	83 x 65 x 135	28	40	123	5.55	0.11	2,460
282311	1,000	83 x 65 x 235	28	61	73	4.93	0.11	1,460





See sterile models on Chapter 1. **Microbiology**

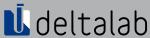

Drum with screen and childproof cap

Made of white colour high density polyethylene.

They embody a handle and a graduated screen which allows the user to see the content.


Childproof, tamper evident, ribbed screw cap with insert plug (see how to use a childproof cap in page 164).

Made with materials suitables for alimentary use.

Drums are supplied non screwed.

Dimensions: 133,8 x 192,5 x 391 mm

code	capacity I	total cap.	Ø int. mouth mm	body weight g	case quantity	case weight	case volume	pallet quantity
202505	5	5,260	34.5	170	16	3.25	0.14	256

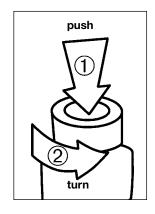



Bottle with screen and childproof cap

Bottle with screen and childproof cap.

Made of white colour high density polyethylene.

It embodies a graduated screen (up to 900 ml) which allows the user to see the content.

Childproof, tamper evident, ribbed screw cap with insert plug (see how to use a childproof cap in page 164).

Made with materials suitables for alimentary use.

Bottles are supplied non screwed.

	code	volume ml	total volume ml	Ø x height mm	body weight g	case quantity	case weight	case volume	pallet quantity
2	02502	1,000	1,120	87.5x226.5	65	74	5.4	0.14	1,184



Cylinder with handle


Made of polypropylene.


The handle provides a secure, solid grip for superior control, specially when handling large volumes.

Pouring is safer and more convenient, as the cylinder is specially designed to ensure a total liquid delivering. 89 x 29 mm hand opening for easy handling. Hexagonal base to prevent tipping and rolling.

Accurate and easy to read graduated double scale on the cylinder. Autoclavable at 121 °C.

code	volume I	graduation ml	height mm	bottom Ø mm	Ø cylinder int. mm	Ø cylinder ext. mm	quantity	weight	volume
191578	1	10	440	146	59	64	1	0.62	0.020
191579	2	20	490	161	81	87	1	0.80	0.020

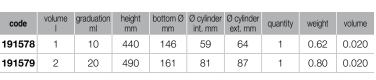

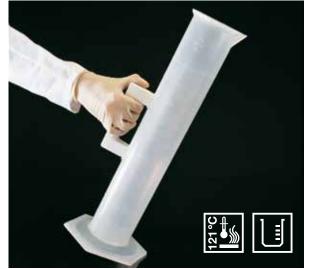





Stackable drums with tamper evident caps

Made of natural colour high density polyethylene. They embody a **handle** and a **tamper evident**, ribbed screw red **cap**, made of HDPE, with internal joint and internal security ring.

Hermetical closure when the cap is totally screwed.

ONU ADR compliant. Drums and caps are supplied non screwed.

code	capacity I	total capacity ml	minimum wall tickness mm	dimensions mm	quantity	weight	volume	pallet quantity
206020	20	20,900	1.20	242 x 298 x 365	1	1.00	0.028	48
206025	25	25,550	1.20	242 x 298 x 443	1	1.20	0.028	48



Wide mouth drums with security cap

Made of high density polyethylene.

Translucent natural color.

Self-sealing green screw cap in polyethylene, with white insert plug. It features a black folding handle.

code	capacity I	Ø external body mm	Ø internal mouth mm	height mm	quantity	weight	volume	pallet quantity
202616	5	169	105	288	1	0.22	0.0065	256





Stackable drum with tamper evident screw cap

Made of high density polyethylene.

It is provided with a sealing screw cap with grip.

This drum has four flat sides suitable for sticking information labels.

They are distinguished by two main characteristics: they can be easily stacked joining the tops and the bottoms together, and they can be palletised thanks to their reduced overall dimensions.

All of them are ONU compliant, meeting all the standards for the transport of hazardous, toxic and noxious products.

code	nominal capacity I	total capacity I	dimensions mm	Ø mouth mm	DIN	quantity	weight	volume	pallet quantity
202605	5	5.5	198 x 148 x 240	44	51	1	0.260	0.007	200
202610	10	11.6	224 x 193 x 305	44	51	1	0.450	0.017	100
202620	20	21.7	245 x 295 x 375	44	51	1	0.900	0.027	48
202625	25	26.9	245 x 295 x 435	44	51	1	1.100	0.030	48
202630	30	34.5	332 x 290 x 428	49	51	1	1.300	0.040	36

Storage bottles

Bottles and caps made of high density polyethylene.

Wide neck bottle (natural colour, translucent) with screw cap (black and blue). It features a double cap: it includes a blue internal cap that works as a plug. Ribbed cap for a better handling.


Bottles include a sturdy aluminium handle for easy carrying.

Not suitable to transport liquids.

code	capacity I	mouth Ø mm	body Ø mm	height mm	unitary weight g	case quantity	case weight	case volume	pallet quantity
202825	5	100	175	275	190	14	4.9	0.14	224
202826	10	100	215	340	380	6	4.5	0.14	96





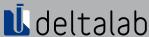



Drums with tap

Made of high density polyeythylene.

Heavy duty drums, wall thickness 2 mm.

Can be used at temperatures below 0 °C.


The tap includes a seal ring to ensure a leakproof seal.

Each drum includes a screw cap and a plug for added security. With two handles for easy transport.

Drums are supplied complete with tap (also available spare tap).

	code	capacity I	Ø neck mm	Ø base mm	height mm	quantity	weight	volume
	19660	5	52.5	165	345	1	0.40	0.008
	19662	10	52.5	210	425	1	0.80	0.012
NEW	19664*	25	79.5	280	565	2	1.4	0.042
NEW	19666*	50	79.5	350	700	2	2.1	0.085
	Spare tap:							
	19375	_	-	_	-	1	0.20	0.001

^{*}Minimum order quantity: 1.

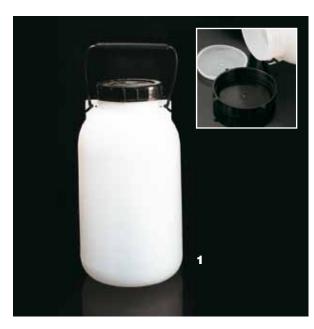

Wide and narrow neck drums

Heavy duty drums made of high density polyeythylene.

Include an insert plug for a positive leakproof seal, and one handle for easy carrying. Tagging points on neck and cap for a tamper proof closure or for labelling.

Wide mouth drums are designed both for liquids and powders, easy to fill and pour.

Wide neck bottles (1)

code	capacity I	Ø neck mm	Ø body mm	height mm	quantity	weight	volume	pallet quantity
191640	5	90	165	305	1	0.45	0.008	200
191642	10	120	210	395	1	0.95	0.017	90

Narrow neck bottles (2)

code	capacity I	Ø neck mm	Ø body mm	height mm	quantity	weight	volume	pallet quantity
191646	10	53	210	420	1	0.95	0.018	90









Wide mouth storage drums with tap and plug

Body and cap made of polyethylene.

Screw cap, green colour, with non hermetical plug.

They embody a tap, and two black folding handles at the top.

Not designed to transport liquids. External mouth diameter: 150 mm.

	code	capacity I	Ø base mm	height mm	quantity	weight kg	volume
Naw	209030	30	345	470	1	1.51	0.059
NEW	209050	50	345	675	1	2.35	0.084

Funnels

Material: polypropylene.

Autoclavable.

code	stem internal Ø mm	stem external Ø mm	bouche mm	quantity	weight	volume
194171	5.7	8.1	80	1	0.027	0.0001
194172	9.1	11.7	100	1	0.022	0.0005
194173	9.6	13.0	120	1	0.032	0.0005
194175	13.9	18.2	180	1	0.080	0.0012
194176	18.5	21.8	220	1	0.160	0.0028
194177	21.7	26.4	260	1	0.20	0.0081
194178	24.7	29.7	300	1	0.34	0.0110









Industrial drums

These drums are produced in high density polyethylene Their particular top, which has a **perfectly stiffened rim obtained** during the moulding phase, prevents the opening from becoming deformed and the lid from becoming detached as a consequence of any drops or violent shocks. They feature a metal washer that ensures a perfect closure where is also suitable to place a tamper evident seal. Thanks to their **airtight seal**, they are **suitable for transport by sea, road or railway** and, thanks to their construction characteristics and qualities; they **fully meet the current ONU ADR standards** for **transport of hazardous toxic and noxious products.** Incinerable (excluding the metal closure).

These drums are produced using selected top-quality raw materials, they feature excellent drop-impact resistance even at low temperature (-20 °C) and stability against corrosion from chemical agents which may cause stress cracking.

Their **internal surface** is **perfectly smooth**, and this makes it easy to fully empty and completely clean the inside from any powder, paste or semi-dense products. Codes **208030** and **208060** have two handles for easy carrying.

code	capacity I	body Ø mm	mouth Ø mm	height mm	quantity	weight kg	pallet quantity
208030	30	320	252	502	1	2.0	48
208060	60	405	327	610	1	3.0	27
208128	120	480	400	800	1	5.3	18



Large deep trays

Made of natural high density polyethylene. Stackable. Include ergonomically designed handles.

code	dimensions mm	height mm	tray volume l	quantity	weight kg
19603	550 x 430	310	50	1	1.80
19604	630 x 440	340	65	1	2.20
19605	720 x 470	360	90	1	4.00






The technical data (dimensions, capacities and resistances) specified in this catalogue are approximated.

The resistances (to chemical products, temperature, centrifugation...) may suffer variations depending in external conditions.

We recommend to test the product under the forecasted conditions of use. Please do not hesitate to ask for samples if necessary.



Description	page	Description	page
		Caps for screw threaded glass tubes	134
A		Caps, general features	126
ABS trays	265	Caps, polyethylene	125-127
Absorbent paper	250	Cardboard waste containers	243
Acetal polymer racks	216	Cardboard storage box	204, 237
Accessories, WILMUT (WILMUT system)	233	Cary Blair. Swab	24
Adhesive labels	246	Cassettes	89-93
Adhesive sealing film for multiwell plates	62	Cassettes, sponge	88
Adhesive tapes for sample identification	246	Cassettes, storage	94-95
Alcohol pads	245	Centrifugation microtubes	136-141, 143-145
Aluminium sealing foil	62	Chlamydia, swab	25
Aluminium racks	221	Cobas Mira cuvette segments	82
Amies. Swab	23	Coliforms, water control	48
Analytical funnels	257	Composed disconnectors / connectors	271
Anticoagulants and preservatives in containers	79	Concentrated deodorant Nilodor	251
Autoclave resistant bags	114, 244	Conical measuring flasks	257
Automatic slide dispenser	106	Conical tubes	34-38, 124, 128,
AYRE, spatula	115		130
		Connectors	270
В		Contact plate rack	29
Bags for autoclave	114, 244	Contact plates	29
Bags	41, 42, 110-114	Container with boric acid	149
Bags, isothermal	116	Container, rectangular	159, 162
Baskets for autoclave	227	Containers for vacuum system	158
Baskets with folding handles	226	Containers with tamper evident	151, 165
Beakers for cell counters	83	Containers	86, 87, 160, 280,
Beakers	273		283, 285
Beakers, graduated	256, 257, 273	Containers, leak proof screw cap	87, 160
Beakers, polyethylene	282	Containers, racks	152, 223, 225
Beckman® microtubes	136	Coplin type staining jar	102
Bed pan	251	Coulter counter cups	82
Bed Scaner for 2D codes W-DATAPAQ 96		Couvette segments, Cobas Mira	82
(WILMUT system)	236	Cover slides	97
Blood group test plate	76	Cryoboxes (WILMUT system)	234
Bottle carrier with folding handlers	226	CRYOINSTANT porous beads	31, 142
Bottle with screen and childproof cap	292	Cryoracks, metal. Space organization	
Bottles for sterilization	244	(WILMUT system)	237
Bottles with standard, narrow neck	287, 288	Cryovials and microtubes, freezing cardboard	
Bottles with standard, wide neck	287	storage box	204
Bottles	45, 116, 244,	Cryovials and microtubes, racks	200-212
	249, 272, 285-	Cryovials	31, 142-145
D" ()	291	Cryovials, acnes	145
Büchner funnels	259	Cryovials, caps	145
Buckets	86, 283, 284	Cryovials, colour coded inserts	144
Burette clamps	260	Cryovials, racks	200-212
		Culture tub cap	132
Comp formials	105	Culture tubes	34, 128
Caps for vials	135	Cups, disposable	157
Canes for cryovials	145	Cups, graduated	164
Cap, culture tub	132	Cuvette rack	212
Caps, plug style (WILMUT system)	233	Cylinder with handle	292
Cappar automatic W CAPPO6	233	Cylinders Cylinderical container with handle	254, 292
Capper, automatic W-CAPP96	005	Cytology by Joh	161
(WILMUT system)	235	Cytology brush	115
Capillary tubes seeling way	235		
Capalfar tubes sealing wax	80 145	December IV/ DECAD (M/II M/I IT avistom)	225
Caps for cryovials		Decapper W-DECAP (WILMUT system)	235 186
Caps for microtubes	50, 51, 140, 141	Daslab tips	100

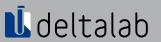

	page
Deltaswab	21
Dessicator discs ce	262
Digralsky spreader. Sterile	33
DNA free swabs, steriles	20, 52
Drum with screen and childproof cap	292
Drums with tap	295
Drums	292-297
DS Barrier, filter tips	186
D3 Barrier, lilter tips	100
E	
E.S.R	74-76
	74-76
E.S.R. citrate tube	
Emergency eyewash bottle	248
Erlenmayer flasks	255
_	
F	
Filter bag, sterile	42
Filter paper	250
Filter tips	170-186
Flask weight ring	263
Flocked swabs	19, 21
Flow indicator	271
Forceps	89, 113
Freezing cardboard storage box for microtubes	
and cryovials	204
Funnel holders	260
Funnels	258-259, 296
Turriois	200 200, 200
G	
Glass dessicator	262
Glass tubes	132-134
	248
Glasses, safety	
Gloves	114, 253, 253
Graduated beakers	256, 257, 273
Graduated cups	164
н	
Haematological stains	80, 99
l la augusta la augusta (agusta (agust	00 00
Haematology, stains	80, 99
Handrop Digital SC	187
Handrop Digital SC	187
Handrop Digital SC Handrop Digital MC Handrop stepper	187 187
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould	187 187 187
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer	187 187 187 88
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for	187 187 187 88 44
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs	187 187 187 88 44 44 247
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for	187 187 187 88 44 44
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles	187 187 187 88 44 44 247
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles	187 187 187 88 44 44 247 20, 52
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles	187 187 187 88 44 44 247 20, 52
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles Imhoff, sedimentation cone Incubation rack	187 187 187 88 44 44 247 20, 52
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles Imhoff, sedimentation cone Incubation rack Industrial drums	187 187 187 88 44 44 247 20, 52 259 208 297
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles Imhoff, sedimentation cone Incubation rack Industrial drums Instrument trays	187 187 187 88 44 44 247 20, 52 259 208 297 226
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles Imhoff, sedimentation cone Incubation rack Industrial drums Instrument trays Isofreeze rack with temperature indicator	187 187 187 88 44 44 247 20, 52 259 208 297 226 207
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles Imhoff, sedimentation cone Incubation rack Industrial drums Instrument trays Isofreeze rack with temperature indicator Isofreeze reversible racks with lid	187 187 187 88 44 44 247 20, 52 259 208 297 226 207 208
Handrop Digital SC Handrop Digital MC Handrop stepper Histology mould Homogenizer Homogenizer, bags for Hospital swabs Human DNA free swabs, steriles Imhoff, sedimentation cone Incubation rack Industrial drums Instrument trays Isofreeze rack with temperature indicator	187 187 187 88 44 44 247 20, 52 259 208 297 226 207

Description	page
J	
Jars for water sampling, sterile	44
Jars with lid	86, 283
Jars, self sealing	285
Jars, wide mouth	86, 87, 280
Jugs	284
K	
Kangaroo bags	112
_	
L	005
Labels W-Labels (WILMUT system)	235
Laboratory holders	270
Laboratory storage cabinet	95
Large deep trays	297
Leak proof screw cap containers	87, 160
Liquid transport modium	89
Liquid transport medium Lithium hengrine tubes	21 72
Lithium heparine tubes Lockers W-LOCK (WILMUT system)	234
Long-handled dippers	46
Long-stem funnels	257
Loops	32-33
Loops	02-00
M	
Macrotips	183-185
Magnetic stirring bars	274-276
Mailing container	76, 113
Manual 2D and 1D bar codes scanner	70, 110
W-READER (WILMUT system)	236
Measuring cylinders	254, 292
Measuring scoops	264
Mega cassettes	93
Metal cryoracks, space organization	
(WILMUT system)	237
Metal tray histology	88
Microincinerator for nichroome loops	32
Microscope slide labels	98
Microtitre plates, sterile	30
Microtubes	50, 51, 136-141
Microtubes, Beckman®	136
Microtubes, centrifugation resistant	136-141, 144,
	145
Microtubes, racks	200-212
Microtubes, safety lock	138
Microtubes, screw thread	50, 51, 140, 141
Microtubes, Vitatron®	136
Mixed racks	206, 210, 211,
	213, 216
Mounting medium for microscopy	98
Multichannel pipette	187
Multiwell sealing mat	64
N N	000
Narrow neck opaque bottles	288
Nichrome loops microincinerator	32
Nilodor, concentrated deodorant	251

A-Z



Description	page	Description	page
Non sterile swabs	17	Round bottom tubes	36, 120, 121, 130
Non-return valve	271	Round jugs with handle	284
		Scaners, bed type and manual	236
P	0.45		
Pads	245	S	0.40.040
Paraffin	88	Safety containers	240-243
Pasteur, glass pipettes	197	Safety glasses	247
Pasteur, polyethylene pipettes	200-193	Safety labeling	240
PCR plates	56-63	Sample containers (PS, PP)	148, 149, 152,
PCR racks	207, 210, 211		153
PCR tubes	53-56	Sample containers	148-158
Petri dish baskets	28, 229	Sample management, software	
Petri dish racks	28	(WILMUT system)	236
Petri dishes	26-30	Sample storage system	64,65
Pickers, W-PICKER and W-MULTIPICKER		Sampling kits: meat or surfaces	110
(WILMUT system)	234	Scintillation vials	88
Pipette and burette rinsing set	266	Screw cap containers with seal and dropper	289
Pipette pump support rack	247	Screw thread microtubes	50, 51, 140, 141
Pipette pumps	40, 195, 247	Self sealing security jars	285
Pipette rack	212	Serological pipettes	39, 194, 196
Pipette stands	261	Serotub, serum glucose	69
Pipettes	189-193	Silicone mat	268
Pipettes, Handrop	187	Silicone tourniquet	68
Pipettes, trays	196, 265	Silicone tubing	268-269
Pipetting bulb	40, 195, 197,	Single use racks	216
	246, 247	Slide labels	98
Pipetting reservoir	182	Slide storage boxes	105, 106
Plastic dessicator	262	Slide tray	103
Polycarbonate racks	200-201	Slidefolder	103
Polyethylene beakers	282	Slides, automatic dispenser	106
Polyethylene bottles	291	Slides, boxes	105, 106
Polyethylene caps	125-127	Slides, Coplin staining jar	102
Polypropylene tubes	123	Slides, cylindrical mailer	105
Polystyrene racks	200	Slides, glass slides	96, 97
Polystyrene tubes	123	Slides, jars and baskets	102, 106, 229
Porous beads CRYOINSTANT	31, 142	Slides, labels	98
Powder funnels	259	Slides, mounting medium	98
PVC antiacid trays	265	Slides, staining racks	102
		Slides, staining	98, 100-102
R		Slides, storage	107
Rack for contact plates	29	Slides, tube mailer	104
Rack for sedimentation cones	259	Sodium thiosulfate Whirl-Pak® bags	46
Rack with temperature indicator	207	Spatula AYRE	115
Racks for bags	229	Spatulas	264
Racks for containers	152, 223, 225	Special techniques	79
Racks for cryovials and microtubes	200-212	Spectrophotometer cuvettes	81
Racks for Petri dishes	28	Sponge for cassettes	88
Racks for tubes	135, 213-224	Spreader, Digralsky	33
Racks, single use	216	Software, sample management	
Re-caps	125	(WILMUT system)	236
Rectangular base for burettes	260	Squared containers	117, 163
Rectangular bottles	284, 290, 291	Stackable drums, tamper evident	293, 294
Rectangular container	159, 162	Stackable Petri dish incubation tray	30
Rectangular waste containers (PP)	243	Staining jar and baskets	100-102
Reinforced silicone tubing	268	Staining racks and jars	100-101
Reticulocyte staining kit	79	Staining racks, slides	102, 229
Ritips®	188	Staining tray	98
Ripette Pro®	188	Stainless steel racks and baskets	220-229
	.00	The state of the s	



Description	page
Stainless steel water still	263
Stains	80, 99
Steel staining basket for blood collection	228
Stepper	187, 188
Sterile bags for STOMACHER® blenders	41
Sterile bottles for water sampling	45
Sterile filter bag	42
Sterile swab in round tube	19
Sterile swabs	16,18
Sterile Whirl-Pak® blender bags	41
Sterile Whirl-Pak® sample bags	111
	245
Sterilization indicator tapes	244
Sterilization, bottles	
Sterilization, indicator tapes	245
Stirring bar assortment pack	277
Stirring bar retriever	274
Stirring bars	274-276
Stopcock, straight bored	271
Storage bottles	294, 295
Storage for laboratory	95
Straw-like pipettes	196
Stuart, swab	24
T Takives, E.S.R. semi micro system	75
Tamper evident containers	151, 165
Techniques, special (Haematology and	
Biochemistry)	79
Test tube baskets	212
Test tube racks, universal	214, 215
Test tube tongs	264
Timer, digital	277
Timer, one-hand	277
Tips for pipettes	170-182
Tips, filter	170-186
Tips, Daslab	186
Tips, DS Barrier	186
Tips, macrotips	183-185
Tongue depressors	250
Transport swabs with media	22-25
Tube, slide tube mailer	104
Tubes Edta dipotassium	72
Tubes Edta tripotassium	72
Tubes for Blood collection	77
Tubes for vacuum system	158
Tubes with clot accelerator and serum	100
	70 71
separator	70, 71
Tubes	34-38, 120-124, 128-134
Tubes, codified 2D datamatrix	
(WILMUT system)	232
Tubes, no codified (WILMUT system)	233
Tubes, conical	34-38, 124, 128, 130
Tubes, Edta dipotassium	73
Tubes, Edta:tripotassium	73
Tubes, glass	132-134
Tubes, lithium heparine	72

Description	page
Tubes, PCR	53-56
Tubes, racks	135, 213-224
Tubes, special pediatric	77
Tubing clamps	269
U	054
Urine collection bottle with lid	251
Urine collection bottle	251 157
Urine container set	107
V	
Vacuum pump	263
Vacuum system, containers and tubes	158
Vaginal speculum	115
Vials, flat bottom screw	135
Vials, scintillation	83
Virus, swab	25
Vitatron® microtubes	136
Volumetric flasks	255
w	
Wash bottles	249, 272, 273
Watch glasses	267
Water control: absence/presence of Coliforms	
(E. Coli)	47
Weighing bottles	267
Weighing dishes	267
Whirl-Pak® bags with flat wires	110
Whirl-Pak® bags	41, 46, 110, 111
NA/In: IN D. L. ® L	113
Whirl-Pak® kangaroo bag	113
Whirl-Pak® stand-up bags Wide mouth drums with tap and plug	110 296
Wide mouth drums, tamper evident	293
Wide neck containers	286, 288-290
W-CAPP10, manual capper (WILMUT system)	235
W-CAPP96, automatic capper (WILMUT system	
W-DECAP, decapper (WILMUT system)	235
W-LABELS, labels (WILMUT system)	235
W-PICKER and W-MULTIPICKER, pickers	
(WILMUT system)	234
W-SEAL, plug style caps (WILMUT system)	233
WILMUT 2D standard (W-2DST)	232
WILMUT 2D Right Position® (W-2DRP)*	232
WILMUT Eco (W-ECO)	233
WILMUT Seroteca (W-SER)	233
WILMUT, 2D datamatrix codified tubes	232
WILMUT, non codified tubes	233
z	
Zip-lock bags	112

A-Z



Code	page	Code	page	Code	page	Code	page
100LG	253	1911108	267	191633	249, 272	19252	196, 265
100LM	253	1911109	267	191634	249, 272	19254	261
100LP	253	191219	266	191637	249, 272	192562	254
100LSP	253	191220	246	191637.04	249, 272	192563	254
102LG	253	191221.1	246	191637.06	249, 272	192564	254
102LM	253	191221.2	246	191637.10	249, 272	192565	254
102LP	253	191221.3	246	191638	249, 272	192566	254
1164	75	191221.4	246	191638.04	249, 272	192567	254
120002	250	191221.5	246	191638.06	249, 272	192571	254
120003	250	191221.6	246	191638.10	249, 272	192572	254
120010	250	191221/1	246	19164	258	192573	254
128030	76	191221/2	246	191640	295	192574	254
1360	75	191221/3	246	191642	295	192575	254
1361	75, 76	191221/4	246	191646	295	192576	254
15001	41	191221/5	246	19165	258	192577	254
15003	41	191221/6	246	19166	259	19262	261
15004	41	191222	245	19167	259	19265	261
15005	42	191223	245	19168	259	19266	260
15006	41, 229	19139	260	19169	259	19267	260
15007	41, 229	191395	263	191693	273	19269	260
15008	42	19140	260	191694	273	19270	253
169001	186	191411	270	19170	259	19274.5	107
169010	186	191412	270	19171	259	19274.7	107
169011	186	191414	270	19172	258	19276	106
169020	186	19145	258	191720	256	19276.1	106
169021	186	19146	258	191721	256	19276.A	105
169100	186	191460	255	191722	256	19276.B	105
169200	186	191461	255	191723	256	19277	106
169300	186	191462	255	191724	256	19277.1	106
191050	245	191463	255	191725	256	19277.A	105
191051	245	191464	255	191726	256	19277.B	105
191052	245	191465	255	191727	256	19278	106
191055	259	19147	258	19173	258	19278.2	106
191056	259	19148	258	19174	258	19278.3	106
191057	264	191484	255	19175	258	19278.A	105
191058	264	191485	255	191822	256	19278.B	105
191059	264	191486	255	191823	256	19279.3	107
191060	264	191487	255	191824	256	19279.5	107
191062	264	191488	255	191825	256	19279.7	107
191063	264	191489	255	191826	256	19280	265
191064	264	19150	258	191827	256	19281	265
191065	264	19152	258	191828	256	19282	265
191066	264	19153	258	191829	256	19283	265
191067	264	191537	273	19194	197	19284	265
191068	264	191578	292	19195	197	19285	265
191069	264	191579	292 287	19200	40, 197, 247	192922.04	88 245
191070 191071	264 264	191597	287	19201	40, 195, 246 265	192928	88
191071	102	191599 191601	287	19212 19213	266	192931 192932	88
1911007	100	191603	287	19215	266	192933	88
191101	100	191610	287	19217	266	192934	88
191104	101	191612	287	19218	266	192935	88
191104	101	191614	287	19219	266	19317	288
191105	101	191616	287	19219	266	19317	288
191107	101	191618	287	19230	262	19319	288
191107	101	_	258			19323	
1911103	267	19162 191623	258	19232 19235	262 262	19324	288 288
		_					
1911104	267	191625	287	192470	263	19325	288
1911105	267	191627	287	192471	263	193251	288
1911106	267	191629	287	192472	263	19326	288
1911107	267	191631	287	192503	82	19340	267



ode	page	Code	page	Code	page	Code	page
9345	267	19513	270	19763	275	200007	191
9351	100	19530	270	19764	275	200008	192
9353	101	19535	271	19782	275	200009	176
9354	101	19538	271	19783	275	200009B	176
9355	100	19539	271	19784	276	200009R	176
9360	102	19540	106	19786	276	200010FR	171
9361	102	19553	262	19787	276	200012	179
9371	271	19554	262	19787.4	276	200012B	179
9374	271	19563	215	19788	276	200012R	179
9375	295	19564	215	19788.10	276	200015	185
9379	271	19565	215	19788.4	276	200015R	185
384	248	19566	215	19793	276	200015RF	185
385	248	19567	215	19794	276	200016	174
386	248	19568	215	19797	276	200016B	174
387	248	19570	215	19806	276	200016R	174
400	277	195701	265	19809	276	200017	178
400.1	277	195702	265	19810	276	200017R	178
408	286	195703	265	19811	276	200018	181
409	286	195705	265	19851	275	200019	178
410	286	195706	265	19852	275	200019 200019R	178
411	286	195707	265	19853	275	2000191	181
411 412	286		265	19854	275		
		195709				200021	181
4171	258, 296	195711	265	19855	275	200022	182
4172	258, 296	195713	265	19856	275	200024	170
4173	258, 296	19575	46	19857	275	200026	175
4175	258, 296	19576	46	1988314	83	200028	184
4176	258, 296	19577	46	1988315	83	200028R	184
4177	258, 296	19580	263	19901	277	200028RF	184
4178	258, 296	19586	251	199014N	174	200029	180
418	271	19591	251	199014RC	174	200029F	180
423	257	19591/T	251	19905	284	200029R	180
424	257	19593	264	19906	284	200029RF	180
425	257	19596	264	19907	284	200032.0	191
434	271	19603	297	199075	183	200034	191
435	271	19604	297	199076RC	177	200037	191
436	271	19605	297	199078N	173	200038	192
437	259	196102	42	199078RC	173	200050	193
438	259	196105	42	199083	177	200060	193
440	259	196110	42	199084RC	170	200062.0	191
441	259	196115	42	19919	274	200070	179
442	259	196302	42	19920	94	200070NR	179
443	259	196305	42	19921	94	200070R	179
445	259	196310	42	19922	105	20007011	174
452	267	196315	42	19923	104	200072	177
454	267	19660	295	19923	104	200074	183
	267	19662	295	19924	104	200075 200075RF	183
455 460							
460	270	19664	295	19930	94	200076	178
461 460	270	19666	295	19931	94	200076R	178
462	270	19698	275	199620	177	200076RF	178
165 10 5	270	19700	226	19970	260	200077	179
167	270	19705	275	19971	270	200077F	179
168	270	19718	273	19972	270	200077FR	179
500	89, 113	19748	275	19973	212	200077R	179
501	89, 113	19749	275	19996	196, 265	200078	173
503	89, 113	19750	275	200000	192	200079	175
504	269	19751	275	200001	192	200079F	175
505	269	19753	275	200002	191	200079FS	175
509	264	19756	275	200003	191	200079R	175
510	270	19757	275	200004	191	200080.1	180
512	270	19758	275	200005	192	200080.1R	180
					191		



Code	page	Code	page	Code	page	Code	page
200082F	180	200397.V	43	202847	87	249010	186
200082NR	180	200400	137	202848	160	249011	186
200082R	180	200400P	137	202911N	285	249020	186
200083	184	200401	137	202912N	285	249021	186
200083R	184	200405	137	202913N	285	249100	186
200100	244	200407	137	202914N	285	249200	186
200102	82	200410	137	202915N	285	249300	186
200103	82	200420	209	202916N	285	250050	243
200121	192	200500	33	202VG	252	260000	242
200200	26	200510	33	202VM	252	27	74
200200.4	26	200VG	252	202VP	252	270045	243
200200.5	27	200VM	252	206020	293	272301	290
200201	27	200VP	252	206025	293	272303	290
200201.B	27	201110	190	208030	297	272304	290
200202	26	201250	190	208060	297	28	76
200202	26	202248	157, 282	208128	297	282300	47
200203	27	202317	157, 282	209030	296	282303	289
200204	29	202502	292	209050	296	282304	289
200208	29	202505	292	210002	193	282305	289
						282305	
200210	26	202605	294	210003	193		289
200214	27	202610	294	210004	193	282310	291
200218	29	202616	293	210005	190	282311	291
200219	27	202620	294	210006	193	282320	45
200310	244	202625	294	210022	193	282321	45
200311	244	202630	294	210023	193	282323	45
200312	244	202801	86, 280	210024	193	282323.BU	45
200318	244	202802	86, 280	210030	190	282324	45
200320	244	202802.1	280	222802	86, 283	282324.BU	45
200321	244	202802B	280	222803	86, 283	282330	45
200322	244	202803	86, 280	222804	86, 283	282331	45
200325	111	202804	281	222805	86, 283	282333	45
200326	111	202808	86, 280	232801	242	282334	45
200329	111	202809	86, 280	232802	242	282340	44
200330	111	202810	86, 280	232806	242	282341	44
200332	111	202811	86, 280	232808	242	282350	44
200333	111	202813	86, 281	232809	242	282351	44
200340	111	202813.10	281	232810	242	29	75
200341	111	202814	86, 280	232811	242	292300	289
200342	41, 111	202814.1	280	240007	240	292303	289
200343	41, 111	202815	281	240011	241	292304	289
200345	111	202815B.O	281	240021	241	292305	289
200346	111	202818	86, 281	240028	241	292306	289
200347	111	202820	281	240035	243	292604	46
200349	111	202820.O	281	240065	243	292818	291
200351	41, 111	202821	86, 281	241010	86, 283	292819	291
200356	110	202821.1	280	241011	86, 283	292819G	290
200357	111	202821.O	281	241013	86, 283	292820	291
200358	111	202823	281	241014	86, 283	292820G	290
200359	111	202823.1	280	241400	240	292821	291
200361	110	202825	294	241500	240	292821G	290
200363	111	202826	294	242000	240	300100	121
200364	111	202828	281	242801	282	300200	18
200365	110	202828G	281	242802	282	300201	18
200372	113	202839	281	242803	282	300202	18
200373	41	202840	160	242804	282	300203	18
200373	41	202841	160	242805	282	300205	17
200374	41	202842	160	242806	282	300203	18
200370 200393.V	43	202843	160	242811	282	300210	18
200393. v 200393P	43			242811		300211.1	
		202844	160		240	_	17
200394.V	43	202845	87	246000	240	300232	17
200396.V	43	202846	87	249001	186	300234	247



Code	page	Code	page	Code	page	Code	page
300236	247	301-10	171	308205	125	327-37RX	171
300237	247	301-11	171	308502	125	327-40	176
300243	17	301200	124	308506	125	330	213
300250	19	301201	124	308602	125	330089	182
300250.1	19	301202	124	308606	125	330129	182
300251	19	301212	124	309202	148	350002	268
300252	19	301213	124	309222	148	350004	268
300252 300252DNA	20, 52	301402	34, 128	309222.O	148	350004	269
300252DNA 300253	19	301403			157	_	269
300253	19	301700	34, 128 121	309316 309316.E	157	350024 350034	269
						_	
300259	19	302000	81	309318	157	350035	269
300260	17	302100	81	309324	157	350037	269
300261	19	302400	81	309324.E	157	350046	269
300261DNA	20, 52	302703	33	309324.E3	157	350048	269
300263	19	302704	33	309402	150	350051	269
300263DNA	20, 52	302723	33	309502	152	350057	269
300265	19	302724	33	309526	152	350059	269
300268	17	302762	32	309552	152	350069	269
300270	17	302764	32	309752	154	350079	269
300280	24	302771	32	309903	156	350115	269
300280.2	24	302772	32	310200	18	350610	269
300281	23	302780	32	310202	18	350612	269
300281/1	23	302792	32	310211.1	18	350613	268
300281T	23	302794	32	310253.1	18	350710	269
300281TC	23	303100	81	316-314.1	188	350810	269
300284	23	303101	81	316-314.2	188	350812	269
300285	23	303102	81	316-315.1	188	350814	269
300287	23	303103	81	316-315.2	188	350815	268
300290	24	304000	81	316-316.1	188	350912	269
300291	24	305002	125	316-316.2	188	350915	269
300292	24	305007	125	316-317.1	188	351014	269
300294	25	305008	125	316-317.2	188	351420	269
300295	24	305302	125	316N	188	351622	269
300295C	24	305306	125	317-314.1	188	352030	269
300297	25	305307	125	317-315.1	188	359517	268
300299	25	305308	125	317-316.1	188	400400	122
300300	120	305502	125	317-317.1	188	400500	122
300302	125	305507	125	317-317.1	188	400700	122
300302	120	305802	125	317-318.1	188	400705	122
300500	120	305805	125	326-19	170	400703	122
300700	121	305806	125	326-37	171	400800	122
						_	122
300702	123	305807	125	327-15	173	400900	
300704	123	305808	125	327-16	180	400S	253
300705	121	305809	125	327-16X	180	401100	122
300800	120	306601	125	327-19	170	401200	124
300800.1	120	307006	125	327-19X	170	401201	124
300800.2	120	307101	125	327-23	172	401204	124
300804	123	307102	125	327-24	172	401301	76, 113
300807	34, 128	307105	125	327-24X	172	401302	76, 113
300808	34, 128	307106	125	327-26	175	401402	34, 128
300900	121	307107	125	327-30	179	401403	34, 128
300903	123	307108	125	327-30R	179	401410	144
300904	123	307109	125	327-30RX	179	401700	122
300907	123	307202	125	327-33	174	402001	244
300908	123	307205	125	327-33R	174	402002	244
301-01	170	307206	125	327-34	174	402005	244
301-01F	170	307207	125	327-34R	174	402010	244
301-01FR	170	308102	125	327-34RX	174	402020	244
					170	400010	050
301-01FRX	170	308106	125	327-36	172	402NG	252
301-01FRX 301-01R	170 170	308106 308107	125 125	327-36 327-36R	172	402NM	252



Code	page	Code	page	Code	page	Code	page
407001	162	409100	256	4092.5N	138	409920	38, 130
407001/T	162	409102.1	145	4092.6N	139	409922	38, 130
407001OBT	162	409103.1	144	4092.6NS	139	410046	165
407003	162	409104.1	145	4092.7N	139	410056	165
407005	117, 163	409105.1	144	409202	149	414141	93
407005.O	117, 163	409106.1	144	409202.B	149	414143	93
407005/T	117, 163	409106.2	145	409222	149	414147	93
407006	117, 163	409107	144	409222.O	149	415141	93
407006/T	117, 163	409107.1	145	409250	256	415143	93
407008	162	409108	144	4094.1A	53	415147	93
407008.O	162	409108.1	145	4094.1AM	53	418003	65
407008OBT	162	409109	144	4094.1N	53	419802	151
408001	161	409109.1	145	4094.1R	53	419826	151
408001.1	161	409110	145	4094.1V	53	429900	37, 131
408001.0	161	409110.1	50, 140	4094.2A	56	429900SP	37, 131
408001-86	161	409110.2	50, 140	4094.2AM	56	429901	37, 131
4080010/T	161	409110.2T	51, 141	4094.2N	56	429910	36, 131
408001OBT	161	409110.3	50, 140	4094.2R	56	429920	36, 131
408002	65	409110.4	50, 140	4094.2V	56	429926	37, 131
408003 408005	65 65	409110.4T 409110.5T	51, 141 51, 141	4094.3N 4094.4N	54 54	429926.10 429926.25	37, 131 37, 131
408025	254	409110.51	145	4094.5BP	55	429926.25	37, 131
408050	254	409110/1	145	4094.5N	53	429927.10	37, 131
408030	254	409110/4	51, 141	409402	150	429927.25	37, 131
408250	254	409111/3	51, 141	409426	150	429930	35, 129
408500	254	409111/4	51, 141	4095.1BP	55	429931	35, 129
408600	159	409111/5	51, 141	4095.1N	54	429940	35, 129
408610	159	409111/6	51, 141	4095.1NP	54	429942	35, 129
408702	158	409111/7	51, 141	4095.2N	54	429945	35, 129
408702.G	158	409111A	50, 140	4095.4N	54	429946	36, 131
408720	158	409111AM	50, 140	4095.5N	56	429950	35, 129
408721	158	409111B	50, 140	4095.6N	54	429951	35, 129
408726	158	409111R	50, 140	4095.7N	54	430100.1	164
408726.G	158	409112/0	51, 141	4095.9N	53	430130	164
408900	158	409112/1	51, 141	409500	256	430130.1	164
408910	158	409112/2	51, 141	409501	153	439905	117, 130
409001	144	409112/4	51, 141	409502	153	439910	117, 130
409002	144	409112/6	51, 141	409502G	153	440131	88
409002.1	144	409113.1	50, 140	409511	156	440132	88
409003	144	409113.2	50, 140	409526	153	440133	88
409003.1	144	409113.3	50, 140	409526.1	153	440139	88
409003.2	144	409113.4	50, 140	409526.O	153	440139P	88
409003A	144	409113/1	31, 142	409526G	153	440141	89
409003AM	144	409113/2	31, 142	409552	153	440142	115
409003R	144	409113/3	31, 142	4096.2N	54	440142.0	115
409003V 409004	144 65	409113/4	31, 142	4096.3N 409602	54 149	440143 440144	89 89
409004	64	409113/6	31, 142 31, 142	409700	154	440147	89
409003 409007.M	50, 140	409115/2	50, 140	409701	155	440149	89
409007.N	50, 140	409115/3	50, 140	409702	155	440150	115
409008	64	409115/4	50, 140	409703	155	440151	115
409008.1	64	409115/6	50, 140	409711	156	440157	115
409008.M	50, 140	409120	145	409726	155	440158	115
409008.N	50, 140	409132	145	409726.G	155	440159	115
409009	64	409133	145	409752	154	441141	91
409010	64	409135	145	409756	154	441142	250
409011	64	4092.1N	137	409756.B	154	441142.1	250
409012	64	4092.1NS	137	409802	150	441142.2	250
409025	256	4092.2N	138	409826	150	441143	91
409050	256	4092.3N	139	409852	150	441145	91
4091.1N	138	4092.3NS	139	409902	156	441242.2	250
						_	



Code	page	Code	page	Code	page	Code	page
44125	93	470002	284	620200	77	808100	80, 99
44140	91	470003	284	620300	77	808101	80, 99
44141	91	470005	284	620400	77	808200	80, 99
44143	91	470006	284	621101	77	813100	133
44145	91	470007	284	621102	77	816100	133
44147	91	481000	254	621610	77	816150	133
44151	92	482000	254	621611	77	816160	133
44240	89	490050	257	621613	77	818150	133
44241	89	490100	257	701	197	820150	133
44245	89	490250	257	702	197	820200	133
44260000	288	490400	257	705000	79	831-312	183
44270000	288	490900	257	705010	79	831-312B	183
44280000	288	491000	256	707094	71	900000	268
44300000	249, 272	492000	256	712	197	900005.1	83
44310000	249, 272	493000	256	713100	134	900007	136
44330000	249, 272	495000	256	716100	134	900008	83
44351	92	53020	289	716125	134	900010	30
44355	92	53125	289	716150	134	900010.1	30
444142	115	53250	289	720150	134	900010.1	30
444150	115	53500	289	721008	187	900011	30
444150	286	554111	90	721008	80	900011.1	30
444242	115	554113	90	7311	80	900012.1	30
44430000	286	554117	90	7401	80	900015	30
44440000	286	554140	90	7411	80	900015.1	30
44450000	286	554141	90	7600	80	900019	82
44460000	286	554143	90	770000	187	900020	83
444602801	164	554147	90	770005	187	900022	82
444602901	164	554148	90	770040	187	900023	82
444602903	164	555141	92	770050	187	900024	82
444603202	164	555251	92	770060	187	900025	104
444603204	164	555258	92	770070	187	900026	82
444603300	164	600300	71	770080	187	900028	104
444603402	164	600400	71	770110	187	900030	39, 194
444608	285	600602	69	770210	187	900031	39, 194
444609	285	600800	70	770220	187	900032	39, 194
444610	285	600801	70	770240	187	900032.1	196
444611	285	601006	74, 75, 76	770310	187	900033	39, 194
444612	285	601102	69	770320	187	900034	39, 194
444613	285	601103	69	770340	187	900036	39, 194
444614	285	601203	69	800000	79	900037	39, 194
44470000	286	601402	73	801000	79	900038	39, 194
445140	91	601403	73	801075	133	900041	39, 194
445141	91	601412	73	801175T	133	900043	39, 194
445143	91	601413	73	801275	133	900050	196
445145	91	601603	73	802000	79	900073	116
445147	91	601612	73	805000	99	900074	116
451241	91	601613	73	805001	99	900076	116
451243	91	601702	73	805010	99	900078	116
451245	91	601802	72	805011	99	900093	59
454110	90	601803	72	805013	80, 99	900095	59
454111	90	601810	72	805020	99	900098	59
454113	90	602002	72	805030	99	900100	83
454117	90	602003	72	805040	99	900101	83
454118	90	611603	73	805050	99	900110	60
454140	90	611604	73	805120	99	900111	60
454141	90	617000	134	805130	99	900111B	61
454143	90	617100	134	805140	99	900112	61
454147	90	617145	134	805150	99	900113B	61
454148	90	617200	134	806030	99	900117	136
470000	284	617300	134	808000	80, 99	900117	136
470001	284	617402	132	808001	80, 99	900118	60
470001	204	017402	132	000001	ou, 99	900122	00



Code	page	Code	page	Code	page	Code	page
900123	59	Α		DS12	95	M 300400	112
900130	39, 194	A100010	96	DS35	95	M 5555	112
900131	39, 194	A100011	96	DS40	94	M 5555B	112
900134N	39, 194	A100013	96	DS42	94	M 6080	112
900136N	39, 194	A100014	96	DS61	95	M 6080B	112
900144N	39, 194	A100018	96	DS62	95	M 70100	112
900146N	39, 194	A100019	96			M 70110B	112
900155	63	A-200	222	F		M 80120	112
900170	63	A-201	222	F-200	222	M 80120B	112
900195	63	A20100	98	F-210	222	M 80160	112
900196	63	A-202	222	F-211	222	M-010.1	206
900197	63	A20250	98	F-240	222	M-010.2	206
900204	135	A20500	98			M-010.5	206
900205	135	A-210	222	G		M-010.6	206
900211	135	A-211	222	G-050	209	M-010.7	206
900212	135	A-212	222	G-100	209	M-011.1	206
900220	135	A-222	222	GS-01	68	M-019AM	210
900225	135	A-230	222	GS-02	68	M-019R	210
900300	62	A-231	222			M-019V	210
900305	62	A-240	222	н		M-020AM	202
900306	63	A-241	222	H-200	227	M-020N	202
900309	63	A-242	222	H-210	227	M-020R	202
900320	62	A-250	222	H-215	227	M-020V	202
900330	62	A-251	222	H-220	227	M-021AM	211
900384	56	A-260	222	H-240	227	M-021R	211
900400	277			H-260	226	M-021V	211
900600	277	В		H-261	226	M-022	208
901016	246	B-200	222	H-290	227	M-024	208
901030	246	B-210	222	H-295	227	M-026	208
901031	246	B-220	222	H-401	228	M-028.5	209
901031.2	246	B-240	222	H-600	28, 229	M-028.6	209
901041	246			H-601	28, 229	M-033AM	201
901075	133	С		H-610	29	M-033R	201
901275	133	C-200	223	H-611	28	M-033V	201
901407	216	C-211	223	H-615	30	M-034AM	201
901410	216	C-254	223			M-034R	201
901531	241	C-256	223	L		M-034V	201
901533	241	CP-30 H	102, 229	I-240	152, 225	M-035AM	201
901538	241	CP-30	102, 229	I-250	152, 225	M-035R	201
901600	98	CP-45 H	102, 229			M-035V	201
910000	251	_		K		M-036A	201
910022	82	D		KIT199083R	177	M-036AM	200
910023	83	D100001	97	KIT200076R	178	M-036B	200
910026	83	D100002	97	_		M-036P	200
911000	251	D100003	97	L	005	M-040	212
913100	133	D100004	97	L-240	225	M-050	212
916100	133	D101818	97	L-250	225	M-052	205
916150	133	D102020	97			M-053	205
917000	135	D102222	97	M 100150	110	M-060	212
918000	135	D102240	97	M 100150	112	M-061AM	200
918150	133	D102424	97	M 100150B	112	M-061NA	200
922000	135	D102432	97	M 110110	112	M-062B	200
928000	135	D102440	97	M 120180	112	M-062R	200
979930	104	D102450	97	M 120180B	112	M-063B	200
000040	103	D102460	97	M 150220	112	M-063R	200
989910	103	D-200	135, 223	M 160220B	112	M-064AM	200
989920		D 000			117		
989920 989945	103	D-230	135, 223	M 180250		M-064V	200
989920 989945 989951	103 98	D-240	135, 223	M 180250B	112	M-100	212
989920 989945 989951 989953	103 98 98	D-240 D-250	135, 223 135, 223	M 180250B M 200300	112 112	M-100 M-110	212 212
989920 989945 989951	103 98	D-240	135, 223	M 180250B	112	M-100	212



INDEX	
Code	page
M-210	182
M-300	228
M-410	205
M-450	205
M-500	205
M-501	205
M-502	205
M-503	205
M-510	205
M-511	205
M-512	205
M-513	205
M-520	213
M-521	213
M-522	213
M-523	213
M-530	213
M-531	213
M-532	213
M-533	213
M-540	210
M-541	210
M-542	210
M-550 M-551	210
M-552	210
M-553	210
M-554	211
M-560A	216
M-561A	216
M-563A	206
M-563NA	206
M-563P	206
M-570	208
M-571	208
M-572	207
M-581	207
M-585	207
M-600	204
M-603	204
M-604	204
M-606	204
M-607	204
M-610	204
M-640	044
	211
M-641	211
M-641	
M-641	211
M-641 N N 160230	112
M-641 N N 160230 N-240	211 112 225
M-641 N N 160230 N-240 N-250	211 112 225 225
M-641 N N 160230 N-240	211 112 225
M-641 N N 160230 N-240 N-250 NT-250	211 112 225 225
M-641 N N 160230 N-240 N-250 NT-250 R	211 112 225 225 225 225
M-641 N N 160230 N-240 N-250 NT-250 R R-200	211 112 225 225 225 225 227
M-641 N N 160230 N-240 N-250 NT-250 R R-200 R-215	211 112 225 225 225 227
M-641 N 160230 N-240 N-250 NT-250 R R-200 R-215 R-220	211 112 225 225 225 227 227
M-641 N 160230 N-240 N-250 NT-250 R R-200 R-215 R-220 R-230	211 112 225 225 225 227 227 227
M-641 N 160230 N-240 N-250 NT-250 R R-200 R-215 R-220 R-230 R-240	211 112 225 225 225 227 227 227 227 227
M-641 N 160230 N-240 N-250 NT-250 R R-200 R-215 R-220 R-230	211 112 225 225 225 227 227 227


Code	page
R-282	224
R-283	224
R-292	224
R-293	224
K-293	224
S	
S-002	102, 229
S-003	102, 229
S-004	102, 229
Т	
T-200	226
T-210	226
U	
U-0537A	204
U-0537N	204
U-0537NA	204
U-105105A	219
U-105105Y	219
U-1145A	217
U-1145AM	217
U-1145N	217
U-1145NA	217
	217
U-1170A	
U-1170AM	217
U-1170N	217
U-15128A	219
U-15128AM	219
U-15128N	219
U-4445A	204
U-4445AM	204
U-4445N	204
U-4445NE	204
U-50128A	216
U-50128N	216
U-50128R	216
U-50128V	216
U-7750A	203
U-7750AM	203
U-7750N	203
U-7780A	218
U-7780N	218
U-7780NA	218
U-7796A	218
U-7796V U-8845A	218
	203
U-8845N	203
U-8850A	203
U-8850N	203
U-8850NA	203
U-9581N	202
U-9950CAM	202
U-9950CN	202
U-9950CNA	202
U-9950CR	202
U-9950CV	202
U-9970A	217
U-9970N	217
U-9970R	217
U-9980A	218

Code	page
U-9980N	
	218
U-9996A	219
U-9996N	219
V	
V 160230	112
w	
	005
W000010	235
W000011	235
W000018	235
W000020	235
W000040	234
W000048	234
W000050	234
W000051	234
W000059	234
W000120	234
W000121	234
W000129	234
W001400	236
W003500	236
W009680	236
W-011	214
W-012	214
W-013	214
W-014	214
W-015	214
W-016	214
W-017	214
W-018	214
W025100	235
W030930	237
W031030	237
W037100	235
W051100	233
W051960	233
W052960	232
W053100	232
W053960	232
W058100	233
W058960	233
W059960	233
W092006	235
W092170	235
W092171	235
W092221	235
W-100	40, 195, 247
W-110	40, 195, 247
W-120	40, 195, 247
W121100	233
W121960	233
W122960	232
W123100	232
W123960	232
W128100	233
W128960	233
W129960	233
W-150	247
W310900	237
W310990	237
44910990	201

Code	page
W350900	237
W350990	237
W370900	237
W370990	237
W700000	234
W700001	234
W700002	234
W700003	234
W700004	234
W990000	233
W990001	233
W990002	233
W990003	233
W990004	233
W990082	233
W992960	233
W999600	233
W999601	233
W999602	233
W999603	233
W999604	233
Z Z-025	221
Z-050	221
Z-075	221
Z-100	221
Z-250	221
Z-300	221
Z-500	221
Z-550	221
Z-600	221
Z-700	221
Z-750	221
Z-800	221



IIdeltalab

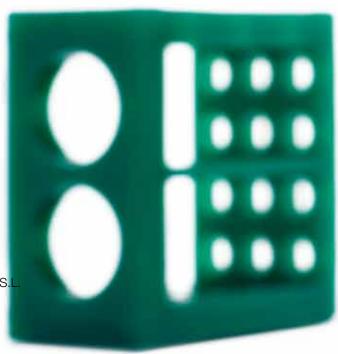

Direction and coordination

Marketing Department, DELTALAB S.I

Photography
DELTALAB S.L. files
Christian Ribas

Artwork and digital image Innovagrafic - Pere Mola

Graphic concept and cover design
Marketing Department, DELTALAB S.L.

Printed by Vanguard Grafic, S.A.

All rights reserved. Both partial or total reproduction of this catalogue without the written permission from DELTALAB S.L. is forbidden

Printed in Spain

November, 2013







All information presented in this catalogue is correct except in

the case of typing errors.

DELTALAB S.L. reserves the right to alter the characteristics of products presented in our catalogue when said changes are beyond control of this company.



Plaza de la Verneda, 1 Pol. Ind. La Llana 08191 Rubí, Barcelona (SPAIN)

t. +34 936 995 000 f +34 936 994 512

www.deltalab.es marketing@deltalab.es